

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH


PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II

County Alleghany Instrument Location Alleghany County Jail
Instrument Serial No. 808890 Sparta, N. C.

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 6 day of September, 2018, the foregoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

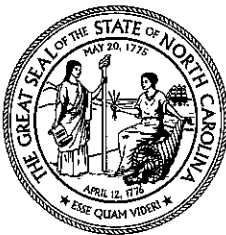

[Signature]
Signature of Certifying Official

657
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

ALLEGHANY COUNTY ALLEGHANY CO JAIL 020

Serial Number: 008890
Test Date: 09/06/2018

Citation Number: M0000000-0
Subject's Name:
PREVENTIVE, MAINTENANCE
Subject's Date of Birth: 11/11/1911
Subject's Sex: Male
Driver's License State: XX
Driver's License Number: NONE

Analyst's Name:
BENFIELD II, KENNETH R
Permit Number: 22067E
Effective:
06/01/2018-06/01/2020

Officer's Name: NONE, NONE
Type of Agency: FTA
Agency: DHHS
Test Type: Breath Test

Lot Number: AG702302
Exp Date: 01/23/2019


Test	g/210L	Time
DIAG	Pass	11:52am
AIR BLK	.00	11:53am
ACCY CHK	.07	11:54am
AIR BLK	.00	11:54am
SUB TEST	.00	11:55am
AIR BLK	.00	11:56am
SUB TEST	.00	11:57am
AIR BLK	.00	11:58am

Reported AC: .00 g/210L



Signature of Chemical Analyst

Court CVR

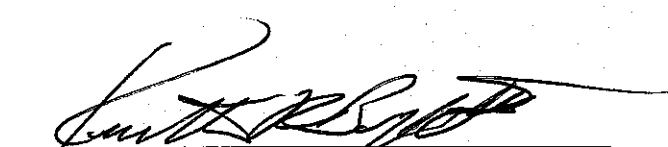

Analyst

Intox EC/IR-II: Preventive Maintenance

ALLEGHANY COUNTY ALLEGHANY CO JAIL 020

Serial Number: 008890 Test Record Number: 721
Test Date: 09/06/2018 Test Time: 11:59am EDT

System Check: Passed

Baseline Tests

Test	Status	Time
IR	Pass	11:59am
FLO	Pass	11:59am
FC	Pass	11:59am

Temperature Tests

Test	Status	Time
FC1	Pass	12:00pm
SRC	Pass	12:00pm
DET	Pass	12:00pm
BAR	Pass	12:00pm
BT	Pass	12:00pm

Blank Tests

Test	Status	Time
AIR	Pass	12:00pm


Printer Tests

Test	Status	Time
PRNT	Pass	12:00pm

CRC Tests

Test	Status	Time
COMP	Pass	12:00pm
CAL	Pass	12:00pm

Preventive Maintenance
Status: Pass

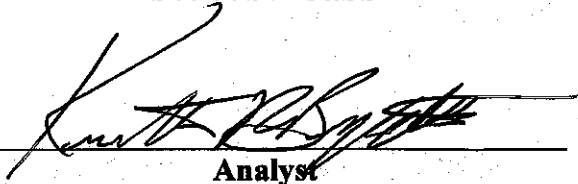
Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH


**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**

County Anson Co. Instrument Location Anson Co. Sheriff's Office
Instrument Serial No. 08139 Wadesboro, NC

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 26th day of September, 20 18 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

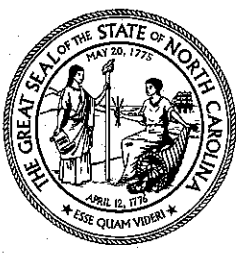

[Signature]
Signature of Certifying Official

654
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

ANSON COUNTY ANSON COUNTY SD 030

Serial Number: 008739
Test Date: 09/26/2018

Citation Number: M0000000-0

Subject's Name:

PREVENTIVE, MAINTENANCE

Subject's Date of Birth: 11/11/1911

Subject's Sex: Male

Driver's License State: XX

Driver's License Number: NONE

Analyst's Name: KEESLER, GRAYHAM C

Permit Number: 07682E

Effective:

12/01/2017-12/01/2019

Officer's Name: NONE, NONE

Type of Agency: FTA

Agency: DHHS

Test Type: Breath Test

Lot Number: AG734101

Exp Date: 12/07/2019

Test	g/210L	Time
DIAG	Pass	12:36pm
AIR BLK	.00	12:37pm
ACCY CHK	.08	12:38pm
AIR BLK	.00	12:39pm
SUB TEST	.00	12:40pm
AIR BLK	.00	12:41pm
SUB TEST	.00	12:42pm
AIR BLK	.00	12:43pm

Reported AC: .00 g/210L


Signature of Chemical Analyst

Court CVR


Analyst

**This form is used when performing Preventive Maintenance procedures
Forensic Tests for Alcohol Branch
Department of Health and Human Services
Rev. 12/2007**

Intox EC/IR-II: Preventive Maintenance

ANSON COUNTY ANSON COUNTY SD 030

Serial Number: 008739 Test Record Number: 321
Test Date: 09/26/2018 Test Time: 12:44pm EDT

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	12:44pm
FLO	Pass	12:44pm
FC	Pass	12:44pm

Temperature Tests

Test	Status	Time
FC1	Pass	12:44pm
SRC	Pass	12:44pm
DET	Pass	12:44pm
BAR	Pass	12:44pm
BT	Pass	12:44pm

Blank Tests

Test	Status	Time
AIR	Pass	12:45pm


Printer Tests

Test	Status	Time
PRNT	Pass	12:45pm

CRC Tests

Test	Status	Time
COMP	Pass	12:45pm
CAL	Pass	12:45pm

Preventive Maintenance
Status: Pass

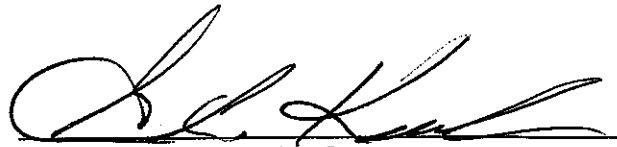

Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH

**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**

County ANSON Co. Instrument Location ANSON Co. Sheriff Office

Instrument Serial No. 008597 Wadesboro

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 26th day of September, 20 18 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.



[Signature]
Signature of Certifying Official

654
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

ANSON COUNTY ANSON COUNTY SD 030

Serial Number: 008597
Test Date: 09/26/2018

Citation Number: M0000000-0
Subject's Name:

PREVENTIVE, MAINTENANCE

Subject's Date of Birth: 11/11/1911

Subject's Sex: Male

Driver's License State: XX

Driver's License Number: NONE

Analyst's Name: KEESLER, GRAYHAM C

Permit Number: 07682E

Effective:

12/01/2017-12/01/2019

Officer's Name: NONE, NONE

Type of Agency: FTA

Agency: DHHS

Test Type: Breath Test

Lot Number: AG716201

Exp Date: 06/11/2019

Test	g/210L	Time
DIAG	Pass	12:39pm
AIR BLK	.00	12:40pm
ACCY CHK	.08	12:41pm
AIR BLK	.00	12:42pm
SUB TEST	.00	12:43pm
AIR BLK	.00	12:44pm
SUB TEST	.00	12:45pm
AIR BLK	.00	12:46pm

Reported AC: .00 g/210L


Signature of Chemical Analyst

Court CVR


Analyst

**This form is used when performing Preventive Maintenance procedures
Forensic Tests for Alcohol Branch
Department of Health and Human Services
Rev. 12/2007**

Intox EC/IR-II: Preventive Maintenance

ANSON COUNTY ANSON COUNTY SD 030

Serial Number: 008597 Test Record Number: 1709
Test Date: 09/26/2018 Test Time: 12:47pm EDT

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	12:47pm
FLO	Pass	12:47pm
FC	Pass	12:47pm

Temperature Tests

Test	Status	Time
FC1	Pass	12:47pm
SRC	Pass	12:47pm
DET	Pass	12:47pm
BAR	Pass	12:47pm
BT	Pass	12:47pm

Blank Tests

Test	Status	Time
AIR	Pass	12:48pm


Printer Tests

Test	Status	Time
PRNT	Pass	12:48pm

CRC Tests

Test	Status	Time
COMP	Pass	12:48pm
CAL	Pass	12:48pm

Preventive Maintenance
Status: Pass

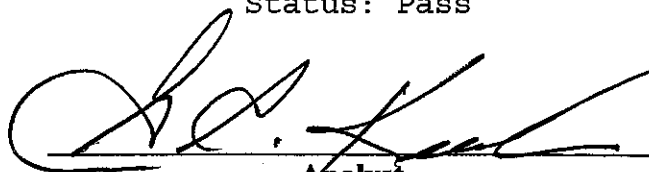

Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH


**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**

County Ashe Instrument Location Ashe County Jail
Instrument Serial No. 008849 Jefferson, N.C.

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 6 day of September, 2018, the foregoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

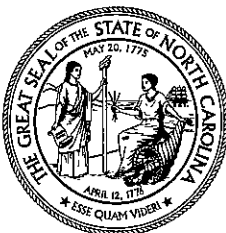

[Signature] 657
Signature of Certifying Official Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

ASHE COUNTY ASHE COUNTY JAIL 040

Serial Number: 008849

Test Date: 09/06/2018

Citation Number: M0000000-0

Subject's Name:

PREVENTIVE, MAINTENANCE

Subject's Date of Birth: 11/11/1911

Subject's Sex: Male

Driver's License State: XX

Driver's License Number: NONE

Analyst's Name:

BENFIELD II, KENNETH R

Permit Number: 22067E

Effective:

06/01/2018-06/01/2020

Officer's Name: NONE, NONE

Type of Agency: FTA

Agency: DHHS


Test Type: Breath Test

Lot Number: AG805801

Exp Date: 02/27/2020

Test	g/210L	Time
DIAG	Pass	10:42am
AIR BLK	.00	10:42am
ACCY CHK	.07	10:43am
AIR BLK	.00	10:44am
SUB TEST	.00	10:44am
AIR BLK	.00	10:45am
SUB TEST	.00	10:47am
AIR BLK	.00	10:48am

Reported AC: .00 g/210L

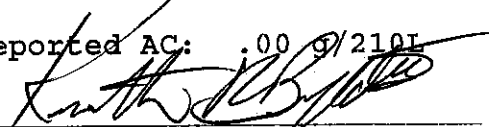
Signature of Chemical Analyst

Court CVR


Analyst

Intox EC/IR-II: Preventive Maintenance

ASHE COUNTY ASHE COUNTY JAIL 040

Serial Number: 008849 Test Record Number: 1176
Test Date: 09/06/2018 Test Time: 10:49am EDT

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	10:49am
FLO	Pass	10:49am
FC	Pass	10:49am

Temperature Tests

Test	Status	Time
FC1	Pass	10:49am
SRC	Pass	10:49am
DET	Pass	10:49am
BAR	Pass	10:49am
BT	Pass	10:49am

Blank Tests

Test	Status	Time
AIR	Pass	10:50am


Printer Tests

Test	Status	Time
PRNT	Pass	10:50am

CRC Tests

Test	Status	Time
COMP	Pass	10:50am
CAL	Pass	10:50am

Preventive Maintenance
Status: Pass

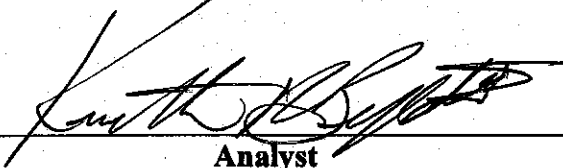

Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH

**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**


County Bertie Instrument Location Bertie Co. S.O.

Instrument Serial No. 008897 222 County Farm Rd., Windsor, NC.

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 11th day of September, 2018 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

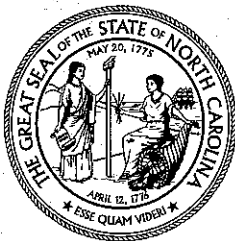

Key [Signature]
Signature of Certifying Official

643
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

BERTIE COUNTY BERTIE CO SO 070

Serial Number: 008897
Test Date: 09/11/2018

Citation Number: M0000000-0
Subject's Name:
PREVENTIVE, MAINTENANCE
Subject's Date of Birth: 11/11/1911
Subject's Sex: Male
Driver's License State: XX
Driver's License Number: NONE


Analyst's Name: GUARD, KELLY G
Permit Number: 12955E
Effective:
06/01/2017-06/01/2019

Officer's Name: NONE, NONE
Type of Agency: FTA
Agency: DHHS
Test Type: Breath Test

Lot Number: AG721401
Exp Date: 08/02/2019


Test	g/210L	Time
DIAG	Pass	9:27am
AIR BLK	.00	9:27am
ACCY CHK	.08	9:28am
AIR BLK	.00	9:29am
SUB TEST	.00	9:30am
AIR BLK	.00	9:31am
SUB TEST	.00	9:32am
AIR BLK	.00	9:33am

Reported AC: .00 g/210L

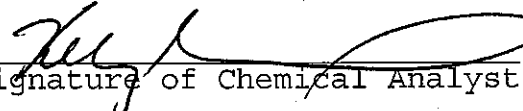

Signature of Chemical Analyst

Court CVR



Analyst

Intox EC/IR-II: Preventive Maintenance

BERTIE COUNTY BERTIE CO SO 070

Serial Number: 008897 Test Record Number: 1151
Test Date: 09/11/2018 Test Time: 9:35am EDT

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	9:35am
FLO	Pass	9:35am
FC	Pass	9:35am

Temperature Tests

Test	Status	Time
FC1	Pass	9:36am
SRC	Pass	9:36am
DET	Pass	9:36am
BAR	Pass	9:36am
BT	Pass	9:36am

Blank Tests

Test	Status	Time
AIR	Pass	9:36am


Printer Tests

Test	Status	Time
PRNT	Pass	9:36am

CRC Tests

Test	Status	Time
COMP	Pass	9:36am
CAL	Pass	9:36am

Preventive Maintenance
Status: Pass

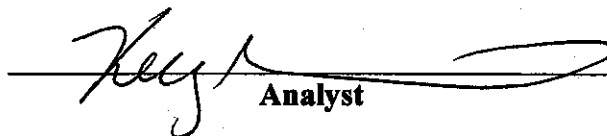
Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH


**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**

County Cabarrus Instrument Location Cabarrus County SD
Instrument Serial No. 008694 30 Corban Ave, Concord

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 11th day of September, 2018 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

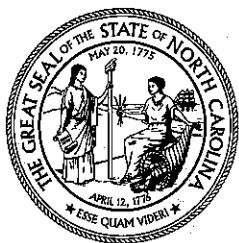

[Handwritten Signature]
Signature of Certifying Official

656
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

CABARRUS CABARRUS COUNTY SD 120

Serial Number: 008694
Test Date: 09/11/2018

Citation Number: M0000000-0
Subject's Name:
PREVENTIVE, MAINTENANCE
Subject's Date of Birth: 11/11/1911
Subject's Sex: Male
Driver's License State: XX
Driver's License Number: NONE

Analyst's Name: HAYS, MARK D
Permit Number: 15924E
Effective:
01/01/2018-01/01/2020

Officer's Name: NONE, NONE
Type of Agency: FTA
Agency: DHHS
Test Type: Breath Test

Lot Number: AG722408
Exp Date: 08/12/2019

Test	g/210L	Time
DIAG	Pass	10:19am
AIR BLK	.00	10:20am
ACCY CHK	.07	10:21am
AIR BLK	.00	10:22am
SUB TEST	.00	10:22am
AIR BLK	.00	10:23am
SUB TEST	.00	10:25am
AIR BLK	.00	10:26am

Reported AC: .00 g/210L



Signature of Chemical Analyst

Court CVR



Analyst

**This form is used when performing Preventive Maintenance procedures
Forensic Tests for Alcohol Branch
Department of Health and Human Services
Rev. 12/2007**

Intox EC/IR-II: Preventive Maintenance

CABARRUS CABARRUS COUNTY SD 120

Serial Number: 008694 Test Record Number: 1401
Test Date: 09/11/2018 Test Time: 10:28am EDT

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	10:28am
FLO	Pass	10:28am
FC	Pass	10:28am

Temperature Tests

Test	Status	Time
FC1	Pass	10:28am
SRC	Pass	10:28am
DET	Pass	10:28am
BAR	Pass	10:28am
BT	Pass	10:28am

Blank Tests

Test	Status	Time
AIR	Pass	10:29am


Printer Tests

Test	Status	Time
PRNT	Pass	10:29am

CRC Tests

Test	Status	Time
COMP	Pass	10:29am
CAL	Pass	10:29am

Preventive Maintenance
Status: Pass



Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH


**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**

County Cabarrus Instrument Location Cabarrus County SD
Instrument Serial No. 008590 30 Corban Ave., Concord

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 11th day of September, 20 18 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

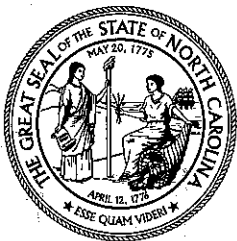

[Handwritten Signature]
Signature of Certifying Official

656

Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

CABARRUS COUNTY CABARRUS COUNTY SD 120

Serial Number: 008590
Test Date: 09/11/2018

Citation Number: M0000000-0
Subject's Name:

PREVENTIVE, MAINTENANCE

Subject's Date of Birth: 11/11/1911

Subject's Sex: Male

Driver's License State: XX

Driver's License Number: NONE

Analyst's Name: HAYS, MARK D

Permit Number: 15924E

Effective:

01/01/2018-01/01/2020

Officer's Name: NONE, NONE

Type of Agency: FTA

Agency: DHHS


Test Type: Breath Test

Lot Number: AG734101

Exp Date: 12/07/2019

Test	g/210L	Time
DIAG	Pass	10:01am
AIR BLK	.00	10:02am
ACCY CHK	.08	10:02am
AIR BLK	.00	10:03am
SUB TEST	.00	10:04am
AIR BLK	.00	10:05am
SUB TEST	.00	10:06am
AIR BLK	.00	10:07am

Reported AC: .00 g/210L

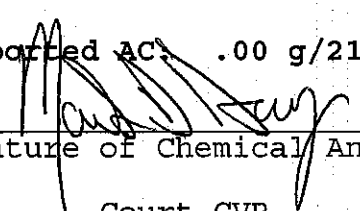

Signature of Chemical Analyst

Court CVR



Analyst

**This form is used when performing Preventive Maintenance procedures
Forensic Tests for Alcohol Branch
Department of Health and Human Services
Rev. 12/2007**

Intox EC/IR-II: Preventive Maintenance

CABARRUS COUNTY CABARRUS COUNTY SD 120

Serial Number: 008590 Test Record Number: 2949
Test Date: 09/11/2018 Test Time: 10:09am EDT

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	10:09am
FLO	Pass	10:09am
FC	Pass	10:09am

Temperature Tests

Test	Status	Time
FC1	Pass	10:10am
SRC	Pass	10:10am
DET	Pass	10:10am
BAR	Pass	10:10am
BT	Pass	10:10am

Blank Tests

Test	Status	Time
AIR	Pass	10:10am


Printer Tests

Test	Status	Time
PRNT	Pass	10:10am

CRC Tests

Test	Status	Time
COMP	Pass	10:10am
CAL	Pass	10:10am

Preventive Maintenance
Status: Pass

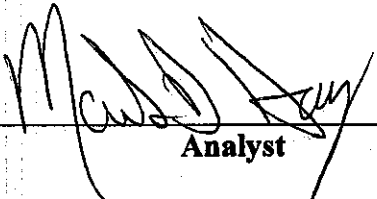

Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH


**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**

County Camden Instrument Location Camden Co. S.O.
Instrument Serial No. 008847 113 Hwy 343, Camden, NC

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 26th day of September, 2018 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

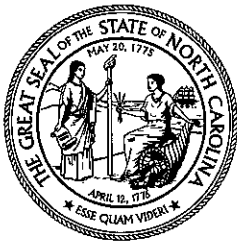

[Signature]
Signature of Certifying Official

643
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

CAMDEN COUNTY CAMDEN CO SO 140

Serial Number: 008847

Test Date: 09/26/2018

Citation Number: M0000000-0

Subject's Name:

PREVENTIVE, MAINTENANCE

Subject's Date of Birth: 11/11/1911

Subject's Sex: Male

Driver's License State: XX

Driver's License Number: NONE

Analyst's Name: GUARD, KELLY G

Permit Number: 12955E

Effective:

06/01/2017-06/01/2019

Officer's Name: NONE, NONE

Type of Agency: FTA

Agency: DHHS


Test Type: Breath Test

Lot Number: AG721401

Exp Date: 08/02/2019


Test	g/210L	Time
DIAG	Pass	11:31am
AIR BLK	.00	11:31am
ACCY CHK	.08	11:32am
AIR BLK	.00	11:33am
SUB TEST	.00	11:33am
AIR BLK	.00	11:34am
SUB TEST	.00	11:36am
AIR BLK	.00	11:37am

Reported AC: .00 g/210L

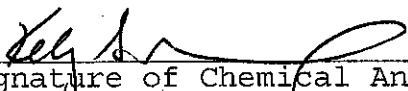

Signature of Chemical Analyst

Court CVR

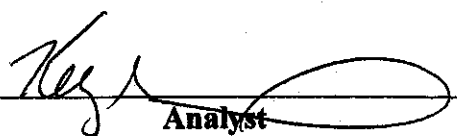

Analyst

Intox EC/IR-II: Preventive Maintenance

CAMDEN COUNTY CAMDEN CO SO 140

Serial Number: 008847 Test Record Number: 655
Test Date: 09/26/2018 Test Time: 11:38am EDT

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	11:38am
FLO	Pass	11:38am
FC	Pass	11:38am

Temperature Tests

Test	Status	Time
FC1	Pass	11:38am
SRC	Pass	11:38am
DET	Pass	11:38am
BAR	Pass	11:38am
BT	Pass	11:38am

Blank Tests

Test	Status	Time
AIR	Pass	11:39am

Printer Tests

Test	Status	Time
PRNT	Pass	11:39am

CRC Tests

Test	Status	Time
COMP	Pass	11:39am
CAL	Pass	11:39am

Preventive Maintenance
Status: Pass



Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH

**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**

County Carteret


Instrument Location Atlantic Beach PD

Instrument Serial No. 008917

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 4 day of September, 2018, the foregoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

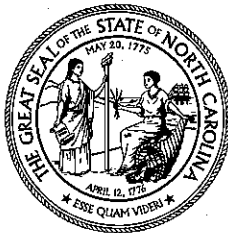

Randy Hall
Signature of Certifying Official

354
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

CARTERET COUNTY ATLANTIC BEACH PD 150

Serial Number: 008917
Test Date: 09/04/2018

Citation Number: M0000000-0
Subject's Name:
PREVENTIVE, MAINTENANCE
Subject's Date of Birth: 11/11/1911
Subject's Sex: Male
Driver's License State: XX
Driver's License Number: NONE


Analyst's Name: HALL, RANDY E
Permit Number: 3462E
Effective:
05/01/2017-05/01/2019

Officer's Name: NONE, NONE
Type of Agency: FTA
Agency: DHHS
Test Type: Breath Test

Lot Number: AG702302
Exp Date: 01/23/2019

Test	g/210L	Time
DIAG	Pass	11:04am
AIR BLK	.00	11:05am
ACCY CHK	.08	11:06am
AIR BLK	.00	11:07am
SUB TEST	.00	11:07am
AIR BLK	.00	11:08am
SUB TEST	.00	11:09am
AIR BLK	.00	11:10am

Reported AC: .00 g/210L

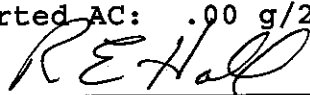

Signature of Chemical Analyst

Court CVR



Analyst

Intox EC/IR-II: Preventive Maintenance

CARTERET COUNTY ATLANTIC BEACH PD 150

Serial Number: 008917 Test Record Number: 755
Test Date: 09/04/2018 Test Time: 11:22am EDT

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	11:23am
FLO	Pass	11:23am
FC	Pass	11:23am

Temperature Tests

Test	Status	Time
FC1	Pass	11:23am
SRC	Pass	11:23am
DET	Pass	11:23am
BAR	Pass	11:23am
BT	Pass	11:23am

Blank Tests

Test	Status	Time
AIR	Pass	11:24am

Printer Tests

Test	Status	Time
PRNT	Pass	11:24am

CRC Tests

Test	Status	Time
COMP	Pass	11:24am
CAL	Pass	11:24am

Preventive Maintenance
Status: Pass



Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH

PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II


County Cumberland Instrument Location BAT Mobile Unit 4

Instrument Serial No. 008929 Fayetteville PO

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 7 day of September, 2018, the foregoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

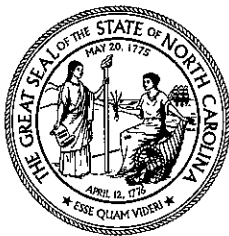

[Signature]
Signature of Certifying Official

660
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

CUMBERLAND COUNTY BAT MOBILE UNIT 4
250

Serial Number: 008929
Test Date: 09/07/2018

Citation Number: M0000000-0
Subject's Name:

PREVENTIVE, MAINTENANCE

Subject's Date of Birth: 11/11/1911

Subject's Sex: Male

Driver's License State: XX

Driver's License Number: NONE

Analyst's Name: SMITH, JASON R

Permit Number: 19145E

Effective:

03/01/2018-03/01/2020

Officer's Name: NONE, NONE

Type of Agency: FTA

Agency: DHHS

Test Type: Breath Test

Lot Number: AG814901

Exp Date: 05/29/2020

Test	g/210L	Time
DIAG	Pass	10:25pm
AIR BLK	.00	10:26pm
ACCY CHK	.08	10:27pm
AIR BLK	.00	10:28pm
SUB TEST	.00	10:28pm
AIR BLK	.00	10:29pm
SUB TEST	.00	10:31pm
AIR BLK	.00	10:32pm

Reported AC: .00 g/210L


Signature of Chemical Analyst

Court CVR


Analyst

**This form is used when performing Preventive Maintenance procedures
Forensic Tests for Alcohol Branch
Department of Health and Human Services
Rev. 12/2007**

Intox EC/IR-II: Preventive Maintenance

CUMBERLAND COUNTY BAT MOBILE UNIT 4 250

Serial Number: 008929 Test Record Number: 1017
Test Date: 09/07/2018 Test Time: 10:36pm EDT

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	10:37pm
FLO	Pass	10:37pm
FC	Pass	10:37pm

Temperature Tests

Test	Status	Time
FC1	Pass	10:37pm
SRC	Pass	10:37pm
DET	Pass	10:37pm
BAR	Pass	10:37pm
BT	Pass	10:37pm

Blank Tests

Test	Status	Time
AIR	Pass	10:38pm


Printer Tests

Test	Status	Time
PRNT	Pass	10:38pm

CRC Tests

Test	Status	Time
COMP	Pass	10:38pm
CAL	Pass	10:38pm

Preventive Maintenance
Status: Pass



Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH


**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**

County Camden Instrument Location BAT Mobile Unit 4
Instrument Serial No. 008601 Payetteville PD

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 7 day of September, 2018, the foregoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

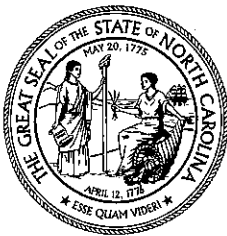

[Signature]
Signature of Certifying Official

660
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

CUMBERLAND COUNTY BAT MOBILE UNIT 4

250

Serial Number: 008601

Test Date: 09/07/2018

Citation Number: M0000000-0

Subject's Name:

PREVENTIVE, MAINTENANCE

Subject's Date of Birth: 11/11/1911

Subject's Sex: Male

Driver's License State: XX

Driver's License Number: NONE

Analyst's Name: SMITH, JASON R

Permit Number: 19145E

Effective:

03/01/2018-03/01/2020

Officer's Name: NONE, NONE

Type of Agency: FTA

Agency: DHHS

Test Type: Breath Test

Lot Number: AG702401

Exp Date: 01/24/2019

Test	g/210L	Time
DIAG	Pass	10:26pm
AIR BLK	.00	10:27pm
ACCY CHK	.08	10:28pm
AIR BLK	.00	10:29pm
SUB TEST	.00	10:30pm
AIR BLK	.00	10:30pm
SUB TEST	.00	10:32pm
AIR BLK	.00	10:33pm

Reported AC: .00 g/210L


Signature of Chemical Analyst

Court CVR


Analyst

**This form is used when performing Preventive Maintenance procedures
Forensic Tests for Alcohol Branch
Department of Health and Human Services
Rev. 12/2007**

Intox EC/IR-II: Preventive Maintenance

CUMBERLAND COUNTY BAT MOBILE UNIT 4 250

Serial Number: 008601 Test Record Number: 1303
Test Date: 09/07/2018 Test Time: 10:34pm EDT

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	10:35pm
FLO	Pass	10:35pm
FC	Pass	10:35pm

Temperature Tests

Test	Status	Time
FC1	Pass	10:35pm
SRC	Pass	10:35pm
DET	Pass	10:35pm
BAR	Pass	10:35pm
BT	Pass	10:35pm

Blank Tests

Test	Status	Time
AIR	Pass	10:35pm


Printer Tests

Test	Status	Time
PRNT	Pass	10:35pm

CRC Tests

Test	Status	Time
COMP	Pass	10:36pm
CAL	Pass	10:36pm

Preventive Maintenance
Status: Pass

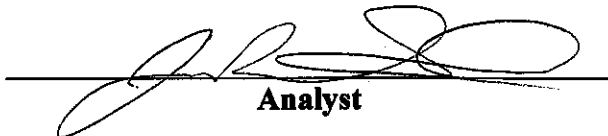
Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH

**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**


County Cumberland Co. Instrument Location FORT BRASS

Instrument Serial No. 008908 PMO

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 10th day of September, 2018 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

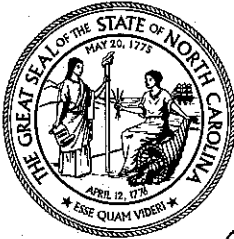

[Signature]
Signature of Certifying Official

654
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

CUMBERLAND COUNTY FT BRAGG LEC 250

Serial Number: 008908
Test Date: 09/10/2018

Citation Number: M0000000-0
Subject's Name:
PREVENTIVE, MAINTENANCE
Subject's Date of Birth: 11/11/1911
Subject's Sex: Male
Driver's License State: XX
Driver's License Number: NONE

Analyst's Name: KEESLER, GRAYHAM C
Permit Number: 07682E
Effective:
12/01/2017-12/01/2019

Officer's Name: NONE, NONE
Type of Agency: FTA
Agency: DHHS
Test Type: Breath Test

Lot Number: AG814902
Exp Date: 05/29/2020

Test	g/210L	Time
DIAG	Pass	2:07pm
AIR BLK	.00	2:07pm
ACCY CHK	.08	2:08pm
AIR BLK	.00	2:09pm
SUB TEST	.00	2:10pm
AIR BLK	.00	2:11pm
SUB TEST	.00	2:12pm
AIR BLK	.00	2:13pm

Reported AC: .00 g/210L


Signature of Chemical Analyst

Court CVR


Analyst

**This form is used when performing Preventive Maintenance procedures
Forensic Tests for Alcohol Branch
Department of Health and Human Services
Rev. 12/2007**

Intox EC/IR-II: Preventive Maintenance

CUMBERLAND COUNTY FT BRAGG LEC 250

Serial Number: 008908 Test Record Number: 1659
Test Date: 09/10/2018 Test Time: 2:14pm EDT

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	2:15pm
FLO	Pass	2:15pm
FC	Pass	2:15pm

Temperature Tests

Test	Status	Time
FC1	Pass	2:15pm
SRC	Pass	2:15pm
DET	Pass	2:15pm
BAR	Pass	2:15pm
BT	Pass	2:15pm

Blank Tests

Test	Status	Time
AIR	Pass	2:16pm


Printer Tests

Test	Status	Time
PRNT	Pass	2:16pm

CRC Tests

Test	Status	Time
COMP	Pass	2:16pm
CAL	Pass	2:16pm

Preventive Maintenance
Status: Pass

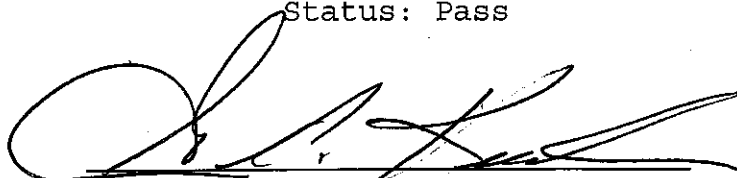

Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH

**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**


County Cumberland Co. Instrument Location FORT BRAGG

Instrument Serial No. 008903 P110

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 10th day of September, 2018 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

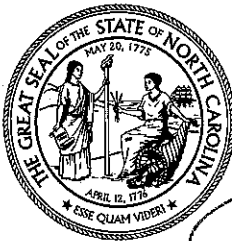

[Signature]
Signature of Certifying Official

654
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

CUMBERLAND COUNTY FORT BRAGG LEC 250

Serial Number: 008903
Test Date: 09/10/2018

Citation Number: M0000000-0
Subject's Name:

PREVENTIVE, MAINTENANCE

Subject's Date of Birth: 11/11/1911

Subject's Sex: Male

Driver's License State: XX

Driver's License Number: NONE

Analyst's Name: KEESLER, GRAYHAM C

Permit Number: 07682E

Effective:

12/01/2017-12/01/2019

Officer's Name: NONE, NONE

Type of Agency: FTA

Agency: DHHS

Test Type: Breath Test

Lot Number: AG814902

Exp Date: 05/29/2020

Test	g/210L	Time
DIAG	Pass	1:59pm
AIR BLK	.00	2:00pm
ACCY CHK	.07	2:01pm
AIR BLK	.00	2:02pm
SUB TEST	.00	2:02pm
AIR BLK	.00	2:03pm
SUB TEST	.00	2:05pm
AIR BLK	.00	2:06pm

Reported AC: .00 g/210L


Signature of Chemical Analyst

Court CVR


Analyst

**This form is used when performing Preventive Maintenance procedures
Forensic Tests for Alcohol Branch
Department of Health and Human Services
Rev. 12/2007**

Intox EC/IR-II: Preventive Maintenance

CUMBERLAND COUNTY FORT BRAGG LEC 250

Serial Number: 008903 Test Record Number: 2235
Test Date: 09/10/2018 Test Time: 2:07pm EDT

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	2:07pm
FLO	Pass	2:07pm
FC	Pass	2:07pm

Temperature Tests

Test	Status	Time
FC1	Pass	2:08pm
SRC	Pass	2:08pm
DET	Pass	2:08pm
BAR	Pass	2:08pm
BT	Pass	2:08pm

Blank Tests

Test	Status	Time
AIR	Pass	2:08pm


Printer Tests

Test	Status	Time
PRNT	Pass	2:08pm

CRC Tests

Test	Status	Time
COMP	Pass	2:08pm
CAL	Pass	2:08pm

Preventive Maintenance
Status: Pass

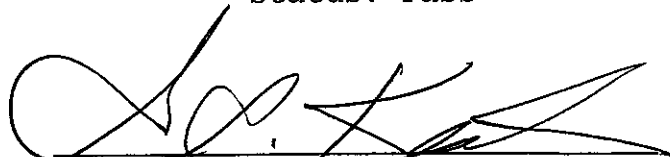

Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH

PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II

County CDWELCOMDE Instrument Location ATM MOUNT UNIT 6

Instrument Serial No. 009570 ROCKY MOUNT

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 29 day of SEPTEMBER, 20 11, the foregoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.



[Signature]

Signature of Certifying Official

663

Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

EDGECOMBE COUNTY BAT MOBILE UNIT 6 320

Serial Number: 008580
Test Date: 09/29/2018

Citation Number: M0000000-0
Subject's Name:
PREVENTIVE, MAINTENANCE
Subject's Date of Birth: 11/11/1911
Subject's Sex: Male
Driver's License State: XX
Driver's License Number: NONE


Analyst's Name: VARNELL, BRYON L
Permit Number: 16896E
Effective:
09/22/2017-09/01/2019

Officer's Name: NONE, NONE
Type of Agency: FTA
Agency: DHHS
Test Type: Breath Test

Lot Number: AG710701
Exp Date: 04/17/2019


Test	g/210L	Time
DIAG	Pass	9:14pm
AIR BLK	.00	9:15pm
ACCY CHK	.07	9:15pm
AIR BLK	.00	9:16pm
SUB TEST	.00	9:17pm
AIR BLK	.00	9:17pm
SUB TEST	.00	9:19pm
AIR BLK	.00	9:20pm

Reported AC: .00 g/210L

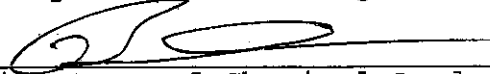

Signature of Chemical Analyst

Court CVR

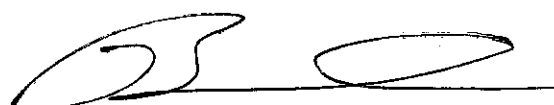

Analyst

Intox EC/IR-II: Preventive Maintenance

EDGECOMBE COUNTY BAT MOBILE UNIT 6 320

Serial Number: 008580 Test Record Number: 2461
Test Date: 09/29/2018 Test Time: 9:26pm EDT

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	9:27pm
FLO	Pass	9:27pm
FC	Pass	9:27pm

Temperature Tests

Test	Status	Time
FC1	Pass	9:27pm
SRC	Pass	9:27pm
DET	Pass	9:27pm
BAR	Pass	9:27pm
BT	Pass	9:27pm

Blank Tests

Test	Status	Time
AIR	Pass	9:27pm

Printer Tests

Test	Status	Time
PRNT	Pass	9:28pm

CRC Tests

Test	Status	Time
COMP	Pass	9:28pm
CAL	Pass	9:28pm

Preventive Maintenance
Status: Pass



Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH

**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**

County FRANKLIN Instrument Location FRANKLINTON PD

Instrument Serial No. 008815 #7 W. MASON ST FRANKLINTON, NC

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 04 day of SEPTEMBER, 20 18 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.



B. D. Small

Signature of Certifying Official

637

Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

FRANKLIN COUNTY FRANKLINTON PD 340

Serial Number: 008815
Test Date: 09/04/2018

Citation Number: M0000000-0
Subject's Name:
PREVENTIVE, MAINTENANCE
Subject's Date of Birth: 11/11/1911
Subject's Sex: Male
Driver's License State: XX
Driver's License Number: NONE

Analyst's Name: SMITH, BRIAN D
Permit Number: 08937E
Effective:
08/01/2017-08/01/2019

Officer's Name: NONE, NONE
Type of Agency: FTA
Agency: DHHS
Test Type: Breath Test

Lot Number: AG814901
Exp Date: 05/29/2020

Test	g/210L	Time
DIAG	Pass	9:16am
AIR BLK	.00	9:16am
ACCY CHK	.07	9:17am
AIR BLK	.00	9:18am
SUB TEST	.00	9:19am
AIR BLK	.00	9:20am
SUB TEST	.00	9:21am
AIR BLK	.00	9:22am

Reported AC: .00 g/210L



Signature of Chemical Analyst

Court CVR



Analyst

Intox EC/IR-II: Preventive Maintenance

FRANKLIN COUNTY FRANKLINTON PD 340

Serial Number: 008815 Test Record Number: 1156
Test Date: 09/04/2018 Test Time: 9:28am EDT

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	9:29am
FLO	Pass	9:29am
FC	Pass	9:29am

Temperature Tests

Test	Status	Time
FC1	Pass	9:29am
SRC	Pass	9:29am
DET	Pass	9:29am
BAR	Pass	9:29am
BT	Pass	9:29am

Blank Tests

Test	Status	Time
AIR	Pass	9:30am

Printer Tests

Test	Status	Time
PRNT	Pass	9:30am

CRC Tests

Test	Status	Time
COMP	Pass	9:30am
CAL	Pass	9:30am

Preventive Maintenance
Status: Pass



Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH

**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**

County Gaston Instrument Location Belmont P.D.

Instrument Serial No. 008733 201 Chronicle Street, Belmont

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 19th day of September, 2018, the foregoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.



Joseph E. [Signature] 650
Signature of Certifying Official Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

GASTON COUNTY BELMONT PD 350

Serial Number: 008733
Test Date: 09/19/2018

Citation Number: M0000000-0

Subject's Name:

PREVENTIVE, MAINTENANCE

Subject's Date of Birth: 11/11/1911

Subject's Sex: Male

Driver's License State: XX

Driver's License Number: NONE

Analyst's Name: HUTCHINSON, JOSEPH E

Permit Number: 19951E

Effective:

08/01/2017-08/01/2019

Officer's Name: NONE, NONE

Type of Agency: FTA

Agency: DHHS


Test Type: Breath Test

Lot Number: AG722408

Exp Date: 08/12/2019

Test	g/210L	Time
DIAG	Pass	11:44am
AIR BLK	.00	11:45am
ACCY CHK	.08	11:45am
AIR BLK	.00	11:46am
SUB TEST	.00	11:47am
AIR BLK	.00	11:48am
SUB TEST	.00	11:49am
AIR BLK	.00	11:51am

Reported AC: .00 g/210L

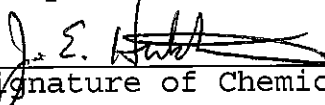

Signature of Chemical Analyst

Court CVR



Analyst

**This form is used when performing Preventive Maintenance procedures
Forensic Tests for Alcohol Branch
Department of Health and Human Services
Rev. 12/2007**

Intox EC/IR-II: Preventive Maintenance

GASTON COUNTY BELMONT PD 350

Serial Number: 008733 Test Record Number: 1111
Test Date: 09/19/2018 Test Time: 11:51am EDT

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	11:52am
FLO	Pass	11:52am
FC	Pass	11:52am

Temperature Tests

Test	Status	Time
FC1	Pass	11:52am
SRC	Pass	11:52am
DET	Pass	11:52am
BAR	Pass	11:52am
BT	Pass	11:52am

Blank Tests

Test	Status	Time
AIR	Pass	11:53am


Printer Tests

Test	Status	Time
PRNT	Pass	11:53am

CRC Tests

Test	Status	Time
COMP	Pass	11:53am
CAL	Pass	11:53am

Preventive Maintenance
Status: Pass


Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH


**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**

County Gaston Instrument Location Gaston County SD
Instrument Serial No. 008684 425 N. Marlette St, Gastonia

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 5th day of September, 2018 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

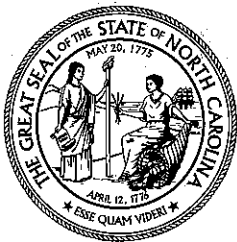

[Signature]
Signature of Certifying Official

656
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

GASTON COUNTY GASTON COUNTY SD 350

Serial Number: 008684
Test Date: 09/05/2018

Citation Number: M0000000-0
Subject's Name:
PREVENTIVE, MAINTENANCE
Subject's Date of Birth: 11/11/1911
Subject's Sex: Male
Driver's License State: XX
Driver's License Number: NONE

Analyst's Name: HAYS, MARK D
Permit Number: 15924E
Effective:
01/01/2018-01/01/2020

Officer's Name: NONE, NONE
Type of Agency: FTA
Agency: DHHS
Test Type: Breath Test

Lot Number: AG721401
Exp Date: 08/02/2019

Test	g/210L	Time
DIAG	Pass	10:21am
AIR BLK	.00	10:22am
ACCY CHK	.08	10:23am
AIR BLK	.00	10:24am
SUB TEST	.00	10:24am
AIR BLK	.00	10:25am
SUB TEST	.00	10:27am
AIR BLK	.00	10:28am

Reported AC:  .00 g/210L

Signature of Chemical Analyst

Court CVR


Analyst

**This form is used when performing Preventive Maintenance procedures
Forensic Tests for Alcohol Branch
Department of Health and Human Services
Rev. 12/2007**

Intox EC/IR-II: Preventive Maintenance

GASTON COUNTY GASTON COUNTY SD 350

Serial Number: 008684 Test Record Number: 4106
Test Date: 09/05/2018 Test Time: 10:29am EDT

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	10:29am
FLO	Pass	10:29am
FC	Pass	10:29am

Temperature Tests

Test	Status	Time
FC1	Pass	10:29am
SRC	Pass	10:29am
DET	Pass	10:29am
BAR	Pass	10:29am
BT	Pass	10:29am

Blank Tests

Test	Status	Time
AIR	Pass	10:30am


Printer Tests

Test	Status	Time
PRNT	Pass	10:30am

CRC Tests

Test	Status	Time
COMP	Pass	10:30am
CAL	Pass	10:30am

Preventive Maintenance
Status: Pass

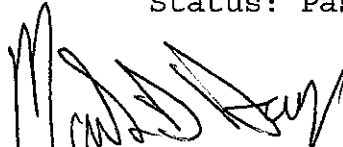

Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH


**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**

County Guilford Instrument Location UNC - Greensboro
Instrument Serial No. 008604 Police Department

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 27 day of September, 2018, the foregoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

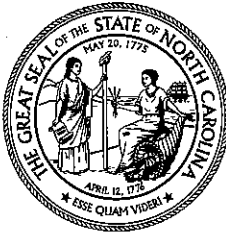

S. Kevin Dean
Signature of Certifying Official

642
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

GUILFORD COUNTY UNC-G POLICE DEPT 400

Serial Number: 008604
Test Date: 09/27/2018

Citation Number: M0000000-0

Subject's Name:

PREVENTIVE, MAINTENANCE

Subject's Date of Birth: 11/11/1911

Subject's Sex: Male

Driver's License State: XX

Driver's License Number: NONE

Analyst's Name: DEAN, L

Permit Number: 11598E

Effective:

04/01/2017-04/01/2019

Officer's Name: NONE, NONE

Type of Agency: FTA

Agency: DHHS


Test Type: Breath Test

Lot Number: AG807101

Exp Date: 03/12/2020

Test	g/210L	Time
DIAG	Pass	10:58am
AIR BLK	.00	10:58am
ACCY CHK	.08	10:59am
AIR BLK	.00	11:00am
SUB TEST	.00	11:01am
AIR BLK	.00	11:02am
SUB TEST	.00	11:03am
AIR BLK	.00	11:04am

Reported AC: .00 g/210L



Signature of Chemical Analyst

Court CVR



Analyst

This form is used when performing Preventive Maintenance procedures

Forensic Tests for Alcohol Branch

Department of Health and Human Services

Rev. 12/2007

Intox EC/IR-II: Preventive Maintenance

GUILFORD COUNTY UNC-G POLICE DEPT 400

Serial Number: 008604 Test Record Number: 1702
Test Date: 09/27/2018 Test Time: 11:05am EDT

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	11:05am
FLO	Pass	11:05am
FC	Pass	11:06am

Temperature Tests

Test	Status	Time
FC1	Pass	11:06am
SRC	Pass	11:06am
DET	Pass	11:06am
BAR	Pass	11:06am
BT	Pass	11:06am

Blank Tests

Test	Status	Time
AIR	Pass	11:06am

Printer Tests

Test	Status	Time
PRNT	Pass	11:06am

CRC Tests

Test	Status	Time
COMP	Pass	11:07am
CAL	Pass	11:07am

Preventive Maintenance
Status: Pass



Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH

**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**


County Guilford Instrument Location Greensboro PD

Instrument Serial No. 008725 100 Police Plaza, Greensboro, N.C.

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 27 day of September, 2018, the foregoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

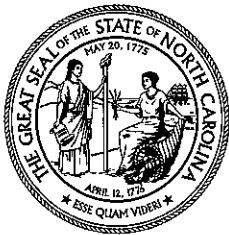

L. Kevin Dean
Signature of Certifying Official

642
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

GUILFORD COUNTY GREENSBORO PD 400

Serial Number: 008725
Test Date: 09/27/2018

Citation Number: M0000000-0
Subject's Name:
PREVENTIVE, MAINTENANCE
Subject's Date of Birth: 11/11/1911
Subject's Sex: Male
Driver's License State: XX
Driver's License Number: NONE


Analyst's Name: DEAN, L
Permit Number: 11598E
Effective:
04/01/2017-04/01/2019

Officer's Name: NONE, NONE
Type of Agency: FTA
Agency: DHHS
Test Type: Breath Test

Lot Number: AG734101
Exp Date: 12/07/2019


Test	g/210L	Time
DIAG	Pass	10:08am
AIR BLK	.00	10:09am
ACCY CHK	.08	10:09am
AIR BLK	.00	10:10am
SUB TEST	.00	10:11am
AIR BLK	.00	10:12am
SUB TEST	.00	10:13am
AIR BLK	.00	10:14am

Reported AC: .00 g/210L



Signature of Chemical Analyst

Court CVR

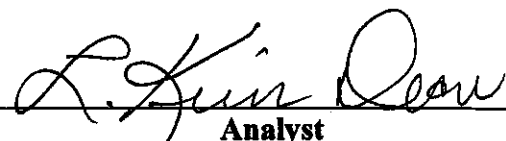

Analyst

Intox EC/IR-II: Preventive Maintenance

GUILFORD COUNTY GREENSBORO PD 400

Serial Number: 008725 Test Record Number: 4281
Test Date: 09/27/2018 Test Time: 10:15am EDT

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	10:16am
FLO	Pass	10:16am
FC	Pass	10:16am

Temperature Tests

Test	Status	Time
FC1	Pass	10:16am
SRC	Pass	10:16am
DET	Pass	10:16am
BAR	Pass	10:16am
BT	Pass	10:16am

Blank Tests

Test	Status	Time
AIR	Pass	10:16am

Printer Tests

Test	Status	Time
PRNT	Pass	10:16am

CRC Tests

Test	Status	Time
COMP	Pass	10:17am
CAL	Pass	10:17am

Preventive Maintenance
Status: Pass



Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH

**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**

County Guilford


Instrument Location Greensboro Jail

Instrument Serial No. 00 8638

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 26 day of September, 20 18, the foregoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

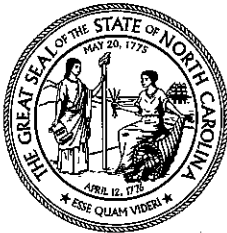

L. Kevin Dean
Signature of Certifying Official

642
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

GUILFORD COUNTY GREENSBORO JAIL 400

Serial Number: 008638
Test Date: 09/26/2018

Citation Number: M0000000-0
Subject's Name:
PREVENTIVE, MAINTENANCE
Subject's Date of Birth: 11/11/1911
Subject's Sex: Male
Driver's License State: XX
Driver's License Number: NONE

Analyst's Name: DEAN, L
Permit Number: 11598E
Effective:
04/01/2017-04/01/2019

Officer's Name: NONE, NONE
Type of Agency: FTA
Agency: DHHS
Test Type: Breath Test

Lot Number: AG807101
Exp Date: 03/12/2020


Test	g/210L	Time
DIAG	Pass	5:01pm
AIR BLK	.00	5:01pm
ACCY CHK	.08	5:02pm
AIR BLK	.00	5:03pm
SUB TEST	.00	5:04pm
AIR BLK	.00	5:05pm
SUB TEST	.00	5:06pm
AIR BLK	.00	5:07pm

Reported AC: .00 g/210L



Signature of Chemical Analyst

Court CVR

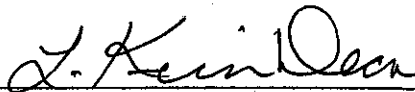

Analyst

**This form is used when performing Preventive Maintenance procedures
Forensic Tests for Alcohol Branch
Department of Health and Human Services
Rev. 12/2007**

Intox EC/IR-II: Preventive Maintenance

GUILFORD COUNTY GREENSBORO JAIL 400

Serial Number: 008638 Test Record Number: 3775
Test Date: 09/26/2018 Test Time: 5:09pm EDT

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	5:09pm
FLO	Pass	5:09pm
FC	Pass	5:09pm

Temperature Tests

Test	Status	Time
FC1	Pass	5:09pm
SRC	Pass	5:09pm
DET	Pass	5:09pm
BAR	Pass	5:09pm
BT	Pass	5:09pm

Blank Tests

Test	Status	Time
AIR	Pass	5:10pm

Printer Tests

Test	Status	Time
PRNT	Pass	5:10pm

CRC Tests

Test	Status	Time
COMP	Pass	5:10pm
CAL	Pass	5:10pm

Preventive Maintenance
Status: Pass



Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH

**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**

County Guilford


Instrument Location Greensboro Jail

Instrument Serial No. 008794

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 26 day of September, 2018, the foregoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

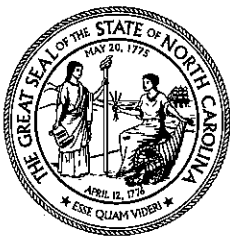

L. Kevin Dean
Signature of Certifying Official

642
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

GUILFORD COUNTY GREENSBORO JAIL 400

Serial Number: 008794

Test Date: 09/26/2018

Citation Number: M0000000-0

Subject's Name:

PREVENTIVE, MAINTENANCE

Subject's Date of Birth: 11/11/1911

Subject's Sex: Male

Driver's License State: XX

Driver's License Number: NONE

Analyst's Name: DEAN, L

Permit Number: 11598E

Effective:

04/01/2017-04/01/2019

Officer's Name: NONE, NONE

Type of Agency: FTA

Agency: DHHS


Test Type: Breath Test

Lot Number: AG807101

Exp Date: 03/12/2020

Test	g/210L	Time
DIAG	Pass	4:38pm
AIR BLK	.00	4:38pm
ACCY CHK	.07	4:39pm
AIR BLK	.00	4:40pm
SUB TEST	.00	4:41pm
AIR BLK	.00	4:41pm
SUB TEST	.00	4:43pm
AIR BLK	.00	4:44pm

Reported AC: .00 g/210L


Signature of Chemical Analyst

Court CVR


Analyst

Intox EC/IR-II: Preventive Maintenance

GUILFORD COUNTY GREENSBORO JAIL 400

Serial Number: 008794 Test Record Number: 6231
Test Date: 09/26/2018 Test Time: 4:45pm EDT

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	4:45pm
FLO	Pass	4:45pm
FC	Pass	4:45pm

Temperature Tests

Test	Status	Time
FC1	Pass	4:45pm
SRC	Pass	4:45pm
DET	Pass	4:45pm
BAR	Pass	4:45pm
BT	Pass	4:45pm

Blank Tests

Test	Status	Time
AIR	Pass	4:46pm


Printer Tests

Test	Status	Time
PRNT	Pass	4:46pm

CRC Tests

Test	Status	Time
COMP	Pass	4:46pm
CAL	Pass	4:46pm

Preventive Maintenance
Status: Pass



Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH

PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II

County Guilford


Instrument Location Greensboro Jail

Instrument Serial No. 008718

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 26 day of September, 2018, the foregoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

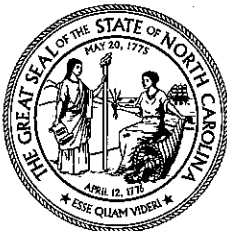

L. Kevin Dean
Signature of Certifying Official

642
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

GUILFORD COUNTY GREENSBORO JAIL 400

Serial Number: 008718
Test Date: 09/26/2018

Citation Number: M0000000-0
Subject's Name:
PREVENTIVE, MAINTENANCE
Subject's Date of Birth: 11/11/1911
Subject's Sex: Male
Driver's License State: XX
Driver's License Number: NONE


Analyst's Name: DEAN, L
Permit Number: 11598E
Effective:
04/01/2017-04/01/2019

Officer's Name: NONE, NONE
Type of Agency: FTA
Agency: DHHS
Test Type: Breath Test

Lot Number: AG807101
Exp Date: 03/12/2020

Test	g/210L	Time
DIAG	Pass	4:08pm
AIR BLK	.00	4:08pm
ACCY CHK	.08	4:09pm
AIR BLK	.00	4:10pm
SUB TEST	.00	4:10pm
AIR BLK	.00	4:11pm
SUB TEST	.00	4:13pm
AIR BLK	.00	4:14pm

Reported AC: .00 g/210L


Signature of Chemical Analyst

Court CVR


Analyst

Intox EC/IR-II: Preventive Maintenance

GUILFORD COUNTY GREENSBORO JAIL 400

Serial Number: 008718 Test Record Number: 1843
Test Date: 09/26/2018 Test Time: 4:15pm EDT

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	4:15pm
FLO	Pass	4:15pm
FC	Pass	4:15pm

Temperature Tests

Test	Status	Time
FC1	Pass	4:15pm
SRC	Pass	4:15pm
DET	Pass	4:15pm
BAR	Pass	4:15pm
BT	Pass	4:15pm

Blank Tests

Test	Status	Time
AIR	Pass	4:16pm


Printer Tests

Test	Status	Time
PRNT	Pass	4:16pm

CRC Tests

Test	Status	Time
COMP	Pass	4:16pm
CAL	Pass	4:16pm

Preventive Maintenance
Status: Pass

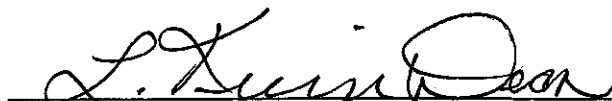

Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH


**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**

County Guilford Instrument Location High Point
Instrument Serial No. 008828 Police Department

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 26 day of September, 2018, the foregoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

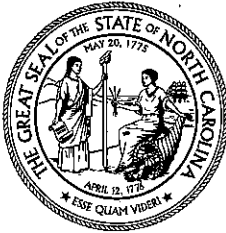

L. Kevin Owen
Signature of Certifying Official

642
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

GUILFORD COUNTY HIGH POINT PD 401

Serial Number: 008828
Test Date: 09/26/2018

Citation Number: M0000000-0
Subject's Name:
PREVENTIVE, MAINTENANCE
Subject's Date of Birth: 11/11/1911
Subject's Sex: Male
Driver's License State: XX
Driver's License Number: NONE


Analyst's Name: DEAN, L
Permit Number: 11598E
Effective:
04/01/2017-04/01/2019

Officer's Name: NONE, NONE
Type of Agency: FTA
Agency: DHHS
Test Type: Breath Test

Lot Number: AG807101
Exp Date: 03/12/2020

Test	g/210L	Time
DIAG	Pass	1:58pm
AIR BLK	.00	1:59pm
ACCY CHK	.08	2:00pm
AIR BLK	.00	2:00pm
SUB TEST	.00	2:01pm
AIR BLK	.00	2:02pm
SUB TEST	.00	2:04pm
AIR BLK	.00	2:05pm

Reported AC: .00 g/210L

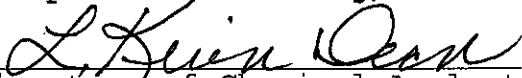
Signature of Chemical Analyst

Court CVR


Analyst

Intox EC/IR-II: Preventive Maintenance

GUILFORD COUNTY HIGH POINT PD 401

Serial Number: 008828 Test Record Number: 2790
Test Date: 09/26/2018 Test Time: 2:06pm EDT

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	2:06pm
FLO	Pass	2:06pm
FC	Pass	2:06pm

Temperature Tests

Test	Status	Time
FC1	Pass	2:07pm
SRC	Pass	2:07pm
DET	Pass	2:07pm
BAR	Pass	2:07pm
BT	Pass	2:07pm

Blank Tests

Test	Status	Time
AIR	Pass	2:07pm

Printer Tests

Test	Status	Time
PRNT	Pass	2:07pm

CRC Tests

Test	Status	Time
COMP	Pass	2:07pm
CAL	Pass	2:07pm

Preventive Maintenance
Status: Pass



Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH

**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**


County HARNETT Instrument Location HARNETT COUNTY

Instrument Serial No. 008850 DETENTION CENTER
LILLINGTON, NC

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 25 day of SEPTEMBER, 2018, the foregoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

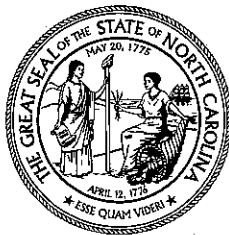

Alan R. Barnes
Signature of Certifying Official

648
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

HARNETT COUNTY DETENTION CENTER 420

Serial Number: 008850
Test Date: 09/25/2018

Citation Number: M0000000-0
Subject's Name:
PREVENTIVE, MAINTENANCE
Subject's Date of Birth: 11/11/1911
Subject's Sex: Male
Driver's License State: XX
Driver's License Number: NONE

Analyst's Name: BARNES, ALVIN R
Permit Number: 15671E
Effective:
07/01/2017-07/01/2019

Officer's Name: NONE, NONE
Type of Agency: FTA
Agency: DHHS
Test Type: Breath Test

Lot Number: AG702302
Exp Date: 01/23/2019

Test	g/210L	Time
DIAG	Pass	6:37pm
AIR BLK	.00	6:38pm
ACCY CHK	.07	6:38pm
AIR BLK	.00	6:39pm
SUB TEST	.00	6:40pm
AIR BLK	.00	6:41pm
SUB TEST	.00	6:42pm
AIR BLK	.00	6:43pm

Reported AC: .00 g/210L

Signature of Chemical Analyst

Court CVR



Analyst

Intox EC/IR-II: Preventive Maintenance

HARNETT COUNTY DETENTION CENTER 420

Serial Number: 008850 Test Record Number: 681
Test Date: 09/25/2018 Test Time: 6:43pm EDT

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	6:44pm
FLO	Pass	6:44pm
FC	Pass	6:44pm

Temperature Tests

Test	Status	Time
FC1	Pass	6:44pm
SRC	Pass	6:44pm
DET	Pass	6:44pm
BAR	Pass	6:44pm
BT	Pass	6:44pm

Blank Tests

Test	Status	Time
AIR	Pass	6:45pm


Printer Tests

Test	Status	Time
PRNT	Pass	6:45pm

CRC Tests

Test	Status	Time
COMP	Pass	6:45pm
CAL	Pass	6:45pm

Preventive Maintenance
Status: Pass

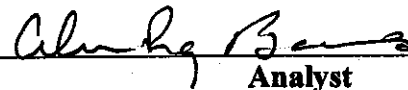
Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH


PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II

County Hertford Instrument Location Ahoskie P.D.
Instrument Serial No. 008848 705 W. Main St., Ahoskie, NC.

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 19th day of September, 2018 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

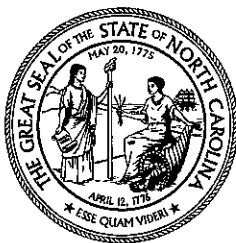

[Signature]
Signature of Certifying Official

643
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

HERTFORD COUNTY AHOSKIE PD 450

Serial Number: 008848
Test Date: 09/19/2018

Citation Number: M0000000-0
Subject's Name:

PREVENTIVE, MAINTENANCE

Subject's Date of Birth: 11/11/1911

Subject's Sex: Male

Driver's License State: XX

Driver's License Number: NONE

Analyst's Name: GUARD, KELLY G

Permit Number: 12955E

Effective:

06/01/2017-06/01/2019

Officer's Name: NONE, NONE

Type of Agency: FTA

Agency: DHHS


Test Type: Breath Test

Lot Number: AG702402

Exp Date: 01/24/2019


Test	g/210L	Time
DIAG	Pass	10:41am
AIR BLK	.00	10:42am
ACCY CHK	.07	10:43am
AIR BLK	.00	10:44am
SUB TEST	.00	10:44am
AIR BLK	.00	10:45am
SUB TEST	.00	10:47am
AIR BLK	.00	10:48am

Reported AC: .00 g/210L

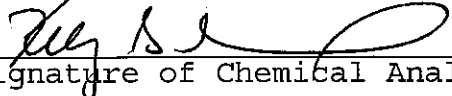

Signature of Chemical Analyst

Court CVR

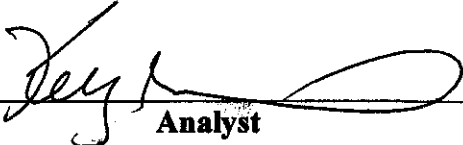

Analyst

**This form is used when performing Preventive Maintenance procedures
Forensic Tests for Alcohol Branch
Department of Health and Human Services
Rev. 12/2007**

Intox EC/IR-II: Preventive Maintenance

HERTFORD COUNTY AHOSKIE PD 450

Serial Number: 008848 Test Record Number: 1359
Test Date: 09/19/2018 Test Time: 10:50am EDT

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	10:50am
FLO	Pass	10:50am
FC	Pass	10:50am

Temperature Tests

Test	Status	Time
FC1	Pass	10:50am
SRC	Pass	10:50am
DET	Pass	10:50am
BAR	Pass	10:50am
BT	Pass	10:50am

Blank Tests

Test	Status	Time
AIR	Pass	10:51am


Printer Tests

Test	Status	Time
PRNT	Pass	10:51am

CRC Tests

Test	Status	Time
COMP	Pass	10:51am
CAL	Pass	10:51am

Preventive Maintenance
Status: Pass

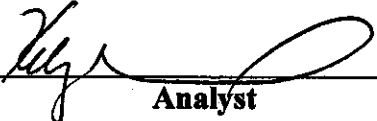

Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH

**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**

County Hoke Co. Instrument Location Hoke Co. Detention CTR

Instrument Serial No. 008855 Raeford, NC

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 12th day of September, 2018 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.



[Signature]
Signature of Certifying Official

654
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

HOKE COUNTY DETENTION CENTER 460

Serial Number: 008855
Test Date: 09/12/2018

Citation Number: M0000000-0
Subject's Name:
PREVENTIVE, MAINTENANCE
Subject's Date of Birth: 11/11/1911
Subject's Sex: Male
Driver's License State: XX
Driver's License Number: NONE


Analyst's Name: KEESLER, GRAYHAM C
Permit Number: 07682E
Effective:
12/01/2017-12/01/2019

Officer's Name: NONE, NONE
Type of Agency: FTA
Agency: DHHS
Test Type: Breath Test

Lot Number: AG710701
Exp Date: 04/17/2019


Test	g/210L	Time
DIAG	Pass	12:08pm
AIR BLK	.00	12:08pm
ACCY CHK	.07	12:09pm
AIR BLK	.00	12:10pm
SUB TEST	.00	12:10pm
AIR BLK	.00	12:11pm
SUB TEST	.00	12:13pm
AIR BLK	.00	12:14pm

Reported AC: .00 g/210L

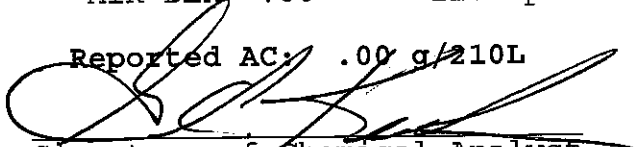

Signature of Chemical Analyst

Court CVR

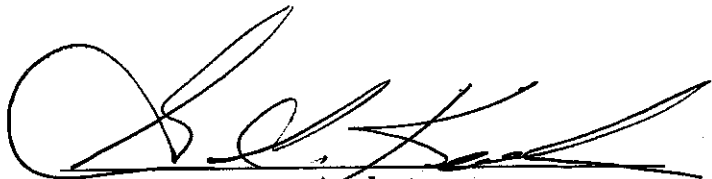

Analyst

Intox EC/IR-II: Preventive Maintenance

HOKE COUNTY DETENTION CENTER 460

Serial Number: 008855 Test Record Number: 1419
Test Date: 09/12/2018 Test Time: 12:15pm EDT

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	12:15pm
FLO	Pass	12:15pm
FC	Pass	12:15pm

Temperature Tests

Test	Status	Time
FC1	Pass	12:15pm
SRC	Pass	12:15pm
DET	Pass	12:15pm
BAR	Pass	12:15pm
BT	Pass	12:15pm

Blank Tests

Test	Status	Time
AIR	Pass	12:16pm


Printer Tests

Test	Status	Time
PRNT	Pass	12:16pm

CRC Tests

Test	Status	Time
COMP	Pass	12:16pm
CAL	Pass	12:16pm

Preventive Maintenance
Status: Pass

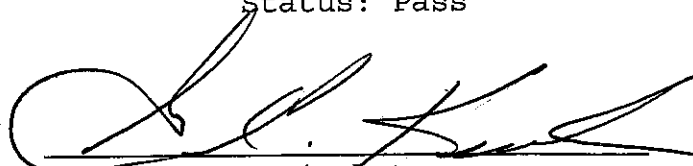

Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH

**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**

County Hoke Co. Instrument Location Hoke Co Detention CTR
Instrument Serial No. 008852 RAE Ford, NC

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 12th day of September, 20 18 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.



[Signature]
Signature of Certifying Official

654
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

HOKE COUNTY DETENTION CENTER 460

Serial Number: 008852
Test Date: 09/12/2018

Citation Number: M0000000-0
Subject's Name:
PREVENTIVE, MAINTENANCE
Subject's Date of Birth: 11/11/1911
Subject's Sex: Male
Driver's License State: XX
Driver's License Number: NONE

Analyst's Name: KEESLER, GRAYHAM C
Permit Number: 07682E
Effective:
12/01/2017-12/01/2019

Officer's Name: NONE, NONE
Type of Agency: FTA
Agency: DHHS
Test Type: Breath Test

Lot Number: AG805802
Exp Date: 02/27/2020

Test	g/210L	Time
DIAG	Pass	12:03pm
AIR BLK	.00	12:03pm
ACCY CHK	.08	12:04pm
AIR BLK	.00	12:05pm
SUB TEST	.00	12:05pm
AIR BLK	.00	12:06pm
SUB TEST	.00	12:08pm
AIR BLK	.00	12:09pm

Reported ACt .00 g/210L


Signature of Chemical Analyst

Court CVR


Analyst

**This form is used when performing Preventive Maintenance procedures
Forensic Tests for Alcohol Branch
Department of Health and Human Services
Rev. 12/2007**

Intox EC/IR-II: Preventive Maintenance

HOKE COUNTY DETENTION CENTER 460

Serial Number: 008852 Test Record Number: 875
Test Date: 09/12/2018 Test Time: 12:09pm EDT

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	12:10pm
FLO	Pass	12:10pm
FC	Pass	12:10pm

Temperature Tests

Test	Status	Time
FC1	Pass	12:10pm
SRC	Pass	12:10pm
DET	Pass	12:10pm
BAR	Pass	12:10pm
BT	Pass	12:10pm

Blank Tests

Test	Status	Time
AIR	Pass	12:11pm

Printer Tests


Test	Status	Time
PRNT	Pass	12:11pm

CRC Tests

Test	Status	Time
COMP	Pass	12:11pm
CAL	Pass	12:11pm

Preventive Maintenance

Status: Pass

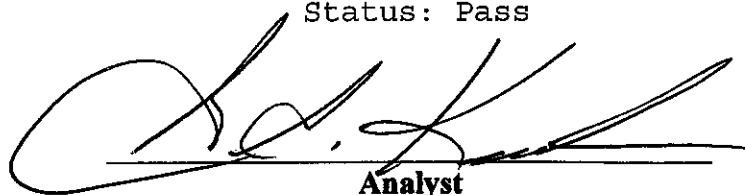

Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH

**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**

County Hyde Instrument Location Hyde Co. S.O.

Instrument Serial No. 008588 1233 MAIN ST., SWAN QUETER, N.C.

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 28th day of SEPTEMBER, 20 18 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.



Linda A. Kessel
Signature of Certifying Official

647
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

HYDE COUNTY HYDE CO SO SWAN QUAR 470

Serial Number: 008588
Test Date: 09/28/2018

Citation Number: M0000000-0
Subject's Name:
PREVENTIVE, MAINTENANCE
Subject's Date of Birth: 11/11/1911
Subject's Sex: Male
Driver's License State: XX
Driver's License Number: NONE


Analyst's Name: KEESLER, LINDA A
Permit Number: 11646E
Effective:
08/01/2017-08/01/2019

Officer's Name: NONE, NONE
Type of Agency: FTA
Agency: DHHS
Test Type: Breath Test

Lot Number: AG702402
Exp Date: 01/24/2019


Test	g/210L	Time
DIAG	Pass	12:45pm
AIR BLK	.00	12:46pm
ACCY CHK	.08	12:46pm
AIR BLK	.00	12:47pm
SUB TEST	.00	12:48pm
AIR BLK	.00	12:49pm
SUB TEST	.00	12:50pm
AIR BLK	.00	12:51pm

Reported AC: .00 g/210L

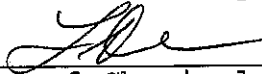

Signature of Chemical Analyst

Court CVR

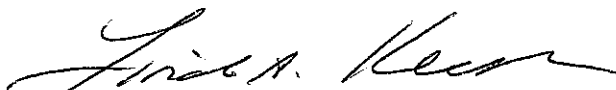

Analyst

Intox EC/IR-II: Preventive Maintenance

HYDE COUNTY HYDE CO SO SWAN QUAR 470

Serial Number: 008588 Test Record Number: 990
Test Date: 09/28/2018 Test Time: 12:52pm EDT

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	12:53pm
FLO	Pass	12:53pm
FC	Pass	12:53pm

Temperature Tests

Test	Status	Time
FC1	Pass	12:53pm
SRC	Pass	12:53pm
DET	Pass	12:53pm
BAR	Pass	12:53pm
BT	Pass	12:53pm

Blank Tests

Test	Status	Time
AIR	Pass	12:54pm


Printer Tests

Test	Status	Time
PRNT	Pass	12:54pm

CRC Tests

Test	Status	Time
COMP	Pass	12:54pm
CAL	Pass	12:54pm

Preventive Maintenance
Status: Pass

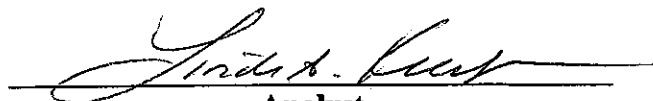
Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH


PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II

County Iredell Instrument Location BAT mobile Unit 2
Instrument Serial No. 608970 Willife

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 2 day of September, 20 10, the foregoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

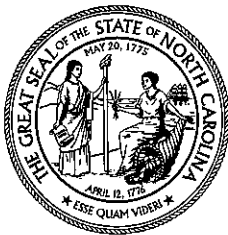

Phil V. Jones

Signature of Certifying Official

656

Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

IREDELL BAT MOBILE UNIT 02 480

Serial Number: 008970

Test Date: 09/02/2018

Citation Number: M0000000-0

Subject's Name:

PREVENTIVE, MAINTENANCE

Subject's Date of Birth: 11/11/1911

Subject's Sex: Male

Driver's License State: XX

Driver's License Number: NONE

Analyst's Name: TOWERY, CHAD V

Permit Number: 26632E

Effective:

06/01/2017-06/01/2019

Officer's Name: NONE, NONE

Type of Agency: FTA

Agency: DHHS


Test Type: Breath Test

Lot Number: AG807101

Exp Date: 03/12/2020

Test	g/210L	Time
DIAG	Pass	2:06pm
AIR BLK	.00	2:07pm
ACCY CHK	.07	2:08pm
AIR BLK	.00	2:09pm
SUB TEST	.00	2:09pm
AIR BLK	.00	2:10pm
SUB TEST	.00	2:12pm
AIR BLK	.00	2:13pm

Reported AC: .00 g/210L

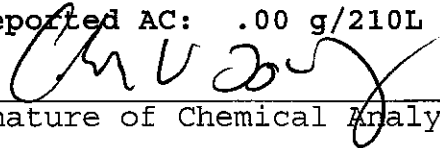

Signature of Chemical Analyst

Court CVR



Analyst

**This form is used when performing Preventive Maintenance procedures
Forensic Tests for Alcohol Branch
Department of Health and Human Services
Rev. 12/2007**

Intox EC/IR-II: Preventive Maintenance

IREDELL BAT MOBILE UNIT 02 480

Serial Number: 008970 Test Record Number: 509
Test Date: 09/02/2018 Test Time: 2:13pm EDT

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	2:14pm
FLO	Pass	2:14pm
FC	Pass	2:14pm

Temperature Tests

Test	Status	Time
FC1	Pass	2:14pm
SRC	Pass	2:14pm
DET	Pass	2:14pm
BAR	Pass	2:14pm
BT	Pass	2:14pm

Blank Tests

Test	Status	Time
AIR	Pass	2:15pm


Printer Tests

Test	Status	Time
PRNT	Pass	2:15pm

CRC Tests

Test	Status	Time
COMP	Pass	2:15pm
CAL	Pass	2:15pm

Preventive Maintenance
Status: Pass

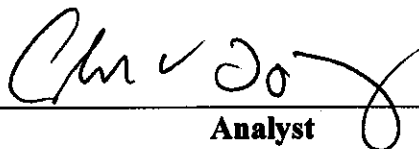

Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH


PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II

County Tredell Instrument Location Bat Mobile Unit 2
Instrument Serial No. 0073 Widdie

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 2 day of September, 2018, the foregoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

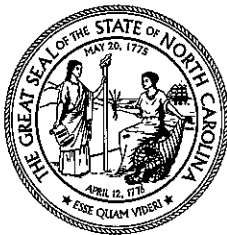

Chris Doy

Signature of Certifying Official

658

Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

IREDELL BAT MOBILE UNIT 02 480

Serial Number: 008973
Test Date: 09/02/2018

Citation Number: M0000000-0
Subject's Name:
PREVENTIVE, MAINTENANCE
Subject's Date of Birth: 11/11/1911
Subject's Sex: Male
Driver's License State: XX
Driver's License Number: NONE


Analyst's Name: TOWERY, CHAD V
Permit Number: 26632E
Effective:
06/01/2017-06/01/2019

Officer's Name: NONE, NONE
Type of Agency: FTA
Agency: DHHS
Test Type: Breath Test

Lot Number: AG807101
Exp Date: 03/12/2020

Test	g/210L	Time
DIAG	Pass	2:08pm
AIR BLK	.00	2:09pm
ACCY CHK	.07	2:09pm
AIR BLK	.00	2:10pm
SUB TEST	.00	2:11pm
AIR BLK	.00	2:12pm
SUB TEST	.00	2:13pm
AIR BLK	.00	2:14pm

Reported AC: .00 g/210L

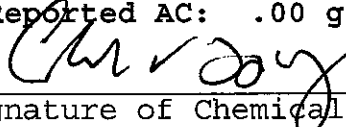

Signature of Chemical Analyst

Court CVR



Analyst

**This form is used when performing Preventive Maintenance procedures
Forensic Tests for Alcohol Branch
Department of Health and Human Services
Rev. 12/2007**

Intox EC/IR-II: Preventive Maintenance

IREDELL BAT MOBILE UNIT 02 480

Serial Number: 008973 Test Record Number: 548
Test Date: 09/02/2018 Test Time: 2:15pm EDT

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	2:15pm
FLO	Pass	2:15pm
FC	Pass	2:15pm

Temperature Tests

Test	Status	Time
FC1	Pass	2:15pm
SRC	Pass	2:15pm
DET	Pass	2:15pm
BAR	Pass	2:15pm
BT	Pass	2:15pm

Blank Tests

Test	Status	Time
AIR	Pass	2:16pm

Printer Tests

Test	Status	Time
PRNT	Pass	2:16pm

CRC Tests

Test	Status	Time
COMP	Pass	2:16pm
CAL	Pass	2:16pm

Preventive Maintenance
Status: Pass



Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH

**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**


County JOHNSTON Instrument Location JOHNSTON CO JAIL

Instrument Serial No. 008810 SMITHFIELD, NC

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 12 day of SEPTEMBER, 2018, the foregoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

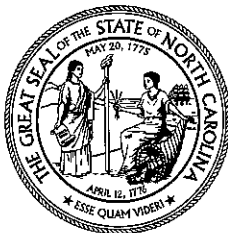

Alan R. Bowers
Signature of Certifying Official

648
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

JOHNSTON COUNTY JOHNSTON CO JAIL 500

Serial Number: 008810
Test Date: 09/12/2018

Citation Number: M0000000-0
Subject's Name:
PREVENTIVE, MAINTENANCE
Subject's Date of Birth: 11/11/1911
Subject's Sex: Male
Driver's License State: XX
Driver's License Number: NONE

Analyst's Name: BARNES, ALVIN R
Permit Number: 15671E
Effective:
07/01/2017-07/01/2019

Officer's Name: NONE, NONE
Type of Agency: FTA
Agency: DHHS
Test Type: Breath Test

Lot Number: AG716201
Exp Date: 06/11/2019

Test	g/210L	Time
DIAG	Pass	7:35am
AIR BLK	.00	7:36am
ACCY CHK	.08	7:36am
AIR BLK	.00	7:37am
SUB TEST	.00	7:38am
AIR BLK	.00	7:39am
SUB TEST	.00	7:40am
AIR BLK	.00	7:41am

Reported AC: .00 g/210L

Signature of Chemical Analyst

Court CVR



Analyst

Intox EC/IR-II: Preventive Maintenance

JOHNSTON COUNTY JOHNSTON CO JAIL 500

Serial Number: 008810 Test Record Number: 3291
Test Date: 09/12/2018 Test Time: 7:42am EDT

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	7:42am
FLO	Pass	7:42am
FC	Pass	7:42am

Temperature Tests

Test	Status	Time
FC1	Pass	7:42am
SRC	Pass	7:42am
DET	Pass	7:42am
BAR	Pass	7:42am
BT	Pass	7:42am

Blank Tests

Test	Status	Time
AIR	Pass	7:43am

Printer Tests

Test	Status	Time
PRNT	Pass	7:43am

CRC Tests

Test	Status	Time
COMP	Pass	7:43am
CAL	Pass	7:43am

Preventive Maintenance
Status: Pass



Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH

PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II

County JOHNSTON Instrument Location JOHNSTON CO JAIL
Instrument Serial No. 008846 SMITHFIELD, NC

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 26 day of SEPTEMBER, 20 18, the foregoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.



Alvin P. Barnes
Signature of Certifying Official

648
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

JOHNSTON COUNTY JOHNSTON CO JAIL 500

Serial Number: 008846
Test Date: 09/26/2018

Citation Number: M0000000-0
Subject's Name:

PREVENTIVE, MAINTENANCE

Subject's Date of Birth: 11/11/1911

Subject's Sex: Male

Driver's License State: XX

Driver's License Number: NONE

Analyst's Name: BARNES, ALVIN R

Permit Number: 15671E

Effective:

07/01/2017-07/01/2019

Officer's Name: NONE, NONE

Type of Agency: FTA

Agency: DHHS

Test Type: Breath Test

Lot Number: AG805802


Exp Date: 02/27/2020

Test	g/210L	Time
DIAG	Pass	5:46pm
AIR BLK	.00	5:47pm
ACCY CHK	.08	5:47pm
AIR BLK	.00	5:48pm
SUB TEST	.00	5:49pm
AIR BLK	.00	5:50pm
SUB TEST	.00	5:51pm
AIR BLK	.00	5:52pm

Reported AC: .00 g/210L

Signature of Chemical Analyst

Court CVR

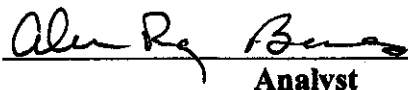
Analyst

Intox EC/IR-II: Preventive Maintenance

JOHNSTON COUNTY JOHNSTON CO JAIL 500

Serial Number: 008846 Test Record Number: 4377

Test Date: 09/26/2018 Test Time: 5:53pm EDT

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	5:53pm
FLO	Pass	5:53pm
FC	Pass	5:53pm

Temperature Tests

Test	Status	Time
FC1	Pass	5:54pm
SRC	Pass	5:54pm
DET	Pass	5:54pm
BAR	Pass	5:54pm
BT	Pass	5:54pm

Blank Tests

Test	Status	Time
AIR	Pass	5:54pm


Printer Tests

Test	Status	Time
PRNT	Pass	5:54pm

CRC Tests

Test	Status	Time
COMP	Pass	5:54pm
CAL	Pass	5:54pm

Preventive Maintenance
Status: Pass

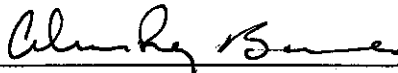

Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH

**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**


County Lenoir Instrument Location Kinston P.D.

Instrument Serial No. 008624 205 E. King St., Kinston, N.C.

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 28th day of September, 2018 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

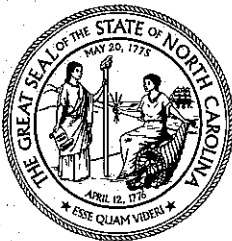

[Signature]
Signature of Certifying Official

643
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

LENOIR COUNTY KINSTON PD 530

Serial Number: 008624
Test Date: 09/28/2018

Citation Number: M0000000-0
Subject's Name:
PREVENTIVE, MAINTENANCE
Subject's Date of Birth: 11/11/1911
Subject's Sex: Male
Driver's License State: XX
Driver's License Number: NONE


Analyst's Name: GUARD, KELLY G
Permit Number: 12955E
Effective:
06/01/2017-06/01/2019

Officer's Name: NONE, NONE
Type of Agency: FTA
Agency: DHHS
Test Type: Breath Test

Lot Number: AG805801
Exp Date: 02/27/2020


Test	g/210L	Time
DIAG	Pass	10:11am
AIR BLK	.00	10:11am
ACCY CHK	.08	10:12am
AIR BLK	.00	10:13am
SUB TEST	.00	10:13am
AIR BLK	.00	10:14am
SUB TEST	.00	10:16am
AIR BLK	.00	10:17am

Reported AC: .00 g/210L

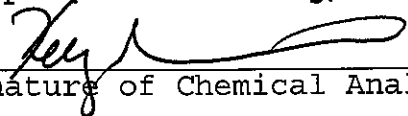

Signature of Chemical Analyst

Court CVR

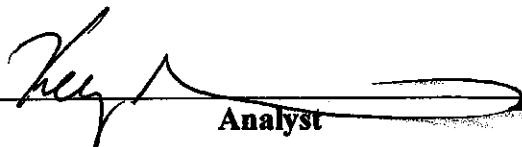

Analyst

Intox EC/IR-II: Preventive Maintenance

LENOIR COUNTY KINSTON PD 530

Serial Number: 008624 Test Record Number: 1720
Test Date: 09/28/2018 Test Time: 10:19am EDT

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	10:19am
FLO	Pass	10:19am
FC	Pass	10:19am

Temperature Tests

Test	Status	Time
FC1	Pass	10:19am
SRC	Pass	10:19am
DET	Pass	10:19am
BAR	Pass	10:19am
BT	Pass	10:19am

Blank Tests

Test	Status	Time
AIR	Pass	10:20am

Printer Tests

Test	Status	Time
PRNT	Pass	10:20am

CRC Tests

Test	Status	Time
COMP	Pass	10:21am
CAL	Pass	10:21am

Preventive Maintenance
Status: Pass



Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH


**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**

County Lenoir Instrument Location Lenoir Co. S.O.
Instrument Serial No. 008639 130 Queen St., Kinston, NC

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 28th day of September, 2018 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

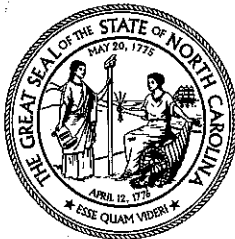

[Signature]
Signature of Certifying Official

643
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

LENOIR COUNTY LENOIR CO SO 530

Serial Number: 008639
Test Date: 09/28/2018

Citation Number: M0000000-0
Subject's Name:

PREVENTIVE, MAINTENANCE

Subject's Date of Birth: 11/11/1911

Subject's Sex: Male

Driver's License State: XX

Driver's License Number: NONE

Analyst's Name: GUARD, KELLY G

Permit Number: 12955E

Effective:

06/01/2017-06/01/2019

Officer's Name: NONE, NONE

Type of Agency: FTA

Agency: DHHS


Test Type: Breath Test

Lot Number: AG734101

Exp Date: 12/07/2019


Test	g/210L	Time
DIAG	Pass	9:32am
AIR BLK	.00	9:33am
ACCY CHK	.07	9:33am
AIR BLK	.00	9:34am
SUB TEST	.00	9:35am
AIR BLK	.00	9:36am
SUB TEST	.00	9:37am
AIR BLK	.00	9:38am

Reported AC: .00 g/210L

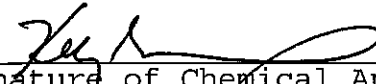

Signature of Chemical Analyst

Court CVR



Analyst

**This form is used when performing Preventive Maintenance procedures
Forensic Tests for Alcohol Branch
Department of Health and Human Services
Rev. 12/2007**

Intox EC/IR-II: Preventive Maintenance

LENOIR COUNTY LENOIR CO SO 530

Serial Number: 008639 Test Record Number: 3210
Test Date: 09/28/2018 Test Time: 9:39am EDT

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	9:40am
FLO	Pass	9:40am
FC	Pass	9:40am

Temperature Tests

Test	Status	Time
FC1	Pass	9:40am
SRC	Pass	9:40am
DET	Pass	9:40am
BAR	Pass	9:40am
BT	Pass	9:40am

Blank Tests

Test	Status	Time
AIR	Pass	9:41am

Printer Tests

Test	Status	Time
PRNT	Pass	9:41am

CRC Tests

Test	Status	Time
COMP	Pass	9:41am
CAL	Pass	9:41am

Preventive Maintenance
Status: Pass


Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH


**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**

County MARTIN Instrument Location MARTIN CO. S.D.
Instrument Serial No. 008912 305 E. MAIN ST., WILKINSON, N.C.

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 25th day of SEPTEMBER, 20 18 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

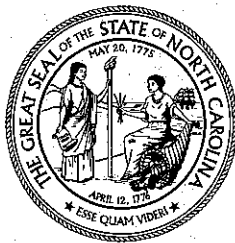

Linda Keel
Signature of Certifying Official

647
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

MARTIN COUNTY SHERIFF'S OFFICE 570

Serial Number: 008912

Test Date: 09/25/2018

Citation Number: M0000000-0

Subject's Name:

PREVENTIVE, MAINTENANCE

Subject's Date of Birth: 11/11/1911

Subject's Sex: Male

Driver's License State: XX

Driver's License Number: NONE

Analyst's Name: KEESLER, LINDA A

Permit Number: 11646E

Effective:

08/01/2017-08/01/2019

Officer's Name: NONE, NONE

Type of Agency: FTA

Agency: DHHS


Test Type: Breath Test

Lot Number: AG814901

Exp Date: 05/29/2020

Test	g/210L	Time
DIAG	Pass	1:22pm
AIR BLK	.00	1:23pm
ACCY CHK	.07	1:24pm
AIR BLK	.00	1:25pm
SUB TEST	.00	1:27pm
AIR BLK	.00	1:28pm
SUB TEST	.00	1:29pm
AIR BLK	.00	1:30pm

Reported AC: .00 g/210L


Signature of Chemical Analyst

Court CVR


Analyst

This form is used when performing Preventive Maintenance procedures
Forensic Tests for Alcohol Branch
Department of Health and Human Services
Rev. 12/2007

Intox EC/IR-II: Preventive Maintenance

MARTIN COUNTY SHERIFF'S OFFICE 570

Serial Number: 008912 Test Record Number: 1408
Test Date: 09/25/2018 Test Time: 1:31pm EDT

System Check: Passed

Baseline Tests

Test	Status	Time
IR	Pass	1:32pm
FLO	Pass	1:32pm
FC	Pass	1:32pm

Temperature Tests

Test	Status	Time
FC1	Pass	1:32pm
SRC	Pass	1:32pm
DET	Pass	1:32pm
BAR	Pass	1:32pm
BT	Pass	1:32pm

Blank Tests

Test	Status	Time
AIR	Pass	1:33pm


Printer Tests

Test	Status	Time
PRNT	Pass	1:33pm

CRC Tests

Test	Status	Time
COMP	Pass	1:33pm
CAL	Pass	1:33pm

Preventive Maintenance
Status: Pass

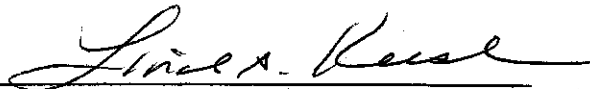

Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH

**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**

County Mecklenburg Instrument Location Mecklenburg County SD
Instrument Serial No. 008690 801 E 4th St, Charlotte

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 4th day of September, 2018 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.



M. J. [Signature]
Signature of Certifying Official

656
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

MECKLENBURG COUNTY SHERIFFS DEPARTMENT
590

Serial Number: 008690
Test Date: 09/04/2018

Citation Number: M0000000-0
Subject's Name:
PREVENTIVE, MAINTENANCE
Subject's Date of Birth: 11/11/1911
Subject's Sex: Male
Driver's License State: XX
Driver's License Number: NONE


Analyst's Name: HAYS, MARK D
Permit Number: 15924E
Effective:
01/01/2018-01/01/2020

Officer's Name: NONE, NONE
Type of Agency: FTA
Agency: DHHS
Test Type: Breath Test

Lot Number: AG807101
Exp Date: 03/12/2020

Test	g/210L	Time
DIAG	Pass	9:18am
AIR BLK	.00	9:19am
ACCY CHK	.08	9:20am
AIR BLK	.00	9:21am
SUB TEST	.00	9:22am
AIR BLK	.00	9:22am
SUB TEST	.00	9:24am
AIR BLK	.00	9:25am

Reported AC: .00 g/210L

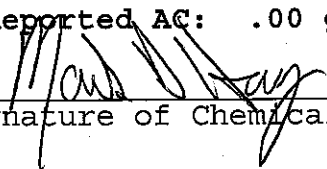

Signature of Chemical Analyst

Court CVR



Analyst

This form is used when performing Preventive Maintenance procedures
Forensic Tests for Alcohol Branch
Department of Health and Human Services
Rev. 12/2007

Intox EC/IR-II: Preventive Maintenance

MECKLENBURG COUNTY SHERIFFS DEPARTMENT 590

Serial Number: 008690 Test Record Number: 5989
Test Date: 09/04/2018 Test Time: 9:31am EDT

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	9:31am
FLO	Pass	9:31am
FC	Pass	9:32am

Temperature Tests

Test	Status	Time
FC1	Pass	9:32am
SRC	Pass	9:32am
DET	Pass	9:32am
BAR	Pass	9:32am
BT	Pass	9:32am

Blank Tests

Test	Status	Time
AIR	Pass	9:32am

Printer Tests

Test	Status	Time
PRNT	Pass	9:32am

CRC Tests

Test	Status	Time
COMP	Pass	9:32am
CAL	Pass	9:32am

Preventive Maintenance
Status: Pass



Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH

**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**

County Mecklenburg Instrument Location Mecklenburg County SD
Instrument Serial No. 008665 801 E 4th St., Charlotte

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 4th day of September, 2018 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.



[Signature]
Signature of Certifying Official

656
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

MECKLENBURG COUNTY SHERIFFS DEPARTMENT
590

Serial Number: 008665
Test Date: 09/04/2018

Citation Number: M0000000-0
Subject's Name:
PREVENTIVE, MAINTENANCE
Subject's Date of Birth: 11/11/1911
Subject's Sex: Male
Driver's License State: XX
Driver's License Number: NONE


Analyst's Name: HAYS, MARK D
Permit Number: 15924E
Effective:
01/01/2018-01/01/2020

Officer's Name: NONE, NONE
Type of Agency: FTA
Agency: DHHS
Test Type: Breath Test

Lot Number: AG807101
Exp Date: 03/12/2020

Test	g/210L	Time
DIAG	Pass	9:16am
AIR BLK	.00	9:17am
ACCY CHK	.07	9:18am
AIR BLK	.00	9:18am
SUB TEST	.00	9:19am
AIR BLK	.00	9:20am
SUB TEST	.00	9:22am
AIR BLK	.00	9:23am

Reported AC: .00 g/210L

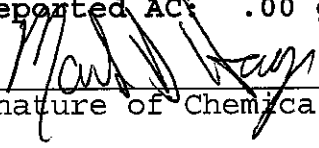

Signature of Chemical Analyst

Court CVR



Analyst

**This form is used when performing Preventive Maintenance procedures
Forensic Tests for Alcohol Branch
Department of Health and Human Services
Rev. 12/2007**

Intox EC/IR-II: Preventive Maintenance

MECKLENBURG COUNTY SHERIFFS DEPARTMENT 590

Serial Number: 008665 Test Record Number: 4580
Test Date: 09/04/2018 Test Time: 9:26am EDT

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	9:26am
FLO	Pass	9:26am
FC	Pass	9:27am

Temperature Tests

Test	Status	Time
FC1	Pass	9:27am
SRC	Pass	9:27am
DET	Pass	9:27am
BAR	Pass	9:27am
BT	Pass	9:27am

Blank Tests

Test	Status	Time
AIR	Pass	9:27am


Printer Tests

Test	Status	Time
PRNT	Pass	9:27am

CRC Tests

Test	Status	Time
COMP	Pass	9:27am
CAL	Pass	9:27am

Preventive Maintenance
Status: Pass

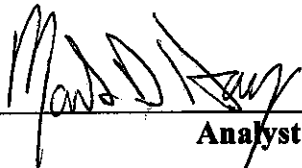

Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH

**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**

County Mecklenburg Instrument Location CMFD-LEC
Instrument Serial No. 008594 601 E. Trade St. Charlotte

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 4th day of September, 2018 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.



M. D. [Signature]
Signature of Certifying Official

656
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

MECKLENBURG COUNTY CMPD LEC 590

Serial Number: 008594
Test Date: 09/04/2018

Citation Number: M0000000-0
Subject's Name:
PREVENTIVE, MAINTENANCE
Subject's Date of Birth: 11/11/1911
Subject's Sex: Male
Driver's License State: XX
Driver's License Number: NONE


Analyst's Name: HAYS, MARK D
Permit Number: 15924E
Effective:
01/01/2018-01/01/2020

Officer's Name: NONE, NONE
Type of Agency: FTA
Agency: DHHS
Test Type: Breath Test

Lot Number: AG814902
Exp Date: 05/29/2020


Test	g/210L	Time
DIAG	Pass	10:18am
AIR BLK	.00	10:19am
ACCY CHK	.07	10:19am
AIR BLK	.00	10:20am
SUB TEST	.00	10:21am
AIR BLK	.00	10:22am
SUB TEST	.00	10:23am
AIR BLK	.00	10:24am

Reported AS: .00 g/210L

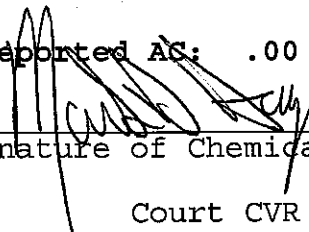

Signature of Chemical Analyst

Court CVR

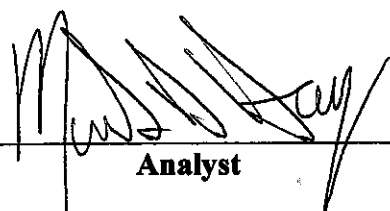

Analyst

**This form is used when performing Preventive Maintenance procedures
Forensic Tests for Alcohol Branch
Department of Health and Human Services
Rev. 12/2007**

Intox EC/IR-II: Preventive Maintenance

MECKLENBURG COUNTY CMPD LEC 590

Serial Number: 008594 Test Record Number: 4170
Test Date: 09/04/2018 Test Time: 10:27am EDT

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	10:28am
FLO	Pass	10:28am
FC	Pass	10:28am

Temperature Tests

Test	Status	Time
FC1	Pass	10:28am
SRC	Pass	10:28am
DET	Pass	10:28am
BAR	Pass	10:28am
BT	Pass	10:28am

Blank Tests

Test	Status	Time
AIR	Pass	10:28am


Printer Tests

Test	Status	Time
PRNT	Pass	10:28am

CRC Tests

Test	Status	Time
COMP	Pass	10:29am
CAL	Pass	10:29am

Preventive Maintenance
Status: Pass

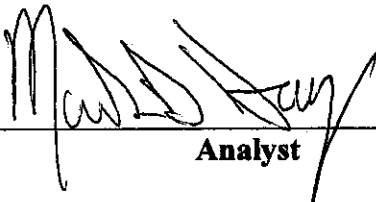

Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH


**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**

County Mecklenburg Instrument Location Pineville PD
Instrument Serial No. 008703 427 Main St., Pineville

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 12th day of September, 2018 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

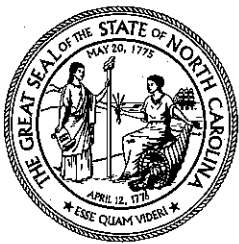

[Signature]
Signature of Certifying Official

656
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

MECKLENBURG COUNTY PINEVILLE PD 590

Serial Number: 008703
Test Date: 09/12/2018

Citation Number: M0000000-0
Subject's Name:

PREVENTIVE, MAINTENANCE

Subject's Date of Birth: 11/11/1911

Subject's Sex: Male

Driver's License State: XX

Driver's License Number: NONE

Analyst's Name: HAYS, MARK D

Permit Number: 15924E

Effective:

01/01/2018-01/01/2020

Officer's Name: NONE, NONE

Type of Agency: FTA

Agency: DHHS

Test Type: Breath Test

Lot Number: AG721401

Exp Date: 08/02/2019

Test	g/210L	Time
DIAG	Pass	11:01am
AIR BLK	.00	11:02am
ACCY CHK	.08	11:03am
AIR BLK	.00	11:04am
SUB TEST	.00	11:04am
AIR BLK	.00	11:05am
SUB TEST	.00	11:07am
AIR BLK	.00	11:08am

Reported AC: .00 g/210L



Signature of Chemical Analyst

Court CVR



Analyst

**This form is used when performing Preventive Maintenance procedures
Forensic Tests for Alcohol Branch
Department of Health and Human Services
Rev. 12/2007**

Intox EC/IR-II: Preventive Maintenance

MECKLENBURG COUNTY PINEVILLE PD 590

Serial Number: 008703 Test Record Number: 5668
Test Date: 09/12/2018 Test Time: 11:09am EDT

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	11:09am
FLO	Pass	11:09am
FC	Pass	11:09am

Temperature Tests

Test	Status	Time
FC1	Pass	11:09am
SRC	Pass	11:09am
DET	Pass	11:09am
BAR	Pass	11:09am
BT	Pass	11:09am

Blank Tests

Test	Status	Time
AIR	Pass	11:10am


Printer Tests

Test	Status	Time
PRNT	Pass	11:10am

CRC Tests

Test	Status	Time
COMP	Pass	11:10am
CAL	Pass	11:10am

Preventive Maintenance
Status: Pass

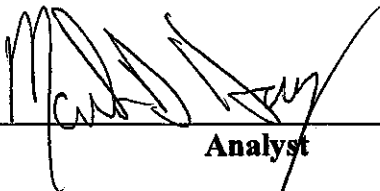

Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH

**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**

County Mitchell Instrument Location Service Pine PD
Instrument Serial No. 007726 Service Pine, NC

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 28 day of September, 2018 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.



[Signature]
Signature of Certifying Official

644
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

MITCHELL COUNTY SPRUCE PINE PD 600

Serial Number: 008726

Test Date: 09/28/2018

Citation Number: M0000000-0

Subject's Name:

PREVENTIVE, MAINTENANCE

Subject's Date of Birth: 11/11/1911

Subject's Sex: Male

Driver's License State: XX

Driver's License Number: NONE

Analyst's Name: BURNETTE, ANTHONY J

Permit Number: 11304E

Effective:

05/01/2017-05/01/2019

Officer's Name: NONE,

Type of Agency: FTA

Agency: DHHS

Test Type: Breath Test

Lot Number: AG716202


Exp Date: 06/11/2019

Test	g/210L	Time
DIAG	Pass	5:26pm
AIR BLK	.00	5:27pm
ACCY CHK	.08	5:27pm
AIR BLK	.00	5:28pm
SUB TEST	.00	5:29pm
AIR BLK	.00	5:30pm
SUB TEST	.00	5:31pm
AIR BLK	.00	5:32pm

Reported AC: .00 g/210L

Signature of Chemical Analyst

Court CVR



Analyst

Intox EC/IR-II: Preventive Maintenance

MITCHELL COUNTY SPRUCE PINE PD 600

Serial Number: 008726 Test Record Number: 946
Test Date: 09/28/2018 Test Time: 5:33pm EDT

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	5:34pm
FLO	Pass	5:34pm
FC	Pass	5:34pm

Temperature Tests

Test	Status	Time
FC1	Pass	5:34pm
SRC	Pass	5:34pm
DET	Pass	5:34pm
BAR	Pass	5:34pm
BT	Pass	5:34pm

Blank Tests

Test	Status	Time
AIR	Pass	5:35pm

Printer Tests

Test	Status	Time
PRNT	Pass	5:35pm

CRC Tests

Test	Status	Time
COMP	Pass	5:35pm
CAL	Pass	5:35pm

Preventive Maintenance
Status: Pass


Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH


**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**


County Montgomery Instrument Location BAT mobile unit 4
Instrument Serial No. 008717

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 2 day of September, 2018 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

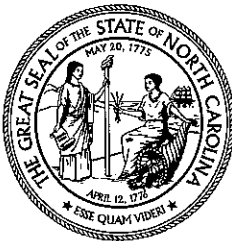

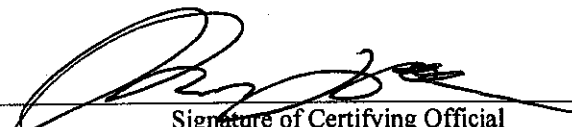
Signature of Certifying Official

864
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

MONTGOMERY COUNTY BAT MOBILE UNIT 4
610

Serial Number: 008717
Test Date: 09/02/2018

Citation Number: M0000000-0
Subject's Name:
PREVENTIVE, MAINTENANCE
Subject's Date of Birth: 11/11/1911
Subject's Sex: Male
Driver's License State: XX
Driver's License Number: NONE


Analyst's Name: TODD, SHANE C
Permit Number: 11391E
Effective:
07/27/2018-07/01/2020

Officer's Name: NONE, NONE
Type of Agency: FTA
Agency: DHHS
Test Type: Breath Test

Lot Number: AG821401
Exp Date: 08/02/2020

Test	g/210L	Time
DIAG	Pass	5:08pm
AIR BLK	.00	5:09pm
ACCY CHK	.07	5:09pm
AIR BLK	.00	5:10pm
SUB TEST	.00	5:11pm
AIR BLK	.00	5:12pm
SUB TEST	.00	5:13pm
AIR BLK	.00	5:14pm

Reported AC: .00 g/210L

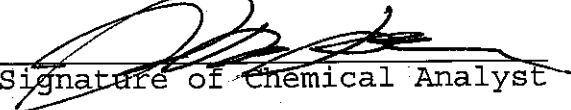
Signature of Chemical Analyst

Court CVR


Analyst

Intox EC/IR-II: Preventive Maintenance

MONTGOMERY COUNTY BAT MOBILE UNIT 4 610

Serial Number: 008717 Test Record Number: 612
Test Date: 09/02/2018 Test Time: 5:16pm EDT

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	5:17pm
FLO	Pass	5:17pm
FC	Pass	5:17pm

Temperature Tests

Test	Status	Time
FC1	Pass	5:17pm
SRC	Pass	5:17pm
DET	Pass	5:17pm
BAR	Pass	5:17pm
BT	Pass	5:17pm

Blank Tests

Test	Status	Time
AIR	Pass	5:18pm

Printer Tests

Test	Status	Time
PRNT	Pass	5:18pm

CRC Tests

Test	Status	Time
COMP	Pass	5:18pm
CAL	Pass	5:18pm

Preventive Maintenance
Status: Pass


Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH

PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II

County Montgomery Instrument Location BAT Mobile Unit 4
Instrument Serial No. 008734

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 2 day of September, 20 18 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.



[Signature]
Signature of Certifying Official

664
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

**MONTGOMERY COUNTY BAT MOBILE UNIT 4
610**

Serial Number: 008734
Test Date: 09/02/2018

Citation Number: M0000000-0
Subject's Name:
PREVENTIVE, MAINTENANCE
Subject's Date of Birth: 11/11/1911
Subject's Sex: Male
Driver's License State: XX
Driver's License Number: NONE

Analyst's Name: TODD, SHANE C
Permit Number: 11391E
Effective:
07/27/2018-07/01/2020

Officer's Name: NONE, NONE
Type of Agency: FTA
Agency: DHHS
Test Type: Breath Test

Lot Number: 702401
Exp Date: 08/02/2020

Test	g/210L	Time
DIAG	Pass	5:11pm
AIR BLK	.00	5:12pm
ACCY CHK	.08	5:13pm
AIR BLK	.00	5:14pm
SUB TEST	.00	5:15pm
AIR BLK	.00	5:16pm
SUB TEST	.00	5:17pm
AIR BLK	.00	5:18pm

Reported AC: .00 g/210L



Signature of Chemical Analyst

Court CVR



Analyst

Intox EC/IR-II: Preventive Maintenance

MONTGOMERY COUNTY BAT MOBILE UNIT 4 610

Serial Number: 008734 Test Record Number: 985
Test Date: 09/02/2018 Test Time: 5:20pm EDT

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	5:21pm
FLO	Pass	5:21pm
FC	Pass	5:21pm

Temperature Tests

Test	Status	Time
FC1	Pass	5:21pm
SRC	Pass	5:21pm
DET	Pass	5:21pm
BAR	Pass	5:21pm
BT	Pass	5:21pm

Blank Tests

Test	Status	Time
AIR	Pass	5:22pm

Printer Tests

Test	Status	Time
PRNT	Pass	5:22pm

CRC Tests

Test	Status	Time
COMP	Pass	5:22pm
CAL	Pass	5:22pm

Preventive Maintenance
Status: Pass



Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH

**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**


County MOORE Instrument Location MOORE CO JAIL

Instrument Serial No. 008735 CARTHAGE, NC

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 27 day of SEPTEMBER, 2018, the foregoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

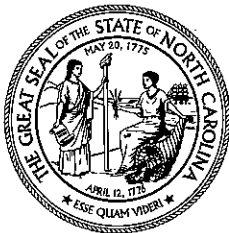

Alan P. Burns

Signature of Certifying Official

648

Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

MOORE COUNTY MOORE CO JAIL 620

Serial Number: 008735
Test Date: 09/27/2018

Citation Number: M0000000-0
Subject's Name:
PREVENTIVE, MAINTENANCE
Subject's Date of Birth: 11/11/1911
Subject's Sex: Male
Driver's License State: XX
Driver's License Number: NONE

Analyst's Name: BARNES, ALVIN R
Permit Number: 15671E
Effective:
07/01/2017-07/01/2019

Officer's Name: NONE, NONE
Type of Agency: FTA
Agency: DHHS
Test Type: Breath Test


Lot Number: AG805802
Exp Date: 02/27/2020

Test	g/210L	Time
DIAG	Pass	10:05am
AIR BLK	.00	10:05am
ACCY CHK	.08	10:06am
AIR BLK	.00	10:07am
SUB TEST	.00	10:07am
AIR BLK	.00	10:08am
SUB TEST	.00	10:10am
AIR BLK	.00	10:10am

Reported AC: .00 g/210L

Signature of Chemical Analyst

Court CVR

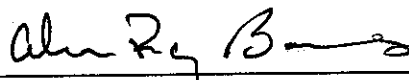

Analyst

Intox EC/IR-II: Preventive Maintenance

MOORE COUNTY MOORE CO JAIL 620

Serial Number: 008735 Test Record Number: 2153
Test Date: 09/27/2018 Test Time: 10:11am EDT

System Check: Passed

Baseline Tests

Test	Status	Time
IR	Pass	10:12am
FLO	Pass	10:12am
FC	Pass	10:12am

Temperature Tests

Test	Status	Time
FC1	Pass	10:12am
SRC	Pass	10:12am
DET	Pass	10:12am
BAR	Pass	10:12am
BT	Pass	10:12am

Blank Tests

Test	Status	Time
AIR	Pass	10:13am


Printer Tests

Test	Status	Time
PRNT	Pass	10:13am

CRC Tests

Test	Status	Time
COMP	Pass	10:13am
CAL	Pass	10:13am

Preventive Maintenance
Status: Pass

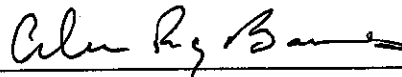

Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH

**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**

County MOORE Instrument Location PINEHURST POLICE DEPT
Instrument Serial No. 008710 PINEHURST, NC

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 27 day of SEPTEMBER, 2018, the foregoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.



Alan R. Bane
Signature of Certifying Official

648
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

MOORE COUNTY PINEHURST PD 620

Serial Number: 008710
Test Date: 09/27/2018

Citation Number: M0000000-0
Subject's Name:

PREVENTIVE, MAINTENANCE

Subject's Date of Birth: 11/11/1911

Subject's Sex: Male

Driver's License State: XX

Driver's License Number: NONE

Analyst's Name: BARNES, ALVIN R

Permit Number: 15671E

Effective:

07/01/2017-07/01/2019

Officer's Name: NONE, NONE

Type of Agency: FTA

Agency: DHHS

Test Type: Breath Test

Lot Number: AG702302


Exp Date: 01/23/2019

Test	g/210L	Time
DIAG	Pass	8:30am
AIR BLK	.00	8:31am
ACCY CHK	.07	8:31am
AIR BLK	.00	8:32am
SUB TEST	.00	8:33am
AIR BLK	.00	8:34am
SUB TEST	.00	8:35am
AIR BLK	.00	8:36am

Reported AC: .00 g/210L

Signature of Chemical Analyst

Court CVR

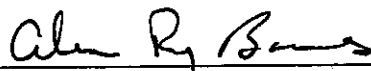
Analyst

Intox EC/IR-II: Preventive Maintenance

MOORE COUNTY PINEHURST PD 620

Serial Number: 008710 Test Record Number: 1536
Test Date: 09/27/2018 Test Time: 8:37am EDT

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	8:37am
FLO	Pass	8:37am
FC	Pass	8:37am

Temperature Tests

Test	Status	Time
FC1	Pass	8:38am
SRC	Pass	8:38am
DET	Pass	8:38am
BAR	Pass	8:38am
BT	Pass	8:38am

Blank Tests

Test	Status	Time
AIR	Pass	8:38am


Printer Tests

Test	Status	Time
PRNT	Pass	8:38am

CRC Tests

Test	Status	Time
COMP	Pass	8:38am
CAL	Pass	8:38am

Preventive Maintenance
Status: Pass

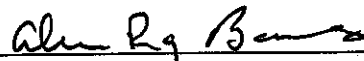

Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH


PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II

County MOORE Instrument Location SOUTHERN PINES PD
Instrument Serial No. 008720 SOUTHERN PINES, NC

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 27 day of SEPTEMBER, 2018, the foregoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

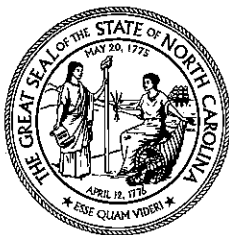

Alan R. Barnes

Signature of Certifying Official

648

Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

MOORE COUNTY SOUTHERN PINES PD 620

Serial Number: 008720
Test Date: 09/27/2018

Citation Number: M0000000-0
Subject's Name:
PREVENTIVE, MAINTENANCE
Subject's Date of Birth: 11/11/1911
Subject's Sex: Male
Driver's License State: XX
Driver's License Number: NONE

Analyst's Name: BARNES, ALVIN R
Permit Number: 15671E
Effective:
07/01/2017-07/01/2019

Officer's Name: NONE, NONE
Type of Agency: FTA
Agency: DHHS
Test Type: Breath Test


Lot Number: AG805802
Exp Date: 02/27/2020

Test	g/210L	Time
DIAG	Pass	7:39am
AIR BLK	.00	7:40am
ACCY CHK	.07	7:40am
AIR BLK	.00	7:41am
SUB TEST	.00	7:42am
AIR BLK	.00	7:43am
SUB TEST	.00	7:44am
AIR BLK	.00	7:45am

Reported AC: .00 g/210L

Signature of Chemical Analyst

Court CVR

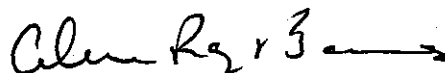

Analyst

This form is used when performing Preventive Maintenance procedures
Forensic Tests for Alcohol Branch
Department of Health and Human Services
Rev. 12/2007

Intox EC/IR-II: Preventive Maintenance

MOORE COUNTY SOUTHERN PINES PD 620

Serial Number: 008720 Test Record Number: 1073
Test Date: 09/27/2018 Test Time: 7:46am EDT

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	7:46am
FLO	Pass	7:46am
FC	Pass	7:46am

Temperature Tests

Test	Status	Time
FC1	Pass	7:46am
SRC	Pass	7:46am
DET	Pass	7:46am
BAR	Pass	7:46am
BT	Pass	7:46am

Blank Tests

Test	Status	Time
AIR	Pass	7:47am

Printer Tests

Test	Status	Time
PRNT	Pass	7:47am

CRC Tests

Test	Status	Time
COMP	Pass	7:47am
CAL	Pass	7:47am

Preventive Maintenance
Status: Pass



Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH

**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**


County New Hanover Instrument Location Mobik PAT unit 4

Instrument Serial No. 008871

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 1 day of September, 20 18 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

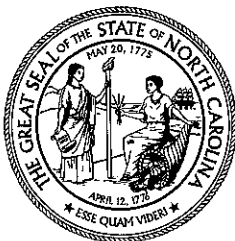

[Signature]
Signature of Certifying Official

664
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

NEW HANOVER COUNTY BAT MOBILE UNIT 4
640

Serial Number: 008871
Test Date: 09/01/2018

Citation Number: M0000000-0
Subject's Name:

PREVENTIVE, MAINTENANCE

Subject's Date of Birth: 11/11/1911

Subject's Sex: Male

Driver's License State: XX

Driver's License Number: NONE

Analyst's Name: TODD, SHANE C

Permit Number: 11391E

Effective:

07/27/2018-07/01/2020

Officer's Name: NONE, NONE

Type of Agency: FTA

Agency: DHHS


Test Type: Breath Test

Lot Number: AG821401

Exp Date: 08/02/2020

Test	g/210L	Time
DIAG	Pass	8:25pm
AIR BLK	.00	8:26pm
ACCY CHK	.08	8:26pm
AIR BLK	.00	8:27pm
SUB TEST	.00	8:28pm
AIR BLK	.00	8:29pm
SUB TEST	.00	8:31pm
AIR BLK	.00	8:32pm

Reported AC: .00 g/210L

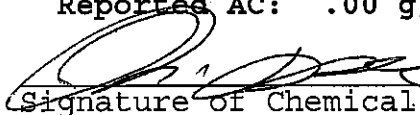
Signature of Chemical Analyst

Court CVR


Analyst

This form is used when performing Preventive Maintenance procedures
Forensic Tests for Alcohol Branch
Department of Health and Human Services
Rev. 12/2007

Intox EC/IR-II: Preventive Maintenance

NEW HANOVER COUNTY BAT MOBILE UNIT 4 640

Serial Number: 008871 Test Record Number: 927
Test Date: 09/01/2018 Test Time: 8:33pm EDT

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	8:34pm
FLO	Pass	8:34pm
FC	Pass	8:34pm

Temperature Tests

Test	Status	Time
FC1	Pass	8:34pm
SRC	Pass	8:34pm
DET	Pass	8:34pm
BAR	Pass	8:34pm
BT	Pass	8:34pm

Blank Tests

Test	Status	Time
AIR	Pass	8:34pm

Printer Tests

Test	Status	Time
PRNT	Pass	8:34pm

CRC Tests

Test	Status	Time
COMP	Pass	8:35pm
CAL	Pass	8:35pm

Preventive Maintenance
Status: Pass


Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH

**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**


County New Hanover Instrument Location BAT mobile unit 4

Instrument Serial No. 008734

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 1 day of September, 2018 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

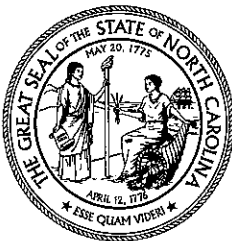


Signature of Certifying Official

664
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

**NEW HANOVER COUNTY BAT MOBILE UNIT 4
640**

Serial Number: 008734
Test Date: 09/01/2018

Citation Number: M0000000-0
Subject's Name:

PREVENTIVE, MAINTENANCE

Subject's Date of Birth: 11/11/1911

Subject's Sex: Male

Driver's License State: XX

Driver's License Number: NONE

Analyst's Name: TODD, SHANE C

Permit Number: 11391E

Effective:

07/27/2018-07/01/2020

Officer's Name: NONE, NONE

Type of Agency: FTA

Agency: DHHS


Test Type: Breath Test

Lot Number: 702401

Exp Date: 08/02/2020

Test	g/210L	Time
DIAG	Pass	8:28pm
AIR BLK	.00	8:29pm
ACCY CHK	.07	8:30pm
AIR BLK	.00	8:30pm
SUB TEST	.00	8:31pm
AIR BLK	.00	8:32pm
SUB TEST	.00	8:34pm
AIR BLK	.00	8:35pm

Reported AC: .00 g/210L



Signature of Chemical Analyst

Court CVR



Analyst

**This form is used when performing Preventive Maintenance procedures
Forensic Tests for Alcohol Branch
Department of Health and Human Services
Rev. 12/2007**

Intox EC/IR-II: Preventive Maintenance

NEW HANOVER COUNTY BAT MOBILE UNIT 4 640

Serial Number: 008734 Test Record Number: 982
Test Date: 09/01/2018 Test Time: 8:37pm EDT

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	8:37pm
FLO	Pass	8:37pm
FC	Pass	8:37pm

Temperature Tests

Test	Status	Time
FC1	Pass	8:37pm
SRC	Pass	8:37pm
DET	Pass	8:37pm
BAR	Pass	8:37pm
BT	Pass	8:37pm

Blank Tests

Test	Status	Time
AIR	Pass	8:38pm

Printer Tests

Test	Status	Time
PRNT	Pass	8:38pm

CRC Tests

Test	Status	Time
COMP	Pass	8:38pm
CAL	Pass	8:38pm

Preventive Maintenance
Status: Pass



Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH

**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**


County New Hanover Instrument Location BAT Mobile Unit 4

Instrument Serial No. 008717

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 1 day of September, 2018 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

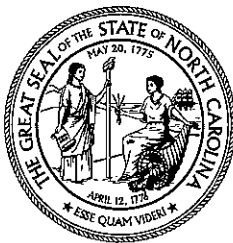

[Signature]
Signature of Certifying Official

664
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

NEW HANOVER COUNTY BAT MOBILE UNIT 4
640

Serial Number: 008717
Test Date: 09/01/2018

Citation Number: M0000000-0
Subject's Name:
PREVENTIVE, MAINTENANCE
Subject's Date of Birth: 11/11/1911
Subject's Sex: Male
Driver's License State: XX
Driver's License Number: NONE


Analyst's Name: TODD, SHANE C
Permit Number: 11391E
Effective:
07/27/2018-07/01/2020

Officer's Name: NONE, NONE
Type of Agency: FTA
Agency: DHHS
Test Type: Breath Test

Lot Number: AG821401
Exp Date: 08/02/2020

Test	g/210L	Time
DIAG	Pass	10:22pm
AIR BLK	.00	10:23pm
ACCY CHK	.08	10:24pm
AIR BLK	.00	10:25pm
SUB TEST	.00	10:25pm
AIR BLK	.00	10:26pm
SUB TEST	.00	10:28pm
AIR BLK	.00	10:29pm

Reported AC: .00 g/210L



Signature of Chemical Analyst

Court CVR



Analyst

Intox EC/IR-II: Preventive Maintenance

NEW HANOVER COUNTY BAT MOBILE UNIT 4 640

Serial Number: 008717 Test Record Number: 606
Test Date: 09/01/2018 Test Time: 10:42pm EDT

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	10:42pm
FLO	Pass	10:42pm
FC	Pass	10:42pm

Temperature Tests

Test	Status	Time
FC1	Pass	10:43pm
SRC	Pass	10:43pm
DET	Pass	10:43pm
BAR	Pass	10:43pm
BT	Pass	10:43pm

Blank Tests

Test	Status	Time
AIR	Pass	10:43pm


Printer Tests

Test	Status	Time
PRNT	Pass	10:43pm

CRC Tests

Test	Status	Time
COMP	Pass	10:43pm
CAL	Pass	10:43pm

Preventive Maintenance
Status: Pass

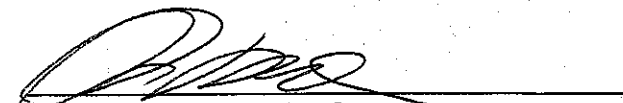
Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH


**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**

County Onslow Instrument Location Mobile BAT unit 4
Instrument Serial No. 008717

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 7 day of September, 2018 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

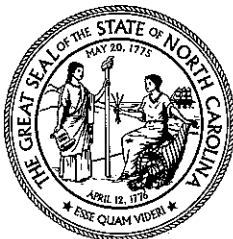

[Signature]
Signature of Certifying Official

664
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

ONslow COUNTY BAT MOBILE UNIT 4 660

Serial Number: 008717
Test Date: 09/07/2018

Citation Number: M0000000-0
Subject's Name:
PREVENTIVE, MAINTENANCE
Subject's Date of Birth: 11/11/1911
Subject's Sex: Male
Driver's License State: XX
Driver's License Number: NONE


Analyst's Name: TODD, SHANE C
Permit Number: 11391E
Effective:
07/27/2018-07/01/2020

Officer's Name: NONE, NONE
Type of Agency: FTA
Agency: DHHS
Test Type: Breath Test


Lot Number: AG821401
Exp Date: 08/02/2020

Test	g/210L	Time
DIAG	Pass	10:33pm
AIR BLK	.00	10:34pm
ACCY CHK	.08	10:35pm
AIR BLK	.00	10:36pm
SUB TEST	.00	10:36pm
AIR BLK	.00	10:37pm
SUB TEST	.00	10:39pm
AIR BLK	.00	10:40pm

Reported AC: .00 g/210L

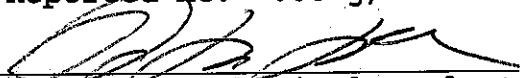
Signature of Chemical Analyst

Court CVR

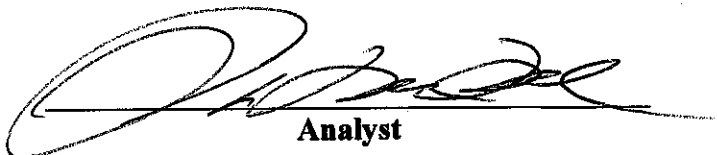
Analyst

Intox EC/IR-II: Preventive Maintenance

ONslow COUNTY BAT MOBILE UNIT 4 660

Serial Number: 008717 Test Record Number: 616
Test Date: 09/07/2018 Test Time: 10:45pm EDT

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	10:46pm
FLO	Pass	10:46pm
FC	Pass	10:46pm

Temperature Tests

Test	Status	Time
FC1	Pass	10:46pm
SRC	Pass	10:46pm
DET	Pass	10:46pm
BAR	Pass	10:46pm
BT	Pass	10:46pm

Blank Tests

Test	Status	Time
AIR	Pass	10:47pm

Printer Tests

Test	Status	Time
PRNT	Pass	10:47pm

CRC Tests

Test	Status	Time
COMP	Pass	10:47pm
CAL	Pass	10:47pm

Preventive Maintenance
Status: Pass


Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH

**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**


County Onslow Instrument Location Mobile BAT Unit 4

Instrument Serial No. 008734

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 7 day of September, 2018 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

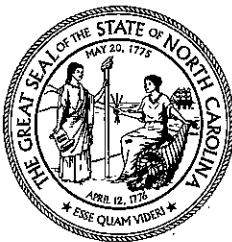

[Signature]
Signature of Certifying Official

664
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

ONslow COUNTY BAT MOBILE UNIT 4 660

Serial Number: 008734
Test Date: 09/07/2018

Citation Number: M0000000-0
Subject's Name:

PREVENTIVE, MAINTENANCE

Subject's Date of Birth: 11/11/1911

Subject's Sex: Male

Driver's License State: XX

Driver's License Number: NONE

Analyst's Name: TODD, SHANE C

Permit Number: 11391E

Effective:

07/27/2018-07/01/2020

Officer's Name: NONE, NONE

Type of Agency: FTA

Agency: DHHS


Test Type: Breath Test

Lot Number: 702401

Exp Date: 08/02/2020

Test	g/210L	Time
DIAG	Pass	9:36pm
AIR BLK	.00	9:37pm
ACCY CHK	.08	9:37pm
AIR BLK	.00	9:38pm
SUB TEST	.00	9:39pm
AIR BLK	.00	9:40pm
SUB TEST	.00	9:41pm
AIR BLK	.00	9:42pm

Reported AC: .00 g/210L

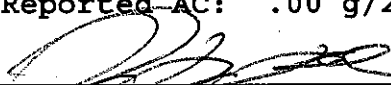
Signature of Chemical Analyst

Court CVR


Analyst

This form is used when performing Preventive Maintenance procedures
Forensic Tests for Alcohol Branch
Department of Health and Human Services
Rev. 12/2007

Intox EC/IR-II: Preventive Maintenance

ONslow COUNTY BAT MOBILE UNIT 4 660

Serial Number: 008734 Test Record Number: 989
Test Date: 09/07/2018 Test Time: 9:44pm EDT

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	9:44pm
FLO	Pass	9:44pm
FC	Pass	9:44pm

Temperature Tests

Test	Status	Time
FC1	Pass	9:44pm
SRC	Pass	9:44pm
DET	Pass	9:44pm
BAR	Pass	9:44pm
BT	Pass	9:44pm

Blank Tests

Test	Status	Time
AIR	Pass	9:45pm


Printer Tests

Test	Status	Time
PRNT	Pass	9:45pm

CRC Tests

Test	Status	Time
COMP	Pass	9:45pm
CAL	Pass	9:45pm

Preventive Maintenance
Status: Pass

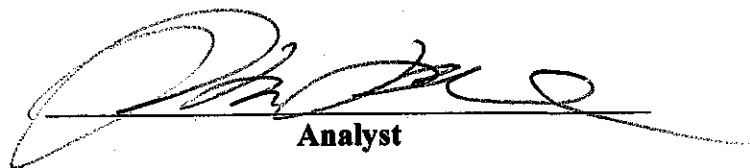

Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH


PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II

County PERQUIMANS Instrument Location Perquimans Co. S.O.
Instrument Serial No. 00 8851 110 Church St., Hartford, N.C.

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 5th day of SEPTEMBER, 2018 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

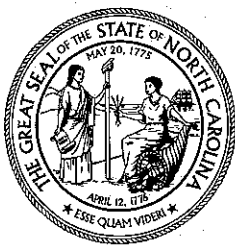

J. Neal
Signature of Certifying Official

647
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

PERQUIMANS COUNTY PERQUIMANS CO SO 710

Serial Number: 008851
Test Date: 09/05/2018

Citation Number: M0000000-0
Subject's Name:
PREVENTIVE, MAINTENANCE
Subject's Date of Birth: 11/11/1911
Subject's Sex: Male
Driver's License State: XX
Driver's License Number: NONE

Analyst's Name: KEESLER, LINDA A
Permit Number: 11646E
Effective:
08/01/2017-08/01/2019

Officer's Name: NONE, NONE
Type of Agency: FTA
Agency: DHHS
Test Type: Breath Test

Lot Number: AG734101
Exp Date: 12/07/2019

Test	g/210L	Time
DIAG	Pass	9:39am
AIR BLK	.00	9:40am
ACCY CHK	.08	9:40am
AIR BLK	.00	9:41am
SUB TEST	.00	9:42am
AIR BLK	.00	9:43am
SUB TEST	.00	9:44am
AIR BLK	.00	9:45am

Reported AC: .00 g/210L



Signature of Chemical Analyst

Court CVR



Analyst

Intox EC/IR-II: Preventive Maintenance

PERQUIMANS COUNTY PERQUIMANS CO SO 710

Serial Number: 008851 Test Record Number: 621
Test Date: 09/05/2018 Test Time: 9:46am EDT

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	9:47am
FLO	Pass	9:47am
FC	Pass	9:47am

Temperature Tests

Test	Status	Time
FC1	Pass	9:47am
SRC	Pass	9:47am
DET	Pass	9:47am
BAR	Pass	9:47am
BT	Pass	9:47am

Blank Tests

Test	Status	Time
AIR	Pass	9:47am

Printer Tests

Test	Status	Time
PRNT	Pass	9:47am

CRC Tests

Test	Status	Time
COMP	Pass	9:48am
CAL	Pass	9:48am

Preventive Maintenance
Status: Pass


Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH

**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**

County Perquimans Instrument Location Perquimans Co. S.O.
Instrument Serial No. 008921 110 Church St., Hertford, N.C.

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 26th day of September, 20 18 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.



Key A. [Signature]
Signature of Certifying Official

643
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

PERQUIMANS COUNTY PERQUIMANS CO SO 710

Serial Number: 008921
Test Date: 09/26/2018

Citation Number: M0000000-0
Subject's Name:
PREVENTIVE, MAINTENANCE
Subject's Date of Birth: 11/11/1911
Subject's Sex: Male
Driver's License State: XX
Driver's License Number: NONE


Analyst's Name: GUARD, KELLY G
Permit Number: 12955E
Effective:
06/01/2017-06/01/2019

Officer's Name: NONE, NONE
Type of Agency: FTA
Agency: DHHS
Test Type: Breath Test


Lot Number: AG734101
Exp Date: 12/07/2019

Test	g/210L	Time
DIAG	Pass	1:23pm
AIR BLK	.00	1:23pm
ACCY CHK	.08	1:24pm
AIR BLK	.00	1:25pm
SUB TEST	.00	1:25pm
AIR BLK	.00	1:26pm
SUB TEST	.00	1:28pm
AIR BLK	.00	1:29pm

Reported AC: .00 g/210L

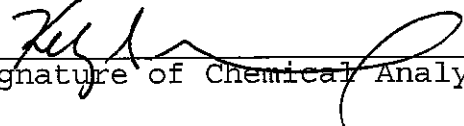

Signature of Chemical Analyst
Court CVR

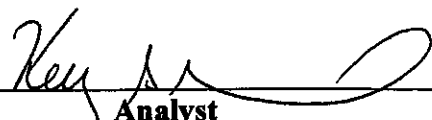

Analyst

Intox EC/IR-II: Preventive Maintenance

PERQUIMANS COUNTY PERQUIMANS CO SO 710

Serial Number: 008921 Test Record Number: 801
Test Date: 09/26/2018 Test Time: 1:32pm EDT

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	1:33pm
FLO	Pass	1:33pm
FC	Pass	1:33pm

Temperature Tests

Test	Status	Time
FC1	Pass	1:33pm
SRC	Pass	1:33pm
DET	Pass	1:33pm
BAR	Pass	1:33pm
BT	Pass	1:33pm

Blank Tests

Test	Status	Time
AIR	Pass	1:34pm


Printer Tests

Test	Status	Time
PRNT	Pass	1:34pm

CRC Tests

Test	Status	Time
COMP	Pass	1:34pm
CAL	Pass	1:34pm

Preventive Maintenance
Status: Pass



Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH

**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**

County PERSON Instrument Location PERSON CO. LEC

Instrument Serial No. 008880 120 COURT ST. ROXBORO, NC

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 05 day of SEPTEMBER, 20 18 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.



Bruce D. Smith
Signature of Certifying Official

637
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

PERSON COUNTY PERSON CO. LEC 720

Serial Number: 008880
Test Date: 09/05/2018

Citation Number: M0000000-0
Subject's Name:
PREVENTIVE, MAINTENANCE
Subject's Date of Birth: 11/11/1911
Subject's Sex: Male
Driver's License State: XX
Driver's License Number: NONE

Analyst's Name: SMITH, BRIAN D
Permit Number: 08937E
Effective:
08/01/2017-08/01/2019

Officer's Name: NONE, NONE
Type of Agency: FTA
Agency: DHHS
Test Type: Breath Test

Lot Number: AG702302
Exp Date: 01/23/2019

Test	g/210L	Time
DIAG	Pass	11:05am
AIR BLK	.00	11:06am
ACCY CHK	.08	11:06am
AIR BLK	.00	11:07am
SUB TEST	.00	11:08am
AIR BLK	.00	11:09am
SUB TEST	.00	11:10am
AIR BLK	.00	11:11am

Reported AC: .00 g/210L


Signature of Chemical Analyst

Court CVR


Analyst

**This form is used when performing Preventive Maintenance procedures
Forensic Tests for Alcohol Branch
Department of Health and Human Services
Rev. 12/2007**

Intox EC/IR-II: Preventive Maintenance

PERSON COUNTY PERSON CO. LEC 720

Serial Number: 008880 Test Record Number: 1360
Test Date: 09/05/2018 Test Time: 11:12am EDT

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	11:12am
FLO	Pass	11:12am
FC	Pass	11:13am

Temperature Tests

Test	Status	Time
FC1	Pass	11:13am
SRC	Pass	11:13am
DET	Pass	11:13am
BAR	Pass	11:13am
BT	Pass	11:13am

Blank Tests

Test	Status	Time
AIR	Pass	11:13am


Printer Tests

Test	Status	Time
PRNT	Pass	11:13am

CRC Tests

Test	Status	Time
COMP	Pass	11:14am
CAL	Pass	11:14am

Preventive Maintenance
Status: Pass

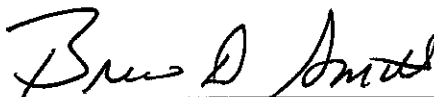

Analyst

**This form is used when performing Preventive Maintenance procedures
Forensic Tests for Alcohol Branch
Department of Health and Human Services
Rev. 12/2007**

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH

**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**

County Polk Instrument Location Polk County LEC

Instrument Serial No. 008927 880 East NC 108, Columbus

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 5th day of September, 20 18 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.



[Signature]
Signature of Certifying Official

656
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

POLK COUNTY POLK COUNTY SD 740

Serial Number: 008927
Test Date: 09/05/2018

Citation Number: M0000000-0
Subject's Name:
PREVENTIVE, MAINTENANCE
Subject's Date of Birth: 11/11/1911
Subject's Sex: Male
Driver's License State: XX
Driver's License Number: NONE

Analyst's Name: HAYS, MARK D
Permit Number: 15924E
Effective:
01/01/2018-01/01/2020

Officer's Name: NONE, NONE
Type of Agency: FTA
Agency: DHHS
Test Type: Breath Test

Lot Number: AG821401
Exp Date: 08/02/2020

Test	g/210L	Time
DIAG	Pass	12:33pm
AIR BLK	.00	12:33pm
ACCY CHK	.08	12:34pm
AIR BLK	.00	12:35pm
SUB TEST	.00	12:36pm
AIR BLK	.00	12:36pm
SUB TEST	.00	12:38pm
AIR BLK	.00	12:39pm

Reported AC: .00 g/210L



Signature of Chemical Analyst

Court CVR



Analyst

**This form is used when performing Preventive Maintenance procedures
Forensic Tests for Alcohol Branch
Department of Health and Human Services
Rev. 12/2007**

Intox EC/IR-II: Preventive Maintenance

POLK COUNTY POLK COUNTY SD 740

Serial Number: 008927 Test Record Number: 539
Test Date: 09/05/2018 Test Time: 12:40pm EDT

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	12:40pm
FLO	Pass	12:40pm
FC	Pass	12:40pm

Temperature Tests

Test	Status	Time
FC1	Pass	12:40pm
SRC	Pass	12:40pm
DET	Pass	12:40pm
BAR	Pass	12:40pm
BT	Pass	12:40pm

Blank Tests

Test	Status	Time
AIR	Pass	12:41pm


Printer Tests

Test	Status	Time
PRNT	Pass	12:41pm

CRC Tests

Test	Status	Time
COMP	Pass	12:41pm
CAL	Pass	12:41pm

Preventive Maintenance
Status: Pass

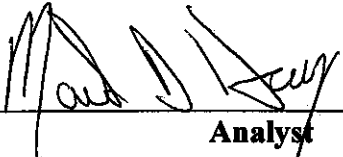

Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH

**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**

County RANDOLPH Instrument Location Archdate
Instrument Serial No. 008791 Police Dept.

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 24 day of September, 2018, the foregoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.



L. Kevin Dean
Signature of Certifying Official

642
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

RANDOLPH COUNTY ARCHDALE PD 750

Serial Number: 008791
Test Date: 09/24/2018

Citation Number: M0000000-0

Subject's Name:

PREVENTIVE, MAINTENANCE

Subject's Date of Birth: 11/11/1911

Subject's Sex: Male

Driver's License State: XX

Driver's License Number: NONE

Analyst's Name: DEAN, L

Permit Number: 11598E

Effective:

04/01/2017-04/01/2019

Officer's Name: NONE, NONE

Type of Agency: FTA

Agency: DHHS


Test Type: Breath Test

Lot Number: AG807101


Exp Date: 03/12/2020

Test	g/210L	Time
DIAG	Pass	12:15pm
AIR BLK	.00	12:15pm
ACCY CHK	.08	12:16pm
AIR BLK	.00	12:17pm
SUB TEST	.00	12:17pm
AIR BLK	.00	12:18pm
SUB TEST	.00	12:20pm
AIR BLK	.00	12:21pm

Reported AC: .00 g/210L


Signature of Chemical Analyst

Court CVR

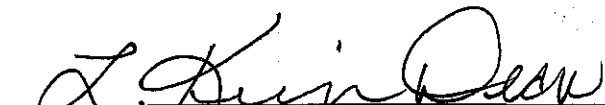
Analyst

**This form is used when performing Preventive Maintenance procedures
Forensic Tests for Alcohol Branch
Department of Health and Human Services
Rev. 12/2007**

Intox EC/IR-II: Preventive Maintenance

RANDOLPH COUNTY ARCHDALE PD 750

Serial Number: 008791 Test Record Number: 1279
Test Date: 09/24/2018 Test Time: 12:23pm EDT

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	12:23pm
FLO	Pass	12:23pm
FC	Pass	12:23pm

Temperature Tests

Test	Status	Time
FC1	Pass	12:24pm
SRC	Pass	12:24pm
DET	Pass	12:24pm
BAR	Pass	12:24pm
BT	Pass	12:24pm

Blank Tests

Test	Status	Time
AIR	Pass	12:24pm


Printer Tests

Test	Status	Time
PRNT	Pass	12:24pm

CRC Tests

Test	Status	Time
COMP	Pass	12:24pm
CAL	Pass	12:24pm

Preventive Maintenance
Status: Pass

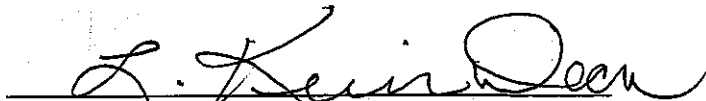

Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH

**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**

County RANDOLPH


Instrument Location LIBERTY POLICE DEPT

Instrument Serial No. 008830 LIBERTY, NC

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 27 day of SEPTEMBER, 2018, the foregoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

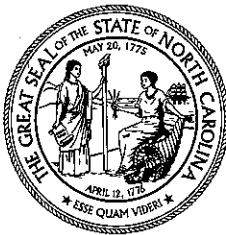

Alan R. Barnes

Signature of Certifying Official

648

Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

RANDOLPH COUNTY LIBERTY PD 750

Serial Number: 008830

Test Date: 09/27/2018

Citation Number: M0000000-0

Subject's Name:

PREVENTIVE, MAINTENANCE

Subject's Date of Birth: 11/11/1911

Subject's Sex: Male

Driver's License State: XX

Driver's License Number: NONE

Analyst's Name: BARNES, ALVIN R

Permit Number: 15671E

Effective:

07/01/2017-07/01/2019

Officer's Name: NONE, NONE

Type of Agency: FTA

Agency: DHHS

Test Type: Breath Test

Lot Number: AG805802

Exp Date: 02/27/2020

Test	g/210L	Time
DIAG	Pass	4:24pm
AIR BLK	.00	4:25pm
ACCY CHK	.08	4:25pm
AIR BLK	.00	4:26pm
SUB TEST	.00	4:26pm
AIR BLK	.00	4:27pm
SUB TEST	.00	4:29pm
AIR BLK	.00	4:30pm

Reported AC: .00 g/210L

Signature of Chemical Analyst

Court CVR



Analyst

Intox EC/IR-II: Preventive Maintenance

RANDOLPH COUNTY LIBERTY PD 750

Serial Number: 008830 Test Record Number: 598
Test Date: 09/27/2018 Test Time: 4:31pm EDT

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	4:31pm
FLO	Pass	4:31pm
FC	Pass	4:31pm

Temperature Tests

Test	Status	Time
FC1	Pass	4:31pm
SRC	Pass	4:31pm
DET	Pass	4:31pm
BAR	Pass	4:31pm
BT	Pass	4:31pm

Blank Tests

Test	Status	Time
AIR	Pass	4:32pm


Printer Tests

Test	Status	Time
PRNT	Pass	4:32pm

CRC Tests

Test	Status	Time
COMP	Pass	4:32pm
CAL	Pass	4:32pm

Preventive Maintenance
Status: Pass

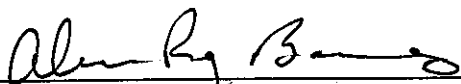

Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH

PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II


County RANDOLPH Instrument Location RANDLEMAN POLICE DEPT.

Instrument Serial No. 008737 RANDLEMAN, NC

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 27 day of SEPTEMBER, 2018, the foregoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

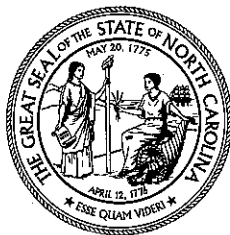

Alan Rg Bams
Signature of Certifying Official

648
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

RANDOLPH COUNTY RANDLEMAN PD 750

Serial Number: 008737
Test Date: 09/27/2018

Citation Number: M0000000-0
Subject's Name:
PREVENTIVE, MAINTENANCE
Subject's Date of Birth: 11/11/1911
Subject's Sex: Male
Driver's License State: XX
Driver's License Number: NONE

Analyst's Name: BARNES, ALVIN R
Permit Number: 15671E
Effective:
07/01/2017-07/01/2019

Officer's Name: NONE, NONE
Type of Agency: FTA
Agency: DHHS
Test Type: Breath Test

Lot Number: AG805802
Exp Date: 02/27/2020

Test	g/210L	Time
DIAG	Pass	2:06pm
AIR BLK	.00	2:06pm
ACCY CHK	.08	2:07pm
AIR BLK	.00	2:08pm
SUB TEST	.00	2:08pm
AIR BLK	.00	2:09pm
SUB TEST	.00	2:11pm
AIR BLK	.00	2:12pm

Reported AC: .00 g/210L

Signature of Chemical Analyst

Court CVR



Analyst

Intox EC/IR-II: Preventive Maintenance

RANDOLPH COUNTY RANDLEMAN PD 750

Serial Number: 008737 Test Record Number: 1055
Test Date: 09/27/2018 Test Time: 2:12pm EDT

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	2:13pm
FLO	Pass	2:13pm
FC	Pass	2:13pm

Temperature Tests

Test	Status	Time
FC1	Pass	2:13pm
SRC	Pass	2:13pm
DET	Pass	2:13pm
BAR	Pass	2:13pm
BT	Pass	2:13pm

Blank Tests

Test	Status	Time
AIR	Pass	2:14pm

Printer Tests

Test	Status	Time
PRNT	Pass	2:14pm

CRC Tests

Test	Status	Time
COMP	Pass	2:14pm
CAL	Pass	2:14pm

Preventive Maintenance
Status: Pass



Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH


**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**

County Richmond Co. Instrument Location Richmond Co. Magistrate office
Instrument Serial No. 008101 Rockingham, NC

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 26th day of September, 2018 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

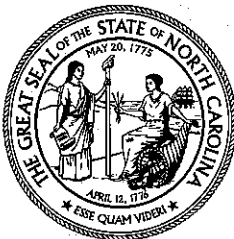

[Signature]
Signature of Certifying Official

654
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

RICHMOND COUNTY RICHMOND CO. MAG OFF
760

Serial Number: 008701
Test Date: 09/26/2018

Citation Number: M0000000-0
Subject's Name:

PREVENTIVE, MAINTENANCE

Subject's Date of Birth: 11/11/1911

Subject's Sex: Male

Driver's License State: XX

Driver's License Number: NONE

Analyst's Name: KEESLER, GRAYHAM C

Permit Number: 07682E

Effective:

12/01/2017-12/01/2019

Officer's Name: NONE, NONE

Type of Agency: FTA

Agency: DHHS

Test Type: Breath Test

Lot Number: AG805802

Exp Date: 02/27/2020

Test	g/210L	Time
DIAG	Pass	3:10pm
AIR BLK	.00	3:11pm
ACCY CHK	.08	3:11pm
AIR BLK	.00	3:12pm
SUB TEST	.00	3:13pm
AIR BLK	.00	3:14pm
SUB TEST	.00	3:16pm
AIR BLK	.00	3:17pm

Reported AC: .00 g/210L

Signature of Chemical Analyst

Court CVR

Analyst

**This form is used when performing Preventive Maintenance procedures
Forensic Tests for Alcohol Branch
Department of Health and Human Services
Rev. 12/2007**

Intox EC/IR-II: Preventive Maintenance

RICHMOND COUNTY RICHMOND CO. MAG OFF 760

Serial Number: 008701 Test Record Number: 1171
Test Date: 09/26/2018 Test Time: 3:18pm EDT

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	3:18pm
FLO	Pass	3:18pm
FC	Pass	3:18pm

Temperature Tests

Test	Status	Time
FC1	Pass	3:18pm
SRC	Pass	3:18pm
DET	Pass	3:18pm
BAR	Pass	3:18pm
BT	Pass	3:18pm

Blank Tests

Test	Status	Time
AIR	Pass	3:19pm


Printer Tests

Test	Status	Time
PRNT	Pass	3:19pm

CRC Tests

Test	Status	Time
COMP	Pass	3:19pm
CAL	Pass	3:19pm

Preventive Maintenance
Status: Pass

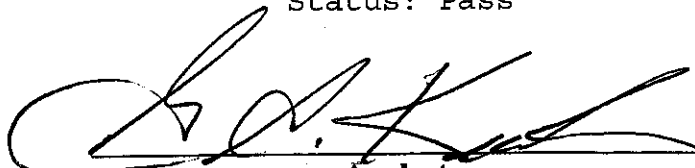
Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH


**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**

County Richmond Co. Instrument Location Richmond Co. Magistrate Office
Instrument Serial No. 008840 Rockingham, NC

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 26th day of September, 2018 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

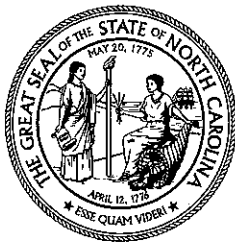

[Signature]
Signature of Certifying Official

654
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

RICHMOND COUNTY MAGISTRATE'S OFFICE
760

Serial Number: 008840
Test Date: 09/26/2018

Citation Number: M0000000-0
Subject's Name:
PREVENTIVE, MAINTENANCE
Subject's Date of Birth: 11/11/1911
Subject's Sex: Male
Driver's License State: XX
Driver's License Number: NONE

Analyst's Name: KEESLER, GRAYHAM C
Permit Number: 07682E
Effective:
12/01/2017-12/01/2019

Officer's Name: NONE, NONE
Type of Agency: FTA
Agency: DHHS
Test Type: Breath Test

Lot Number: AG805802
Exp Date: 02/27/2020

Test	g/210L	Time
DIAG	Pass	3:17pm
AIR BLK	.00	3:18pm
ACCY CHK	.07	3:19pm
AIR BLK	.00	3:20pm
SUB TEST	.00	3:21pm
AIR BLK	.00	3:22pm
SUB TEST	.00	3:23pm
AIR BLK	.00	3:24pm

Reported AG: .00 g/210L


Signature of Chemical Analyst

Court CVR


Analyst

**This form is used when performing Preventive Maintenance procedures
Forensic Tests for Alcohol Branch
Department of Health and Human Services
Rev. 12/2007**

Intox EC/IR-II: Preventive Maintenance

RICHMOND COUNTY MAGISTRATE'S OFFICE 760

Serial Number: 008840 Test Record Number: 2174
Test Date: 09/26/2018 Test Time: 3:26pm EDT

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	3:26pm
FLO	Pass	3:26pm
FC	Pass	3:26pm

Temperature Tests

Test	Status	Time
FC1	Pass	3:26pm
SRC	Pass	3:26pm
DET	Pass	3:26pm
BAR	Pass	3:26pm
BT	Pass	3:26pm

Blank Tests

Test	Status	Time
AIR	Pass	3:27pm


Printer Tests

Test	Status	Time
PRNT	Pass	3:27pm

CRC Tests

Test	Status	Time
COMP	Pass	3:27pm
CAL	Pass	3:27pm

Preventive Maintenance
Status: Pass

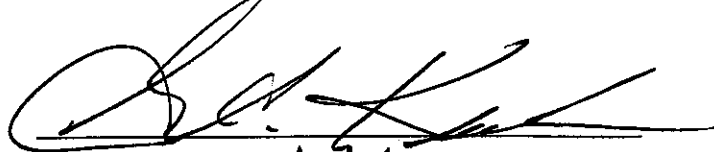

Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH

**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**


County Robeson Instrument Location Mobile BAT unit 4

Instrument Serial No. 008734

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 8 day of September, 2018 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

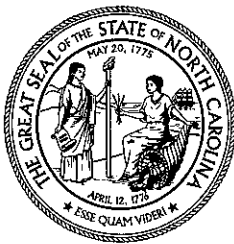

[Signature]
Signature of Certifying Official

664
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

ROBESON COUNTY BAT MOBILE UNIT 4 770

Serial Number: 008734
Test Date: 09/08/2018

Citation Number: M0000000-0
Subject's Name:

PRVENTIVE, MAINTENANCE

Subject's Date of Birth: 11/11/1911

Subject's Sex: Male

Driver's License State: XX

Driver's License Number: NONE

Analyst's Name: TODD, SHANE C

Permit Number: 11391E

Effective:

07/27/2018-07/01/2020

Officer's Name: NONE, NONE

Type of Agency: FTA

Agency: DHHS


Test Type: Breath Test

Lot Number: 702401

Exp Date: 08/02/2020

Test	g/210L	Time
DIAG	Pass	8:00pm
AIR BLK	.00	8:01pm
ACCY CHK	.08	8:01pm
AIR BLK	.00	8:02pm
SUB TEST	.00	8:03pm
AIR BLK	.00	8:04pm
SUB TEST	.00	8:06pm
AIR BLK	.00	8:07pm

Reported AC: .00 g/210L

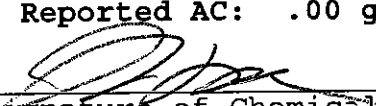

Signature of Chemical Analyst

Court CVR



Analyst

This form is used when performing Preventive Maintenance procedures
Forensic Tests for Alcohol Branch
Department of Health and Human Services
Rev. 12/2007

Intox EC/IR-II: Preventive Maintenance

ROBESON COUNTY BAT MOBILE UNIT 4 770

Serial Number: 008734 Test Record Number: 991
Test Date: 09/08/2018 Test Time: 8:13pm EDT

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	8:13pm
FLO	Pass	8:13pm
FC	Pass	8:13pm

Temperature Tests

Test	Status	Time
FC1	Pass	8:13pm
SRC	Pass	8:13pm
DET	Pass	8:13pm
BAR	Pass	8:13pm
BT	Pass	8:13pm

Blank Tests

Test	Status	Time
AIR	Pass	8:14pm

Printer Tests

Test	Status	Time
PRNT	Pass	8:14pm

CRC Tests

Test	Status	Time
COMP	Pass	8:14pm
CAL	Pass	8:14pm

Preventive Maintenance
Status: Pass



Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH

**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**


County Robeson Instrument Location Mobile BAT unit 4

Instrument Serial No. 008717

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 8 day of September, 2018 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

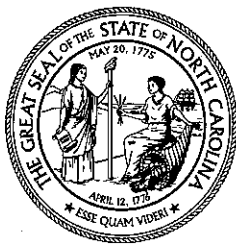

[Signature]
Signature of Certifying Official

664
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

ROBESON COUNTY BAT MOBILE UNIT 4 770

Serial Number: 008717
Test Date: 09/08/2018

Citation Number: M0000000-0
Subject's Name:
PREVENTIVE, MAINTENANCE
Subject's Date of Birth: 11/11/1911
Subject's Sex: Male
Driver's License State: XX
Driver's License Number: NONE


Analyst's Name: TODD, SHANE C
Permit Number: 11391E
Effective:
07/27/2018-07/01/2020

Officer's Name: NONE, NONE
Type of Agency: FTA
Agency: DHHS
Test Type: Breath Test

Lot Number: AG821401
Exp Date: 08/02/2020

Test	g/210L	Time
DIAG	Pass	8:02pm
AIR BLK	.00	8:03pm
ACCY CHK	.08	8:04pm
AIR BLK	.00	8:04pm
SUB TEST	.00	8:05pm
AIR BLK	.00	8:06pm
SUB TEST	.00	8:08pm
AIR BLK	.00	8:09pm

Reported AC: .00 g/210L

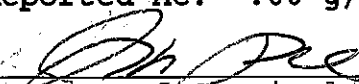

Signature of Chemical Analyst

Court CVR



Analyst

Intox EC/IR-II: Preventive Maintenance

ROBESON COUNTY BAT MOBILE UNIT 4 770

Serial Number: 008717 Test Record Number: 619
Test Date: 09/08/2018 Test Time: 8:14pm EDT

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	8:14pm
FLO	Pass	8:14pm
FC	Pass	8:15pm

Temperature Tests

Test	Status	Time
FC1	Pass	8:15pm
SRC	Pass	8:15pm
DET	Pass	8:15pm
BAR	Pass	8:15pm
BT	Pass	8:15pm

Blank Tests

Test	Status	Time
AIR	Pass	8:15pm

Printer Tests

Test	Status	Time
PRNT	Pass	8:15pm

CRC Tests

Test	Status	Time
COMP	Pass	8:15pm
CAL	Pass	8:15pm

Preventive Maintenance
Status: Pass



Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH


PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II

County SAMPSON Instrument Location SAMPSON COUNTY
Instrument Serial No. 008825 SHERIFF'S OFFICE

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 25 day of SEPTEMBER, 2018, the foregoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

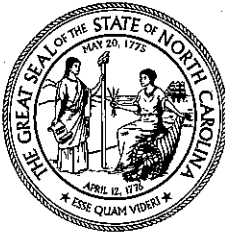

Al R. B...
Signature of Certifying Official

648
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

SAMPSON COUNTY SAMPSON CO SD 810

Serial Number: 008825

Test Date: 09/25/2018

Citation Number: M0000000-0

Subject's Name:

PREVENTIVE, MAINTENANCE

Subject's Date of Birth: 11/11/1911

Subject's Sex: Male

Driver's License State: XX

Driver's License Number: NONE

Analyst's Name: BARNES, ALVIN R

Permit Number: 15671E

Effective:

07/01/2017-07/01/2019

Officer's Name: NONE, NONE

Type of Agency: FTA

Agency: DHHS

Test Type: Breath Test

Lot Number: AG702401


Exp Date: 01/24/2019

Test	g/210L	Time
DIAG	Pass	3:48pm
AIR BLK	.00	3:49pm
ACCY CHK	.07	3:49pm
AIR BLK	.00	3:50pm
SUB TEST	.00	3:51pm
AIR BLK	.00	3:52pm
SUB TEST	.00	3:53pm
AIR BLK	.00	3:54pm

Reported AC: .00 g/210L

Signature of Chemical Analyst

Court CVR

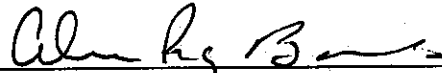

Analyst

Intox EC/IR-II: Preventive Maintenance

SAMPSON COUNTY SAMPSON CO SD 810

Serial Number: 008825 Test Record Number: 2544
Test Date: 09/25/2018 Test Time: 3:56pm EDT

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	3:56pm
FLO	Pass	3:56pm
FC	Pass	3:56pm

Temperature Tests

Test	Status	Time
FC1	Pass	3:56pm
SRC	Pass	3:56pm
DET	Pass	3:56pm
BAR	Pass	3:56pm
BT	Pass	3:56pm

Blank Tests

Test	Status	Time
AIR	Pass	3:57pm


Printer Tests

Test	Status	Time
PRNT	Pass	3:57pm

CRC Tests

Test	Status	Time
COMP	Pass	3:57pm
CAL	Pass	3:57pm

Preventive Maintenance
Status: Pass

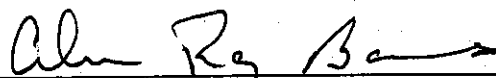

Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH

**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**

County SAMPSON Instrument Location SAMPSON COUNTY

Instrument Serial No. 008877 SHERIFF'S OFFICE

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 25 day of SEPTEMBER, 2018, the foregoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.



Alan R. Bunn
Signature of Certifying Official

648
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

SAMPSON COUNTY SAMPSON COUNTY SD 810

Serial Number: 008877
Test Date: 09/25/2018

Citation Number: M0000000-0
Subject's Name:
PREVENTIVE, MAINTENANCE
Subject's Date of Birth: 11/11/1911
Subject's Sex: Male
Driver's License State: XX
Driver's License Number: NONE

Analyst's Name: BARNES, ALVIN R
Permit Number: 15671E
Effective:
07/01/2017-07/01/2019

Officer's Name: NONE, NONE
Type of Agency: FTA
Agency: DHHS
Test Type: Breath Test


Lot Number: AG805802
Exp Date: 02/27/2020

Test	g/210L	Time
DIAG	Pass	3:41pm
AIR BLK	.00	3:42pm
ACCY CHK	.08	3:42pm
AIR BLK	.00	3:43pm
SUB TEST	.00	3:44pm
AIR BLK	.00	3:45pm
SUB TEST	.00	3:46pm
AIR BLK	.00	3:47pm

Reported AC: .00 g/210L

Signature of Chemical Analyst

Court CVR

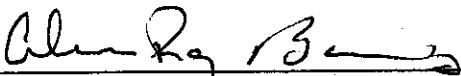

Analyst

Intox EC/IR-II: Preventive Maintenance

SAMPSON COUNTY SAMPSON COUNTY SD 810

Serial Number: 008877 Test Record Number: 2919
Test Date: 09/25/2018 Test Time: 3:51pm EDT

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	3:52pm
FLO	Pass	3:52pm
FC	Pass	3:52pm

Temperature Tests

Test	Status	Time
FC1	Pass	3:52pm
SRC	Pass	3:52pm
DET	Pass	3:52pm
BAR	Pass	3:52pm
BT	Pass	3:52pm

Blank Tests

Test	Status	Time
AIR	Pass	3:53pm

Printer Tests

Test	Status	Time
PRNT	Pass	3:53pm

CRC Tests

Test	Status	Time
COMP	Pass	3:53pm
CAL	Pass	3:53pm

Preventive Maintenance
Status: Pass



Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH


**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**

County Surry Instrument Location Mount Airy
Instrument Serial No. 008943 Police Department

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 26 day of September, 2018, the foregoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

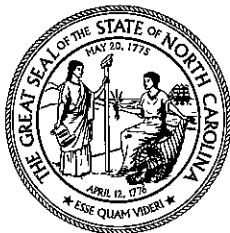

Kent D. Smith
Signature of Certifying Official

657
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

SURRY COUNTY MOUNT AIRY PD 850

Serial Number: 008943
Test Date: 09/26/2018

Citation Number: M0000000-0
Subject's Name:
PREVENTIVE, MAINTENANCE
Subject's Date of Birth: 11/11/1911
Subject's Sex: Male
Driver's License State: XX
Driver's License Number: NONE


Analyst's Name:
BENFIELD II, KENNETH R
Permit Number: 22067E
Effective:
06/01/2018-06/01/2020

Officer's Name: NONE, NONE
Type of Agency: FTA
Agency: DHHS
Test Type: Breath Test

Lot Number: AG702302
Exp Date: 01/23/2019


Test	g/210L	Time
DIAG	Pass	4:22pm
AIR BLK	.00	4:23pm
ACCY CHK	.08	4:23pm
AIR BLK	.00	4:24pm
SUB TEST	.00	4:25pm
AIR BLK	.00	4:26pm
SUB TEST	.00	4:27pm
AIR BLK	.00	4:28pm

Reported AC: .00 g/210L

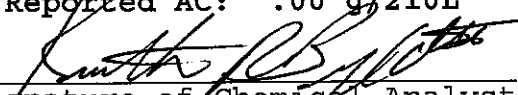

Signature of Chemical Analyst

Court CVR

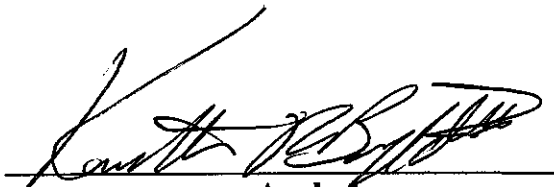

Analyst

Intox EC/IR-II: Preventive Maintenance

SURRY COUNTY MOUNT AIRY PD 850

Serial Number: 008943 Test Record Number: 2041
Test Date: 09/26/2018 Test Time: 4:29pm EDT

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	4:29pm
FLO	Pass	4:29pm
FC	Pass	4:29pm

Temperature Tests

Test	Status	Time
FC1	Pass	4:29pm
SRC	Pass	4:29pm
DET	Pass	4:29pm
BAR	Pass	4:29pm
BT	Pass	4:29pm

Blank Tests

Test	Status	Time
AIR	Pass	4:30pm


Printer Tests

Test	Status	Time
PRNT	Pass	4:30pm

CRC Tests

Test	Status	Time
COMP	Pass	4:30pm
CAL	Pass	4:30pm

Preventive Maintenance
Status: Pass

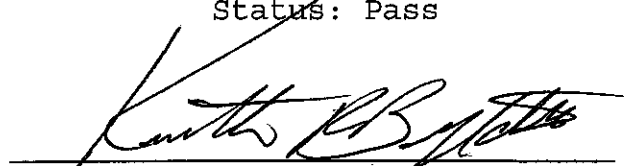

Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH


PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II

County Surry Instrument Location Pilot Mountain
Instrument Serial No. 008938 Police Department

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 26 day of September, 2018, the foregoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

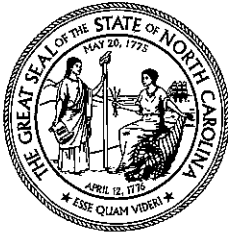

[Signature]
Signature of Certifying Official

657
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

SURRY COUNTY PILOT MOUNTAIN PD 850

Serial Number: 008938
Test Date: 09/26/2018

Citation Number: M0000000-0
Subject's Name:
PREVENTIVE, MAINTENANCE
Subject's Date of Birth: 11/11/1911
Subject's Sex: Male
Driver's License State: XX
Driver's License Number: NONE


Analyst's Name:
BENFIELD II, KENNETH R
Permit Number: 22067E
Effective:
06/01/2018-06/01/2020

Officer's Name: NONE, NONE
Type of Agency: FTA
Agency: DHHS
Test Type: Breath Test

Lot Number: AG734101
Exp Date: 12/07/2019


Test	g/210L	Time
DIAG	Pass	2:40pm
AIR BLK	.00	2:41pm
ACCY CHK	.08	2:41pm
AIR BLK	.00	2:42pm
SUB TEST	.00	2:43pm
AIR BLK	.00	2:44pm
SUB TEST	.00	2:45pm
AIR BLK	.00	2:46pm

Reported AC: .00 g/210L

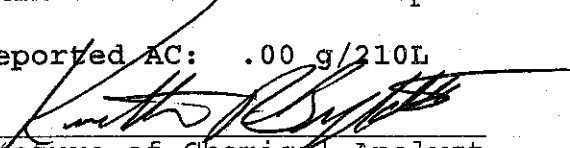

Signature of Chemical Analyst

Court CVR

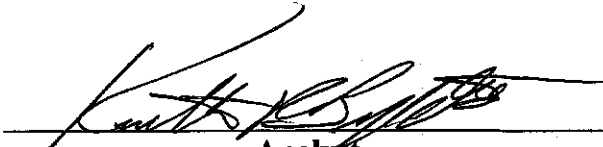

Analyst

Intox EC/IR-II: Preventive Maintenance

SURRY COUNTY PILOT MOUNTAIN PD 850

Serial Number: 008938 Test Record Number: 620
Test Date: 09/26/2018 Test Time: 2:47pm EDT

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	2:47pm
FLO	Pass	2:47pm
FC	Pass	2:47pm

Temperature Tests

Test	Status	Time
FC1	Pass	2:47pm
SRC	Pass	2:47pm
DET	Pass	2:47pm
BAR	Pass	2:47pm
BT	Pass	2:47pm

Blank Tests

Test	Status	Time
AIR	Pass	2:48pm


Printer Tests

Test	Status	Time
PRNT	Pass	2:48pm

CRC Tests

Test	Status	Time
COMP	Pass	2:48pm
CAL	Pass	2:48pm

Preventive Maintenance
Status: Pass

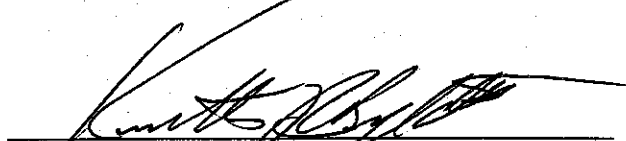

Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH


**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**

County Surry Instrument Location Surry County Jail
Instrument Serial No. 008934 Dobson, N.C.

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 26 day of September, 2018, the foregoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

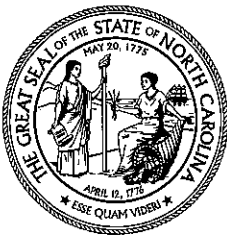

[Signature]
Signature of Certifying Official

657

Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

SURRY COUNTY SURRY CO JAIL 850

Serial Number: 008934
Test Date: 09/26/2018

Citation Number: M0000000-0
Subject's Name:
PREVENTIVE, MAINTENANCE
Subject's Date of Birth: 11/11/1911
Subject's Sex: Male
Driver's License State: XX
Driver's License Number: NONE


Analyst's Name:
BENFIELD II, KENNETH R
Permit Number: 22067E
Effective:
06/01/2018-06/01/2020

Officer's Name: NONE, NONE
Type of Agency: FTA
Agency: DHHS
Test Type: Breath Test

Lot Number: AG702302
Exp Date: 01/23/2019


Test	g/210L	Time
DIAG	Pass	3:34pm
AIR BLK	.00	3:35pm
ACCY CHK	.07	3:36pm
AIR BLK	.00	3:36pm
SUB TEST	.00	3:37pm
AIR BLK	.00	3:38pm
SUB TEST	.00	3:39pm
AIR BLK	.00	3:40pm

Reported AC: .00 g/210L

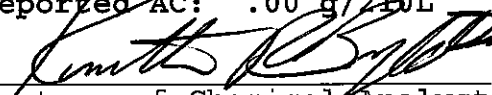

Signature of Chemical Analyst

Court CVR



Analyst

Intox EC/IR-II: Preventive Maintenance

SURRY COUNTY SURRY CO JAIL 850

Serial Number: 008934 Test Record Number: 1945
Test Date: 09/26/2018 Test Time: 3:41pm EDT

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	3:42pm
FLO	Pass	3:42pm
FC	Pass	3:42pm

Temperature Tests

Test	Status	Time
FC1	Pass	3:42pm
SRC	Pass	3:42pm
DET	Pass	3:42pm
BAR	Pass	3:42pm
BT	Pass	3:42pm

Blank Tests

Test	Status	Time
AIR	Pass	3:42pm


Printer Tests

Test	Status	Time
PRNT	Pass	3:42pm

CRC Tests

Test	Status	Time
COMP	Pass	3:43pm
CAL	Pass	3:43pm

Preventive Maintenance
Status: Pass

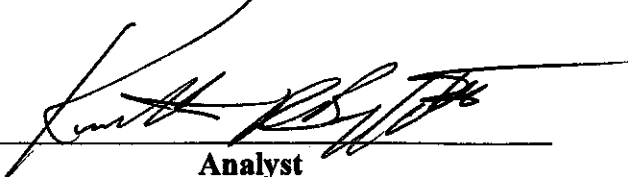

Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH

**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**

County Swain

Instrument Location Cherokee Tribal Det.


Instrument Serial No. 008782

Cherokee, NC

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 25 day of September, 2018 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

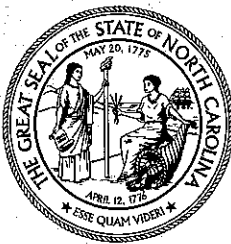

David R. Cuthbert

Signature of Certifying Official

635

Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

SWAIN COUNTY CHEROKEE DETENTION 860

Serial Number: 008782
Test Date: 09/25/2018

Citation Number: M0000000-0
Subject's Name:
PREVENTIVE, MAINTENANCE
Subject's Date of Birth: 11/11/1911
Subject's Sex: Male
Driver's License State: XX
Driver's License Number: NONE

Analyst's Name: CUTLER, DANIEL R
Permit Number: 08457E
Effective:
09/01/2017-09/01/2019

Officer's Name: NONE,
Type of Agency: FTA
Agency: DHHS
Test Type: Breath Test


Lot Number: AG807101
Exp Date: 03/12/2020

Test	g/210L	Time
DIAG	Pass	9:32am
AIR BLK	.00	9:33am
ACCY CHK	.08	9:34am
AIR BLK	.00	9:34am
SUB TEST	.00	9:35am
AIR BLK	.00	9:36am
SUB TEST	.00	9:37am
AIR BLK	.00	9:38am

Reported AC: .00 g/210L

Signature of Chemical Analyst

Court CVR

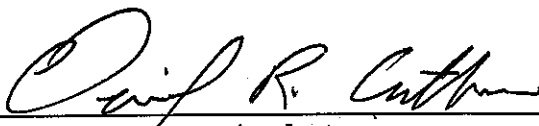
Analyst

Intox EC/IR-II: Preventive Maintenance

SWAIN COUNTY CHEROKEE DETENTION 860

Serial Number: 008782 Test Record Number: 1067
Test Date: 09/25/2018 Test Time: 9:39am EDT

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	9:40am
FLO	Pass	9:40am
FC	Pass	9:40am

Temperature Tests

Test	Status	Time
FC1	Pass	9:40am
SRC	Pass	9:40am
DET	Pass	9:40am
BAR	Pass	9:40am
BT	Pass	9:40am

Blank Tests

Test	Status	Time
AIR	Pass	9:40am


Printer Tests

Test	Status	Time
PRNT	Pass	9:40am

CRC Tests

Test	Status	Time
COMP	Pass	9:41am
CAL	Pass	9:41am

Preventive Maintenance
Status: Pass

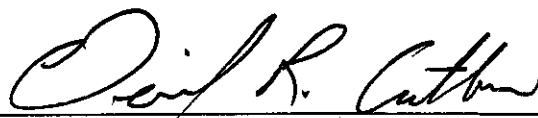

Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH


PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II

County Tyrrell Instrument Location Tyrrell Co. S.O.
Instrument Serial No. 008902 412 MAIN ST., Columbia, N.C.

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 24th day of September, 20 18 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

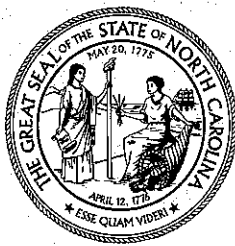

Jim A. Keene
Signature of Certifying Official

647
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

TYRRELL COUNTY SHERIFF'S OFFICE 880

Serial Number: 008902
Test Date: 09/24/2018

Citation Number: M0000000-0
Subject's Name:
PREVENTIVE, MAINTENANCE
Subject's Date of Birth: 11/11/1911
Subject's Sex: Male
Driver's License State: XX
Driver's License Number: NONE


Analyst's Name: KEESLER, LINDA A
Permit Number: 11646E
Effective:
08/01/2017-08/01/2019

Officer's Name: NONE, NONE
Type of Agency: FTA
Agency: DHHS
Test Type: Breath Test

Lot Number: AG734102
Exp Date: 12/07/2019

Test	g/210L	Time
DIAG	Pass	1:06pm
AIR BLK	.00	1:07pm
ACCY CHK	.07	1:07pm
AIR BLK	.00	1:08pm
SUB TEST	.00	1:09pm
AIR BLK	.00	1:10pm
SUB TEST	.00	1:11pm
AIR BLK	.00	1:12pm

Reported AC: .00 g/210L

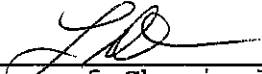

Signature of Chemical Analyst

Court CVR



Analyst

Intox EC/IR-II: Preventive Maintenance

TYRRELL COUNTY SHERIFF'S OFFICE 880

Serial Number: 008902 Test Record Number: 774
Test Date: 09/24/2018 Test Time: 1:13pm EDT

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	1:13pm
FLO	Pass	1:13pm
FC	Pass	1:13pm

Temperature Tests

Test	Status	Time
FC1	Pass	1:14pm
SRC	Pass	1:14pm
DET	Pass	1:14pm
BAR	Pass	1:14pm
BT	Pass	1:14pm

Blank Tests

Test	Status	Time
AIR	Pass	1:14pm


Printer Tests

Test	Status	Time
PRNT	Pass	1:14pm

CRC Tests

Test	Status	Time
COMP	Pass	1:14pm
CAL	Pass	1:14pm

Preventive Maintenance
Status: Pass

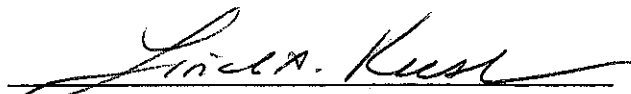
Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH

PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II


County VAANCE Instrument Location NAT MONTIE UNIT 6

Instrument Serial No. 009584 HENDERSON

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 21 day of SEPTEMBER, 2010, the foregoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

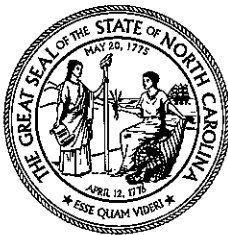

[Signature]
Signature of Certifying Official

663
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

VANCE COUNTY BAT MOBILE UNIT 6 900

Serial Number: 008584

Test Date: 09/28/2018

Citation Number: M0000000-0

Subject's Name:

PREVENTIVE, MAINTENANCE

Subject's Date of Birth: 11/11/1911

Subject's Sex: Male

Driver's License State: XX

Driver's License Number: NONE

Analyst's Name: VARNELL, BRYON L

Permit Number: 16896E

Effective:

09/22/2017-09/01/2019

Officer's Name: NONE, NONE

Type of Agency: FTA

Agency: DHHS


Test Type: Breath Test

Lot Number: AG807101

Exp Date: 03/12/2020

Test	g/210L	Time
DIAG	Pass	9:15pm
AIR BLK	.00	9:16pm
ACCY CHK	.07	9:16pm
AIR BLK	.00	9:17pm
SUB TEST	.00	9:18pm
AIR BLK	.00	9:19pm
SUB TEST	.00	9:20pm
AIR BLK	.00	9:21pm

Reported AC: .00 g/210L

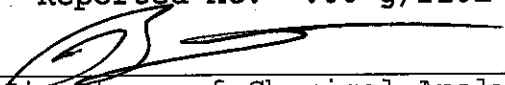
Signature of Chemical Analyst

Court CVR


Analyst

**This form is used when performing Preventive Maintenance procedures
Forensic Tests for Alcohol Branch
Department of Health and Human Services
Rev. 12/2007**

Intox EC/IR-II: Preventive Maintenance

VANCE COUNTY BAT MOBILE UNIT 6 900

Serial Number: 008584 Test Record Number: 2231

Test Date: 09/28/2018 Test Time: 9:29pm EDT

System Check: Passed

Baseline Tests

Test	Status	Time
IR	Pass	9:29pm
FLO	Pass	9:29pm
FC	Pass	9:29pm

Temperature Tests

Test	Status	Time
FC1	Pass	9:29pm
SRC	Pass	9:29pm
DET	Pass	9:29pm
BAR	Pass	9:29pm
BT	Pass	9:29pm

Blank Tests

Test	Status	Time
AIR	Pass	9:30pm

Printer Tests

Test	Status	Time
PRNT	Pass	9:30pm

CRC Tests

Test	Status	Time
COMP	Pass	9:30pm
CAL	Pass	9:30pm

Preventive Maintenance
Status: Pass



Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH

PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II


County VANCE Instrument Location 11th MONTGOMERY UNIT 6

Instrument Serial No. 008779 HENDERSON

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 24 day of SEPTEMBER, 2014, the foregoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

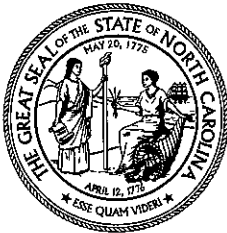

[Signature]
Signature of Certifying Official

667
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

VANCE COUNTY BAT MOBILE UNIT 6 900

Serial Number: 008779
Test Date: 09/28/2018

Citation Number: M0000000-0
Subject's Name:
PREVENTIVE, MAINTENANCE
Subject's Date of Birth: 11/11/1911
Subject's Sex: Male
Driver's License State: XX
Driver's License Number: NONE


Analyst's Name: VARNELL, BRYON L
Permit Number: 16896E
Effective:
09/22/2017-09/01/2019

Officer's Name: NONE, NONE
Type of Agency: FTA
Agency: DHHS
Test Type: Breath Test


Lot Number: AG716201
Exp Date: 06/11/2019

Test	g/210L	Time
DIAG	Pass	9:14pm
AIR BLK	.00	9:15pm
ACCY CHK	.07	9:15pm
AIR BLK	.00	9:16pm
SUB TEST	.00	9:17pm
AIR BLK	.00	9:18pm
SUB TEST	.00	9:19pm
AIR BLK	.00	9:20pm

Reported AC: .00 g/210L


Signature of Chemical Analyst

Court CVR

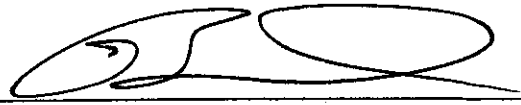

Analyst

**This form is used when performing Preventive Maintenance procedures
Forensic Tests for Alcohol Branch
Department of Health and Human Services
Rev. 12/2007**

Intox EC/IR-II: Preventive Maintenance

VANCE COUNTY BAT MOBILE UNIT 6 900

Serial Number: 008779 Test Record Number: 3510
Test Date: 09/28/2018 Test Time: 9:23pm EDT

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	9:24pm
FLO	Pass	9:24pm
FC	Pass	9:24pm

Temperature Tests

Test	Status	Time
FC1	Pass	9:24pm
SRC	Pass	9:24pm
DET	Pass	9:24pm
BAR	Pass	9:24pm
BT	Pass	9:24pm

Blank Tests

Test	Status	Time
AIR	Pass	9:25pm


Printer Tests

Test	Status	Time
PRNT	Pass	9:25pm

CRC Tests

Test	Status	Time
COMP	Pass	9:25pm
CAL	Pass	9:25pm

Preventive Maintenance
Status: Pass

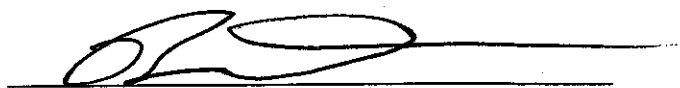

Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH


PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II

County VANCE Instrument Location SAF MONZIE UNIT 6
Instrument Serial No. 008637 HENDERSON

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 28 day of SEPTEMBER, 20 11, the foregoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

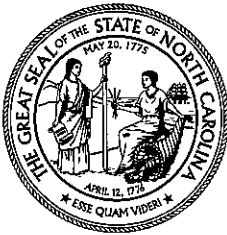

[Signature]

Signature of Certifying Official

6663

Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

VANCE COUNTY BAT MOBILE UNIT 6 900

Serial Number: 008637

Test Date: 09/28/2018

Citation Number: M0000000-0

Subject's Name:

PREVENTIVE, MAINTENANCE

Subject's Date of Birth: 11/11/1911

Subject's Sex: Male

Driver's License State: XX

Driver's License Number: NONE

Analyst's Name: VARNELL, BRYON L

Permit Number: 16896E

Effective:

09/22/2017-09/01/2019

Officer's Name: NONE, NONE

Type of Agency: FTA

Agency: DHHS


Test Type: Breath Test

Lot Number: AG716201

Exp Date: 06/11/2019

Test	g/210L	Time
DIAG	Pass	9:12pm
AIR BLK	.00	9:13pm
ACCY CHK	.07	9:13pm
AIR BLK	.00	9:14pm
SUB TEST	.00	9:15pm
AIR BLK	.00	9:16pm
SUB TEST	.00	9:17pm
AIR BLK	.00	9:18pm

Reported AC: .00 g/210L

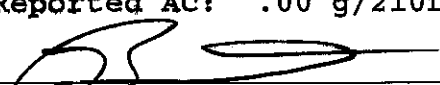

Signature of Chemical Analyst

Court CVR



Analyst

**This form is used when performing Preventive Maintenance procedures
Forensic Tests for Alcohol Branch
Department of Health and Human Services
Rev. 12/2007**

Intox EC/IR-II: Preventive Maintenance

VANCE COUNTY BAT MOBILE UNIT 6 900

Serial Number: 008637 Test Record Number: 2946
Test Date: 09/28/2018 Test Time: 9:22pm EDT

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	9:23pm
FLO	Pass	9:23pm
FC	Pass	9:23pm

Temperature Tests

Test	Status	Time
FC1	Pass	9:23pm
SRC	Pass	9:23pm
DET	Pass	9:23pm
BAR	Pass	9:23pm
BT	Pass	9:23pm

Blank Tests

Test	Status	Time
AIR	Pass	9:23pm

Printer Tests

Test	Status	Time
PRNT	Pass	9:24pm

CRC Tests

Test	Status	Time
COMP	Pass	9:24pm
CAL	Pass	9:24pm

Preventive Maintenance
Status: Pass



Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH


**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**

County Washington Instrument Location Washington Co. S.O.
Instrument Serial No. 008829 Adams St., Plymouth, NC.

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 27th day of September, 2018 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

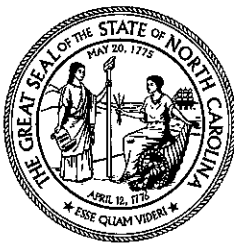

Keey A. D. P.
Signature of Certifying Official

643
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

WASHINGTON COUNTY SHERIFF'S OFFICE 930

Serial Number: 008829
Test Date: 09/27/2018

Citation Number: M0000000-0
Subject's Name:

PREVENTIVE, MAINTENANCE

Subject's Date of Birth: 11/11/1911

Subject's Sex: Male

Driver's License State: XX

Driver's License Number: NONE

Analyst's Name: GUARD, KELLY G

Permit Number: 12955E

Effective:

06/01/2017-06/01/2019

Officer's Name: NONE, NONE

Type of Agency: FTA

Agency: DHHS

Test Type: Breath Test

Lot Number: AG710701

Exp Date: 04/17/2019


Test	g/210L	Time
DIAG	Pass	10:22am
AIR BLK	.00	10:23am
ACCY CHK	.07	10:23am
AIR BLK	.00	10:24am
SUB TEST	.00	10:25am
AIR BLK	.00	10:26am
SUB TEST	.00	10:27am
AIR BLK	.00	10:28am

Reported AC: .00 g/210L



Signature of Chemical Analyst

Court CVR



Analyst

This form is used when performing Preventive Maintenance procedures
Forensic Tests for Alcohol Branch
Department of Health and Human Services
Rev. 12/2007

Intox EC/IR-II: Preventive Maintenance

WASHINGTON COUNTY SHERIFF'S OFFICE 930

Serial Number: 008829 Test Record Number: 873
Test Date: 09/27/2018 Test Time: 10:29am EDT

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	10:29am
FLO	Pass	10:29am
FC	Pass	10:30am

Temperature Tests

Test	Status	Time
FC1	Pass	10:30am
SRC	Pass	10:30am
DET	Pass	10:30am
BAR	Pass	10:30am
BT	Pass	10:30am

Blank Tests

Test	Status	Time
AIR	Pass	10:30am

Printer Tests

Test	Status	Time
PRNT	Pass	10:30am

CRC Tests

Test	Status	Time
COMP	Pass	10:31am
CAL	Pass	10:31am

Preventive Maintenance
Status: Pass



Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH


**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**

County Washington Instrument Location Washington Co. S.O.
Instrument Serial No. 008829 Adams St., Plymouth, N.C.

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 12th day of September, 2018 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

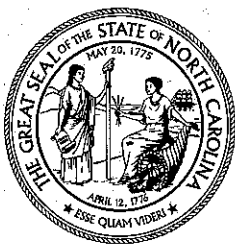

[Signature]
Signature of Certifying Official

643

Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

WASHINGTON COUNTY SHERIFF'S OFFICE 930

Serial Number: 008829
Test Date: 09/12/2018

Citation Number: M0000000-0
Subject's Name:

PREVENTIVE, MAINTENANCE

Subject's Date of Birth: 11/11/1911

Subject's Sex: Male

Driver's License State: XX

Driver's License Number: NONE

Analyst's Name: GUARD, KELLY G

Permit Number: 12955E

Effective:

06/01/2017-06/01/2019

Officer's Name: NONE, NONE

Type of Agency: FTA

Agency: DHHS


Test Type: Breath Test

Lot Number: AG710701

Exp Date: 04/17/2019


Test	g/210L	Time
DIAG	Pass	10:20am
AIR BLK	.00	10:21am
ACCY CHK	.07	10:22am
AIR BLK	.00	10:23am
SUB TEST	.00	10:24am
AIR BLK	.00	10:24am
SUB TEST	.00	10:26am
AIR BLK	.00	10:27am

Reported AC: .00 g/210L

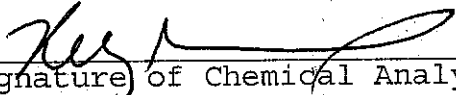

Signature of Chemical Analyst

Court CVR



Analyst

Intox EC/IR-II: Preventive Maintenance

WASHINGTON COUNTY SHERIFF'S OFFICE 930

Serial Number: 008829 Test Record Number: 869
Test Date: 09/12/2018 Test Time: 10:28am EDT

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	10:28am
FLO	Pass	10:28am
FC	Pass	10:28am

Temperature Tests

Test	Status	Time
FC1	Pass	10:28am
SRC	Pass	10:28am
DET	Pass	10:28am
BAR	Pass	10:28am
BT	Pass	10:28am

Blank Tests

Test	Status	Time
AIR	Pass	10:29am

Printer Tests

Test	Status	Time
PRNT	Pass	10:29am

CRC Tests

Test	Status	Time
COMP	Pass	10:29am
CAL	Pass	10:29am

Preventive Maintenance
Status: Pass



Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH

PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II


County Wayne Instrument Location Wayne Co. Detention CTR.

Instrument Serial No. 008879 207 E. CHESTNUT ST., GOLDSBORO, N.C.

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 26th day of SEPTEMBER, 20 18 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

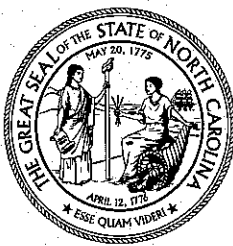

Jintha A. Kease
Signature of Certifying Official

647
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

WAYNE COUNTY WAYNE CO DETENTION 950

Serial Number: 008879
Test Date: 09/26/2018

Citation Number: M0000000-0
Subject's Name:
PREVENTIVE, MAINTENANCE
Subject's Date of Birth: 11/11/1911
Subject's Sex: Male
Driver's License State: XX
Driver's License Number: NONE

Analyst's Name: KEESLER, LINDA A
Permit Number: 11646E
Effective:
08/01/2017-08/01/2019

Officer's Name: NONE, NONE
Type of Agency: FTA
Agency: DHHS
Test Type: Breath Test

Lot Number: AG814901
Exp Date: 05/29/2020


Test	g/210L	Time
DIAG	Pass	1:32pm
AIR BLK	.00	1:33pm
ACCY CHK	.08	1:34pm
AIR BLK	.00	1:35pm
SUB TEST	.00	1:35pm
AIR BLK	.00	1:36pm
SUB TEST	.00	1:38pm
AIR BLK	.00	1:39pm

Reported AC: .00 g/210L



Signature of Chemical Analyst

Court CVR

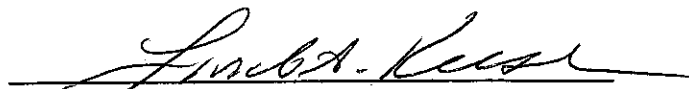

Analyst

Intox EC/IR-II: Preventive Maintenance

WAYNE COUNTY WAYNE CO DETENTION 950

Serial Number: 008879 Test Record Number: 1072
Test Date: 09/26/2018 Test Time: 1:39pm EDT

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	1:40pm
FLO	Pass	1:40pm
FC	Pass	1:40pm

Temperature Tests

Test	Status	Time
FC1	Pass	1:40pm
SRC	Pass	1:40pm
DET	Pass	1:40pm
BAR	Pass	1:40pm
BT	Pass	1:40pm

Blank Tests

Test	Status	Time
AIR	Pass	1:40pm

Printer Tests

Test	Status	Time
PRNT	Pass	1:41pm

CRC Tests

Test	Status	Time
COMP	Pass	1:41pm
CAL	Pass	1:41pm

Preventive Maintenance
Status: Pass


Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH

**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**

County Watauga Instrument Location Watauga Co Jail
Instrument Serial No. 008715 Boone, NC

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 29 day of September, 2018 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.



[Signature]
Signature of Certifying Official

649
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

WATAUGA COUNTY WATAUGA JAIL 940

Serial Number: 008715

Test Date: 09/29/2018

Citation Number: M0000000-0

Subject's Name:

PREVENTIVE, MAINTENANCE

Subject's Date of Birth: 11/11/1911

Subject's Sex: Male

Driver's License State: XX

Driver's License Number: NONE

Analyst's Name: BURNETTE, ANTHONY J

Permit Number: 11304E

Effective:

05/01/2017-05/01/2019

Officer's Name: NONE,

Type of Agency: FTA

Agency: DHHS

Test Type: Breath Test

Lot Number: AG722408

Exp Date: 08/12/2019

Test	g/210L	Time
DIAG	Pass	12:45pm
AIR BLK	.00	12:46pm
ACCY CHK	.08	12:47pm
AIR BLK	.00	12:48pm
SUB TEST	.00	12:49pm
AIR BLK	.00	12:50pm
SUB TEST	.00	12:51pm
AIR BLK	.00	12:52pm

Reported AC: .00 g/210L

Signature of Chemical Analyst

Court CVR



Analyst

**This form is used when performing Preventive Maintenance procedures
Forensic Tests for Alcohol Branch
Department of Health and Human Services
Rev. 12/2007**

Intox EC/IR-II: Preventive Maintenance

WATAUGA COUNTY WATAUGA JAIL 940

Serial Number: 008715 Test Record Number: 2165
Test Date: 09/29/2018 Test Time: 12:53pm EDT

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	12:53pm
FLO	Pass	12:53pm
FC	Pass	12:53pm

Temperature Tests

Test	Status	Time
FC1	Pass	12:53pm
SRC	Pass	12:53pm
DET	Pass	12:53pm
BAR	Pass	12:53pm
BT	Pass	12:53pm

Blank Tests

Test	Status	Time
AIR	Pass	12:54pm

Printer Tests

Test	Status	Time
PRNT	Pass	12:54pm

CRC Tests

Test	Status	Time
COMP	Pass	12:54pm
CAL	Pass	12:54pm

Preventive Maintenance
Status: Pass


Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH

PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II

County Wilson

Instrument Location Wilson Co. Detention CTR.


Instrument Serial No. 008627

100 E. Greene ST., Wilson, N.C.

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 20th day of September, 2018 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

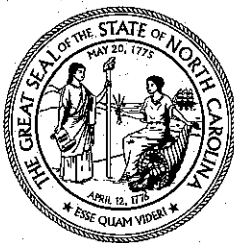

Shirley A. Keen
Signature of Certifying Official

647
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

WILSON COUNTY DETENTION CENTER 970

Serial Number: 008627
Test Date: 09/20/2018

Citation Number: M0000000-0
Subject's Name:
PREVENTIVE, MAINTENANCE
Subject's Date of Birth: 11/11/1911
Subject's Sex: Male
Driver's License State: XX
Driver's License Number: NONE

Analyst's Name: KEESLER, LINDA A
Permit Number: 11646
Effective:
08/01/2017-08/01/2019

Officer's Name: NONE, NONE
Type of Agency: FTA
Agency: DHHS
Test Type: Breath Test

Lot Number: AG721401
Exp Date: 08/02/2019

Test	g/210L	Time
DIAG	Pass	12:34pm
AIR BLK	.00	12:35pm
ACCY CHK	.07	12:36pm
AIR BLK	.00	12:37pm
SUB TEST	.00	12:37pm
AIR BLK	.00	12:38pm
SUB TEST	.00	12:40pm
AIR BLK	.00	12:40pm

Reported AC: 00 g/210L



Signature of Chemical Analyst

Court CVR



Analyst

Intox EC/IR-II: Preventive Maintenance

WILSON COUNTY DETENTION CENTER 970

Serial Number: 008627 Test Record Number: 2249
Test Date: 09/20/2018 Test Time: 12:42pm EDT

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	12:42pm
FLO	Pass	12:42pm
FC	Pass	12:43pm

Temperature Tests

Test	Status	Time
FC1	Pass	12:43pm
SRC	Pass	12:43pm
DET	Pass	12:43pm
BAR	Pass	12:43pm
BT	Pass	12:43pm

Blank Tests

Test	Status	Time
AIR	Pass	12:43pm

Printer Tests

Test	Status	Time
PRNT	Pass	12:43pm

CRC Tests

Test	Status	Time
COMP	Pass	12:43pm
CAL	Pass	12:43pm

Preventive Maintenance
Status: Pass


Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH

**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**


County Wilson Instrument Location Wilson Co. Detention Ctr.

Instrument Serial No. 008652 100 E. Green ST., Wilson, N.C.

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 20th day of September, 2018 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

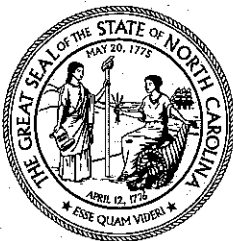

Linba Kue
Signature of Certifying Official

647
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

WILSON COUNTY DETENTION CENTER 970

Serial Number: 008652

Test Date: 09/20/2018

Citation Number: M0000000-0

Subject's Name:

PREVENTIVE, MAINTENANCE

Subject's Date of Birth: 11/11/1911

Subject's Sex: Male

Driver's License State: XX

Driver's License Number: NONE

Analyst's Name: KEESLER, LINDA A

Permit Number: 11646E

Effective:

08/01/2017-08/01/2019

Officer's Name: NONE, NONE

Type of Agency: FTA

Agency: DHHS

Test Type: Breath Test

Lot Number: AG814901

Exp Date: 05/29/2020

Test	g/210L	Time
DIAG	Pass	12:36pm
AIR BLK	.00	12:37pm
ACCY CHK	.08	12:37pm
AIR BLK	.00	12:38pm
SUB TEST	.00	12:39pm
AIR BLK	.00	12:40pm
SUB TEST	.00	12:42pm
AIR BLK	.00	12:43pm

Reported AC: .00 g/210L



Signature of Chemical Analyst

Court CVR



Analyst

**This form is used when performing Preventive Maintenance procedures
Forensic Tests for Alcohol Branch
Department of Health and Human Services
Rev. 12/2007**

Intox EC/IR-II: Preventive Maintenance

WILSON COUNTY DETENTION CENTER 970

Serial Number: 008652 Test Record Number: 3163
Test Date: 09/20/2018 Test Time: 12:44pm EDT

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	12:44pm
FLO	Pass	12:44pm
FC	Pass	12:44pm

Temperature Tests

Test	Status	Time
FC1	Pass	12:44pm
SRC	Pass	12:44pm
DET	Pass	12:44pm
BAR	Pass	12:44pm
BT	Pass	12:44pm

Blank Tests

Test	Status	Time
AIR	Pass	12:45pm

Printer Tests

Test	Status	Time
PRNT	Pass	12:45pm

CRC Tests

Test	Status	Time
COMP	Pass	12:45pm
CAL	Pass	12:45pm

Preventive Maintenance
Status: Pass


Analyst