

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH

**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**


County Macon Instrument Location Macon Co. Jail

Instrument Serial No. 008618 Franklin, NC

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 12 day of Nov., 2008 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

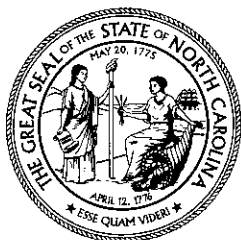

Paul R. Little
Signature of Certifying Official

635
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

MACON COUNTY MACON COUNTY JAIL 550

Serial Number: 008618
Test Date: 11/12/2008

Citation Number: M0000000-0
Subject's Name:
PREVENTIVE, MAINTENANCE
Subject's Date of Birth: 11/11/1911
Subject's Sex: Male
Driver's License State: XX
Driver's License Number: NONE

Analyst's Name: CUTLER, DANIEL R
Permit Number: 08457E
Effective:
12/01/2007-12/01/2009

Officer's Name: NONE,
Type of Agency: FTA
Agency: DHHS
Test Type: Breath Test


Lot Number: AG722701
Exp Date: 08/14/2009

Test	g/210L	Time
DIAG	Pass	11:38am
AIR BLK	.00	11:39am
ACCY CHK	.07	11:39am
AIR BLK	.00	11:40am
SUB TEST	.00	11:41am
AIR BLK	.00	11:42am
SUB TEST	.00	11:43am
AIR BLK	.00	11:44am

Reported AC: .00 g/210L

Signature of Chemical Analyst

Court CVR



Analyst

Intox EC/IR-II: Preventive Maintenance

MACON COUNTY MACON COUNTY JAIL 550

Serial Number: 008618 Test Record Number: 467
Test Date: 11/12/2008 Test Time: 11:45am EST

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	11:45am
FLO	Pass	11:45am
FC	Pass	11:45am

Temperature Tests

Test	Status	Time
FC1	Pass	11:45am
SRC	Pass	11:45am
DET	Pass	11:45am
BAR	Pass	11:45am
BT	Pass	11:45am

Blank Tests

Test	Status	Time
AIR	Pass	11:46am

Printer Tests

Test	Status	Time
PRNT	Pass	11:46am

CRC Tests

Test	Status	Time
COMP	Pass	11:46am
CAL	Pass	11:46am

Preventive Maintenance
Status: Pass



Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH

**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**

County Macon

Instrument Location Macon Co. Jail


Instrument Serial No. 1108789

Franklin, NC

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 12 day of Nov., 20 08 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

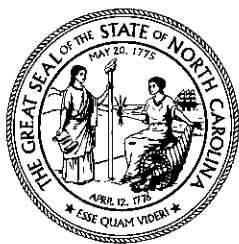

Paul R. [Signature]
Signature of Certifying Official

635
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

MACON COUNTY MACON COUNTY JAIL 550

Serial Number: 008789
Test Date: 11/12/2008

Citation Number: M0000000-0

Subject's Name:

PREVENTIVE, MAINTENANCE

Subject's Date of Birth: 11/11/1911

Subject's Sex: Male

Driver's License State: XX

Driver's License Number: NONE

Analyst's Name: CUTLER, DANIEL R

Permit Number: 08457E

Effective:

12/01/2007-12/01/2009

Officer's Name: NONE,

Type of Agency: FTA

Agency: DHHS

Test Type: Breath Test

Lot Number: AG722701


Exp Date: 08/14/2009

Test	g/210L	Time
DIAG	Pass	11:39am
AIR BLK	.00	11:40am
ACCY CHK	.07	11:40am
AIR BLK	.00	11:41am
SUB TEST	.00	11:42am
AIR BLK	.00	11:43am
SUB TEST	.00	11:44am
AIR BLK	.00	11:45am

Reported AC: .00 g/210L

Signature of Chemical Analyst

Court CVR

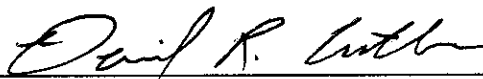

Analyst

Intox EC/IR-II: Preventive Maintenance

MACON COUNTY MACON COUNTY JAIL 550

Serial Number: 008789 Test Record Number: 69
Test Date: 11/12/2008 Test Time: 11:45am EST

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	11:46am
FLO	Pass	11:46am
FC	Pass	11:46am

Temperature Tests

Test	Status	Time
FC1	Pass	11:46am
SRC	Pass	11:46am
DET	Pass	11:46am
BAR	Pass	11:46am
BT	Pass	11:46am

Blank Tests

Test	Status	Time
AIR	Pass	11:47am

Printer Tests

Test	Status	Time
PRNT	Pass	11:47am

CRC Tests

Test	Status	Time
COMP	Pass	11:47am
CAL	Pass	11:47am

Preventive Maintenance
Status: Pass



Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH

**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**


County Macon Instrument Location Macon Co. Magistrate

Instrument Serial No. 008795 81 Maple St. Highlands, NC

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 12 day of Nov., 20 08 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

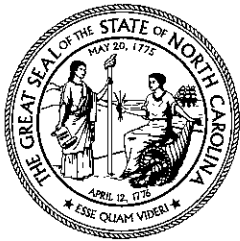

Edward R. Cuth

Signature of Certifying Official

635

Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

MACON COUNTY MACON CO MAGISTRATE 550

Serial Number: 008795
Test Date: 11/12/2008

Citation Number: M0000000-0
Subject's Name:
PREVENTIVE, MAINTENANCE
Subject's Date of Birth: 11/11/1911
Subject's Sex: Male
Driver's License State: XX
Driver's License Number: NONE

Analyst's Name: CUTLER, DANIEL R
Permit Number: 08457E
Effective:
12/01/2007-12/01/2009

Officer's Name: NONE,
Type of Agency: FTA
Agency: DHHS
Test Type: Breath Test

Lot Number: AG722802
Exp Date: 08/15/2009

Test	g/210L	Time
DIAG	Pass	1:36pm
AIR BLK	.00	1:36pm
ACCY CHK	.07	1:37pm
AIR BLK	.00	1:38pm
SUB TEST	.00	1:38pm
AIR BLK	.00	1:39pm
SUB TEST	.00	1:41pm
AIR BLK	.00	1:42pm

Reported AC: .00 g/210L

Signature of Chemical Analyst

Court CVR



Analyst

Intox EC/IR-II: Preventive Maintenance

MACON COUNTY MACON CO MAGISTRATE 550

Serial Number: 008795 Test Record Number: 77
Test Date: 11/12/2008 Test Time: 1:42pm EST

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	1:42pm
FLO	Pass	1:42pm
FC	Pass	1:42pm

Temperature Tests

Test	Status	Time
FC1	Pass	1:43pm
SRC	Pass	1:43pm
DET	Pass	1:43pm
BAR	Pass	1:43pm
BT	Pass	1:43pm

Blank Tests

Test	Status	Time
AIR	Pass	1:43pm

Printer Tests

Test	Status	Time
PRNT	Pass	1:43pm

CRC Tests

Test	Status	Time
COMP	Pass	1:43pm
CAL	Pass	1:43pm

Preventive Maintenance
Status: Pass



Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH

**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**


County Pitt Instrument Location Pitt Co. Detention Center

Instrument Serial No. 008668 124 Detention Dr, Greenville, N.C.

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 25th day of November, 20 08 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

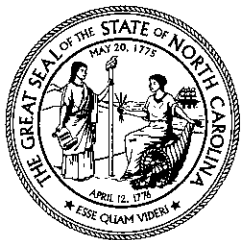

[Signature]
Signature of Certifying Official

643
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

PITT COUNTY PITT CO DETENTION 730

Serial Number: 008668

Test Date: 11/25/2008

Citation Number: M0000000-0

Subject's Name:

PREVENTIVE, MAINTENANCE

Subject's Date of Birth: 11/11/1911

Subject's Sex: Male

Driver's License State: XX

Driver's License Number: NONE

Analyst's Name: GUARD, KELLY G

Permit Number: 12955E

Effective:

12/01/2007-12/01/2009

Officer's Name: NONE, NONE

Type of Agency: FTA

Agency: DHHS

Test Type: Breath Test

Lot Number: AG722602

Exp Date: 08/13/2009


Test	g/210L	Time
DIAG	Pass	1:28pm
AIR BLK	.00	1:29pm
ACCY CHK	.07	1:29pm
AIR BLK	.00	1:30pm
SUB TEST	.00	1:30pm
AIR BLK	.00	1:32pm
SUB TEST	.00	1:33pm
AIR BLK	.00	1:34pm

Reported AC: .00 g/210L



Signature of Chemical Analyst

Court CVR

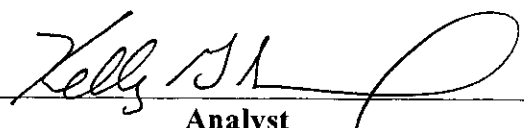

Analyst

Intox EC/IR-II: Preventive Maintenance

PITT COUNTY PITT CO DETENTION 730

Serial Number: 008668 Test Record Number: 669
Test Date: 11/25/2008 Test Time: 1:35pm EST

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	1:36pm
FLO	Pass	1:36pm
FC	Pass	1:36pm

Temperature Tests

Test	Status	Time
FC1	Pass	1:36pm
SRC	Pass	1:36pm
DET	Pass	1:36pm
BAR	Pass	1:36pm
BT	Pass	1:36pm

Blank Tests

Test	Status	Time
AIR	Pass	1:37pm

Printer Tests

Test	Status	Time
PRNT	Pass	1:37pm

CRC Tests

Test	Status	Time
COMP	Pass	1:37pm
CAL	Pass	1:37pm

Preventive Maintenance
Status: Pass



Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH

**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**


County Pitt Instrument Location Pitt Co. Detention Center

Instrument Serial No. 008662 124 Detention Dr., Greenville, N.C.

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 25th day of November, 2008 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

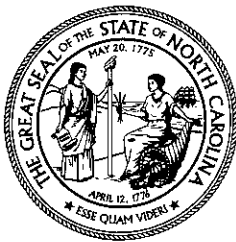

Kelly M. O
Signature of Certifying Official

643
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

PITT COUNTY PITT CO DETENTION 730

Serial Number: 008662

Test Date: 11/25/2008

Citation Number: M0000000-0

Subject's Name:

PREVENTIVE, MAINTENANCE

Subject's Date of Birth: 11/11/1911

Subject's Sex: Male

Driver's License State: XX

Driver's License Number: NONE

Analyst's Name: GUARD, KELLY G

Permit Number: 12955E

Effective:

12/01/2007-12/01/2009

Officer's Name: NONE, NONE

Type of Agency: FTA

Agency: DHHS


Test Type: Breath Test

Lot Number: AG722602

Exp Date: 08/13/2009


Test	g/210L	Time
DIAG	Pass	1:26pm
AIR BLK	.00	1:27pm
ACCY CHK	.07	1:28pm
AIR BLK	.00	1:29pm
SUB TEST	.00	1:29pm
AIR BLK	.00	1:30pm
SUB TEST	.00	1:31pm
AIR BLK	.00	1:32pm

Reported AC: .00 g/210L



Signature of Chemical Analyst

Court CVR

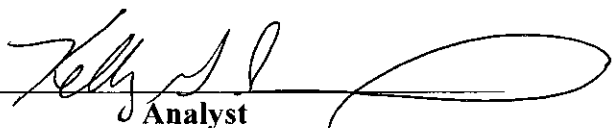

Analyst

Intox EC/IR-II: Preventive Maintenance

PITT COUNTY PITT CO DETENTION 730

Serial Number: 008662 Test Record Number: 269
Test Date: 11/25/2008 Test Time: 1:34pm EST

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	1:35pm
FLO	Pass	1:35pm
FC	Pass	1:35pm

Temperature Tests

Test	Status	Time
FC1	Pass	1:35pm
SRC	Pass	1:35pm
DET	Pass	1:35pm
BAR	Pass	1:35pm
BT	Pass	1:35pm

Blank Tests

Test	Status	Time
AIR	Pass	1:36pm


Printer Tests

Test	Status	Time
PRNT	Pass	1:36pm

CRC Tests

Test	Status	Time
COMP	Pass	1:36pm
CAL	Pass	1:36pm

Preventive Maintenance
Status: Pass

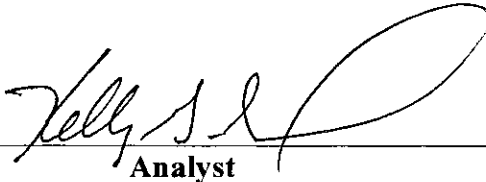

Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH

**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**


County Pitt Instrument Location Ayden Police Dept

Instrument Serial No. 0086666 4144 West Ave, Ayden, NC 28513

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 25th day of November, 2008 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

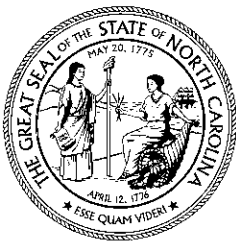

Kelly M. O
Signature of Certifying Official

643
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

PITT AYDEN PD 730

Serial Number: 008666
Test Date: 11/25/2008

Citation Number: M0000000-0
Subject's Name:
PREVENTIVE, MAINTENANCE
Subject's Date of Birth: 11/11/1911
Subject's Sex: Male
Driver's License State: XX
Driver's License Number: NONE


Analyst's Name: GUARD, KELLY G
Permit Number: 12955E
Effective:
12/01/2007-12/01/2009

Officer's Name: NONE, NONE
Type of Agency: FTA
Agency: DHHS
Test Type: Breath Test

Lot Number: AG722602
Exp Date: 08/13/2009

Test	g/210L	Time
DIAG	Pass	9:51am
AIR BLK	.00	9:52am
ACCY CHK	.07	9:52am
AIR BLK	.00	9:53am
SUB TEST	.00	9:54am
AIR BLK	.00	9:55am
SUB TEST	.00	9:56am
AIR BLK	.00	9:57am

Reported AC: .00 g/210L



Signature of Chemical Analyst

Court CVR



Analyst

Intox EC/IR-II: Preventive Maintenance

PITT AYDEN PD 730

Serial Number: 008666 Test Record Number: 261
Test Date: 11/25/2008 Test Time: 9:58am EST

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	9:59am
FLO	Pass	9:59am
FC	Pass	9:59am

Temperature Tests

Test	Status	Time
FC1	Pass	9:59am
SRC	Pass	9:59am
DET	Pass	9:59am
BAR	Pass	9:59am
BT	Pass	9:59am

Blank Tests

Test	Status	Time
AIR	Pass	9:59am


Printer Tests

Test	Status	Time
PRNT	Pass	10:00am

CRC Tests

Test	Status	Time
COMP	Pass	10:00am
CAL	Pass	10:00am

Preventive Maintenance
Status: Pass

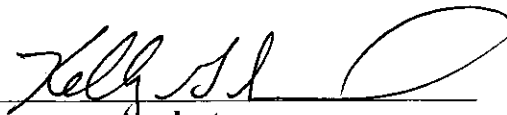

Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH

**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**


County Hyde Instrument Location Hyde Co. Sheriff's Office

Instrument Serial No. 008801 1223 Main St., Swan Quarter, N.C.

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 10th day of November, 2008 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

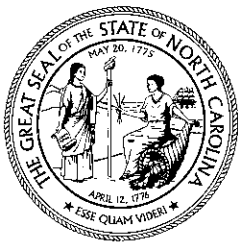

Kelly D. [Signature]
Signature of Certifying Official

1643
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

HYDE COUNTY HYDE CO SO SWAN QUAR 470

Serial Number: 008801
Test Date: 11/10/2008

Citation Number: M0000000-0

Subject's Name:

PREVENTIVE, MAINTENANCE

Subject's Date of Birth: 11/11/1911

Subject's Sex: Male

Driver's License State: XX

Driver's License Number: NONE

Analyst's Name: GUARD, KELLY G

Permit Number: 12955E

Effective:

12/01/2007-12/01/2009

Officer's Name: NONE, NONE

Type of Agency: FTA

Agency: DHHS


Test Type: Breath Test

Lot Number: AG722701

Exp Date: 08/14/2009

Test	g/210L	Time
DIAG	Pass	10:21am
AIR BLK	.00	10:22am
ACCY CHK	.08	10:23am
AIR BLK	.00	10:24am
SUB TEST	.00	10:24am
AIR BLK	.00	10:25am
SUB TEST	.00	10:27am
AIR BLK	.00	10:28am

Reported AC: .00 g/210L



Signature of Chemical Analyst

Court CVR



Analyst

**This form is used when performing Preventive Maintenance procedures
Forensic Tests for Alcohol Branch
Department of Health and Human Services
Rev. 12/2007**

Intox EC/IR-II: Preventive Maintenance

HYDE COUNTY HYDE CO SO SWAN QUAR 470

Serial Number: 008801 Test Record Number: 83
Test Date: 11/10/2008 Test Time: 10:29am EST

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	10:29am
FLO	Pass	10:29am
FC	Pass	10:29am

Temperature Tests

Test	Status	Time
FC1	Pass	10:30am
SRC	Pass	10:30am
DET	Pass	10:30am
BAR	Pass	10:30am
BT	Pass	10:30am

Blank Tests

Test	Status	Time
AIR	Pass	10:30am

Printer Tests

Test	Status	Time
PRNT	Pass	10:30am

CRC Tests

Test	Status	Time
COMP	Pass	10:30am
CAL	Pass	10:30am

Preventive Maintenance
Status: Pass



Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH


**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**

County Dare Instrument Location Dare Co. Detention Center
Instrument Serial No. 008804 1044 Driftwood Dr, Manteo, N.C.

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 4th day of November, 2008 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

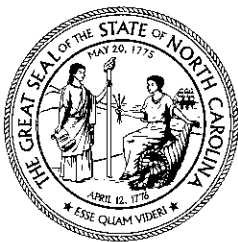

Kelly D. J.
Signature of Certifying Official

643
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

DARE COUNTY DARE CO DETENTION CE 270

Serial Number: 008804
Test Date: 11/04/2008

Citation Number: M0000000-0
Subject's Name:
PREVENTIVE, MAINTENANCE
Subject's Date of Birth: 11/11/1911
Subject's Sex: Male
Driver's License State: XX
Driver's License Number: NONE


Analyst's Name: GUARD, KELLY G
Permit Number: 12955E
Effective:
12/01/2007-12/01/2009

Officer's Name: NONE, NONE
Type of Agency: FTA
Agency: DHHS
Test Type: Breath Test

Lot Number: AG722802
Exp Date: 08/15/2009

Test	g/210L	Time
DIAG	Pass	2:28pm
AIR BLK	.00	2:29pm
ACCY CHK	.07	2:30pm
AIR BLK	.00	2:31pm
SUB TEST	.00	2:31pm
AIR BLK	.00	2:32pm
SUB TEST	.00	2:34pm
AIR BLK	.00	2:35pm

Reported AC: .00 g/210L



Signature of Chemical Analyst

Court CVR



Analyst

Intox EC/IR-II: Preventive Maintenance

DARE COUNTY DARE CO DETENTION CE 270

Serial Number: 008804 Test Record Number: 236
Test Date: 11/04/2008 Test Time: 2:36pm EST

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	2:36pm
FLO	Pass	2:36pm
FC	Pass	2:36pm

Temperature Tests

Test	Status	Time
FC1	Pass	2:36pm
SRC	Pass	2:36pm
DET	Pass	2:36pm
BAR	Pass	2:36pm
BT	Pass	2:36pm

Blank Tests

Test	Status	Time
AIR	Pass	2:37pm

Printer Tests

Test	Status	Time
PRNT	Pass	2:37pm

CRC Tests

Test	Status	Time
COMP	Pass	2:37pm
CAL	Pass	2:37pm

Preventive Maintenance
Status: Pass


Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH

**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**


County Dare Instrument Location Dare Co. Detention

Instrument Serial No. 008783 1044 Driftwood Dr., Manteo, N.C.

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 4th day of November, 2008 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

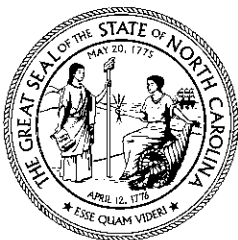

Kelly A. L. P.
Signature of Certifying Official

643
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

DARE COUNTY DARE CO DETENTION CE 270

Serial Number: 008783
Test Date: 11/04/2008

Citation Number: M0000000-0
Subject's Name:
PREVENTIVE, MAINTENANCE
Subject's Date of Birth: 11/11/1911
Subject's Sex: Male
Driver's License State: XX
Driver's License Number: NONE


Analyst's Name: GUARD, KELLY G
Permit Number: 12955E
Effective:
12/01/2007-12/01/2009

Officer's Name: NONE, NONE
Type of Agency: FTA
Agency: DHHS
Test Type: Breath Test

Lot Number: AG722701
Exp Date: 08/14/2009

Test	g/210L	Time
DIAG	Pass	2:27pm
AIR BLK	.00	2:28pm
ACCY CHK	.07	2:29pm
AIR BLK	.00	2:30pm
SUB TEST	.00	2:30pm
AIR BLK	.00	2:31pm
SUB TEST	.00	2:32pm
AIR BLK	.00	2:33pm

Reported AC: .00 g/210L

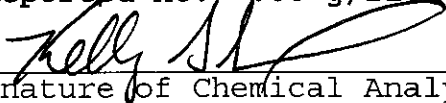

Signature of Chemical Analyst

Court CVR



Analyst

Intox EC/IR-II: Preventive Maintenance

DARE COUNTY DARE CO DETENTION CE 270

Serial Number: 008783 Test Record Number: 83
Test Date: 11/04/2008 Test Time: 2:35pm EST

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	2:35pm
FLO	Pass	2:35pm
FC	Pass	2:35pm

Temperature Tests

Test	Status	Time
FC1	Pass	2:35pm
SRC	Pass	2:35pm
DET	Pass	2:35pm
BAR	Pass	2:35pm
BT	Pass	2:35pm

Blank Tests

Test	Status	Time
AIR	Pass	2:36pm


Printer Tests

Test	Status	Time
PRNT	Pass	2:36pm

CRC Tests

Test	Status	Time
COMP	Pass	2:36pm
CAL	Pass	2:36pm

Preventive Maintenance
Status: Pass

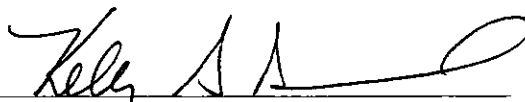

Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH

**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**

County Dare


Instrument Location Kill Devil Hills P.D.

Instrument Serial No. 008844 102, Town Hall Dr, Kill Devil Hills, N.C.

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 4th day of November, 2008 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

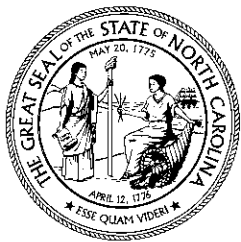

Kelly A. L. O.
Signature of Certifying Official

643
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

DARE COUNTY KILL DEVIL HILLS PD 270

Serial Number: 008844
Test Date: 11/04/2008

Citation Number: M0000000-0
Subject's Name:
PREVENTIVE, MAINTENANCE
Subject's Date of Birth: 11/11/1911
Subject's Sex: Male
Driver's License State: XX
Driver's License Number: NONE


Analyst's Name: GUARD, KELLY G
Permit Number: 12955E
Effective:
12/01/2007-12/01/2009

Officer's Name: NONE, NONE
Type of Agency: FTA
Agency: DHHS
Test Type: Breath Test

Lot Number: AG722701
Exp Date: 08/14/2009


Test	g/210L	Time
DIAG	Pass	1:46pm
AIR BLK	.00	1:47pm
ACCY CHK	.07	1:48pm
AIR BLK	.00	1:48pm
SUB TEST	.00	1:49pm
AIR BLK	.00	1:50pm
SUB TEST	.00	1:51pm
AIR BLK	.00	1:52pm

Reported AC: .00 g/210L

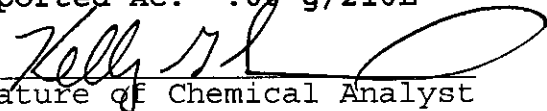

Signature of Chemical Analyst

Court CVR

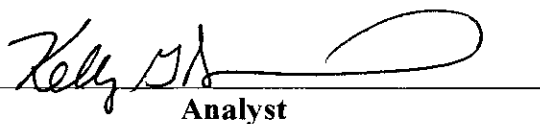

Analyst

Intox EC/IR-II: Preventive Maintenance

DARE COUNTY KILL DEVIL HILLS PD 270

Serial Number: 008844 Test Record Number: 216
Test Date: 11/04/2008 Test Time: 1:54pm EST

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	1:54pm
FLO	Pass	1:54pm
FC	Pass	1:54pm

Temperature Tests

Test	Status	Time
FC1	Pass	1:54pm
SRC	Pass	1:54pm
DET	Pass	1:54pm
BAR	Pass	1:54pm
BT	Pass	1:54pm

Blank Tests

Test	Status	Time
AIR	Pass	1:55pm


Printer Tests

Test	Status	Time
PRNT	Pass	1:55pm

CRC Tests

Test	Status	Time
COMP	Pass	1:55pm
CAL	Pass	1:55pm

Preventive Maintenance
Status: Pass



Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH

**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**


County Caldwell Instrument Location Caldwell Co. Jail

Instrument Serial No. 108803 Verat, NC

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 26 day of November, 2008 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

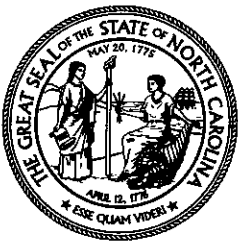

[Signature]
Signature of Certifying Official

649
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

CALDWELL COUNTY CALDWELL COUNTY JAIL
130

Serial Number: 008803
Test Date: 11/26/2008

Citation Number: M0000000-0
Subject's Name:
PREVENTIVE, MAINTENANCE
Subject's Date of Birth: 11/11/1911
Subject's Sex: Male
Driver's License State: XX
Driver's License Number: NONE

Analyst's Name: BURNETTE, ANTHONY J
Permit Number: 11304E
Effective:
12/01/2007-12/01/2009

Officer's Name: NONE,
Type of Agency: FTA
Agency: DHHS
Test Type: Breath Test

Lot Number: AG722502
Exp Date: 08/12/2009

Test	g/210L	Time
DIAG	Pass	4:56pm
AIR BLK	.00	4:57pm
ACCY CHK	.08	4:57pm
AIR BLK	.00	4:58pm
SUB TEST	.00	4:59pm
AIR BLK	.00	4:59pm
SUB TEST	.00	5:01pm
AIR BLK	.00	5:02pm

Reported AC: .00 g/210L

Signature of Chemical Analyst

Court CVR



Analyst

**This form is used when performing Preventive Maintenance procedures
Forensic Tests for Alcohol Branch
Department of Health and Human Services
Rev. 12/2007**

Intox EC/IR-II: Preventive Maintenance

CALDWELL COUNTY CALDWELL COUNTY JAIL 130

Serial Number: 008803 Test Record Number: 111
Test Date: 11/26/2008 Test Time: 5:03pm EST

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	5:03pm
FLO	Pass	5:03pm
FC	Pass	5:03pm

Temperature Tests

Test	Status	Time
FC1	Pass	5:03pm
SRC	Pass	5:03pm
DET	Pass	5:03pm
BAR	Pass	5:03pm
BT	Pass	5:03pm

Blank Tests

Test	Status	Time
AIR	Pass	5:04pm

Printer Tests

Test	Status	Time
PRNT	Pass	5:04pm

CRC Tests

Test	Status	Time
COMP	Pass	5:04pm
CAL	Pass	5:04pm

Preventive Maintenance
Status: Pass


Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH

**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**


County Caldwell Instrument Location Caldwell Co Jail

Instrument Serial No. 108719 66001, MS

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 26 day of November, 2008 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

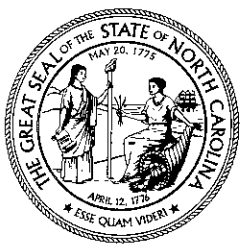

[Signature]
Signature of Certifying Official

649
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

CALDWELL COUNTY CALDWELL COUNTY JAIL
130

Serial Number: 008719
Test Date: 11/26/2008

Citation Number: M0000000-0
Subject's Name:
PREVENTIVE, MAINTENANCE
Subject's Date of Birth: 11/11/1911
Subject's Sex: Male
Driver's License State: XX
Driver's License Number: NONE

Analyst's Name: BURNETTE, ANTHONY J
Permit Number: 11304E
Effective:
12/01/2007-12/01/2009

Officer's Name: NONE,
Type of Agency: FTA
Agency: DHHS
Test Type: Breath Test

Lot Number: AG722503
Exp Date: 08/12/2009

Test	g/210L	Time
DIAG	Pass	4:50pm
AIR BLK	.00	4:51pm
ACCY CHK	.07	4:51pm
AIR BLK	.00	4:52pm
SUB TEST	.00	4:53pm
AIR BLK	.00	4:54pm
SUB TEST	.00	4:55pm
AIR BLK	.00	4:56pm

Reported AC: .00 g/210L

Signature of Chemical Analyst

Court CVR


Analyst

Intox EC/IR-II: Preventive Maintenance

CALDWELL COUNTY CALDWELL COUNTY JAIL 130

Serial Number: 008719 Test Record Number: 270
Test Date: 11/26/2008 Test Time: 4:56pm EST

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	4:57pm
FLO	Pass	4:57pm
FC	Pass	4:57pm

Temperature Tests

Test	Status	Time
FC1	Pass	4:57pm
SRC	Pass	4:57pm
DET	Pass	4:57pm
BAR	Pass	4:57pm
BT	Pass	4:57pm

Blank Tests

Test	Status	Time
AIR	Pass	4:58pm

Printer Tests

Test	Status	Time
PRNT	Pass	4:58pm

CRC Tests

Test	Status	Time
COMP	Pass	4:58pm
CAL	Pass	4:58pm

Preventive Maintenance
Status: Pass


Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH

**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**


County Madison Instrument Location Moss Hill PD

Instrument Serial No. 008594 Moss Hill, NC

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 12 day of November, 2008 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

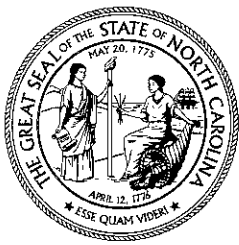

[Signature]
Signature of Certifying Official

649
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

MADISON COUNTY MARS HILL PD 560

Serial Number: 008599
Test Date: 11/12/2008

Citation Number: M0000000-0
Subject's Name:

PREVENTIVE, MAINTENANCE

Subject's Date of Birth: 11/11/1911

Subject's Sex: Male

Driver's License State: XX

Driver's License Number: NONE

Analyst's Name: BURNETTE, ANTHONY J

Permit Number: 11304E

Effective:

12/01/2007-12/01/2009

Officer's Name: NONE,

Type of Agency: FTA

Agency: DHHS

Test Type: Breath Test

Lot Number: AG722502


Exp Date: 08/12/2009

Test	g/210L	Time
DIAG	Pass	6:17pm
AIR BLK	.00	6:18pm
ACCY CHK	.07	6:19pm
AIR BLK	.00	6:19pm
SUB TEST	.00	6:20pm
AIR BLK	.00	6:21pm
SUB TEST	.00	6:22pm
AIR BLK	.00	6:23pm

Reported AC: .00 g/210L

Signature of Chemical Analyst

Court CVR

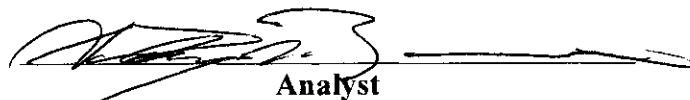
Analyst

Intox EC/IR-II: Preventive Maintenance

MADISON COUNTY MARS HILL PD 560

Serial Number: 008599 Test Record Number: 292
Test Date: 11/12/2008 Test Time: 6:24pm EST

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	6:25pm
FLO	Pass	6:25pm
FC	Pass	6:25pm

Temperature Tests

Test	Status	Time
FC1	Pass	6:25pm
SRC	Pass	6:25pm
DET	Pass	6:25pm
BAR	Pass	6:25pm
BT	Pass	6:25pm

Blank Tests

Test	Status	Time
AIR	Pass	6:26pm

Printer Tests

Test	Status	Time
PRNT	Pass	6:26pm

CRC Tests

Test	Status	Time
COMP	Pass	6:26pm
CAL	Pass	6:26pm

Preventive Maintenance
Status: Pass


Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH

**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**


County Madison Instrument Location Mars Hill P.D.

Instrument Serial No. 008582 Mars Hill, NC

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 12 day of November, 2008 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

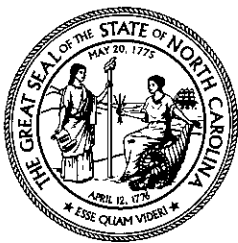

[Signature]
Signature of Certifying Official

649
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

MADISON COUNTY MARS HILL PD 560

Serial Number: 008582
Test Date: 11/12/2008

Citation Number: M0000000-0
Subject's Name:
PREVENTIVE, MAINTENANCE
Subject's Date of Birth: 11/11/1911
Subject's Sex: Male
Driver's License State: XX
Driver's License Number: NONE

Analyst's Name: BURNETTE, ANTHONY J
Permit Number: 11304E
Effective:
12/01/2007-12/01/2009

Officer's Name: NONE,
Type of Agency: FTA
Agency: DHHS
Test Type: Breath Test

Lot Number: AG723301
Exp Date: 08/20/2009

Test	g/210L	Time
DIAG	Pass	6:16pm
AIR BLK	.00	6:17pm
ACCY CHK	.07	6:18pm
AIR BLK	.00	6:19pm
SUB TEST	.00	6:20pm
AIR BLK	.00	6:21pm
SUB TEST	.00	6:22pm
AIR BLK	.00	6:23pm

Reported AC: .00 g/210L

Signature of Chemical Analyst

Court CVR


Analyst

Intox EC/IR-II: Preventive Maintenance

MADISON COUNTY MARS HILL PD 560

Serial Number: 008582 Test Record Number: 490
Test Date: 11/12/2008 Test Time: 6:24pm EST

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	6:24pm
FLO	Pass	6:24pm
FC	Pass	6:24pm

Temperature Tests

Test	Status	Time
FC1	Pass	6:24pm
SRC	Pass	6:24pm
DET	Pass	6:24pm
BAR	Pass	6:24pm
BT	Pass	6:24pm

Blank Tests

Test	Status	Time
AIR	Pass	6:25pm

Printer Tests

Test	Status	Time
PRNT	Pass	6:25pm

CRC Tests

Test	Status	Time
COMP	Pass	6:25pm
CAL	Pass	6:25pm

Preventive Maintenance
Status: Pass


Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH

**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**


County Forsyth Instrument Location Kernersville P. D.

Instrument Serial No. 008650

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 25 day of November, 2008 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

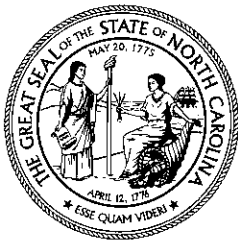

[Signature]
Signature of Certifying Official

682
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

FORSYTH COUNTY KERNERSVILLE PD 330

Serial Number: 008650

Test Date: 11/25/2008

Citation Number: M0000000-0

Subject's Name:

PREVENTIVE, MAINTENANCE

Subject's Date of Birth: 11/11/1911

Subject's Sex: Male

Driver's License State: XX

Driver's License Number: NONE

Analyst's Name: WEAVER, GEORGE A

Permit Number: 09442E

Effective:

12/01/2007-12/01/2009

Officer's Name: NONE, NONE

Type of Agency: FTA

Agency: DHHS


Test Type: Breath Test

Lot Number: AG722501-07

Exp Date: 08/12/2009

Test	g/210L	Time
DIAG	Pass	11:00am
AIR BLK	.00	11:01am
ACCY CHK	.07	11:01am
AIR BLK	.00	11:02am
SUB TEST	.00	11:03am
AIR BLK	.00	11:04am
SUB TEST	.00	11:06am
AIR BLK	.00	11:07am

Reported AC: .00 g/210L

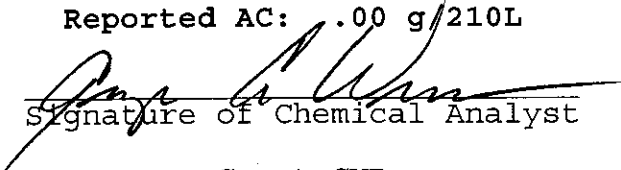
Signature of Chemical Analyst

Court CVR


Analyst

Intox EC/IR-II: Preventive Maintenance

FORSYTH COUNTY KERNERSVILLE PD 330

Serial Number: 008650 Test Record Number: 262

Test Date: 11/25/2008 Test Time: 11:08am EST

System Check: Passed

Baseline Tests

Test	Status	Time
IR	Pass	11:08am
FLO	Pass	11:08am
FC	Pass	11:08am

Temperature Tests

Test	Status	Time
FC1	Pass	11:08am
SRC	Pass	11:08am
DET	Pass	11:08am
BAR	Pass	11:08am
BT	Pass	11:08am

Blank Tests

Test	Status	Time
AIR	Pass	11:09am


Printer Tests

Test	Status	Time
PRNT	Pass	11:09am

CRC Tests

Test	Status	Time
COMP	Pass	11:09am
CAL	Pass	11:09am

Preventive Maintenance
Status: Pass


Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH

**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**


County Ashe Instrument Location Ashe Co Jail

Instrument Serial No. 008349

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 24 day of November, 20 08 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

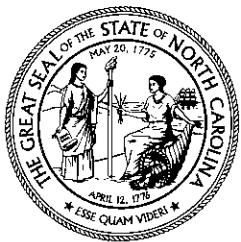

[Signature]
Signature of Certifying Official

632
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

ASHE COUNTY ASHE COUNTY JAIL 040

Serial Number: 008849
Test Date: 11/24/2008

Citation Number: M0000000-0
Subject's Name:

PREVENTIVE, MAINTENANCE

Subject's Date of Birth: 11/11/1911

Subject's Sex: Male

Driver's License State: XX

Driver's License Number: NONE

Analyst's Name: WEAVER, GEORGE A

Permit Number: 09442E

Effective:

12/01/2007-12/01/2009

Officer's Name: NONE, NONE

Type of Agency: FTA

Agency: DHHS


Test Type: Breath Test

Lot Number: AG816302

Exp Date: 06/11/2010

Test	g/210L	Time
DIAG	Pass	12:14pm
AIR BLK	.00	12:15pm
ACCY CHK	.07	12:16pm
AIR BLK	.00	12:16pm
SUB TEST	.00	12:17pm
AIR BLK	.00	12:18pm
SUB TEST	.00	12:20pm
AIR BLK	.00	12:21pm

Reported AC: .00 g/210L

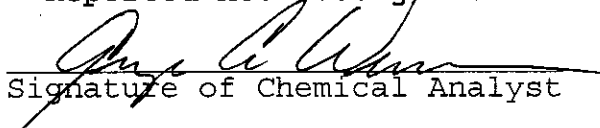
Signature of Chemical Analyst

Court CVR


Analyst

**This form is used when performing Preventive Maintenance procedures
Forensic Tests for Alcohol Branch
Department of Health and Human Services
Rev. 12/2007**

Intox EC/IR-II: Preventive Maintenance

ASHE COUNTY ASHE COUNTY JAIL 040

Serial Number: 008849 Test Record Number: 129
Test Date: 11/24/2008 Test Time: 12:22pm EST

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	12:22pm
FLO	Pass	12:22pm
FC	Pass	12:22pm

Temperature Tests

Test	Status	Time
FC1	Pass	12:23pm
SRC	Pass	12:23pm
DET	Pass	12:23pm
BAR	Pass	12:23pm
BT	Pass	12:23pm

Blank Tests

Test	Status	Time
AIR	Pass	12:23pm


Printer Tests

Test	Status	Time
PRNT	Pass	12:23pm

CRC Tests

Test	Status	Time
COMP	Pass	12:23pm
CAL	Pass	12:23pm

Preventive Maintenance
Status: Pass

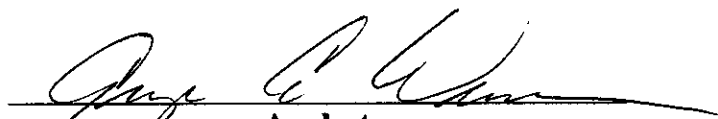
Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH


**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**

County Allegany Instrument Location Allegany Co. Jail
Instrument Serial No. 008890

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 24 day of November, 2008 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

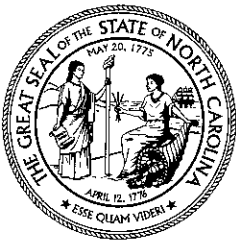

[Signature]
Signature of Certifying Official

632
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

ALLEGHANY COUNTY ALLEGHANY CO JAIL 020

Serial Number: 008890
Test Date: 11/24/2008

Citation Number: M0000000-0
Subject's Name:
PREVENTIVE, MAINTENANCE
Subject's Date of Birth: 11/11/1911
Subject's Sex: Male
Driver's License State: XX
Driver's License Number: NONE

Analyst's Name: WEAVER, GEORGE A
Permit Number: 09442E
Effective:
12/01/2007-12/01/2009

Officer's Name: NONE, NONE
Type of Agency: FTA
Agency: DHHS
Test Type: Breath Test

Lot Number: AG816303
Exp Date: 06/11/2010

Test	g/210L	Time
DIAG	Pass	1:30pm
AIR BLK	.00	1:31pm
ACCY CHK	.07	1:31pm
AIR BLK	.00	1:32pm
SUB TEST	.00	1:33pm
AIR BLK	.00	1:33pm
SUB TEST	.00	1:35pm
AIR BLK	.00	1:36pm

Reported AC: .00/ g/210L


Signature of Chemical Analyst

Court CVR


Analyst

**This form is used when performing Preventive Maintenance procedures
Forensic Tests for Alcohol Branch
Department of Health and Human Services
Rev. 12/2007**

Intox EC/IR-II: Preventive Maintenance

ALLEGHANY COUNTY ALLEGHANY CO JAIL 020

Serial Number: 008890 Test Record Number: 70
Test Date: 11/24/2008 Test Time: 1:37pm EST

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	1:37pm
FLO	Pass	1:37pm
FC	Pass	1:37pm

Temperature Tests

Test	Status	Time
FC1	Pass	1:37pm
SRC	Pass	1:37pm
DET	Pass	1:37pm
BAR	Pass	1:37pm
BT	Pass	1:37pm

Blank Tests

Test	Status	Time
AIR	Pass	1:38pm


Printer Tests

Test	Status	Time
PRNT	Pass	1:38pm

CRC Tests

Test	Status	Time
COMP	Pass	1:38pm
CAL	Pass	1:38pm

Preventive Maintenance
Status: Pass

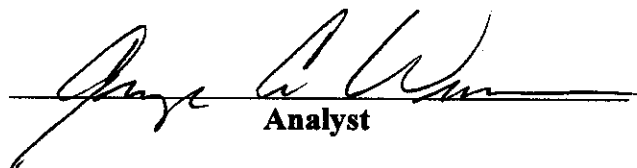
Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH

**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**


County Wilkes Instrument Location N. Wilkesboro P.O.

Instrument Serial No. 008862

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 21 day of November, 2008 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

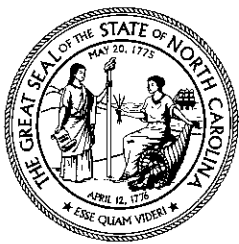

[Signature]
Signature of Certifying Official

Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

WILKES COUNTY NORTH WILKESBORO PD 960

Serial Number: 008862
Test Date: 11/21/2008

Citation Number: M0000000-0
Subject's Name:
PREVENTIVE, MAINTENANCE
Subject's Date of Birth: 11/11/1911
Subject's Sex: Male
Driver's License State: XX
Driver's License Number: NONE


Analyst's Name: WEAVER, GEORGE A
Permit Number: 09442E
Effective:
12/01/2007-12/01/2009

Officer's Name: NONE, NONE
Type of Agency: FTA
Agency: DHHS
Test Type: Breath Test

Lot Number: AG816303
Exp Date: 06/11/2010

Test	g/210L	Time
DIAG	Pass	2:45pm
AIR BLK	.00	2:46pm
ACCY CHK	.07	2:46pm
AIR BLK	.00	2:47pm
SUB TEST	.00	2:48pm
AIR BLK	.00	2:49pm
SUB TEST	.00	2:50pm
AIR BLK	.00	2:51pm

Reported AC: .00 g/210L

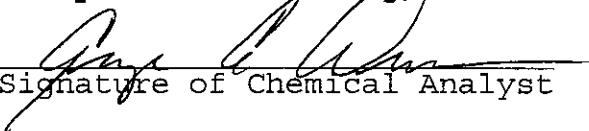
Signature of Chemical Analyst

Court CVR


Analyst

Intox EC/IR-II: Preventive Maintenance

WILKES COUNTY NORTH WILKESBORO PD 960

Serial Number: 008862 Test Record Number: 57
Test Date: 11/21/2008 Test Time: 2:53pm EST

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	2:53pm
FLO	Pass	2:53pm
FC	Pass	2:54pm

Temperature Tests

Test	Status	Time
FC1	Pass	2:54pm
SRC	Pass	2:54pm
DET	Pass	2:54pm
BAR	Pass	2:54pm
BT	Pass	2:54pm

Blank Tests

Test	Status	Time
AIR	Pass	2:54pm

Printer Tests

Test	Status	Time
PRNT	Pass	2:54pm

CRC Tests

Test	Status	Time
COMP	Pass	2:54pm
CAL	Pass	2:54pm

Preventive Maintenance
Status: Pass


Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH

**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**

County Wilkes


Instrument Location Wilkes Co Court House

Instrument Serial No. 008843

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 21 day of November, 2008 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

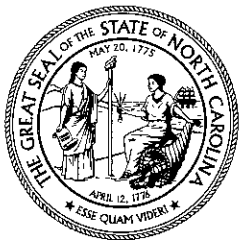

[Signature]
Signature of Certifying Official

632
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

WILKES COUNTY WILKES CO COURTHOUSE 960

Serial Number: 008843
Test Date: 11/21/2008

Citation Number: M0000000-0
Subject's Name:
PREVENTIVE, MAINTENANCE
Subject's Date of Birth: 11/11/1911
Subject's Sex: Male
Driver's License State: XX
Driver's License Number: NONE


Analyst's Name: WEAVER, GEORGE A
Permit Number: 09442E
Effective:
12/01/2007-12/01/2009

Officer's Name: NONE, NONE
Type of Agency: FTA
Agency: DHHS
Test Type: Breath Test

Lot Number: AG816303
Exp Date: 06/11/2010

Test	g/210L	Time
DIAG	Pass	1:55pm
AIR BLK	.00	1:56pm
ACCY CHK	.07	1:57pm
AIR BLK	.00	1:58pm
SUB TEST	.00	1:58pm
AIR BLK	.00	1:59pm
SUB TEST	.00	2:01pm
AIR BLK	.00	2:02pm

Reported AC: .00 g/210L

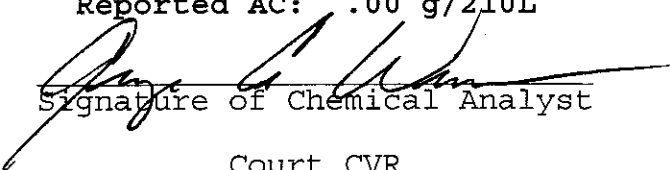
Signature of Chemical Analyst

Court CVR


Analyst

**This form is used when performing Preventive Maintenance procedures
Forensic Tests for Alcohol Branch
Department of Health and Human Services
Rev. 12/2007**

Intox EC/IR-II: Subject Test

WILKES COUNTY WILKES CO COURTHOUSE 960

Serial Number: 008843
Test Date: 11/21/2008

Citation Number: M0000000-0
Subject's Name:

PREVENTIVE, MAINTENANCE

Subject's Date of Birth: 11/11/1911

Subject's Sex: Male

Driver's License State: XX

Driver's License Number: NONE

Analyst's Name: WEAVER, GEORGE A

Permit Number: 09442E

Effective:

12/01/2007-12/01/2009

Officer's Name: NONE, NONE

Type of Agency: FTA


Agency: DHHS

Test Type: Breath Test

Lot Number: AG816303
Exp Date: 06/11/2010

Test	g/210L	Time
DIAG	Pass	1:55pm
AIR BLK	.00	1:56pm
ACCY CHK	.07	1:57pm
AIR BLK	.00	1:58pm
SUB TEST	.00	1:58pm
AIR BLK	.00	1:59pm
SUB TEST	.00	2:01pm
AIR BLK	.00	2:02pm

Reported AC: .00 g/210L

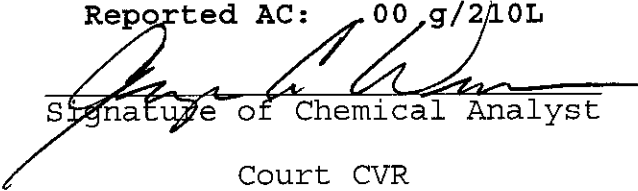
Signature of Chemical Analyst

Court CVR


Analyst

**This form is used when performing Preventive Maintenance procedures
Forensic Tests for Alcohol Branch
Department of Health and Human Services
Rev. 12/2007**

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH

**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**

County Clatsop


Instrument Location BAT Mobile Unit 4

Instrument Serial No. 008734

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 21st day of November, 2008 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

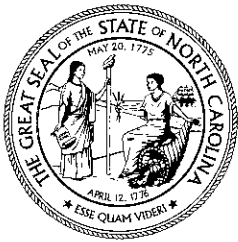

[Signature]
Signature of Certifying Official

657
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Preventive Maintenance

CHATHAM COUNTY BAT MOBILE UNIT 4 180

Serial Number: 008734 Test Record Number: 167
Test Date: 11/21/2008 Test Time: 10:10pm EST

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	10:10pm
FLO	Pass	10:10pm
FC	Pass	10:11pm

Temperature Tests

Test	Status	Time
FC1	Pass	10:11pm
SRC	Pass	10:11pm
DET	Pass	10:11pm
BAR	Pass	10:11pm
BT	Pass	10:11pm

Blank Tests

Test	Status	Time
AIR	Pass	10:11pm


Printer Tests

Test	Status	Time
PRNT	Pass	10:11pm

CRC Tests

Test	Status	Time
COMP	Pass	10:11pm
CAL	Pass	10:11pm

Preventive Maintenance
Status: Pass

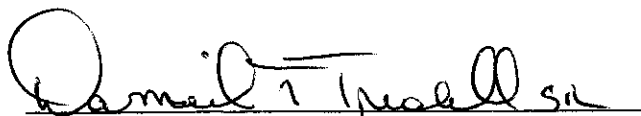
Analyst

Intox EC/IR-II: Subject Test

CHATHAM COUNTY BAT MOBILE UNIT 4 180

Serial Number: 008734
Test Date: 11/21/2008

Citation Number: M0000000-0
Subject's Name:
PREVENTIVE, MAINTENANCE
Subject's Date of Birth: 11/11/1911
Subject's Sex: Male
Driver's License State: XX
Driver's License Number: NONE


Analyst's Name: TRUDELL, SR., DANIEL T
Permit Number: 21535E
Effective:
12/01/2007-12/01/2009

Officer's Name: NONE, NONE
Type of Agency: FTA
Agency: DHHS
Test Type: Breath Test


Lot Number: AG722501
Exp Date: 08/12/2009

Test	g/210L	Time
DIAG	Pass	10:02pm
AIR BLK	.00	10:03pm
ACCY CHK	.07	10:03pm
AIR BLK	.00	10:04pm
SUB TEST	.00	10:05pm
AIR BLK	.00	10:06pm
SUB TEST	.00	10:07pm
AIR BLK	.00	10:08pm

Reported AC: .00 g/210L

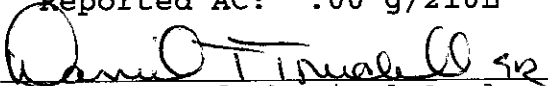
Signature of Chemical Analyst

Court CVR

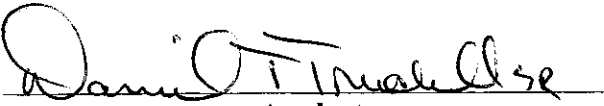
Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH

PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II


County MECKLENBURG Instrument Location BAT MOBILE UNIT 3

Instrument Serial No. 008616 CHARLOTTE, NC

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 20 day of NOVEMBER, 2008 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

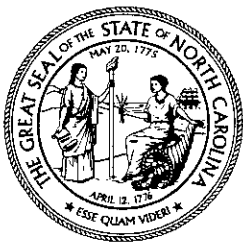

Alvin R. Bunch
Signature of Certifying Official

648
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

~~Intox EC/IR-II: Subject Test~~

MECKLENBURG COUNTY BAT MOBILE UNIT 3
590

Serial Number: 008616
Test Date: 11/20/2008

Citation Number: M0000000-0
Subject's Name:
PREVENTIVE, MAINTENANCE
Subject's Date of Birth: 11/11/1911
Subject's Sex: Male
Driver's License State: XX
Driver's License Number: NONE

Analyst's Name: BARNES, ALVIN R
Permit Number: 15671E
Effective:
12/01/2007-12/01/2009

Officer's Name: NONE, NONE
Type of Agency: FTA
Agency: DHHS
Test Type: Breath Test

Lot Number: AG722601
Exp Date: 08/13/2009

Test	g/210L	Time
DIAG	Pass	6:46pm
AIR BLK	.00	6:47pm
ACCY CHK	.07	6:47pm
AIR BLK	.00	6:48pm
SUB TEST	.00	6:49pm
AIR BLK	.00	6:50pm
SUB TEST	.00	6:51pm
AIR BLK	.00	6:52pm

Reported AC: .00 g/210L

Signature of Chemical Analyst

Court CVR


Analyst

This form is used when performing Preventive Maintenance procedures
Forensic Tests for Alcohol Branch
Department of Health and Human Services
Rev. 12/2007

Intox EC/IR-II: Preventive Maintenance

MECKLENBURG COUNTY BAT MOBILE UNIT 3 590

Serial Number: 008616 Test Record Number: 356
Test Date: 11/20/2008 Test Time: 6:53pm EST

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	6:54pm
FLO	Pass	6:54pm
FC	Pass	6:54pm

Temperature Tests

Test	Status	Time
FC1	Pass	6:54pm
SRC	Pass	6:54pm
DET	Pass	6:54pm
BAR	Pass	6:54pm
BT	Pass	6:54pm

Blank Tests

Test	Status	Time
AIR	Pass	6:55pm


Printer Tests

Test	Status	Time
PRNT	Pass	6:55pm

CRC Tests

Test	Status	Time
COMP	Pass	6:55pm
CAL	Pass	6:55pm

Preventive Maintenance
Status: Pass

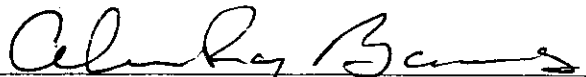
Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH

PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II


County MOORE Instrument Location BAT MOBILE UNIT 3

Instrument Serial No. 008616 SOUTHERN PINES, NC

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 22 day of NOVEMBER, 2008 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

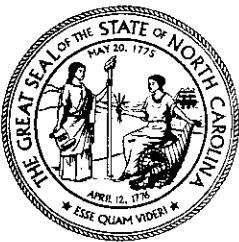

Alvin R. Benz
Signature of Certifying Official

648
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

~~Intox EC/IR-II: Subject Test~~

~~MOORE COUNTY BAT MOBILE UNIT 3 620~~

Serial Number: 008616
Test Date: 11/22/2008

Citation Number: M0000000-0
Subject's Name:

PREVENTIVE, MAINTENANCE

Subject's Date of Birth: 11/11/1911

Subject's Sex: Male

Driver's License State: XX

Driver's License Number: NONE

Analyst's Name: BARNES, ALVIN R

Permit Number: 15671E

Effective:

12/01/2007-12/01/2009

Officer's Name: NONE, NONE

Type of Agency: FTA

Agency: DHHS

Test Type: Breath Test

Lot Number: AG722601


Exp Date: 08/13/2009

Test	g/210L	Time
DIAG	Pass	8:03pm
AIR BLK	.00	8:04pm
ACCY CHK	.07	8:04pm
AIR BLK	.00	8:05pm
SUB TEST	.00	8:06pm
AIR BLK	.00	8:06pm
SUB TEST	.00	8:08pm
AIR BLK	.00	8:09pm

Reported AC: .00 g/210L

Signature of Chemical Analyst

Court CVR

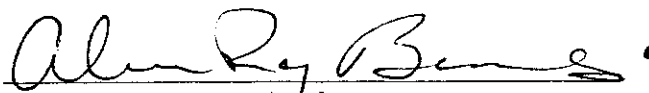
Analyst

Intox-EC/IR-II: Preventive Maintenance

MOORE COUNTY BAT MOBILE UNIT 3 620

Serial Number: 008616 Test Record Number: 362
Test Date: 11/22/2008 Test Time: 8:10pm EST

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	8:10pm
FLO	Pass	8:10pm
FC	Pass	8:10pm

Temperature Tests

Test	Status	Time
FC1	Pass	8:10pm
SRC	Pass	8:10pm
DET	Pass	8:10pm
BAR	Pass	8:10pm
BT	Pass	8:10pm

Blank Tests

Test	Status	Time
AIR	Pass	8:11pm


Printer Tests

Test	Status	Time
PRNT	Pass	8:11pm

CRC Tests

Test	Status	Time
COMP	Pass	8:11pm
CAL	Pass	8:11pm

Preventive Maintenance
Status: Pass

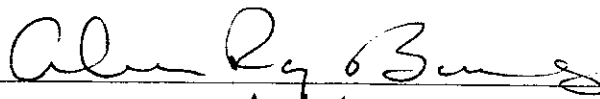
Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH


**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**

County Cleveland Instrument Location Kings Mtn. P. D.
Instrument Serial No. 008900 112 S. Piedmont Ave, Kings Mtn.
704 - 734 - 0444

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 13 day of November, 2008 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

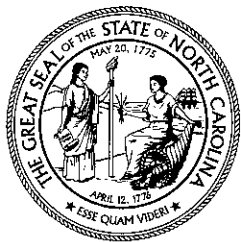

Bobby L. Willis
Signature of Certifying Official

557
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

CLEVELAND COUNTY KINGS MOUNTAIN PD 220

Serial Number: 008900
Test Date: 11/13/2008

Citation Number: M0000000-0
Subject's Name:
PREVENTIVE, MAINTENANCE
Subject's Date of Birth: 11/11/1911
Subject's Sex: Male
Driver's License State: XX
Driver's License Number: NONE

Analyst's Name: WILLIS, BOBBY D
Permit Number: 08010E
Effective:
12/01/2007-12/01/2009

Officer's Name: NONE, NONE
Type of Agency: FTA
Agency: DHHS
Test Type: Breath Test

Lot Number: AG816303
Exp Date: 06/11/2010

Test	g/210L	Time
DIAG	Pass	12:59pm
AIR BLK	.00	1:00pm
ACCY CHK	.07	1:00pm
AIR BLK	.00	1:01pm
SUB TEST	.00	1:02pm
AIR BLK	.00	1:03pm
SUB TEST	.00	1:04pm
AIR BLK	.00	1:05pm

Reported AC: .00 g/210L

Bobby D. Willis
Signature of Chemical Analyst

Court CVR.

Bobby D. Willis
Analyst

Intox EC/IR-II: Preventive Maintenance

CLEVELAND COUNTY KINGS MOUNTAIN PD 220

Serial Number: 008900 Test Record Number: 87
Test Date: 11/13/2008 Test Time: 1:07pm EST

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	1:07pm
FLO	Pass	1:07pm
FC	Pass	1:07pm

Temperature Tests

Test	Status	Time
FC1	Pass	1:07pm
SRC	Pass	1:07pm
DET	Pass	1:07pm
BAR	Pass	1:07pm
BT	Pass	1:07pm

Blank Tests

Test	Status	Time
AIR	Pass	1:08pm

Printer Tests

Test	Status	Time
PRNT	Pass	1:08pm

CRC Tests

Test	Status	Time
COMP	Pass	1:08pm
CAL	Pass	1:08pm

Preventive Maintenance
Status: Pass

Bobby D. Willis
Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH

**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**


County Cleveland Instrument Location Cleveland Co. S.D.

Instrument Serial No. 008893 100 Justice Pl., Shelby
704-484-4888

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 13 day of November, 2008 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

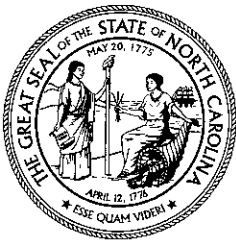

Betsy L. Willis
Signature of Certifying Official

5571
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

CLEVELAND COUNTY CLEVELAND COUNTY SD
220

Serial Number: 008893
Test Date: 11/13/2008

Citation Number: M0000000-0

Subject's Name:

PREVENTIVE, MAINTENANCE

Subject's Date of Birth: 11/11/1911

Subject's Sex: Male

Driver's License State: XX

Driver's License Number: NONE

Analyst's Name: WILLIS, BOBBY D

Permit Number: 08010E

Effective:

12/01/2007-12/01/2009

Officer's Name: NONE, NONE

Type of Agency: FTA

Agency: DHHS

Test Type: Breath Test

Lot Number: AG816303

Exp Date: 06/11/2010

Test	g/210L	Time
DIAG	Pass	10:33am
AIR BLK	.00	10:34am
ACCY CHK	.07	10:35am
AIR BLK	.00	10:36am
SUB TEST	.00	10:36am
AIR BLK	.00	10:37am
SUB TEST	.00	10:39am
AIR BLK	.00	10:40am

Reported AC: .00 g/210L

Bobby D. Willis
Signature of Chemical Analyst

Court CVR

Bobby D. Willis
Analyst

Intox EC/IR-II: Preventive Maintenance

CLEVELAND COUNTY CLEVELAND COUNTY SD 220

Serial Number: 008893 Test Record Number: 165
Test Date: 11/13/2008 Test Time: 10:41am EST

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	10:42am
FLO	Pass	10:42am
FC	Pass	10:42am

Temperature Tests

Test	Status	Time
FC1	Pass	10:42am
SRC	Pass	10:42am
DET	Pass	10:42am
BAR	Pass	10:42am
BT	Pass	10:42am

Blank Tests

Test	Status	Time
AIR	Pass	10:42am


Printer Tests

Test	Status	Time
PRNT	Pass	10:42am

CRC Tests

Test	Status	Time
COMP	Pass	10:43am
CAL	Pass	10:43am

Preventive Maintenance
Status: Pass

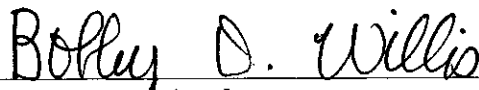

Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH

**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**


County Cleveland Instrument Location Cleveland Co. S.D.

Instrument Serial No. 008887 100 Justice Pl., Shelby
704-484-4888

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 13 day of November, 2008 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

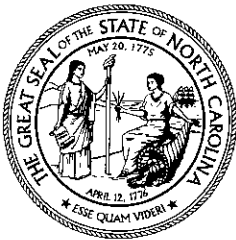

Bobby D. Wilco
Signature of Certifying Official

557
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

CLEVELAND COUNTY CLEVELAND COUNTY SD
220

Serial Number: 008887
Test Date: 11/13/2008

Citation Number: M0000000-0
Subject's Name:
PREVENTIVE, MAINTENANCE
Subject's Date of Birth: 11/11/1911
Subject's Sex: Male
Driver's License State: XX
Driver's License Number: NONE

Analyst's Name: WILLIS, BOBBY D
Permit Number: 08010E
Effective:
12/01/2007-12/01/2009

Officer's Name: NONE, NONE
Type of Agency: FTA
Agency: DHHS
Test Type: Breath Test

Lot Number: AG816303
Exp Date: 06/11/2010

Test	g/210L	Time
DIAG	Pass	9:54am
AIR BLK	.00	9:55am
ACCY CHK	.07	9:56am
AIR BLK	.00	9:56am
SUB TEST	.00	9:57am
AIR BLK	.00	9:58am
SUB TEST	.00	9:59am
AIR BLK	.00	10:00am

Reported AC: .00 g/210L

Bobby D. Willis
Signature of Chemical Analyst

Court CVR

Bobby D. Willis
Analyst

Intox EC/IR-II: Preventive Maintenance

CLEVELAND COUNTY CLEVELAND COUNTY SD 220

Serial Number: 008887 Test Record Number: 112
Test Date: 11/13/2008 Test Time: 10:02am EST

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	10:02am
FLO	Pass	10:02am
FC	Pass	10:02am

Temperature Tests

Test	Status	Time
FC1	Pass	10:02am
SRC	Pass	10:02am
DET	Pass	10:02am
BAR	Pass	10:02am
BT	Pass	10:02am

Blank Tests

Test	Status	Time
AIR	Pass	10:03am

Printer Tests

Test	Status	Time
PRNT	Pass	10:03am

CRC Tests

Test	Status	Time
COMP	Pass	10:03am
CAL	Pass	10:03am

Preventive Maintenance
Status: Pass

Bolley D. Willis
Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH


**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**

County Rutherford Instrument Location Forest City P. D.
Instrument Serial No. 008889 187 S. Church St., Forest City
828-245-5555

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 12 day of November, 20 08 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

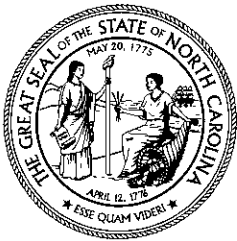

Bobby E. Willis
Signature of Certifying Official

557
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

RUTHERFORD COUNTY FOREST CITY PD 800

Serial Number: 008889
Test Date: 11/12/2008

Citation Number: M0000000-0
Subject's Name:
PREVENTIVE, MAINTENANCE
Subject's Date of Birth: 11/11/1911
Subject's Sex: Male
Driver's License State: XX
Driver's License Number: NONE

Analyst's Name: WILLIS, BOBBY D
Permit Number: 08010E
Effective:
12/01/2007-12/01/2009

Officer's Name: NONE, NONE
Type of Agency: FTA
Agency: DHHS
Test Type: Breath Test

Lot Number: AG816303
Exp Date: 06/11/2010

Test	g/210L	Time
DIAG	Pass	12:43pm
AIR BLK	.00	12:44pm
ACCY CHK	.07	12:44pm
AIR BLK	.00	12:45pm
SUB TEST	.00	12:46pm
AIR BLK	.00	12:47pm
SUB TEST	.00	12:49pm
AIR BLK	.00	12:50pm

Reported AC: .00 g/210L

Bobby D. Willis
Signature of Chemical Analyst

Court CVR

Bobby D. Willis
Analyst

Intox EC/IR-II: Preventive Maintenance

RUTHERFORD COUNTY FOREST CITY PD 800

Serial Number: 008889 Test Record Number: 103
Test Date: 11/12/2008 Test Time: 12:51pm EST

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	12:51pm
FLO	Pass	12:51pm
FC	Pass	12:51pm

Temperature Tests

Test	Status	Time
FC1	Pass	12:51pm
SRC	Pass	12:51pm
DET	Pass	12:51pm
BAR	Pass	12:51pm
BT	Pass	12:51pm

Blank Tests

Test	Status	Time
AIR	Pass	12:52pm

Printer Tests

Test	Status	Time
PRNT	Pass	12:52pm

CRC Tests

Test	Status	Time
COMP	Pass	12:52pm
CAL	Pass	12:52pm

Preventive Maintenance
Status: Pass

Bofley D. Willis
Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH

**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**


County Rutherford Instrument Location Rutherford Co. S.D.

Instrument Serial No. 008914 400 N. Washington St. Rutherfordton
828-627-6747

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 12 day of November, 20 08 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

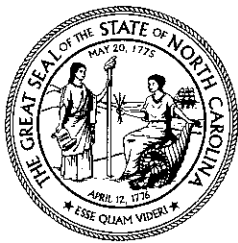

Billy L. Willes
Signature of Certifying Official

557
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Preventive Maintenance

RUTHERFORD COUNTY RUTHERFORD COUNTY SD 800

Serial Number: 008914 Test Record Number: 158
Test Date: 11/12/2008 Test Time: 11:13am EST

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	11:14am
FLO	Pass	11:14am
FC	Pass	11:14am

Temperature Tests

Test	Status	Time
FC1	Pass	11:14am
SRC	Pass	11:14am
DET	Pass	11:14am
BAR	Pass	11:14am
BT	Pass	11:14am

Blank Tests

Test	Status	Time
AIR	Pass	11:15am


Printer Tests

Test	Status	Time
PRNT	Pass	11:15am

CRC Tests

Test	Status	Time
COMP	Pass	11:15am
CAL	Pass	11:15am

Preventive Maintenance
Status: Pass

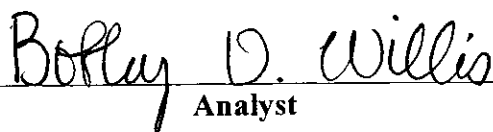
Analyst

Intox EC/IR-II: Subject Test

RUTHERFORD COUNTY RUTHERFORD COUNTY SD
800

Serial Number: 008914
Test Date: 11/12/2008

Citation Number: M0000000-0
Subject's Name:
PREVENTIVE, MAINTENANCE
Subject's Date of Birth: 11/11/1911
Subject's Sex: Male
Driver's License State: XX
Driver's License Number: NONE

Analyst's Name: WILLIS, BOBBY D
Permit Number: 08010E
Effective:
12/01/2007-12/01/2009

Officer's Name: NONE, NONE
Type of Agency: FTA
Agency: DHHS
Test Type: Breath Test

Lot Number: AG816303
Exp Date: 06/11/2010

Test	g/210L	Time
DIAG	Pass	11:03am
AIR BLK	.00	11:04am
ACCY CHK	.07	11:04am
AIR BLK	.00	11:05am
SUB TEST	.00	11:06am
AIR BLK	.00	11:07am
SUB TEST	.00	11:08am
AIR BLK	.00	11:09am

Reported AC: .00 g/210L

Bobby D. Willis
Signature of Chemical Analyst

Court CVR

Bobby D. Willis
Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH

**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**


County Polk Instrument Location Polk Co. S.D.

Instrument Serial No. 008833 46 Ward St., Columbus
828-874-3001

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 10 day of November, 2008 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

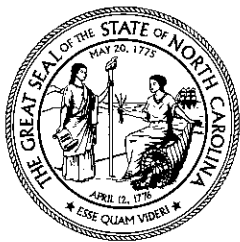

Bobby E. Willis
Signature of Certifying Official

557
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

POLK COUNTY POLK COUNTY SD 740

Serial Number: 008832
Test Date: 11/10/2008

Citation Number: M0000000-0
Subject's Name:

PREVENTIVE, MAINTENANCE

Subject's Date of Birth: 11/11/1911

Subject's Sex: Male

Driver's License State: XX

Driver's License Number: NONE

Analyst's Name: WILLIS, BOBBY D

Permit Number: 08010E

Effective:

12/01/2007-12/01/2009

Officer's Name: NONE, NONE

Type of Agency: FTA

Agency: DHHS


Test Type: Breath Test

Lot Number: AG816303

Exp Date: 06/11/2010


Test	g/210L	Time
DIAG	Pass	10:11am
AIR BLK	.00	10:12am
ACCY CHK	.07	10:13am
AIR BLK	.00	10:14am
SUB TEST	.00	10:15am
AIR BLK	.00	10:16am
SUB TEST	.00	10:17am
AIR BLK	.00	10:18am

Reported AC: .00 g/210L

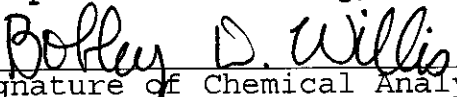

Signature of Chemical Analyst

Court CVR



Analyst

Intox EC/IR-II: Preventive Maintenance

POLK COUNTY POLK COUNTY SD 740

Serial Number: 008832 Test Record Number: 83
Test Date: 11/10/2008 Test Time: 10:21am EST

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	10:21am
FLO	Pass	10:21am
FC	Pass	10:22am

Temperature Tests

Test	Status	Time
FC1	Pass	10:22am
SRC	Pass	10:22am
DET	Pass	10:22am
BAR	Pass	10:22am
BT	Pass	10:22am

Blank Tests

Test	Status	Time
AIR	Pass	10:22am


Printer Tests

Test	Status	Time
PRNT	Pass	10:22am

CRC Tests

Test	Status	Time
COMP	Pass	10:23am
CAL	Pass	10:23am

Preventive Maintenance
Status: Pass

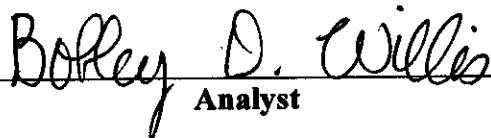
Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH


**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**

County RANDOLPH Instrument Location LIBERTY POLICE
Instrument Serial No. 008830 Dept, LIBERTY NC

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 25 day of NOVEMBER, 2008 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

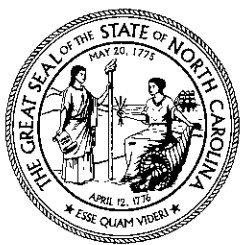

[Signature]
Signature of Certifying Official

371
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

RANDOLPH LIBERTY POLICE DEPT 750

Serial Number: 008830
Test Date: 11/25/2008

Citation Number: M0000000-0
Subject's Name:
PREVENTIVE, MAINTENANCE
Subject's Date of Birth: 11/11/1911
Subject's Sex: Male
Driver's License State: XX
Driver's License Number: NONE


Analyst's Name: RUSSELL, LARRY H
Permit Number: 06108E
Effective:
12/01/2007-12/01/2009

Officer's Name: NONE, NONE
Type of Agency: FTA
Agency: DHHS
Test Type: Breath Test

Lot Number: AG722601
Exp Date: 08/13/2009


Test	g/210L	Time
DIAG	Pass	2:43pm
AIR BLK	.00	2:44pm
ACCY CHK	.07	2:44pm
AIR BLK	.00	2:45pm
SUB TEST	.00	2:46pm
AIR BLK	.00	2:47pm
SUB TEST	.00	2:48pm
AIR BLK	.00	2:49pm

Reported AC: .00 g/210L

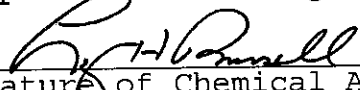

Signature of Chemical Analyst

Court CVR



Analyst

Intox EC/IR-II: Preventive Maintenance

RANDOLPH LIBERTY POLICE DEPT 750

Serial Number: 008830 Test Record Number: 113
Test Date: 11/25/2008 Test Time: 2:50pm EST

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	2:51pm
FLO	Pass	2:51pm
FC	Pass	2:51pm

Temperature Tests

Test	Status	Time
FC1	Pass	2:51pm
SRC	Pass	2:51pm
DET	Pass	2:51pm
BAR	Pass	2:51pm
BT	Pass	2:51pm

Blank Tests

Test	Status	Time
AIR	Pass	2:51pm


Printer Tests

Test	Status	Time
PRNT	Pass	2:51pm

CRC Tests

Test	Status	Time
COMP	Pass	2:52pm
CAL	Pass	2:52pm

Preventive Maintenance
Status: Pass

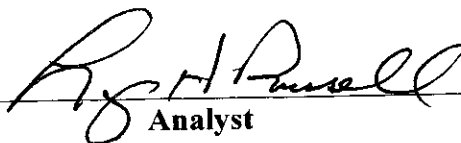
Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH

**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**


County RANDOLPH Instrument Location RANDOLPH Co. Jail

Instrument Serial No. 008899 ASHEBORO NC

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 25 day of NOVEMBER, 20 08 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

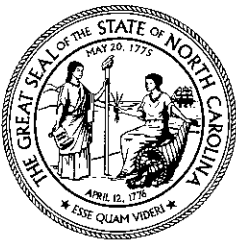

[Signature]
Signature of Certifying Official

371
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

RANDOLPH RANDOLPH CO. JAIL 750

Serial Number: 008899
Test Date: 11/25/2008

Citation Number: M0000000-0

Subject's Name:

PREVENTIVE, MAINTENANCE

Subject's Date of Birth: 11/11/1911

Subject's Sex: Male

Driver's License State: XX

Driver's License Number: NONE

Analyst's Name: RUSSELL, LARRY H

Permit Number: 06108E

Effective:

12/01/2007-12/01/2009

Officer's Name: NONE, NONE

Type of Agency: FTA

Agency: DHHS


Test Type: Breath Test

Lot Number: AG723402

Exp Date: 08/21/2009


Test	g/210L	Time
DIAG	Pass	12:44pm
AIR BLK	.00	12:45pm
ACCY CHK	.07	12:45pm
AIR BLK	.00	12:46pm
SUB TEST	.00	12:47pm
AIR BLK	.00	12:48pm
SUB TEST	.00	12:50pm
AIR BLK	.00	12:51pm

Reported AC: .00 g/210L

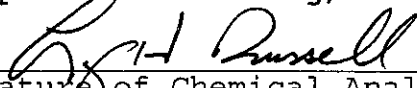

Signature of Chemical Analyst

Court CVR.

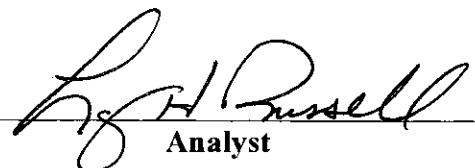

Analyst

Intox EC/IR-II: Preventive Maintenance

RANDOLPH RANDOLPH CO. JAIL 750

Serial Number: 008899 Test Record Number: 297
Test Date: 11/25/2008 Test Time: 12:52pm EST

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	12:53pm
FLO	Pass	12:53pm
FC	Pass	12:53pm

Temperature Tests

Test	Status	Time
FC1	Pass	12:53pm
SRC	Pass	12:53pm
DET	Pass	12:53pm
BAR	Pass	12:53pm
BT	Pass	12:53pm

Blank Tests

Test	Status	Time
AIR	Pass	12:54pm


Printer Tests

Test	Status	Time
PRNT	Pass	12:54pm

CRC Tests

Test	Status	Time
COMP	Pass	12:54pm
CAL	Pass	12:54pm

Preventive Maintenance
Status: Pass

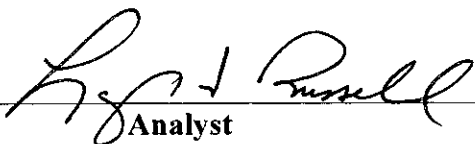

Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH


**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**

County RANDOLPH Instrument Location RANDLEMAN POLICE
Instrument Serial No. 008737 DEPT, RANDLEMAN NC

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 25 day of NOVEMBER, 2008 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

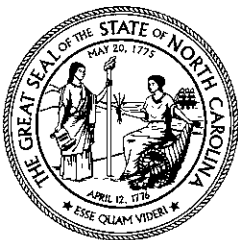

[Signature]
Signature of Certifying Official

371
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

RANDOLPH RANDLEMAN PD 750

Serial Number: 008737

Test Date: 11/25/2008

Citation Number: M0000000-0

Subject's Name:

PREVENTIVE, MAINTENANCE

Subject's Date of Birth: 11/11/1911

Subject's Sex: Male

Driver's License State: XX

Driver's License Number: NONE

Analyst's Name: RUSSELL, LARRY H

Permit Number: 06108E

Effective:

12/01/2007-12/01/2009

Officer's Name: NONE, NONE

Type of Agency: FTA

Agency: DHHS


Test Type: Breath Test

Lot Number: AG723402

Exp Date: 08/21/2009


Test	g/210L	Time
DIAG	Pass	11:42am
AIR BLK	.00	11:43am
ACCY CHK	.07	11:44am
AIR BLK	.00	11:44am
SUB TEST	.00	11:45am
AIR BLK	.00	11:46am
SUB TEST	.00	11:47am
AIR BLK	.00	11:48am

Reported AC: .00 g/210L

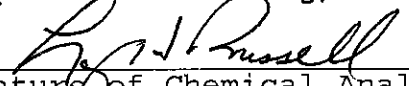

Signature of Chemical Analyst

Court CVR



Analyst

Intox EC/IR-II: Preventive Maintenance

RANDOLPH RANDLEMAN PD 750

Serial Number: 008737 Test Record Number: 136
Test Date: 11/25/2008 Test Time: 11:50am EST

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	11:50am
FLO	Pass	11:50am
FC	Pass	11:50am

Temperature Tests

Test	Status	Time
FC1	Pass	11:50am
SRC	Pass	11:50am
DET	Pass	11:50am
BAR	Pass	11:50am
BT	Pass	11:50am

Blank Tests

Test	Status	Time
AIR	Pass	11:51am


Printer Tests

Test	Status	Time
PRNT	Pass	11:51am

CRC Tests

Test	Status	Time
COMP	Pass	11:51am
CAL	Pass	11:51am

Preventive Maintenance
Status: Pass

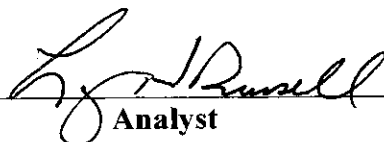
Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH

**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**


County ROBESON Instrument Location ST. PAULS POLICE DEPT.

Instrument Serial No. 008814 ST. PAULS, NC

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 21 day of NOVEMBER, 20 08 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

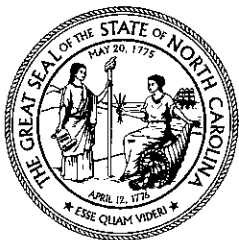

[Signature]
Signature of Certifying Official

371
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

ROBESON ST. PAULS PD. 770

Serial Number: 008814
Test Date: 11/21/2008

Citation Number: M0000000-0

Subject's Name:

PREVENTIVE, MAINTENANCE

Subject's Date of Birth: 11/11/1911

Subject's Sex: Male

Driver's License State: XX

Driver's License Number: NONE

Analyst's Name: RUSSELL, LARRY H

Permit Number: 06108E

Effective:

12/01/2007-12/01/2009

Officer's Name: NONE, NONE

Type of Agency: FTA

Agency: DHHS


Test Type: Breath Test

Lot Number: AG814101

Exp Date: 05/20/2010


Test	g/210L	Time
DIAG	Pass	10:35am
AIR BLK	.00	10:35am
ACCY CHK	.07	10:36am
AIR BLK	.00	10:37am
SUB TEST	.00	10:37am
AIR BLK	.00	10:38am
SUB TEST	.00	10:40am
AIR BLK	.00	10:41am

Reported AC: .00 g/210L

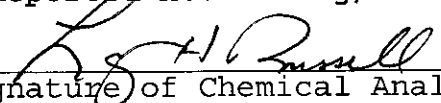

Signature of Chemical Analyst

Court CVR

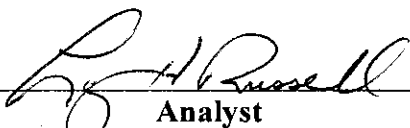

Analyst

Intox EC/IR-II: Preventive Maintenance

ROBESON ST. PAULS PD. 770

Serial Number: 008814 Test Record Number: 55
Test Date: 11/21/2008 Test Time: 10:44am EST

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	10:44am
FLO	Pass	10:44am
FC	Pass	10:44am

Temperature Tests

Test	Status	Time
FC1	Pass	10:44am
SRC	Pass	10:44am
DET	Pass	10:44am
BAR	Pass	10:44am
BT	Pass	10:44am

Blank Tests

Test	Status	Time
AIR	Pass	10:45am


Printer Tests

Test	Status	Time
PRNT	Pass	10:45am

CRC Tests

Test	Status	Time
COMP	Pass	10:45am
CAL	Pass	10:45am

Preventive Maintenance
Status: Pass

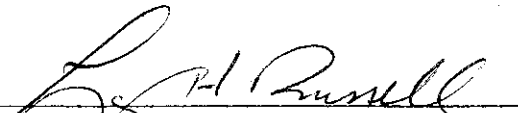

Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH

**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**

County RANDOLPH Instrument Location RANDOLPH Co. Jail

Instrument Serial No. 008860 Asheboro NC

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 25 day of NOVEMBER, 20 08 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.



[Signature]
Signature of Certifying Official

371
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

RANDOLPH RANDOLPH CO. JAIL 750

Serial Number: 008860

Test Date: 11/25/2008

Citation Number: M0000000-0

Subject's Name:

PREVENTIVE, MAINTENANCE

Subject's Date of Birth: 11/11/1911

Subject's Sex: Male

Driver's License State: XX

Driver's License Number: NONE

Analyst's Name: RUSSELL, LARRY H

Permit Number: 06108E

Effective:

12/01/2007-12/01/2009

Officer's Name: NONE, NONE

Type of Agency: FTA

Agency: DHHS


Test Type: Breath Test

Lot Number: AG722601

Exp Date: 08/13/2009


Test	g/210L	Time
DIAG	Pass	1:06pm
AIR BLK	.00	1:07pm
ACCY CHK	.07	1:07pm
AIR BLK	.00	1:08pm
SUB TEST	.00	1:09pm
AIR BLK	.00	1:10pm
SUB TEST	.00	1:11pm
AIR BLK	.00	1:12pm

Reported AC: .00 g/210L

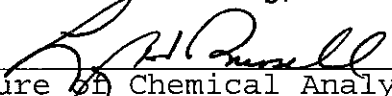

Signature of Chemical Analyst

Court CVR

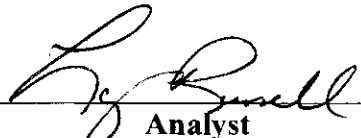

Analyst

Intox EC/IR-II: Preventive Maintenance

RANDOLPH RANDOLPH CO. JAIL 750

Serial Number: 008860 Test Record Number: 142
Test Date: 11/25/2008 Test Time: 1:13pm EST

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	1:14pm
FLO	Pass	1:14pm
FC	Pass	1:14pm

Temperature Tests

Test	Status	Time
FC1	Pass	1:14pm
SRC	Pass	1:14pm
DET	Pass	1:14pm
BAR	Pass	1:14pm
BT	Pass	1:14pm

Blank Tests

Test	Status	Time
AIR	Pass	1:15pm


Printer Tests

Test	Status	Time
PRNT	Pass	1:15pm

CRC Tests

Test	Status	Time
COMP	Pass	1:15pm
CAL	Pass	1:15pm

Preventive Maintenance
Status: Pass

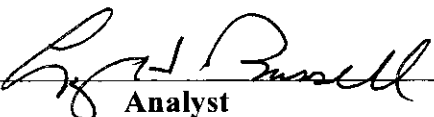
Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH

**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**


County ROBESON Instrument Location Lumberton,

Instrument Serial No. 008805 LEC.

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 25 day of November, 2008 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

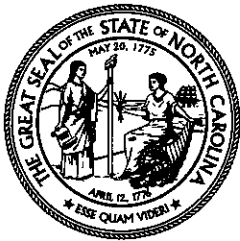

[Handwritten Signature]

Signature of Certifying Official

578

Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

ROBESON LUMBERTON, LEC 770

Serial Number: 008805
Test Date: 11/25/2008

Citation Number: M0000000-0
Subject's Name:
PREVENTIVE, MAINTENANCE
Subject's Date of Birth: 11/11/1911
Subject's Sex: Male
Driver's License State: XX
Driver's License Number: NONE


Analyst's Name: SIMMONS, PAUL T
Permit Number: 08619E
Effective:
12/01/2007-12/01/2009

Officer's Name: NONE, NONE
Type of Agency: FTA
Agency: DHHS
Test Type: Breath Test


Lot Number: AG814002
Exp Date: 05/19/2010

Test	g/210L	Time
DIAG	Pass	3:41pm
AIR BLK	.00	3:42pm
ACCY CHK	.07	3:42pm
AIR BLK	.00	3:43pm
SUB TEST	.00	3:44pm
AIR BLK	.00	3:45pm
SUB TEST	.00	3:46pm
AIR BLK	.00	3:47pm

Reported AC: .00 g/210L

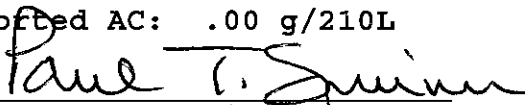
Signature of Chemical Analyst

Court CVR

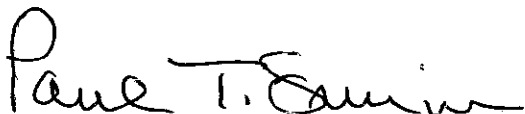
Analyst

Intox EC/IR-II: Preventive Maintenance

ROBESON LUMBERTON, LEC 770

Serial Number: 008805 Test Record Number: 218
Test Date: 11/25/2008 Test Time: 3:49pm EST

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	3:49pm
FLO	Pass	3:49pm
FC	Pass	3:49pm

Temperature Tests

Test	Status	Time
FC1	Pass	3:50pm
SRC	Pass	3:50pm
DET	Pass	3:50pm
BAR	Pass	3:50pm
BT	Pass	3:50pm

Blank Tests

Test	Status	Time
AIR	Pass	3:50pm


Printer Tests

Test	Status	Time
PRNT	Pass	3:50pm

CRC Tests

Test	Status	Time
COMP	Pass	3:50pm
CAL	Pass	3:50pm

Preventive Maintenance
Status: Pass

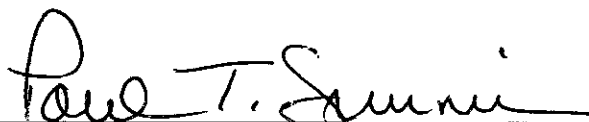
Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH


**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**

County Robeson Instrument Location Lumberton,
Instrument Serial No. 008836 LEC,

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 25 day of November, 2008 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

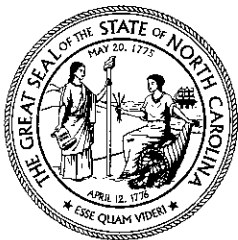

Paul T. [Signature]
Signature of Certifying Official

578
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

ROBESON LUMBERTON, LEC 770

Serial Number: 008836

Test Date: 11/25/2008

Citation Number: M0000000-0

Subject's Name:

PREVENTIVE, MAINTENANCE

Subject's Date of Birth: 11/11/1911

Subject's Sex: Male

Driver's License State: XX

Driver's License Number: NONE

Analyst's Name: SIMMONS, PAUL T

Permit Number: 08619E

Effective:

12/01/2007-12/01/2009

Officer's Name: NONE, NONE

Type of Agency: FTA

Agency: DHHS

Test Type: Breath Test

Lot Number: AG814002

Exp Date: 05/19/2010


Test	g/210L	Time
DIAG	Pass	4:01pm
AIR BLK	.00	4:01pm
ACCY CHK	.07	4:02pm
AIR BLK	.00	4:03pm
SUB TEST	.00	4:03pm
AIR BLK	.00	4:04pm
SUB TEST	.00	4:06pm
AIR BLK	.00	4:07pm

Reported AC: .00 g/210L



Signature of Chemical Analyst

Court CVR

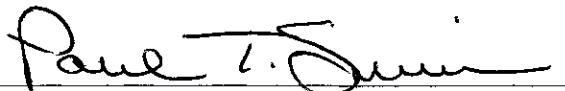

Analyst

Intox EC/IR-II: Preventive Maintenance

ROBESON LUMBERTON, LEC 770

Serial Number: 008836 Test Record Number: 166
Test Date: 11/25/2008 Test Time: 4:09pm EST

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	4:09pm
FLO	Pass	4:09pm
FC	Pass	4:09pm

Temperature Tests

Test	Status	Time
FC1	Pass	4:09pm
SRC	Pass	4:09pm
DET	Pass	4:09pm
BAR	Pass	4:09pm
BT	Pass	4:09pm

Blank Tests

Test	Status	Time
AIR	Pass	4:10pm


Printer Tests

Test	Status	Time
PRNT	Pass	4:10pm

CRC Tests

Test	Status	Time
COMP	Pass	4:10pm
CAL	Pass	4:10pm

Preventive Maintenance
Status: Pass

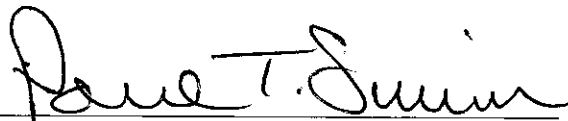

Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH


**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**

County ROBESON Instrument Location PEMBROKE
Instrument Serial No. 008837 POLICE DEPT.

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 25 day of NOVEMBER, 2008 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

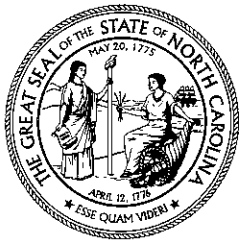

[Signature]
Signature of Certifying Official

578
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

ROBESON PEMBROKE, PD. 770

Serial Number: 008837

Test Date: 11/25/2008

Citation Number: M0000000-0

Subject's Name:

PREVENTIVE, MAINTENANCE

Subject's Date of Birth: 11/11/1911

Subject's Sex: Male

Driver's License State: XX

Driver's License Number: NONE

Analyst's Name: SIMMONS, PAUL T

Permit Number: 08619E

Effective:

12/01/2007-12/01/2009

Officer's Name: NONE, NONE

Type of Agency: FTA

Agency: DHS

Test Type: Breath Test

Lot Number: AG814002

Exp Date: 05/19/2010

Test	g/210L	Time
DIAG	Pass	2:42pm
AIR BLK	.00	2:42pm
ACCY CHK	.07	2:43pm
AIR BLK	.00	2:44pm
SUB TEST	.00	2:44pm
AIR BLK	.00	2:45pm
SUB TEST	.00	2:47pm
AIR BLK	.00	2:48pm

Reported AC: .00 g/210L



Signature of Chemical Analyst

Court CVR.



Analyst

Intox EC/IR-II: Preventive Maintenance

ROBESON PEMBROKE, PD. 770

Serial Number: 008837 Test Record Number: 61
Test Date: 11/25/2008 Test Time: 2:51pm EST

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	2:51pm
FLO	Pass	2:51pm
FC	Pass	2:51pm

Temperature Tests

Test	Status	Time
FC1	Pass	2:51pm
SRC	Pass	2:51pm
DET	Pass	2:51pm
BAR	Pass	2:51pm
BT	Pass	2:51pm

Blank Tests

Test	Status	Time
AIR	Pass	2:52pm


Printer Tests

Test	Status	Time
PRNT	Pass	2:52pm

CRC Tests

Test	Status	Time
COMP	Pass	2:52pm
CAL	Pass	2:52pm

Preventive Maintenance
Status: Pass

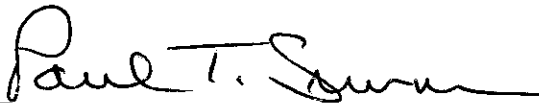

Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH


**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**

County Scotland Instrument Location Scotland
Instrument Serial No. 008850 County Sheriff's Dept.

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 25 day of NOVEMBER, 2008 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

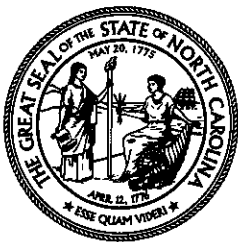

Paul T. Sumner
Signature of Certifying Official

578
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

SCOTLAND SCOTLAND CO. SD. 820

Serial Number: 008850
Test Date: 11/25/2008

Citation Number: M0000000-0
Subject's Name:
PREVENTIVE, MAINTENANCE
Subject's Date of Birth: 11/11/1911
Subject's Sex: Male
Driver's License State: XX
Driver's License Number: NONE

Analyst's Name: SIMMONS, PAUL T
Permit Number: 08619E
Effective:
12/01/2007-12/01/2009

Officer's Name: NONE, NONE
Type of Agency: FTA
Agency: DHHS
Test Type: Breath Test


Lot Number: AG814002
Exp Date: 05/19/2010

Test	g/210L	Time
DIAG	Pass	1:22pm
AIR BLK	.00	1:22pm
ACCY CHK	.07	1:23pm
AIR BLK	.00	1:24pm
SUB TEST	.00	1:24pm
AIR BLK	.00	1:25pm
SUB TEST	.00	1:27pm
AIR BLK	.00	1:28pm

Reported AC: .00 g/210L

Signature of Chemical Analyst

Court CVR


Analyst

Intox EC/IR-II: Preventive Maintenance

SCOTLAND SCOTLAND CO. SD. 820

Serial Number: 008850 Test Record Number: 84
Test Date: 11/25/2008 Test Time: 1:29pm EST

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	1:30pm
FLO	Pass	1:30pm
FC	Pass	1:30pm

Temperature Tests

Test	Status	Time
FC1	Pass	1:30pm
SRC	Pass	1:30pm
DET	Pass	1:30pm
BAR	Pass	1:30pm
BT	Pass	1:30pm

Blank Tests

Test	Status	Time
AIR	Pass	1:30pm


Printer Tests

Test	Status	Time
PRNT	Pass	1:30pm

CRC Tests

Test	Status	Time
COMP	Pass	1:31pm
CAL	Pass	1:31pm

Preventive Maintenance
Status: Pass

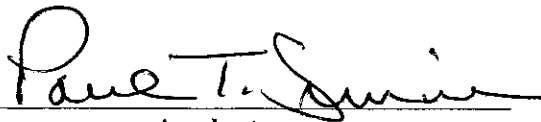

Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH


**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**

County ROBESON Instrument Location RED SPRINGS
Instrument Serial No. 008859 POLICE DEPT

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 25 day of NOVEMBER, 2008 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

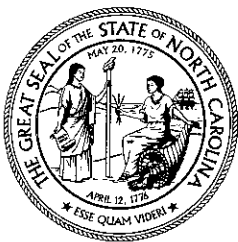

Paul T. Summons
Signature of Certifying Official

578
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

ROBESON RED SPRINGS PD 770

Serial Number: 008857
Test Date: 11/25/2008

Citation Number: M0000000-0
Subject's Name:
PREVENTIVE, MAINTENANCE
Subject's Date of Birth: 11/11/1911
Subject's Sex: Male
Driver's License State: XX
Driver's License Number: NONE


Analyst's Name: SIMMONS, PAUL T
Permit Number: 08619E
Effective:
12/01/2007-12/01/2009

Officer's Name: NONE, NONE
Type of Agency: FTA
Agency: DHHS
Test Type: Breath Test


Lot Number: AG814101
Exp Date: 05/20/2010

Test	g/210L	Time
DIAG	Pass	10:46am
AIR BLK	.00	10:47am
ACCY CHK	.07	10:47am
AIR BLK	.00	10:48am
SUB TEST	.00	10:49am
AIR BLK	.00	10:50am
SUB TEST	.00	10:51am
AIR BLK	.00	10:52am

Reported AC: .00 g/210L

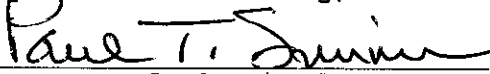
Signature of Chemical Analyst

Court CVR

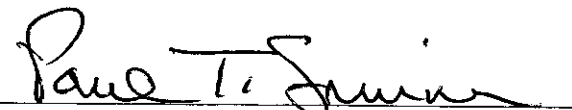
Analyst

Intox EC/IR-II: Preventive Maintenance

ROBESON RED SPRINGS PD 770

Serial Number: 008857 Test Record Number: 52
Test Date: 11/25/2008 Test Time: 10:53am EST

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	10:54am
FLO	Pass	10:54am
FC	Pass	10:54am

Temperature Tests

Test	Status	Time
FC1	Pass	10:54am
SRC	Pass	10:54am
DET	Pass	10:54am
BAR	Pass	10:54am
BT	Pass	10:54am

Blank Tests

Test	Status	Time
AIR	Pass	10:55am


Printer Tests

Test	Status	Time
PRNT	Pass	10:55am

CRC Tests

Test	Status	Time
COMP	Pass	10:55am
CAL	Pass	10:55am

Preventive Maintenance
Status: Pass



Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH


**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**

County SCOTLAND Instrument Location LAURINBURG
Instrument Serial No. 008834 POLICE DEPT.

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 24 day of NOVEMBER, 2008 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

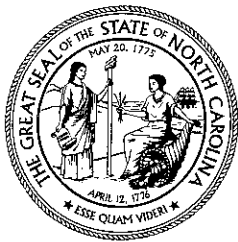

Paul T. Sumner
Signature of Certifying Official

578
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

SCOTLAND LAURINBURG PD 820

Serial Number: 008834
Test Date: 11/24/2008

Citation Number: M0000000-0
Subject's Name:
PREVENTIVE, MAINTENANCE
Subject's Date of Birth: 11/11/1911
Subject's Sex: Male
Driver's License State: XX
Driver's License Number: NONE


Analyst's Name: SIMMONS, PAUL T
Permit Number: 08619E
Effective:
12/01/2007-12/01/2009

Officer's Name: NONE, NONE
Type of Agency: FTA
Agency: DHHS
Test Type: Breath Test

Lot Number: AG814002
Exp Date: 05/19/2010

Test	g/210L	Time
DIAG	Pass	9:56am
AIR BLK	.00	9:56am
ACCY CHK	.07	9:57am
AIR BLK	.00	9:58am
SUB TEST	.00	9:59am
AIR BLK	.00	10:00am
SUB TEST	.00	10:01am
AIR BLK	.00	10:02am

Reported AC: .00 g/210L

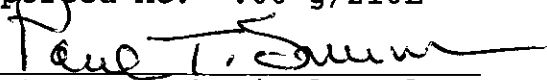

Signature of Chemical Analyst

Court CVR



Analyst

Intox EC/IR-II: Preventive Maintenance

SCOTLAND LAURINBURG PD 820

Serial Number: 008834 Test Record Number: 80
Test Date: 11/24/2008 Test Time: 10:04am EST

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	10:04am
FLO	Pass	10:04am
FC	Pass	10:04am

Temperature Tests

Test	Status	Time
FC1	Pass	10:04am
SRC	Pass	10:04am
DET	Pass	10:04am
BAR	Pass	10:04am
BT	Pass	10:04am

Blank Tests

Test	Status	Time
AIR	Pass	10:05am


Printer Tests

Test	Status	Time
PRNT	Pass	10:05am

CRC Tests

Test	Status	Time
COMP	Pass	10:05am
CAL	Pass	10:05am

Preventive Maintenance
Status: Pass

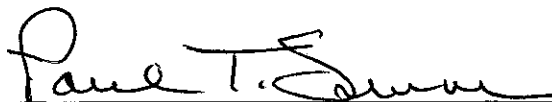

Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH


**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**

County HOKE Instrument Location HOKE COUNTY JAIL
Instrument Serial No. 008855 RAEFORD, NORTH CAROLINA

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 19 day of NOVEMBER, 2008 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

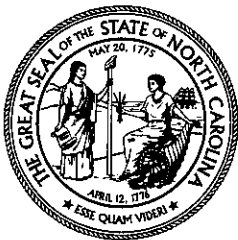

[Signature]
Signature of Certifying Official

572
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

HOKE HOKE CO. SD. 460

Serial Number: 008855
Test Date: 11/19/2008

Citation Number: M0000000-0
Subject's Name:
PREVENTIVE, MAINTENANCE
Subject's Date of Birth: 11/11/1911
Subject's Sex: Male
Driver's License State: XX
Driver's License Number: NONE


Analyst's Name: SIMMONS, PAUL T
Permit Number: 08619E
Effective:
12/01/2007-12/01/2009

Officer's Name: NONE, NONE
Type of Agency: FTA
Agency: DHHS
Test Type: Breath Test

Lot Number: AG814101
Exp Date: 05/20/2010


Test	g/210L	Time
DIAG	Pass	9:40am
AIR BLK	.00	9:41am
ACCY CHK	.07	9:42am
AIR BLK	.00	9:43am
SUB TEST	.00	9:43am
AIR BLK	.00	9:44am
SUB TEST	.00	9:46am
AIR BLK	.00	9:46am

Reported AC: .00 g/210L

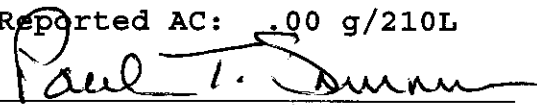

Signature of Chemical Analyst

Court CVR .



Analyst

Intox EC/IR-II: Preventive Maintenance

HOKE HOKE CO. SD. 460

Serial Number: 008855 Test Record Number: 140
Test Date: 11/19/2008 Test Time: 9:47am EST

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	9:48am
FLO	Pass	9:48am
FC	Pass	9:48am

Temperature Tests

Test	Status	Time
FC1	Pass	9:48am
SRC	Pass	9:48am
DET	Pass	9:48am
BAR	Pass	9:48am
BT	Pass	9:48am

Blank Tests

Test	Status	Time
AIR	Pass	9:49am

Printer Tests

Test	Status	Time
PRNT	Pass	9:49am

CRC Tests

Test	Status	Time
COMP	Pass	9:49am
CAL	Pass	9:49am

Preventive Maintenance
Status: Pass



Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH


**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**

County CUMBERLAND Instrument Location CUMBERLAND COUNTY
Instrument Serial No. 008614 Detention Center

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 21 day of November, 2008 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

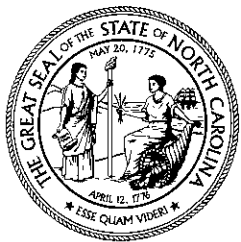

Paul T. Sumner
Signature of Certifying Official

578
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

CUMBERLAND COUNTY DETENTION CTR 250

Serial Number: 008614

Test Date: 11/21/2008

Citation Number: M0000000-0

Subject's Name:

PREVENTIVE, MAINTENANCE

Subject's Date of Birth: 11/11/1911

Subject's Sex: Male

Driver's License State: XX

Driver's License Number: NONE

Analyst's Name: SIMMONS, PAUL T

Permit Number: 08619E

Effective:

12/01/2007-12/01/2009

Officer's Name: NONE, NONE

Type of Agency: FTA

Agency: DHHS


Test Type: Breath Test

Lot Number: AG722702


Exp Date: 08/14/2009

Test	g/210L	Time
DIAG	Pass	5:27pm
AIR BLK	.00	5:28pm
ACCY CHK	.07	5:28pm
AIR BLK	.00	5:29pm
SUB TEST	.00	5:30pm
AIR BLK	.00	5:31pm
SUB TEST	.00	5:33pm
AIR BLK	.00	5:33pm

Reported AC: .00 g/210L

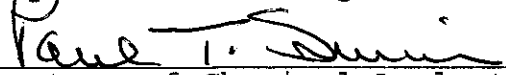
Signature of Chemical Analyst

Court CVR


Analyst

Intox EC/IR-II: Preventive Maintenance

CUMBERLAND COUNTY DETENTION CTR 250

Serial Number: 008614 Test Record Number: 549
Test Date: 11/21/2008 Test Time: 5:34pm EST

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	5:35pm
FLO	Pass	5:35pm
FC	Pass	5:35pm

Temperature Tests

Test	Status	Time
FC1	Pass	5:35pm
SRC	Pass	5:35pm
DET	Pass	5:35pm
BAR	Pass	5:35pm
BT	Pass	5:35pm

Blank Tests

Test	Status	Time
AIR	Pass	5:36pm

Printer Tests

Test	Status	Time
PRNT	Pass	5:36pm

CRC Tests

Test	Status	Time
COMP	Pass	5:36pm
CAL	Pass	5:36pm

Preventive Maintenance
Status: Pass



Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH


**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**

County Cumberland Instrument Location Cumberland County
Instrument Serial No. 008672 Detention Center

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 21 day of November, 2008 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

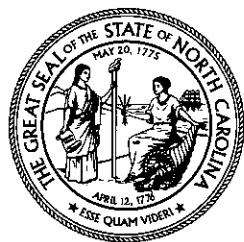

Paul Thomas
Signature of Certifying Official

578
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

CUMBERLAND COUNTY DETENTION CTR 250

Serial Number: 008672
Test Date: 11/21/2008

Citation Number: M0000000-0

Subject's Name:

PREVENTIVE, MAINTENANCE

Subject's Date of Birth: 11/11/1911

Subject's Sex: Male

Driver's License State: XX

Driver's License Number: NONE

Analyst's Name: SIMMONS, PAUL T

Permit Number: 08619E

Effective:

12/01/2007-12/01/2009

Officer's Name: NONE, NONE

Type of Agency: FTA

Agency: DHHS


Test Type: Breath Test

Lot Number: AG722602


Exp Date: 08/13/2009

Test	g/210L	Time
DIAG	Pass	4:49pm
AIR BLK	.00	4:50pm
ACCY CHK	.08	4:51pm
AIR BLK	.00	4:51pm
SUB TEST	.00	4:52pm
AIR BLK	.00	4:53pm
SUB TEST	.00	4:54pm
AIR BLK	.00	4:55pm

Reported AC: .00 g/210L

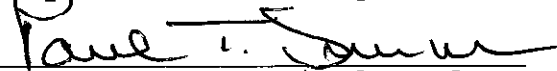
Signature of Chemical Analyst

Court CVR

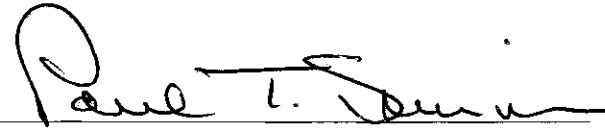
Analyst

Intox EC/IR-II: Preventive Maintenance

CUMBERLAND COUNTY DETENTION CTR 250

Serial Number: 008672 Test Record Number: 840
Test Date: 11/21/2008 Test Time: 4:58pm EST

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	4:58pm
FLO	Pass	4:58pm
FC	Pass	4:58pm

Temperature Tests

Test	Status	Time
FC1	Pass	4:58pm
SRC	Pass	4:58pm
DET	Pass	4:58pm
BAR	Pass	4:58pm
BT	Pass	4:58pm

Blank Tests

Test	Status	Time
AIR	Pass	4:59pm


Printer Tests

Test	Status	Time
PRNT	Pass	4:59pm

CRC Tests

Test	Status	Time
COMP	Pass	4:59pm
CAL	Pass	4:59pm

Preventive Maintenance
Status: Pass

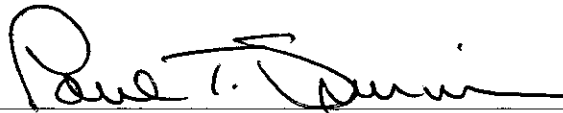

Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH


**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**

County CUMBERLAND Instrument Location CUMBERLAND COUNTY
Instrument Serial No. 008633 DETENTION CENTER

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 21 day of November, 2008 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

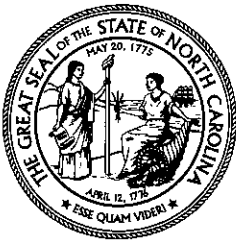

Paul T. Cannon 578
Signature of Certifying Official Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

CUMBERLAND COUNTY DETENTION CTR 250

Serial Number: 008633
Test Date: 11/21/2008

Citation Number: M0000000-0
Subject's Name:
PREVENTIVE, MAINTENANCE
Subject's Date of Birth: 11/11/1911
Subject's Sex: Male
Driver's License State: XX
Driver's License Number: NONE


Analyst's Name: SIMMONS, PAUL T
Permit Number: 08619E
Effective:
12/01/2007-12/01/2009

Officer's Name: NONE, NONE
Type of Agency: FTA
Agency: DHHS
Test Type: Breath Test


Lot Number: AG722802
Exp Date: 08/15/2009

Test	g/210L	Time
DIAG	Pass	4:22pm
AIR BLK	.00	4:23pm
ACCY CHK	.08	4:24pm
AIR BLK	.00	4:24pm
SUB TEST	.00	4:25pm
AIR BLK	.00	4:26pm
SUB TEST	.00	4:27pm
AIR BLK	.00	4:28pm

Reported AC: .00 g/210L


Signature of Chemical Analyst

Court CVR


Analyst

Intox EC/IR-II: Preventive Maintenance

CUMBERLAND COUNTY DETENTION CTR 250

Serial Number: 008633 Test Record Number: 392
Test Date: 11/21/2008 Test Time: 4:30pm EST

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	4:31pm
FLO	Pass	4:31pm
FC	Pass	4:31pm

Temperature Tests

Test	Status	Time
FC1	Pass	4:31pm
SRC	Pass	4:31pm
DET	Pass	4:31pm
BAR	Pass	4:31pm
BT	Pass	4:31pm

Blank Tests

Test	Status	Time
AIR	Pass	4:32pm


Printer Tests

Test	Status	Time
PRNT	Pass	4:32pm

CRC Tests

Test	Status	Time
COMP	Pass	4:32pm
CAL	Pass	4:32pm

Preventive Maintenance
Status: Pass

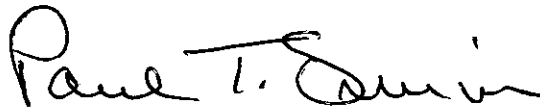

Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH


**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**

County Cumberland Instrument Location Cumberland County
Instrument Serial No. 008632 Detention Center

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 21 day of November, 2008 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

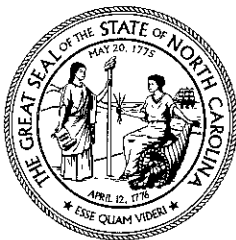

Paul T. Sumner
Signature of Certifying Official

578
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

CUMBERLAND COUNTY DETENTION CTR 250

Serial Number: 008632

Test Date: 11/21/2008

Citation Number: M0000000-0

Subject's Name:

PREVENTIVE, MAINTENANCE

Subject's Date of Birth: 11/11/1911

Subject's Sex: Male

Driver's License State: XX

Driver's License Number: NONE

Analyst's Name: SIMMONS, PAUL T

Permit Number: 08619E

Effective:

12/01/2007-12/01/2009

Officer's Name: NONE, NONE

Type of Agency: FTA

Agency: DHHS

Test Type: Breath Test

Lot Number: AG816303

Exp Date: 11/06/2010

Test	g/210L	Time
DIAG	Pass	3:52pm
AIR BLK	.00	3:52pm
ACCY CHK	.08	3:53pm
AIR BLK	.00	3:54pm
SUB TEST	.00	3:54pm
AIR BLK	.00	3:55pm
SUB TEST	.00	3:57pm
AIR BLK	.00	3:58pm

Reported AC: .00 g/210L



Signature of Chemical Analyst

Court CVR



Analyst

This form is used when performing Preventive Maintenance procedures
Forensic Tests for Alcohol Branch
Department of Health and Human Services
Rev. 12/2007

Intox EC/IR-II: Preventive Maintenance

CUMBERLAND COUNTY DETENTION CTR 250

Serial Number: 008632 Test Record Number: 437
Test Date: 11/21/2008 Test Time: 3:59pm EST

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	3:59pm
FLO	Pass	3:59pm
FC	Pass	3:59pm

Temperature Tests

Test	Status	Time
FC1	Pass	4:00pm
SRC	Pass	4:00pm
DET	Pass	4:00pm
BAR	Pass	4:00pm
BT	Pass	4:00pm

Blank Tests

Test	Status	Time
AIR	Pass	4:00pm


Printer Tests

Test	Status	Time
PRNT	Pass	4:00pm

CRC Tests

Test	Status	Time
COMP	Pass	4:00pm
CAL	Pass	4:00pm

Preventive Maintenance
Status: Pass



Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH

**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**


County CHOWAN Instrument Location CHOWAN CO. PUBLIC SAFETY

Instrument Serial No. 008586 CRA, 305 W. FREEMASON ST.,
EDENTON, N.C.

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 13th day of NOVEMBER, 20 08 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

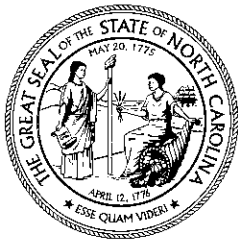

[Signature]
Signature of Certifying Official

647
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

CHOWAN COUNTY PUBLIC SAFETY CENTER 200

Serial Number: 008586
Test Date: 11/13/2008

Citation Number: M0000000-0
Subject's Name:
PREVENTIVE, MAINTENANCE
Subject's Date of Birth: 11/11/1911
Subject's Sex: Male
Driver's License State: XX
Driver's License Number: NONE


Analyst's Name: KEESLER, LINDA A
Permit Number: 11646E
Effective:
12/01/2007-12/01/2009

Officer's Name: NONE, NONE
Type of Agency: FTA
Agency: DHHS
Test Type: Breath Test

Lot Number: AG809301
Exp Date: 04/02/2010


Test	g/210L	Time
DIAG	Pass	1:06pm
AIR BLK	.00	1:07pm
ACCY CHK	.07	1:08pm
AIR BLK	.00	1:09pm
SUB TEST	.00	1:09pm
AIR BLK	.00	1:10pm
SUB TEST	.00	1:12pm
AIR BLK	.00	1:13pm

Reported AC: .00 g/210L

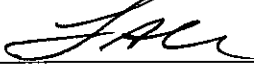

Signature of Chemical Analyst

Court CVR

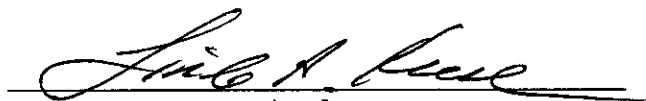

Analyst

Intox EC/IR-II: Preventive Maintenance

CHOWAN COUNTY PUBLIC SAFETY CENTER 200

Serial Number: 008586 Test Record Number: 262
Test Date: 11/13/2008 Test Time: 1:14pm EST

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	1:14pm
FLO	Pass	1:14pm
FC	Pass	1:14pm

Temperature Tests

Test	Status	Time
FC1	Pass	1:15pm
SRC	Pass	1:15pm
DET	Pass	1:15pm
BAR	Pass	1:15pm
BT	Pass	1:15pm

Blank Tests

Test	Status	Time
AIR	Pass	1:15pm


Printer Tests

Test	Status	Time
PRNT	Pass	1:15pm

CRC Tests

Test	Status	Time
COMP	Pass	1:16pm
CAL	Pass	1:16pm

Preventive Maintenance
Status: Pass

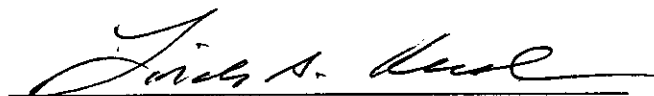

Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH

**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**


County Wake Instrument Location BAT Music & West

Instrument Serial No. 008788 Holly Springs

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 22nd day of November, 2008 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

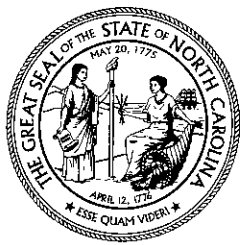

[Signature]
Signature of Certifying Official

636
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Preventive Maintenance

WAKE COUNTY BAT MOBILE UNIT 5 910

Serial Number: 008788 Test Record Number: 132
Test Date: 11/22/2008 Test Time: 11:10pm EST

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	11:11pm
FLO	Pass	11:11pm
FC	Pass	11:11pm

Temperature Tests

Test	Status	Time
FC1	Pass	11:11pm
SRC	Pass	11:11pm
DET	Pass	11:11pm
BAR	Pass	11:11pm
BT	Pass	11:11pm

Blank Tests

Test	Status	Time
AIR	Pass	11:11pm


Printer Tests

Test	Status	Time
PRNT	Pass	11:12pm

CRC Tests

Test	Status	Time
COMP	Pass	11:12pm
CAL	Pass	11:12pm

Preventive Maintenance
Status: Pass

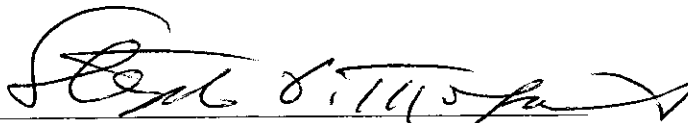
Analyst

Intox EC/IR-II: Subject Test

WAKE COUNTY BAT MOBILE UNIT 5 910

Serial Number: 008788
Test Date: 11/22/2008

Citation Number: M0000000-0
Subject's Name:
PREVENTIVE, MAINTENANCE
Subject's Date of Birth: 11/11/1911
Subject's Sex: Male
Driver's License State: XX
Driver's License Number: NONE

Analyst's Name: MORGART, STEPHEN G
Permit Number: 09372E
Effective:
12/01/2007-12/01/2009

Officer's Name: NONE, NONE
Type of Agency: FTA
Agency: DHHS
Test Type: Breath Test

Lot Number: AG722802
Exp Date: 08/15/2009

Test	g/210L	Time
DIAG	Pass	11:02pm
AIR BLK	.00	11:03pm
ACCY CHK	.07	11:04pm
AIR BLK	.00	11:05pm
SUB TEST	.00	11:05pm
AIR BLK	.00	11:06pm
SUB TEST	.00	11:08pm
AIR BLK	.00	11:09pm

Reported AC: .00 g/210L


Signature of Chemical Analyst

Court CVR


Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH


**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**

County Wake Instrument Location BAT Mobile Unit #5
Instrument Serial No. 008600 Holly Springs

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 22nd day of November, 2008 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

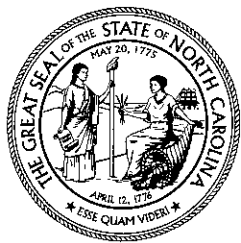

Frank C. Tilton
Signature of Certifying Official

636
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Preventive Maintenance

WAKE COUNTY BAT MOBILE UNIT 5 910

Serial Number: 008600 Test Record Number: 415
Test Date: 11/22/2008 Test Time: 10:48pm EST

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	10:48pm
FLO	Pass	10:48pm
FC	Pass	10:48pm

Temperature Tests

Test	Status	Time
FC1	Pass	10:48pm
SRC	Pass	10:48pm
DET	Pass	10:48pm
BAR	Pass	10:48pm
BT	Pass	10:48pm

Blank Tests

Test	Status	Time
AIR	Pass	10:49pm

Printer Tests

Test	Status	Time
PRNT	Pass	10:49pm

CRC Tests

Test	Status	Time
COMP	Pass	10:49pm
CAL	Pass	10:49pm

Preventive Maintenance
Status: Pass



Analyst

Intox EC/IR-II: Subject Test

WAKE COUNTY BAT MOBILE UNIT 5 910

Serial Number: 008600
Test Date: 11/22/2008

Citation Number: M0000000-0

Subject's Name:

PREVENTIVE, MAINTENANCE

Subject's Date of Birth: 11/11/1911

Subject's Sex: Male

Driver's License State: XX

Driver's License Number: NONE

Analyst's Name: MORGART, STEPHEN G

Permit Number: 09372E

Effective:

12/01/2007-12/01/2009

Officer's Name: NONE, NONE

Type of Agency: FTA

Agency: DHHS

Test Type: Breath Test

Lot Number: AG723401

Exp Date: 08/21/2009

Test	g/210L	Time
DIAG	Pass	10:39pm
AIR BLK	.00	10:39pm
ACCY CHK	.07	10:40pm
AIR BLK	.00	10:41pm
SUB TEST	.00	10:42pm
AIR BLK	.00	10:42pm
SUB TEST	.00	10:44pm
AIR BLK	.00	10:45pm

Reported AC: .00 g/210L

Signature of Chemical Analyst

Court CVR



Analyst

**This form is used when performing Preventive Maintenance procedures
Forensic Tests for Alcohol Branch
Department of Health and Human Services
Rev. 12/2007**

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH

**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**


County Wake Instrument Location BSI Mobile Unit #5

Instrument Serial No. 008698 Holly Springs

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 22nd day of November, 2008 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

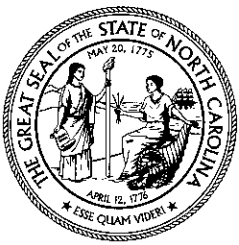

[Signature]
Signature of Certifying Official

636
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Preventive Maintenance

WAKE COUNTY BAT MOBILE UNIT 5 910

Serial Number: 008698 Test Record Number: 261
Test Date: 11/22/2008 Test Time: 10:49pm EST

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	10:50pm
FLO	Pass	10:50pm
FC	Pass	10:50pm

Temperature Tests

Test	Status	Time
FC1	Pass	10:50pm
SRC	Pass	10:50pm
DET	Pass	10:50pm
BAR	Pass	10:50pm
BT	Pass	10:50pm

Blank Tests

Test	Status	Time
AIR	Pass	10:51pm


Printer Tests

Test	Status	Time
PRNT	Pass	10:51pm

CRC Tests

Test	Status	Time
COMP	Pass	10:51pm
CAL	Pass	10:51pm

Preventive Maintenance
Status: Pass

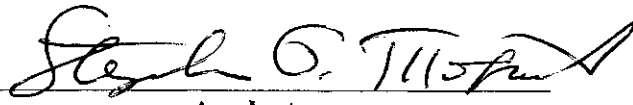

Analyst

Intox EC/IR-II: Subject Test

WAKE COUNTY BAT MOBILE UNIT 5 910

Serial Number: 008698
Test Date: 11/22/2008

Citation Number: M0000000-0
Subject's Name: NONE, NONE
Subject's Date of Birth: 11/11/1911
Subject's Sex: Male
Driver's License State: XX
Driver's License Number: NONE


Analyst's Name: MORGART, STEPHEN G
Permit Number: 09372E
Effective:
12/01/2007-12/01/2009

Officer's Name: NONE, NONE
Type of Agency: FTA
Agency: DHHS
Test Type: Breath Test

Lot Number: AG723401
Exp Date: 08/21/2009


Test	g/210L	Time
DIAG	Pass	10:41pm
AIR BLK	.00	10:42pm
ACCY CHK	.07	10:43pm
AIR BLK	.00	10:44pm
SUB TEST	.00	10:45pm
AIR BLK	.00	10:46pm
SUB TEST	.00	10:47pm
AIR BLK	.00	10:48pm

Reported AC: .00 g/210L

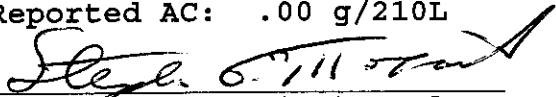

Signature of Chemical Analyst

Court CVR

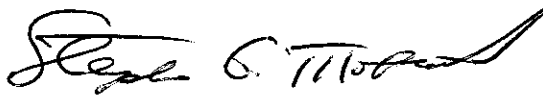

Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH


PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II

County Wake Instrument Location BAF Mobile Unit #5
Instrument Serial No. 608698 Camy

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 21st day of November, 2008 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

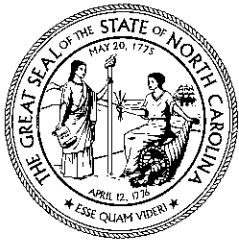

[Signature]

Signature of Certifying Official

636

Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Preventive Maintenance

WAKE COUNTY BAT MOBILE UNIT 5 910

Serial Number: 008698 Test Record Number: 257
Test Date: 11/21/2008 Test Time: 11:49pm EST

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	11:50pm
FLO	Pass	11:50pm
FC	Pass	11:50pm

Temperature Tests

Test	Status	Time
FC1	Pass	11:50pm
SRC	Pass	11:50pm
DET	Pass	11:50pm
BAR	Pass	11:50pm
BT	Pass	11:50pm

Blank Tests

Test	Status	Time
AIR	Pass	11:51pm


Printer Tests

Test	Status	Time
PRNT	Pass	11:51pm

CRC Tests

Test	Status	Time
COMP	Pass	11:51pm
CAL	Pass	11:51pm

Preventive Maintenance
Status: Pass

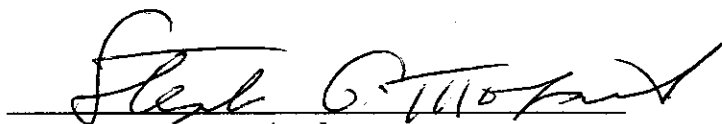

Analyst

Intox EC/IR-II: Subject Test.

WAKE COUNTY BAT MOBILE UNIT 5 910

Serial Number: 008698
Test Date: 11/21/2008

Citation Number: M0000000-0
Subject's Name:
PREVENTIVE, MAINTENANCE
Subject's Date of Birth: 11/11/1911
Subject's Sex: Male
Driver's License State: XX
Driver's License Number: NONE


Analyst's Name: MORGART, STEPHEN G
Permit Number: 09372E
Effective:
12/01/2007-12/01/2009

Officer's Name: NONE, NONE
Type of Agency: FTA
Agency: DHHS
Test Type: Breath Test


Lot Number: AG723401
Exp Date: 08/21/2009

Test	g/210L	Time
DIAG	Pass	11:41pm
AIR BLK	.00	11:42pm
ACCY CHK	.07	11:43pm
AIR BLK	.00	11:44pm
SUB TEST	.00	11:44pm
AIR BLK	.00	11:45pm
SUB TEST	.00	11:47pm
AIR BLK	.00	11:48pm

Reported AC: .00 g/210L

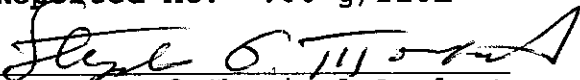
Signature of Chemical Analyst

Court CVR

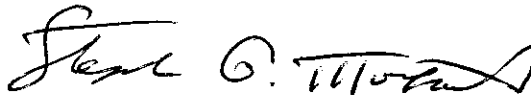

Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH

**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**


County Watauga Instrument Location BAT Mobile Unit #5

Instrument Serial No. 005600 Cary

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 21ST day of NOVEMBER, 2008 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

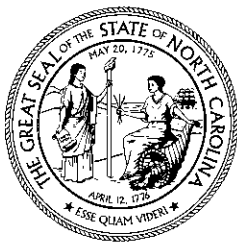

Steph G. Thompson

Signature of Certifying Official

636

Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Preventive Maintenance

WAKE COUNTY BAT MOBILE UNIT 5 910

Serial Number: 008600 Test Record Number: 412
Test Date: 11/21/2008 Test Time: 11:48pm EST

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	11:49pm
FLO	Pass	11:49pm
FC	Pass	11:49pm

Temperature Tests

Test	Status	Time
FC1	Pass	11:49pm
SRC	Pass	11:49pm
DET	Pass	11:49pm
BAR	Pass	11:49pm
BT	Pass	11:49pm

Blank Tests

Test	Status	Time
AIR	Pass	11:49pm

Printer Tests

Test	Status	Time
PRNT	Pass	11:49pm

CRC Tests

Test	Status	Time
COMP	Pass	11:49pm
CAL	Pass	11:49pm

Preventive Maintenance
Status: Pass



Analyst

Intox EC/IR-II: Subject Test.

WAKE COUNTY BAT MOBILE UNIT 5 910

Serial Number: 008600
Test Date: 11/21/2008

Citation Number: M0000000-0
Subject's Name:

PREVENTIVE, MAINTENANCE

Subject's Date of Birth: 11/11/1911

Subject's Sex: Male

Driver's License State: XX

Driver's License Number: NONE

Analyst's Name: MORGART, STEPHEN G

Permit Number: 09372E

Effective:

12/01/2007-12/01/2009

Officer's Name: NONE, NONE

Type of Agency: FTA

Agency: DHHS


Test Type: Breath Test

Lot Number: AG723401


Exp Date: 08/21/2009

Test	g/210L	Time
DIAG	Pass	11:39pm
AIR BLK	.00	11:40pm
ACCY CHK	.07	11:40pm
AIR BLK	.00	11:41pm
SUB TEST	.00	11:42pm
AIR BLK	.00	11:43pm
SUB TEST	.00	11:45pm
AIR BLK	.00	11:46pm

Reported AC: .00 g/210L

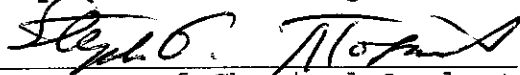
Signature of Chemical Analyst

Court CVR

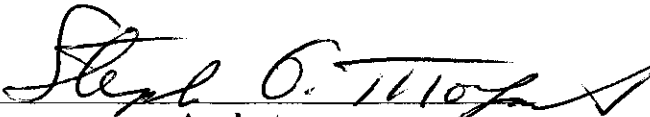
Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH


**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**

County Wake Instrument Location Bot mobile units
Instrument Serial No. 008788 Camy

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 21st day of November, 2008 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

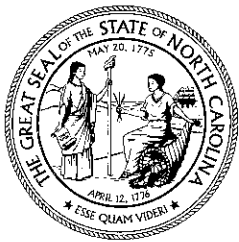

[Signature]

Signature of Certifying Official

636

Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Preventive Maintenance

WAKE COUNTY BAT MOBILE UNIT 5 910

Serial Number: 008788 Test Record Number: 130
Test Date: 11/21/2008 Test Time: 11:57pm EST

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	11:58pm
FLO	Pass	11:58pm
FC	Pass	11:58pm

Temperature Tests

Test	Status	Time
FC1	Pass	11:58pm
SRC	Pass	11:58pm
DET	Pass	11:58pm
BAR	Pass	11:58pm
BT	Pass	11:58pm

Blank Tests

Test	Status	Time
AIR	Pass	11:58pm


Printer Tests

Test	Status	Time
PRNT	Pass	11:59pm

CRC Tests

Test	Status	Time
COMP	Pass	11:59pm
CAL	Pass	11:59pm

Preventive Maintenance
Status: Pass

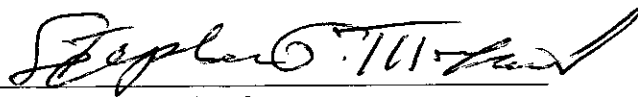

Analyst

Intox EC/IR-II: Subject Test

WAKE COUNTY BAT MOBILE UNIT 5 910

Serial Number: 008788

Test Date: 11/21/2008

Citation Number: M0000000-0

Subject's Name:

PREVENTIVE, MAINTENANCE

Subject's Date of Birth: 11/11/1911

Subject's Sex: Male

Driver's License State: XX

Driver's License Number: NONE

Analyst's Name: MORGART, STEPHEN G

Permit Number: 09372E

Effective:

12/01/2007-12/01/2009

Officer's Name: NONE, NONE

Type of Agency: FTA

Agency: DHHS

Test Type: Breath Test

Lot Number: AG722802


Exp Date: 08/15/2009

Test	g/210L	Time
DIAG	Pass	11:47pm
AIR BLK	.00	11:48pm
ACCY CHK	.07	11:48pm
AIR BLK	.00	11:49pm
SUB TEST	.00	11:50pm
AIR BLK	.00	11:51pm
SUB TEST	.00	11:52pm
AIR BLK	.00	11:53pm

Reported AC: .00 g/210L


Signature of Chemical Analyst

Court CVR


Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH

**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**

County Swain


Instrument Location MOUNT AIRY P.D.

Instrument Serial No. 008865

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 17 day of NOVEMBER, 2008 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

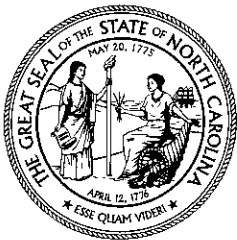

[Signature]
Signature of Certifying Official

632
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

SURRY COUNTY MOUNT AIRY PD 850

Serial Number: 008865
Test Date: 11/17/2008

Citation Number: M0000000-0
Subject's Name:
PREVENTIVE, MAINTENANCE
Subject's Date of Birth: 11/11/1911
Subject's Sex: Male
Driver's License State: XX
Driver's License Number: NONE


Analyst's Name: WEAVER, GEORGE A
Permit Number: 09442E
Effective:
12/01/2007-12/01/2009

Officer's Name: NONE, NONE
Type of Agency: FTA
Agency: DHHS
Test Type: Breath Test

Lot Number: AG825401
Exp Date: 09/10/2010

Test	g/210L	Time
DIAG	Pass	3:09pm
AIR BLK	.00	3:11pm
ACCY CHK	.07	3:11pm
AIR BLK	.00	3:12pm
SUB TEST	.00	3:13pm
AIR BLK	.00	3:14pm
SUB TEST	.00	3:15pm
AIR BLK	.00	3:16pm

Reported AC: .007g/210L

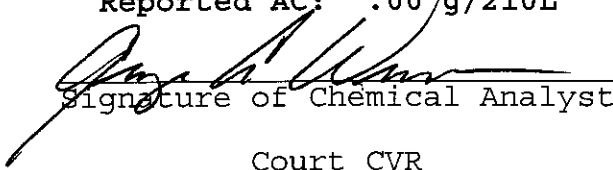
Signature of Chemical Analyst

Court CVR


Analyst

**This form is used when performing Preventive Maintenance procedures
Forensic Tests for Alcohol Branch
Department of Health and Human Services
Rev. 12/2007**

Intox EC/IR-II: Preventive Maintenance

SURRY COUNTY MOUNT AIRY PD 850

Serial Number: 008865 Test Record Number: 48
Test Date: 11/17/2008 Test Time: 3:19pm EST

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	3:19pm
FLO	Pass	3:19pm
FC	Pass	3:19pm

Temperature Tests

Test	Status	Time
FC1	Pass	3:19pm
SRC	Pass	3:19pm
DET	Pass	3:19pm
BAR	Pass	3:19pm
BT	Pass	3:19pm

Blank Tests

Test	Status	Time
AIR	Pass	3:20pm


Printer Tests

Test	Status	Time
PRNT	Pass	3:20pm

CRC Tests

Test	Status	Time
COMP	Pass	3:20pm
CAL	Pass	3:20pm

Preventive Maintenance
Status: Pass

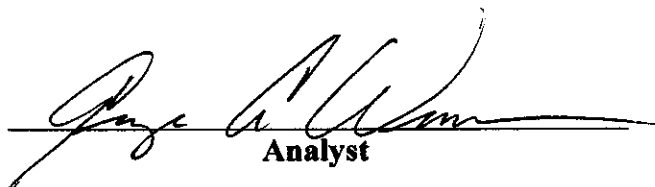
Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH

**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**


County WAKE Instrument Location KINGSDALE PS

Instrument Serial No. 008838 979 STEELE SQUARE CT. KINGSDALE, NC

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 18 day of NOVEMBER, 2008 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

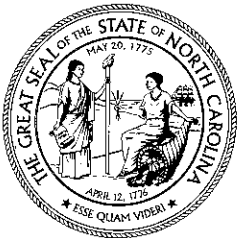

[Handwritten Signature]
Signature of Certifying Official

650
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

WAKE COUNTY KNIGHTDALE PS 910

Serial Number: 008838

Test Date: 11/18/2008

Citation Number: M0000000-0

Subject's Name:

PREVENTIVE, MAINTENANCE

Subject's Date of Birth: 11/11/1911

Subject's Sex: Male

Driver's License State: XX

Driver's License Number: NONE

Analyst's Name:

QUARANTELLA, NICHOLAS J

Permit Number: 21536E

Effective:

01/01/2008-01/01/2010

Officer's Name: NONE, NONE

Type of Agency: FTA

Agency: DHHS

Test Type: Breath Test

Lot Number: AG814101

Exp Date: 05/20/2010


Test	g/210L	Time
DIAG	Pass	11:30am
AIR BLK	.00	11:31am
ACCY CHK	.07	11:31am
AIR BLK	.00	11:32am
SUB TEST	.00	11:33am
AIR BLK	.00	11:33am
SUB TEST	.00	11:36am
AIR BLK	.00	11:37am

Reported AC: .00 g/210L



Signature of Chemical Analyst

Court CVR

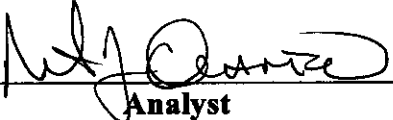

Analyst

Intox EC/IR-II: Preventive Maintenance

WAKE COUNTY KNIGHTDALE PS 910

Serial Number: 008838 Test Record Number: 89
Test Date: 11/18/2008 Test Time: 11:38am EST

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	11:39am
FLO	Pass	11:39am
FC	Pass	11:39am

Temperature Tests

Test	Status	Time
FC1	Pass	11:39am
SRC	Pass	11:39am
DET	Pass	11:39am
BAR	Pass	11:39am
BT	Pass	11:39am

Blank Tests

Test	Status	Time
AIR	Pass	11:39am


Printer Tests

Test	Status	Time
PRNT	Pass	11:39am

CRC Tests

Test	Status	Time
COMP	Pass	11:40am
CAL	Pass	11:40am

Preventive Maintenance
Status: Pass

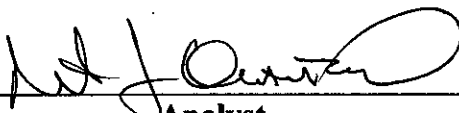

Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH

**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**


County HALIFAX Instrument Location ROANOKE RAPIDS PD

Instrument Serial No. 008656 1040 ROANOKE AVE, ROANOKE RAPIDS, NC

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 14 day of NOVEMBER, 2008 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

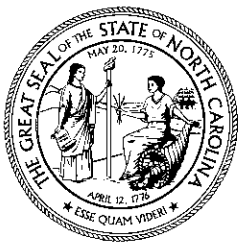

[Signature]
Signature of Certifying Official

652
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Preventive Maintenance

HALIFAX CO. ROANOKE RAPIDS PD 410

Serial Number: 008656 Test Record Number: 229
Test Date: 11/14/2008 Test Time: 2:59pm EST

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	2:59pm
FLO	Pass	2:59pm
FC	Pass	3:00pm

Temperature Tests

Test	Status	Time
FC1	Pass	3:00pm
SRC	Pass	3:00pm
DET	Pass	3:00pm
BAR	Pass	3:00pm
BT	Pass	3:00pm

Blank Tests

Test	Status	Time
AIR	Pass	3:00pm

Printer Tests

Test	Status	Time
PRNT	Pass	3:00pm

CRC Tests

Test	Status	Time
COMP	Pass	3:01pm
CAL	Pass	3:01pm

Preventive Maintenance
Status: Pass



Analyst

Intox EC/IR-II: Subject Test

HALIFAX CO. ROANOKE RAPIDS PD 410

Serial Number: 008656
Test Date: 11/14/2008

Citation Number: M0000000-0
Subject's Name:
PREVENTIVE, MAINTENANCE
Subject's Date of Birth: 11/11/1911
Subject's Sex: Male
Driver's License State: XX
Driver's License Number: NONE


Analyst's Name:
QUARANTELLA, NICHOLAS J
Permit Number: 21536E
Effective:
01/01/2008-01/01/2010

Officer's Name: NONE, NONE
Type of Agency: FTA
Agency: DHHS
Test Type: Breath Test

Lot Number: AG722702
Exp Date: 08/14/2009

Test	g/210L	Time
DIAG	Pass	2:51pm
AIR BLK	.00	2:52pm
ACCY CHK	.07	2:53pm
AIR BLK	.00	2:54pm
SUB TEST	.00	2:54pm
AIR BLK	.00	2:55pm
SUB TEST	.00	2:56pm
AIR BLK	.00	2:57pm

Reported AC: .00 g/210L

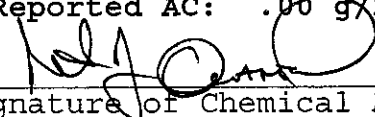

Signature of Chemical Analyst

Court CVR



Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH

**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**


County HALIFAX Instrument Location ROANOKE RAPIDS, PD

Instrument Serial No. 008635 1070 ROANOKE AVE, ROANOKE RAPIDS, NC

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 14 day of NOVEMBER, 2008 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

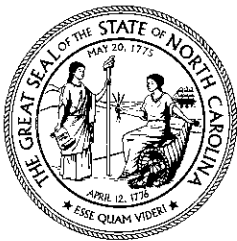

[Signature]
Signature of Certifying Official

652
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

HALIFAX CO ROANOKE RAPIDS PD 410

Serial Number: 008635
Test Date: 11/14/2008

Citation Number: M0000000-0
Subject's Name:
PREVENTIVE, MAINTENANCE
Subject's Date of Birth: 11/11/1911
Subject's Sex: Male
Driver's License State: XX
Driver's License Number: NONE


Analyst's Name:
QUARANTELLA, NICHOLAS J
Permit Number: 21536E
Effective:
01/01/2008-01/01/2010

Officer's Name: NONE, NONE
Type of Agency: FTA
Agency: DHHS
Test Type: Breath Test

Lot Number: AG722503
Exp Date: 08/12/2009

Test	g/210L	Time
DIAG	Pass	2:50pm
AIR BLK	.00	2:51pm
ACCY CHK	.07	2:52pm
AIR BLK	.00	2:53pm
SUB TEST	.00	2:54pm
AIR BLK	.00	2:54pm
SUB TEST	.00	2:56pm
AIR BLK	.00	2:57pm

Reported AC: .00 g/210L



Signature of Chemical Analyst

Court CVR



Analyst

Intox EC/IR-II: Preventive Maintenance

HALIFAX CO ROANOKE RAPIDS PD 410

Serial Number: 008635 Test Record Number: 276
Test Date: 11/14/2008 Test Time: 2:58pm EST

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	2:59pm
FLO	Pass	2:59pm
FC	Pass	2:59pm

Temperature Tests

Test	Status	Time
FC1	Pass	2:59pm
SRC	Pass	2:59pm
DET	Pass	2:59pm
BAR	Pass	2:59pm
BT	Pass	2:59pm

Blank Tests

Test	Status	Time
AIR	Pass	2:59pm

Printer Tests

Test	Status	Time
PRNT	Pass	2:59pm

CRC Tests

Test	Status	Time
COMP	Pass	3:00pm
CAL	Pass	3:00pm

Preventive Maintenance
Status: Pass



Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH

**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**


County NORTHAMPTON Instrument Location NORTHAMPTON

Instrument Serial No. 008688 105 W. JEFFERSON ST. JACKSON, NC

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 14 day of NOVEMBER, 2008 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

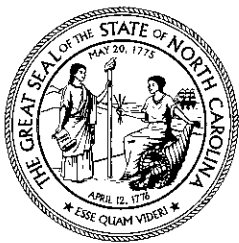

[Signature]
Signature of Certifying Official

652
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

NORTHAMPTON COUNTY SHERIFFS DEPARTMENT
650

Serial Number: 008688
Test Date: 11/14/2008

Citation Number: M0000000-0
Subject's Name:
PREVENTIVE, MAINTENANCE
Subject's Date of Birth: 11/11/1911
Subject's Sex: Male
Driver's License State: XX
Driver's License Number: NONE

Analyst's Name:
QUARANTELLA, NICHOLAS J
Permit Number: 21536E
Effective:
01/01/2008-01/01/2010

Officer's Name: NONE, NONE
Type of Agency: FTA
Agency: DHHS
Test Type: Breath Test

Lot Number: AG723401
Exp Date: 08/21/2009

Test	g/210L	Time
DIAG	Pass	1:41pm
AIR BLK	.00	1:42pm
ACCY CHK	.07	1:42pm
AIR BLK	.00	1:43pm
SUB TEST	.00	1:43pm
AIR BLK	.00	1:44pm
SUB TEST	.00	1:46pm
AIR BLK	.00	1:47pm

Reported AC: .00 g/210L



Signature of Chemical Analyst

Court CVR



Analyst

Intox EC/IR-II: Preventive Maintenance

NORTHAMPTON COUNTY SHERIFFS DEPARTMENT 650

Serial Number: 008688 Test Record Number: 387
Test Date: 11/14/2008 Test Time: 1:48pm EST

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	1:48pm
FLO	Pass	1:48pm
FC	Pass	1:48pm

Temperature Tests

Test	Status	Time
FC1	Pass	1:48pm
SRC	Pass	1:48pm
DET	Pass	1:48pm
BAR	Pass	1:48pm
BT	Pass	1:48pm

Blank Tests

Test	Status	Time
AIR	Pass	1:49pm

Printer Tests

Test	Status	Time
PRNT	Pass	1:49pm

CRC Tests

Test	Status	Time
COMP	Pass	1:49pm
CAL	Pass	1:49pm

Preventive Maintenance
Status: Pass



Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH

**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**


County NORTHAMPTON Instrument Location NORTHAMPTON CO. SHERIFFS DEPT

Instrument Serial No. 008607 105 W. JEFFERSON ST, JACKSON, NC

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 14 day of NOVEMBER, 2008 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

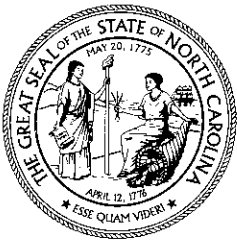

[Signature]
Signature of Certifying Official

652
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

NORTHAMPTON COUNTY SHERIFFS DEPARTMENT
650

Serial Number: 008607
Test Date: 11/14/2008

Citation Number: M0000000-0
Subject's Name:

PREVENTIVE, MAINTENANCE

Subject's Date of Birth: 11/11/1911

Subject's Sex: Male

Driver's License State: XX

Driver's License Number: NONE

Analyst's Name:

QUARANTELLA, NICHOLAS J

Permit Number: 21536E

Effective:

01/01/2008-01/01/2010

Officer's Name: NONE, NONE

Type of Agency: FTA

Agency: DHHS

Test Type: Breath Test

Lot Number: AG723401

Exp Date: 08/21/2009

Test	g/210L	Time
DIAG	Pass	1:34pm
AIR BLK	.00	1:35pm
ACCY CHK	.07	1:36pm
AIR BLK	.00	1:37pm
SUB TEST	.00	1:37pm
AIR BLK	.00	1:38pm
SUB TEST	.00	1:39pm
AIR BLK	.00	1:40pm

Reported AC: .00 g/210L



Signature of Chemical Analyst

Court CVR



Analyst

Intox EC/IR-II: Preventive Maintenance

NORTHAMPTON COUNTY SHERIFFS DEPARTMENT 650

Serial Number: 008607 Test Record Number: 414
Test Date: 11/14/2008 Test Time: 1:42pm EST

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	1:42pm
FLO	Pass	1:42pm
FC	Pass	1:42pm

Temperature Tests

Test	Status	Time
FC1	Pass	1:42pm
SRC	Pass	1:42pm
DET	Pass	1:42pm
BAR	Pass	1:42pm
BT	Pass	1:42pm

Blank Tests

Test	Status	Time
AIR	Pass	1:43pm

Printer Tests

Test	Status	Time
PRNT	Pass	1:43pm

CRC Tests

Test	Status	Time
COMP	Pass	1:43pm
CAL	Pass	1:43pm

Preventive Maintenance
Status: Pass



Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH

**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**

County HALIFAX

Instrument Location HALIFAX CO. SHERIFF DEPT


Instrument Serial No. 008695

FERRILL LN HALIFAX, NC

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 14 day of NOVEMBER, 2008 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

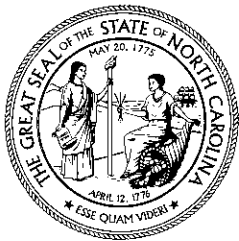

[Handwritten Signature]
Signature of Certifying Official

652
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

HALIFAX CO. HALIFAX CO. SD 410

Serial Number: 008695
Test Date: 11/14/2008

Citation Number: M0000000-0
Subject's Name:
PREVENTIVE, MAINTENANCE
Subject's Date of Birth: 11/11/1911
Subject's Sex: Male
Driver's License State: XX
Driver's License Number: NONE


Analyst's Name:
QUARANTELLA, NICHOLAS J
Permit Number: 21536E
Effective:
01/01/2008-01/01/2010

Officer's Name: NONE, NONE
Type of Agency: FTA
Agency: DHHS
Test Type: Breath Test

Lot Number: AG722503
Exp Date: 08/12/2009

Test	g/210L	Time
DIAG	Pass	12:26pm
AIR BLK	.00	12:27pm
ACCY CHK	.08	12:27pm
AIR BLK	.00	12:28pm
SUB TEST	.00	12:29pm
AIR BLK	.00	12:30pm
SUB TEST	.00	12:31pm
AIR BLK	.00	12:32pm

Reported AC: .00 g/210L

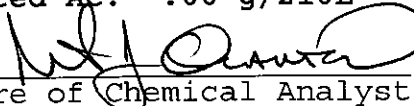

Signature of Chemical Analyst

Court CVR



Analyst

Intox EC/IR-II: Preventive Maintenance

HALIFAX CO. HALIFAX CO. SD 410

Serial Number: 008695 Test Record Number: 218
Test Date: 11/14/2008 Test Time: 12:33pm EST

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	12:34pm
FLO	Pass	12:34pm
FC	Pass	12:34pm

Temperature Tests

Test	Status	Time
FC1	Pass	12:34pm
SRC	Pass	12:34pm
DET	Pass	12:34pm
BAR	Pass	12:34pm
BT	Pass	12:34pm

Blank Tests

Test	Status	Time
AIR	Pass	12:35pm

Printer Tests

Test	Status	Time
PRNT	Pass	12:35pm

CRC Tests

Test	Status	Time
COMP	Pass	12:35pm
CAL	Pass	12:35pm

Preventive Maintenance
Status: Pass



Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH

**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**


County Wake Instrument Location BAT Mobile Unit #5

Instrument Serial No. 008788 Apex

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 15th day of NOVEMBER, 2008 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

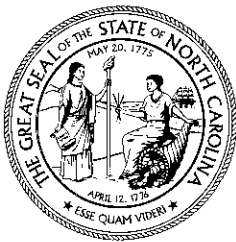

[Handwritten Signature]
Signature of Certifying Official

636
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Preventive Maintenance

WAKE COUNTY BAT MOBILE UNIT 5 910

Serial Number: 008788 Test Record Number: 127
Test Date: 11/15/2008 Test Time: 11:24pm EST

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	11:24pm
FLO	Pass	11:24pm
FC	Pass	11:24pm

Temperature Tests

Test	Status	Time
FC1	Pass	11:24pm
SRC	Pass	11:24pm
DET	Pass	11:24pm
BAR	Pass	11:24pm
BT	Pass	11:24pm

Blank Tests

Test	Status	Time
AIR	Pass	11:25pm


Printer Tests

Test	Status	Time
PRNT	Pass	11:25pm

CRC Tests

Test	Status	Time
COMP	Pass	11:25pm
CAL	Pass	11:25pm

Preventive Maintenance
Status: Pass

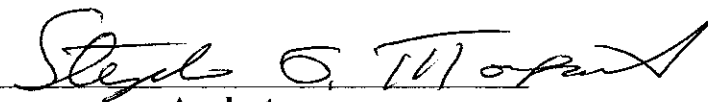

Analyst

Intox EC/IR-II: Subject Test

WAKE COUNTY BAT MOBILE UNIT 5 910

Serial Number: 008788
Test Date: 11/15/2008

Citation Number: M0000000-0

Subject's Name:

PREVENTIVE, MAINTENANCE

Subject's Date of Birth: 11/11/1911

Subject's Sex: Male

Driver's License State: XX

Driver's License Number: NONE

Analyst's Name: MORGART, STEPHEN G

Permit Number: 09372E

Effective:

12/01/2007-12/01/2009

Officer's Name: NONE, NONE

Type of Agency: FTA

Agency: DHHS


Test Type: Breath Test

Lot Number: AG722802


Exp Date: 08/15/2009

Test	g/210L	Time
DIAG	Pass	11:15pm
AIR BLK	.000	11:16pm
ACCY CHK	.080	11:16pm
AIR BLK	.000	11:17pm
SUB TEST	.000	11:18pm
AIR BLK	.000	11:19pm
SUB TEST	.000	11:21pm
AIR BLK	.000	11:22pm

Reported AC: .00 g/210L

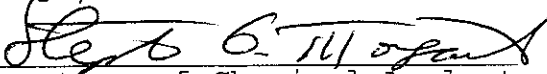
Signature of Chemical Analyst

Court CVR

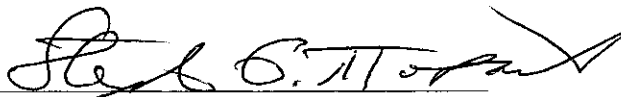
Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH

**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**


County Wake Instrument Location BATMOBILE UNIT #5

Instrument Serial No. 008698 APEX

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 15TH day of NOVEMBER, 2008 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

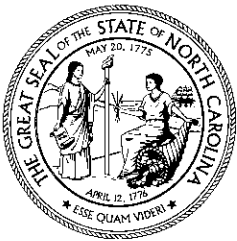

Sept. C. T. T. T. T.
Signature of Certifying Official

636
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Preventive Maintenance

WAKE COUNTY BAT MOBILE UNIT 5 910

Serial Number: 008698 Test Record Number: 251
Test Date: 11/15/2008 Test Time: 11:24pm EST

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	11:24pm
FLO	Pass	11:24pm
FC	Pass	11:24pm

Temperature Tests

Test	Status	Time
FC1	Pass	11:24pm
SRC	Pass	11:24pm
DET	Pass	11:24pm
BAR	Pass	11:24pm
BT	Pass	11:24pm

Blank Tests

Test	Status	Time
AIR	Pass	11:25pm


Printer Tests

Test	Status	Time
PRNT	Pass	11:25pm

CRC Tests

Test	Status	Time
COMP	Pass	11:25pm
CAL	Pass	11:25pm

Preventive Maintenance
Status: Pass

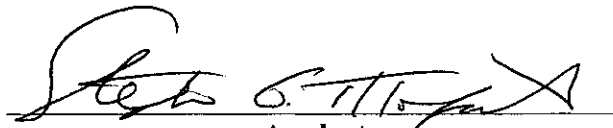
Analyst

Intox EC/IR-II: Subject Test

WAKE COUNTY BAT MOBILE UNIT 5 910

Serial Number: 008698
Test Date: 11/15/2008

Citation Number: M0000000-0
Subject's Name:
PREVENTIVE, MAINTENANCE
Subject's Date of Birth: 11/11/1911
Subject's Sex: Male
Driver's License State: XX
Driver's License Number: NONE


Analyst's Name: MORGART, STEPHEN G
Permit Number: 09372E
Effective:
12/01/2007-12/01/2009

Officer's Name: MORGART, STEPHEN G
Type of Agency: FTA
Agency: DHHS
Test Type: Breath Test


Lot Number: AG723401
Exp Date: 08/21/2009

Test	g/210L	Time
DIAG	Pass	11:11pm
AIR BLK	.00	11:13pm
ACCY CHK	.07	11:13pm
AIR BLK	.00	11:14pm
SUB TEST	.00	11:15pm
AIR BLK	.00	11:16pm
SUB TEST	.00	11:17pm
AIR BLK	.00	11:18pm

Reported AC: .00 g/210L

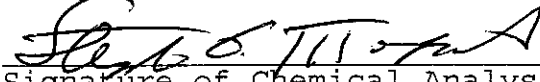
Signature of Chemical Analyst

Court CVR.

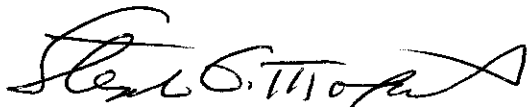

Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH


**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**

County Wake Instrument Location Bot Mobile Unit #5
Instrument Serial No. 008600 Apex

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 15th day of NOVEMBER, 2008 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

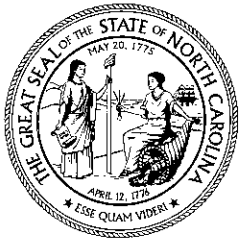

Steph T. H. [Signature]

Signature of Certifying Official

636

Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Preventive Maintenance

WAKE COUNTY BAT MOBILE UNIT 5 910

Serial Number: 008600 Test Record Number: 408
Test Date: 11/15/2008 Test Time: 11:18pm EST

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	11:18pm
FLO	Pass	11:18pm
FC	Pass	11:18pm

Temperature Tests

Test	Status	Time
FC1	Pass	11:18pm
SRC	Pass	11:18pm
DET	Pass	11:18pm
BAR	Pass	11:18pm
BT	Pass	11:18pm

Blank Tests

Test	Status	Time
AIR	Pass	11:19pm


Printer Tests

Test	Status	Time
PRNT	Pass	11:19pm

CRC Tests

Test	Status	Time
COMP	Pass	11:19pm
CAL	Pass	11:19pm

Preventive Maintenance
Status: Pass

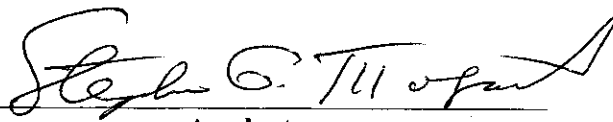

Analyst

Intox EC/IR-II: Subject Test

WAKE COUNTY BAT MOBILE UNIT 5 910

Serial Number: 008600
Test Date: 11/15/2008

Citation Number: M0000000-0
Subject's Name:
PREVENTIVE, MAINTENANCE
Subject's Date of Birth: 11/11/1911
Subject's Sex: Male
Driver's License State: XX
Driver's License Number: NONE


Analyst's Name: MORGART, STEPHEN G
Permit Number: 09372E
Effective:
12/01/2007-12/01/2009

Officer's Name: NONE, NONE
Type of Agency: FTA
Agency: DHHS
Test Type: Breath Test


Lot Number: AG723401
Exp Date: 08/21/2009

Test	g/210L	Time
DIAG	Pass	11:06pm
AIR BLK	.00	11:07pm
ACCY CHK	.07	11:08pm
AIR BLK	.00	11:09pm
SUB TEST	.00	11:09pm
AIR BLK	.00	11:10pm
SUB TEST	.00	11:12pm
AIR BLK	.00	11:13pm

Reported AC: .00 g/210L

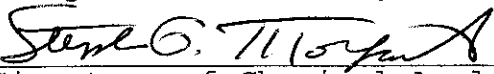
Signature of Chemical Analyst

Court CVR

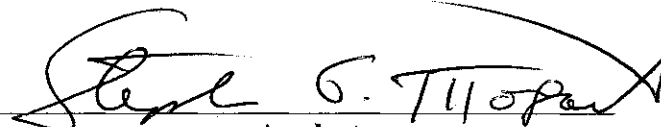
Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH

**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**


County DURHAM Instrument Location DURHAM CO. JAIL

Instrument Serial No. 008891 217 S. MANGUM ST. DURHAM, NC

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 12 day of NOVEMBER, 2008 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

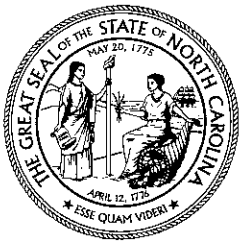

Brian D. Smith
Signature of Certifying Official

637
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

DURHAM COUNTY DURHAM COUNTY JAIL 310

Serial Number: 008891
Test Date: 11/12/2008

Citation Number: M0000000-0
Subject's Name:
PREVENTIVE, MAINTENANCE
Subject's Date of Birth: 11/11/1911
Subject's Sex: Male
Driver's License State: XX
Driver's License Number: NONE

Analyst's Name: SMITH, BRIAN D
Permit Number: 08937E
Effective:
12/01/2007-12/01/2009

Officer's Name: NONE, NONE
Type of Agency: FTA
Agency: DHHS
Test Type: Breath Test

Lot Number: AG814002
Exp Date: 05/19/2010


Test	g/210L	Time
DIAG	Pass	2:17pm
AIR BLK	.00	2:17pm
ACCY CHK	.07	2:18pm
AIR BLK	.00	2:19pm
SUB TEST	.00	2:19pm
AIR BLK	.00	2:20pm
SUB TEST	.00	2:22pm
AIR BLK	.00	2:23pm

Reported AC: .00 g/210L



Signature of Chemical Analyst

Court CVR



Analyst

Intox EC/IR-II: Preventive Maintenance

DURHAM COUNTY DURHAM COUNTY JAIL 310

Serial Number: 008891 Test Record Number: 135
Test Date: 11/12/2008 Test Time: 2:25pm EST

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	2:25pm
FLO	Pass	2:25pm
FC	Pass	2:25pm

Temperature Tests

Test	Status	Time
FC1	Pass	2:26pm
SRC	Pass	2:26pm
DET	Pass	2:26pm
BAR	Pass	2:26pm
BT	Pass	2:26pm

Blank Tests

Test	Status	Time
AIR	Pass	2:26pm


Printer Tests

Test	Status	Time
PRNT	Pass	2:26pm

CRC Tests

Test	Status	Time
COMP	Pass	2:26pm
CAL	Pass	2:26pm

Preventive Maintenance
Status: Pass

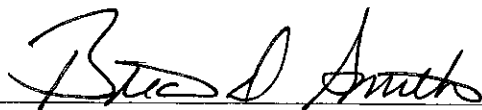

Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH

**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**


County DURHAM Instrument Location DURHAM Co. JAIL

Instrument Serial No. 008878 217 S. MANGUM ST. DURHAM, NC

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 12 day of NOVEMBER, 2008 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

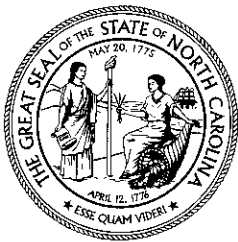

Brian D. Smith
Signature of Certifying Official

637
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

DURHAM COUNTY DURHAM COUNTY JAIL 310

Serial Number: 008878
Test Date: 11/12/2008

Citation Number: M0000000-0
Subject's Name:
PREVENTIVE, MAINTENANCE
Subject's Date of Birth: 11/11/1911
Subject's Sex: Male
Driver's License State: XX
Driver's License Number: NONE

Analyst's Name: SMITH, BRIAN D
Permit Number: 08937E
Effective:
12/01/2007-12/01/2009

Officer's Name: NONE, NONE
Type of Agency: FTA
Agency: DHHS
Test Type: Breath Test

Lot Number: AG814002
Exp Date: 05/19/2010

Test	g/210L	Time
DIAG	Pass	2:16pm
AIR BLK	.00	2:16pm
ACCY CHK	.07	2:17pm
AIR BLK	.00	2:18pm
SUB TEST	.00	2:18pm
AIR BLK	.00	2:19pm
SUB TEST	.00	2:20pm
AIR BLK	.00	2:21pm

Reported AC: .00, g/210L


Signature of Chemical Analyst

Court CVR


Analyst

Intox EC/IR-II: Preventive Maintenance

DURHAM COUNTY DURHAM COUNTY JAIL 310

Serial Number: 008878 Test Record Number: 156
Test Date: 11/12/2008 Test Time: 2:22pm EST

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	2:23pm
FLO	Pass	2:23pm
FC	Pass	2:23pm

Temperature Tests

Test	Status	Time
FC1	Pass	2:23pm
SRC	Pass	2:23pm
DET	Pass	2:23pm
BAR	Pass	2:23pm
BT	Pass	2:23pm

Blank Tests

Test	Status	Time
AIR	Pass	2:23pm


Printer Tests

Test	Status	Time
PRNT	Pass	2:23pm

CRC Tests

Test	Status	Time
COMP	Pass	2:24pm
CAL	Pass	2:24pm

Preventive Maintenance
Status: Pass

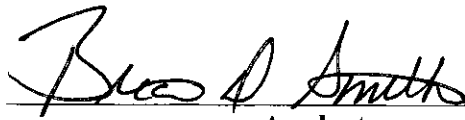

Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH

**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**


County DURHAM Instrument Location DURHAM CO. JAIL

Instrument Serial No. 008859 217 S. MANGUM ST. DURHAM, NC

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 12 day of NOVEMBER, 2008 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

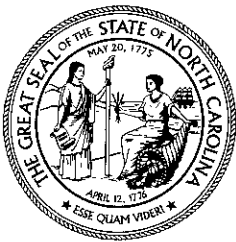

Brian D. Smith
Signature of Certifying Official

637
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

DURHAM COUNTY DURHAM COUNTY JAIL 310

Serial Number: 008859
Test Date: 11/12/2008

Citation Number: M0000000-0
Subject's Name:
PREVENTIVE, MAINTENANCE
Subject's Date of Birth: 11/11/1911
Subject's Sex: Male
Driver's License State: XX
Driver's License Number: NONE


Analyst's Name: SMITH, BRIAN D
Permit Number: 08937E
Effective:
12/01/2007-12/01/2009

Officer's Name: NONE, NONE
Type of Agency: FTA
Agency: DHHS
Test Type: Breath Test

Lot Number: AG814002
Exp Date: 05/19/2010


Test	g/210L	Time
DIAG	Pass	2:14pm
AIR BLK	.00	2:15pm
ACCY CHK	.07	2:16pm
AIR BLK	.00	2:16pm
SUB TEST	.00	2:17pm
AIR BLK	.00	2:18pm
SUB TEST	.00	2:20pm
AIR BLK	.00	2:21pm

Reported AC: .00 g/210L

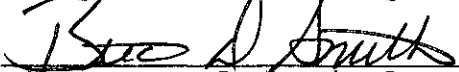

Signature of Chemical Analyst

Court CVR

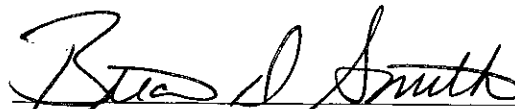

Analyst

Intox EC/IR-II: Preventive Maintenance

DURHAM COUNTY DURHAM COUNTY JAIL 310

Serial Number: 008859 Test Record Number: 126
Test Date: 11/12/2008 Test Time: 2:21pm EST

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	2:22pm
FLO	Pass	2:22pm
FC	Pass	2:22pm

Temperature Tests

Test	Status	Time
FC1	Pass	2:22pm
SRC	Pass	2:22pm
DET	Pass	2:22pm
BAR	Pass	2:22pm
BT	Pass	2:22pm

Blank Tests

Test	Status	Time
AIR	Pass	2:23pm

Printer Tests

Test	Status	Time
PRNT	Pass	2:23pm

CRC Tests

Test	Status	Time
COMP	Pass	2:23pm
CAL	Pass	2:23pm

Preventive Maintenance
Status: Pass



Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH

**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**


County WARREN Instrument Location NORLINA PD

Instrument Serial No. 008651 101 MAIN ST. NORLINA, NC

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 04 day of NOVEMBER, 2008 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

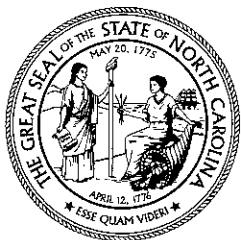

Bud A Smith
Signature of Certifying Official

637
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

WARREN COUNTY NORLINA POLICE DEPT. 920

Serial Number: 008651
Test Date: 11/04/2008

Citation Number: M0000000-0
Subject's Name:
PREVENTIVE, MAINTENANCE
Subject's Date of Birth: 11/11/1911
Subject's Sex: Male
Driver's License State: XX
Driver's License Number: NONE

Analyst's Name: SMITH, BRIAN D
Permit Number: 08937E
Effective:
12/01/2007-12/01/2009

Officer's Name: NONE, NONE
Type of Agency: FTA
Agency: DHHS
Test Type: Breath Test

Lot Number: AG723401
Exp Date: 08/21/2009

Test	g/210L	Time
DIAG	Pass	1:59pm
AIR BLK	.00	2:00pm
ACCY CHK	.07	2:01pm
AIR BLK	.00	2:02pm
SUB TEST	.00	2:03pm
AIR BLK	.00	2:04pm
SUB TEST	.00	2:05pm
AIR BLK	.00	2:06pm

Reported AC: .00 g/210L



Signature of Chemical Analyst

Court CVR



Analyst

Intox EC/IR-II: Preventive Maintenance

WARREN COUNTY NORLINA POLICE DEPT. 920

Serial Number: 008651 Test Record Number: 342
Test Date: 11/04/2008 Test Time: 2:20pm EST

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	2:20pm
FLO	Pass	2:20pm
FC	Pass	2:20pm

Temperature Tests

Test	Status	Time
FC1	Pass	2:21pm
SRC	Pass	2:21pm
DET	Pass	2:21pm
BAR	Pass	2:21pm
BT	Pass	2:21pm

Blank Tests

Test	Status	Time
AIR	Pass	2:21pm

Printer Tests

Test	Status	Time
PRNT	Pass	2:21pm

CRC Tests

Test	Status	Time
COMP	Pass	2:21pm
CAL	Pass	2:21pm

Preventive Maintenance
Status: Pass



Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH

**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**


County Randolph Instrument Location Arckdale Police

Instrument Serial No. 008791 Department

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 12 day of November, 20 08 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

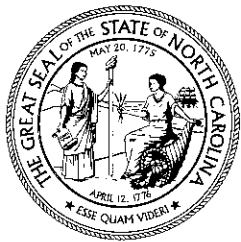

[Signature]
Signature of Certifying Official

642
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

RANDOLPH COUNTY ARCHDALE PD 750

Serial Number: 008791
Test Date: 11/12/2008

Citation Number: M0000000-0
Subject's Name:

PREVENTIVE, MAINTENANCE

Subject's Date of Birth: 11/11/1911

Subject's Sex: Male

Driver's License State: XX

Driver's License Number: NONE

Analyst's Name: DEAN, LARRY K

Permit Number: 11598E

Effective:

12/01/2007-12/01/2009

Officer's Name: NONE, NONE

Type of Agency: FTA

Agency: DHHS


Test Type: Breath Test

Lot Number: AG722601


Exp Date: 08/13/2009

Test	g/210L	Time
DIAG	Pass	11:39am
AIR BLK	.00	11:40am
ACCY CHK	.07	11:40am
AIR BLK	.00	11:41am
SUB TEST	.00	11:42am
AIR BLK	.00	11:43am
SUB TEST	.00	11:44am
AIR BLK	.00	11:45am

Reported AC: .00 g/210L

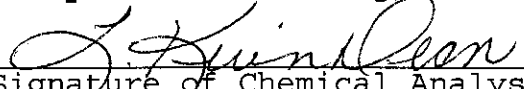
Signature of Chemical Analyst

Court CVR

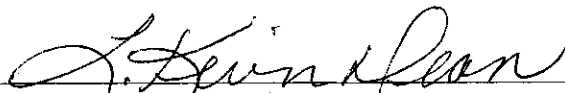
Analyst

Intox EC/IR-II: Preventive Maintenance

RANDOLPH COUNTY ARCHDALE PD 750

Serial Number: 008791 Test Record Number: 149
Test Date: 11/12/2008 Test Time: 11:46am EST

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	11:47am
FLO	Pass	11:47am
FC	Pass	11:47am

Temperature Tests

Test	Status	Time
FC1	Pass	11:47am
SRC	Pass	11:47am
DET	Pass	11:47am
BAR	Pass	11:47am
BT	Pass	11:47am

Blank Tests

Test	Status	Time
AIR	Pass	11:47am


Printer Tests

Test	Status	Time
PRNT	Pass	11:47am

CRC Tests

Test	Status	Time
COMP	Pass	11:48am
CAL	Pass	11:48am

Preventive Maintenance
Status: Pass

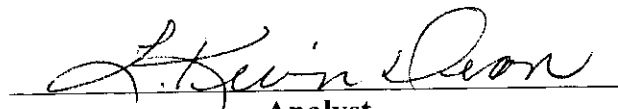
Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH

**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**


County Rockingham Instrument Location Rockingham Co Jail

Instrument Serial No. 008796

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 10 day of November, 2008 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

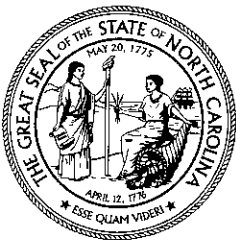

[Handwritten Signature]
Signature of Certifying Official

632
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

ROCKINGHAM COUNTY ROCKINGHAM CO JAIL
780

Serial Number: 008796
Test Date: 11/10/2008

Citation Number: M0000000-0
Subject's Name:
PREVENTIVE, MAINTENANCE
Subject's Date of Birth: 11/11/1911
Subject's Sex: Male
Driver's License State: XX
Driver's License Number: NONE


Analyst's Name: WEAVER, GEORGE A
Permit Number: 09442E
Effective:
12/01/2007-12/01/2009

Officer's Name: NONE, NONE
Type of Agency: FTA
Agency: DHHS
Test Type: Breath Test


Lot Number: AG722601
Exp Date: 08/13/2009

Test	g/210L	Time
DIAG	Pass	4:02pm
AIR BLK	.00	4:03pm
ACCY CHK	.07	4:04pm
AIR BLK	.00	4:05pm
SUB TEST	.00	4:05pm
AIR BLK	.00	4:06pm
SUB TEST	.00	4:08pm
AIR BLK	.00	4:09pm

Reported AC: .00 g/210L

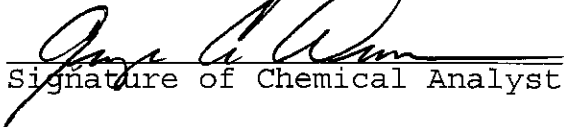
Signature of Chemical Analyst

Court CVR

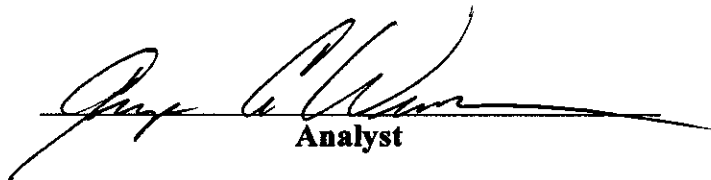
Analyst

Intox EC/IR-II: Preventive Maintenance

ROCKINGHAM COUNTY ROCKINGHAM CO JAIL 780

Serial Number: 008796 Test Record Number: 208
Test Date: 11/10/2008 Test Time: 4:10pm EST

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	4:10pm
FLO	Pass	4:10pm
FC	Pass	4:10pm

Temperature Tests

Test	Status	Time
FC1	Pass	4:10pm
SRC	Pass	4:10pm
DET	Pass	4:10pm
BAR	Pass	4:10pm
BT	Pass	4:10pm

Blank Tests

Test	Status	Time
AIR	Pass	4:11pm


Printer Tests

Test	Status	Time
PRNT	Pass	4:11pm

CRC Tests

Test	Status	Time
COMP	Pass	4:11pm
CAL	Pass	4:11pm

Preventive Maintenance
Status: Pass



Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH

**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**


County Rockingham Instrument Location Reidsville P.A.

Instrument Serial No. 008784

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 10 day of NOVEMBER, 2008 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

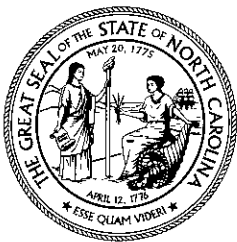

[Signature]
Signature of Certifying Official

632
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

ROCKINGHAM COUNTY REIDSVILLE PD 780

Serial Number: 008784
Test Date: 11/10/2008

Citation Number: M0000000-0

Subject's Name:

PREVENTIVE, MAINTENANCE

Subject's Date of Birth: 11/11/1911

Subject's Sex: Male

Driver's License State: XX

Driver's License Number: NONE

Analyst's Name: WEAVER, GEORGE A

Permit Number: 09442E

Effective:

12/01/2007-12/01/2009

Officer's Name: NONE, NONE

Type of Agency: FTA

Agency: DHHS


Test Type: Breath Test

Lot Number: AG722601-07


Exp Date: 08/13/2009

Test	g/210L	Time
DIAG	Pass	3:18pm
AIR BLK	.00	3:19pm
ACCY CHK	.07	3:20pm
AIR BLK	.00	3:21pm
SUB TEST	.00	3:22pm
AIR BLK	.00	3:23pm
SUB TEST	.00	3:26pm
AIR BLK	.00	3:27pm

Reported AC: .00 g/210L

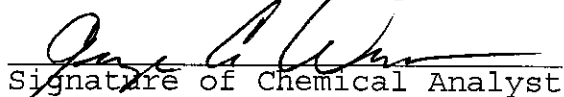
Signature of Chemical Analyst

Court CVR

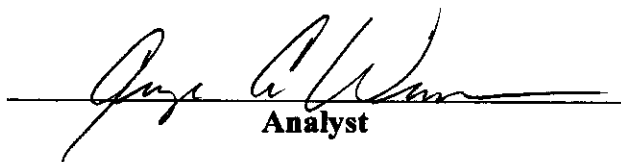
Analyst

Intox EC/IR-II: Preventive Maintenance

ROCKINGHAM COUNTY REIDSVILLE PD 780

Serial Number: 008784 Test Record Number: 130
Test Date: 11/10/2008 Test Time: 3:27pm EST

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	3:28pm
FLO	Pass	3:28pm
FC	Pass	3:28pm

Temperature Tests

Test	Status	Time
FC1	Pass	3:28pm
SRC	Pass	3:28pm
DET	Pass	3:28pm
BAR	Pass	3:28pm
BT	Pass	3:28pm

Blank Tests

Test	Status	Time
AIR	Pass	3:29pm

Printer Tests

Test	Status	Time
PRNT	Pass	3:29pm

CRC Tests

Test	Status	Time
COMP	Pass	3:29pm
CAL	Pass	3:29pm

Preventive Maintenance
Status: Pass


Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH

**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**


County Rockingham Instrument Location MADISON P.D.

Instrument Serial No. 008638

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 10 day of NOVEMBER, 2008 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

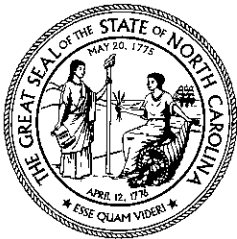

[Signature]
Signature of Certifying Official

632
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

ROCKINGHAM COUNTY MADISON PD 780

Serial Number: 008638
Test Date: 11/10/2008

Citation Number: M0000000-0

Subject's Name:

PREVENTIVE, MAINTENANCE

Subject's Date of Birth: 11/11/1911

Subject's Sex: Male

Driver's License State: XX

Driver's License Number: NONE

Analyst's Name: WEAVER, GEORGE A

Permit Number: 09442E

Effective:

12/01/2007-12/01/2009

Officer's Name: NONE, NONE

Type of Agency: FTA

Agency: DHHS


Test Type: Breath Test

Lot Number: AG722601-18

Exp Date: 08/13/2009

Test	g/210L	Time
DIAG	Pass	8:15pm
AIR BLK	.00	8:16pm
ACCY CHK	.08	8:17pm
AIR BLK	.00	8:18pm
SUB TEST	.00	8:18pm
AIR BLK	.00	8:19pm
SUB TEST	.00	8:21pm
AIR BLK	.00	8:22pm

Reported AC: 00 g/210L

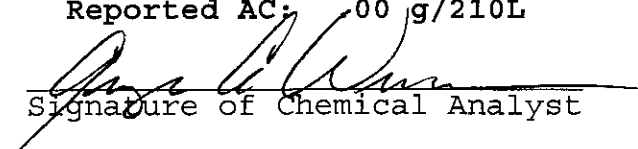
Signature of Chemical Analyst

Court CVR


Analyst

**This form is used when performing Preventive Maintenance procedures
Forensic Tests for Alcohol Branch
Department of Health and Human Services
Rev. 12/2007**

Intox EC/IR-II: Preventive Maintenance

ROCKINGHAM COUNTY MADISON PD 780

Serial Number: 008638 Test Record Number: 298
Test Date: 11/10/2008 Test Time: 8:23pm EST

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	8:24pm
FLO	Pass	8:24pm
FC	Pass	8:24pm

Temperature Tests

Test	Status	Time
FC1	Pass	8:24pm
SRC	Pass	8:24pm
DET	Pass	8:24pm
BAR	Pass	8:24pm
BT	Pass	8:24pm

Blank Tests

Test	Status	Time
AIR	Pass	8:25pm

Printer Tests

Test	Status	Time
PRNT	Pass	8:25pm

CRC Tests

Test	Status	Time
COMP	Pass	8:25pm
CAL	Pass	8:25pm

Preventive Maintenance
Status: Pass


Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH

**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**


County Rockingham Instrument Location Cotter P.D.

Instrument Serial No. 008636

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 10 day of November, 2008 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

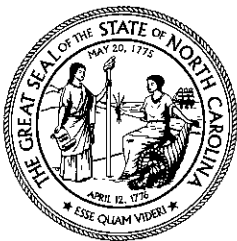

[Signature]
Signature of Certifying Official

632
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

ROCKINGHAM COUNTY EDEN PD 780

Serial Number: 008636
Test Date: 11/10/2008

Citation Number: M0000000-0
Subject's Name:
PREVENTIVE, MAINTENANCE
Subject's Date of Birth: 11/11/1911
Subject's Sex: Male
Driver's License State: XX
Driver's License Number: NONE


Analyst's Name: WEAVER, GEORGE A
Permit Number: 09442E
Effective:
12/01/2007-12/01/2009

Officer's Name: NONE, NONE
Type of Agency: FTA
Agency: DHHS
Test Type: Breath Test

Lot Number: AG722701-07
Exp Date: 08/14/2009

Test	g/210L	Time
DIAG	Pass	2:36pm
AIR BLK	.00	2:37pm
ACCY CHK	.07	2:37pm
AIR BLK	.00	2:38pm
SUB TEST	.00	2:39pm
AIR BLK	.00	2:40pm
SUB TEST	.00	2:41pm
AIR BLK	.00	2:42pm

Reported AC: .00 g/210L

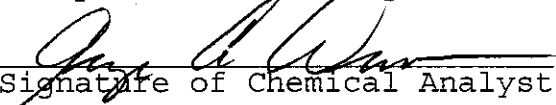
Signature of Chemical Analyst

Court CVR


Analyst

Intox EC/IR-II: Preventive Maintenance

ROCKINGHAM COUNTY EDEN PD 780

Serial Number: 008636 Test Record Number: 494
Test Date: 11/10/2008 Test Time: 2:43pm EST

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	2:44pm
FLO	Pass	2:44pm
FC	Pass	2:44pm

Temperature Tests

Test	Status	Time
FC1	Pass	2:44pm
SRC	Pass	2:44pm
DET	Pass	2:44pm
BAR	Pass	2:44pm
BT	Pass	2:44pm

Blank Tests

Test	Status	Time
AIR	Pass	2:44pm

Printer Tests

Test	Status	Time
PRNT	Pass	2:44pm

CRC Tests

Test	Status	Time
COMP	Pass	2:45pm
CAL	Pass	2:45pm

Preventive Maintenance
Status: Pass


Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH

**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**


County Granville Instrument Location BFI Mobile Unit 4

Instrument Serial No. 008717

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 1st day of November, 2008 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

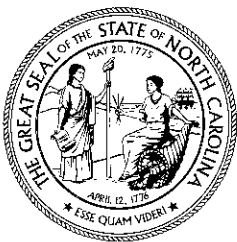

[Signature]
Signature of Certifying Official

657
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Preventive Maintenance

GRANVILLE COUNTY BAT MOBILE UNIT 4 380

Serial Number: 008717 Test Record Number: 115
Test Date: 11/01/2008 Test Time: 9:46pm EDT

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	9:46pm
FLO	Pass	9:46pm
FC	Pass	9:46pm

Temperature Tests

Test	Status	Time
FC1	Pass	9:46pm
SRC	Pass	9:46pm
DET	Pass	9:46pm
BAR	Pass	9:46pm
BT	Pass	9:46pm

Blank Tests

Test	Status	Time
AIR	Pass	9:47pm


Printer Tests

Test	Status	Time
PRNT	Pass	9:47pm

CRC Tests

Test	Status	Time
COMP	Pass	9:47pm
CAL	Pass	9:47pm

Preventive Maintenance
Status: Pass

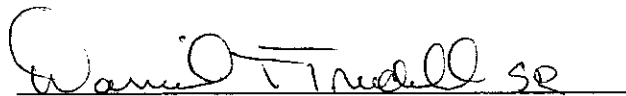
Analyst

Intox EC/IR-II: Subject Test

GRANVILLE COUNTY BAT MOBILE UNIT 4 380

Serial Number: 008717
Test Date: 11/01/2008

Citation Number: M0000000-0

Subject's Name:

PREVENTIVE, MAINTENANCE

Subject's Date of Birth: 11/11/1911

Subject's Sex: Male

Driver's License State: XX

Driver's License Number: NONE

Analyst's Name: TRUDELL, SR., DANIEL T

Permit Number: 21535E

Effective:

12/01/2007-12/01/2009

Officer's Name: NONE, NONE

Type of Agency: FTA

Agency: DHHS


Test Type: Breath Test

Lot Number: AG723401


Exp Date: 08/21/2009

Test	g/210L	Time
DIAG	Pass	9:38pm
AIR BLK	.00	9:39pm
ACCY CHK	.07	9:39pm
AIR BLK	.00	9:40pm
SUB TEST	.00	9:41pm
AIR BLK	.00	9:42pm
SUB TEST	.00	9:43pm
AIR BLK	.00	9:44pm

Reported AC: .00 g/210L

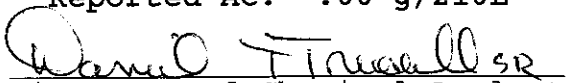
Signature of Chemical Analyst

Court CVR

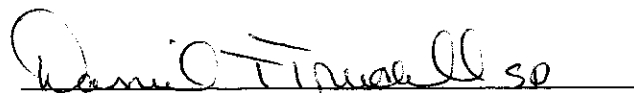
Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH


PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II

County GASTON Instrument Location BAT MOBILE UNIT 3
Instrument Serial No. 008707 BELMONT, NC

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 02 day of NOVEMBER, 2008 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

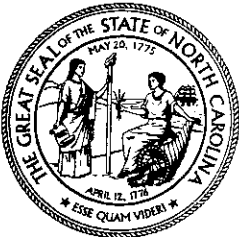

Alvin R. Burns
Signature of Certifying Official

648
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

GASTON COUNTY BAT MOBILE UNIT 3 350

Serial Number: 008707

Test Date: 11/02/2008

Citation Number: M0000000-0

Subject's Name:

PREVENTIVE, MAINTENANCE

Subject's Date of Birth: 11/11/1911

Subject's Sex: Male

Driver's License State: XX

Driver's License Number: NONE

Analyst's Name: BARNES, ALVIN R

Permit Number: 15671E

Effective:

12/01/2007-12/01/2009

Officer's Name: NONE, NONE

Type of Agency: FTA

Agency: DHHS

Test Type: Breath Test

Lot Number: AG722501


Exp Date: 08/12/2009

Test	g/210L	Time
DIAG	Pass	7:26pm
AIR BLK	.00	7:27pm
ACCY CHK	.07	7:28pm
AIR BLK	.00	7:29pm
SUB TEST	.00	7:29pm
AIR BLK	.00	7:30pm
SUB TEST	.00	7:32pm
AIR BLK	.00	7:32pm

Reported AC: .00 g/210L

Signature of Chemical Analyst

Court CVR

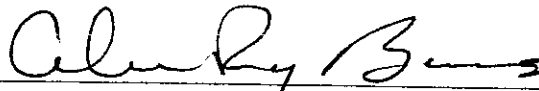
Analyst

Intox EC/IR-II: Preventive Maintenance

GASTON COUNTY BAT MOBILE UNIT 3 350

Serial Number: 008707 Test Record Number: 206
Test Date: 11/02/2008 Test Time: 7:34pm EST

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	7:34pm
FLO	Pass	7:34pm
FC	Pass	7:34pm

Temperature Tests

Test	Status	Time
FC1	Pass	7:34pm
SRC	Pass	7:34pm
DET	Pass	7:34pm
BAR	Pass	7:34pm
BT	Pass	7:34pm

Blank Tests

Test	Status	Time
AIR	Pass	7:35pm


Printer Tests

Test	Status	Time
PRNT	Pass	7:35pm

CRC Tests

Test	Status	Time
COMP	Pass	7:35pm
CAL	Pass	7:35pm

Preventive Maintenance
Status: Pass

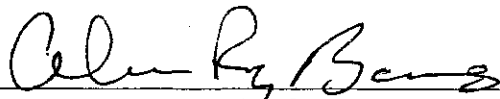
Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH


**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**

County GASTON Instrument Location BAT MOBILE UNIT 3
Instrument Serial No. 008647 BELMONT, NC

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 02 day of NOVEMBER, 2008 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

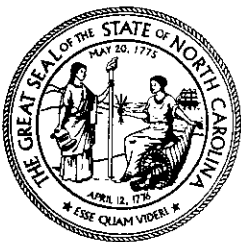

Alan R. B...
Signature of Certifying Official

6048
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

GASTON COUNTY BAT MOBILE UNIT 3 350

Serial Number: 008647

Test Date: 11/02/2008

Citation Number: M0000000-0

Subject's Name:

PREVENTIVE, MAINTENANCE

Subject's Date of Birth: 11/11/1911

Subject's Sex: Male

Driver's License State: XX

Driver's License Number: NONE

Analyst's Name: BARNES, ALVIN R

Permit Number: 15671E

Effective:

12/01/2007-12/01/2009

Officer's Name: NONE, NONE

Type of Agency: FTA

Agency: DHHS

Test Type: Breath Test

Lot Number: AG722602


Exp Date: 08/13/2009

Test	g/210L	Time
DIAG	Pass	11:14pm
AIR BLK	.00	11:15pm
ACCY CHK	.07	11:16pm
AIR BLK	.00	11:17pm
SUB TEST	.00	11:17pm
AIR BLK	.00	11:18pm
SUB TEST	.00	11:20pm
AIR BLK	.00	11:21pm

Reported AC: .00 g/210L

Signature of Chemical Analyst

Court CVR

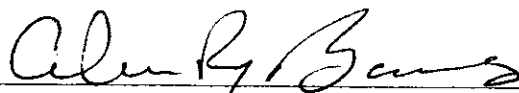
Analyst

Intox EC/IR-II: Preventive Maintenance

GASTON COUNTY BAT MOBILE UNIT 3 350

Serial Number: 008647 Test Record Number: 376
Test Date: 11/02/2008 Test Time: 11:22pm EST

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	11:22pm
FLO	Pass	11:22pm
FC	Pass	11:22pm

Temperature Tests

Test	Status	Time
FC1	Pass	11:22pm
SRC	Pass	11:22pm
DET	Pass	11:22pm
BAR	Pass	11:22pm
BT	Pass	11:22pm

Blank Tests

Test	Status	Time
AIR	Pass	11:23pm


Printer Tests

Test	Status	Time
PRNT	Pass	11:23pm

CRC Tests

Test	Status	Time
COMP	Pass	11:23pm
CAL	Pass	11:23pm

Preventive Maintenance
Status: Pass

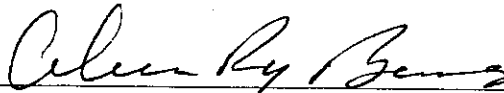

Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH

**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**

County FORSYTH Instrument Location BATMOBILE UNIT 3
Instrument Serial No. 008616 KEBDERSVILLE, NC

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 01 day of NOVEMBER, 20 08 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.



Alan R. Benson
Signature of Certifying Official

648
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

FORSYTH COUNTY BAT MOBILE UNIT 3 330

Serial Number: 008616
Test Date: 11/01/2008

Citation Number: M0000000-0
Subject's Name:

PREVENTIVE, MAINTENANCE

Subject's Date of Birth: 11/11/1911

Subject's Sex: Male

Driver's License State: XX

Driver's License Number: NONE

Analyst's Name: BARNES, ALVIN R

Permit Number: 15671E

Effective:

12/01/2007-12/01/2009

Officer's Name: NONE, NONE

Type of Agency: FTA

Agency: DHHS

Test Type: Breath Test

Lot Number: AG722601

Exp Date: 08/13/2009

Test	g/210L	Time
DIAG	Pass	10:35pm
AIR BLK	.00	10:35pm
ACCY CHK	.07	10:36pm
AIR BLK	.00	10:37pm
SUB TEST	.00	10:37pm
AIR BLK	.00	10:38pm
SUB TEST	.00	10:40pm
AIR BLK	.00	10:40pm

Reported AC: .00 g/210L

Signature of Chemical Analyst

Court CVR


Analyst

Intox EC/IR-II: Preventive Maintenance

FORSYTH COUNTY BAT MOBILE UNIT 3 330

Serial Number: 008616 Test Record Number: 350
Test Date: 11/01/2008 Test Time: 10:41pm EDT

System Check: Passed

Baseline Tests

Test	Status	Time
IR	Pass	10:42pm
FLO	Pass	10:42pm
FC	Pass	10:42pm

Temperature Tests

Test	Status	Time
FC1	Pass	10:42pm
SRC	Pass	10:42pm
DET	Pass	10:42pm
BAR	Pass	10:42pm
BT	Pass	10:42pm

Blank Tests

Test	Status	Time
AIR	Pass	10:42pm


Printer Tests

Test	Status	Time
PRNT	Pass	10:42pm

CRC Tests

Test	Status	Time
COMP	Pass	10:43pm
CAL	Pass	10:43pm

Preventive Maintenance
Status: Pass

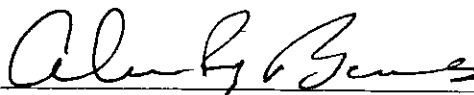

Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH

**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**


County Granville Instrument Location BAT mobile unit 4

Instrument Serial No. 008734

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 1st day of November, 2008 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

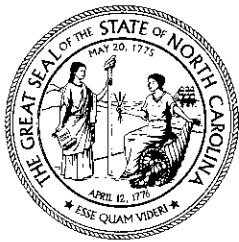

[Signature]
Signature of Certifying Official

651
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Preventive Maintenance

GRANVILLE COUNTY BAT MOBILE UNIT 4 380

Serial Number: 008734 Test Record Number: 161

Test Date: 11/01/2008 Test Time: 9:28pm EDT

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	9:29pm
FLO	Pass	9:29pm
FC	Pass	9:29pm

Temperature Tests

Test	Status	Time
FC1	Pass	9:29pm
SRC	Pass	9:29pm
DET	Pass	9:29pm
BAR	Pass	9:29pm
BT	Pass	9:29pm

Blank Tests

Test	Status	Time
AIR	Pass	9:30pm


Printer Tests

Test	Status	Time
PRNT	Pass	9:30pm

CRC Tests

Test	Status	Time
COMP	Pass	9:30pm
CAL	Pass	9:30pm

Preventive Maintenance
Status: Pass

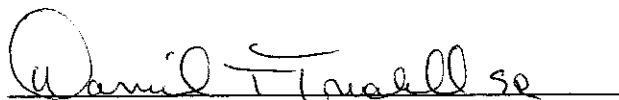

Analyst

Intox EC/IR-II: Subject Test

GRANVILLE COUNTY BAT MOBILE UNIT 4 380

Serial Number: 008734
Test Date: 11/01/2008

Citation Number: M0000000-0
Subject's Name:
PREVENTIVE, MAINTENANCE
Subject's Date of Birth: 11/11/1911
Subject's Sex: Male
Driver's License State: XX
Driver's License Number: NONE


Analyst's Name: TRUDELL, SR., DANIEL T
Permit Number: 21535E
Effective:
12/01/2007-12/01/2009

Officer's Name: NONE, NONE
Type of Agency: FTA
Agency: DHHS
Test Type: Breath Test


Lot Number: AG722501
Exp Date: 08/12/2009

Test	g/210L	Time
DIAG	Pass	9:20pm
AIR BLK	.00	9:21pm
ACCY CHK	.07	9:22pm
AIR BLK	.00	9:23pm
SUB TEST	.00	9:23pm
AIR BLK	.00	9:24pm
SUB TEST	.00	9:26pm
AIR BLK	.00	9:27pm

Reported AC: .00 g/210L

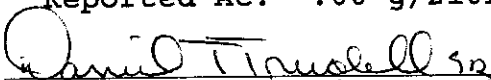
Signature of Chemical Analyst

Court CVR

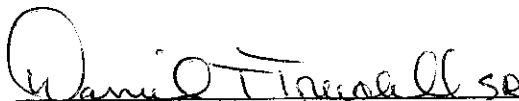
Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH


**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**

County Wake Instrument Location BAT mobile unit #5
Instrument Serial No. 008600 Wendell

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 1st day of November, 2008 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

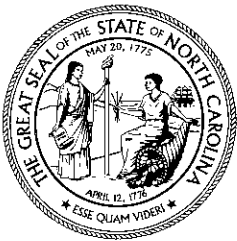

Joseph C. Tilton 636
Signature of Certifying Official Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Preventive Maintenance

WAKE COUNTY BAT MOBILE UNIT 5 910

Serial Number: 008600 Test Record Number: 404
Test Date: 11/01/2008 Test Time: 9:47pm EDT

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	9:47pm
FLO	Pass	9:47pm
FC	Pass	9:47pm

Temperature Tests

Test	Status	Time
FC1	Pass	9:47pm
SRC	Pass	9:47pm
DET	Pass	9:47pm
BAR	Pass	9:47pm
BT	Pass	9:47pm

Blank Tests

Test	Status	Time
AIR	Pass	9:48pm


Printer Tests

Test	Status	Time
PRNT	Pass	9:48pm

CRC Tests

Test	Status	Time
COMP	Pass	9:48pm
CAL	Pass	9:48pm

Preventive Maintenance
Status: Pass

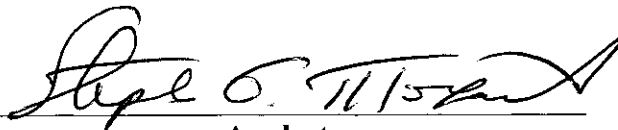
Analyst

Intox EC/IR-II: Subject Test

WAKE COUNTY BAT MOBILE UNIT 5 910

Serial Number: 008600

Test Date: 11/01/2008

Citation Number: M0000000-0

Subject's Name:

PREVENTIVE, MAINTENANCE

Subject's Date of Birth: 11/11/1911

Subject's Sex: Male

Driver's License State: XX

Driver's License Number: NONE

Analyst's Name: MORGART, STEPHEN G

Permit Number: 09372E

Effective:

12/01/2007-12/01/2009

Officer's Name: NONE, NONE

Type of Agency: FTA

Agency: DHHS


Test Type: Breath Test

Lot Number: AG723401


Exp Date: 08/21/2009

Test	g/210L	Time
DIAG	Pass	9:38pm
AIR BLK	.00	9:39pm
ACCY CHK	.07	9:40pm
AIR BLK	.00	9:41pm
SUB TEST	.00	9:41pm
AIR BLK	.00	9:43pm
SUB TEST	.00	9:44pm
AIR BLK	.00	9:45pm

Reported AC: .00 g/210L

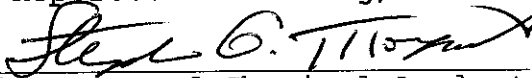
Signature of Chemical Analyst

Court CVR

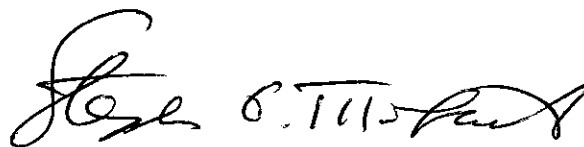

Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH


**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**

County Wake Instrument Location BAT MOBILE UNIT #5
Instrument Serial No. 008698 Wenell

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 1ST day of NOVEMBER, 2009 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

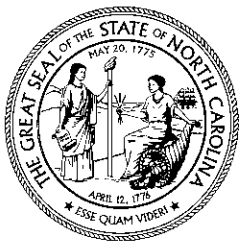

[Signature]
Signature of Certifying Official

636
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Preventive Maintenance

WAKE COUNTY BAT MOBILE UNIT 5 910

Serial Number: 008698 Test Record Number: 246
Test Date: 11/01/2008 Test Time: 10:40pm EDT

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	10:40pm
FLO	Pass	10:40pm
FC	Pass	10:40pm

Temperature Tests

Test	Status	Time
FC1	Pass	10:40pm
SRC	Pass	10:40pm
DET	Pass	10:40pm
BAR	Pass	10:40pm
BT	Pass	10:40pm

Blank Tests

Test	Status	Time
AIR	Pass	10:41pm


Printer Tests

Test	Status	Time
PRNT	Pass	10:41pm

CRC Tests

Test	Status	Time
COMP	Pass	10:41pm
CAL	Pass	10:41pm

Preventive Maintenance
Status: Pass

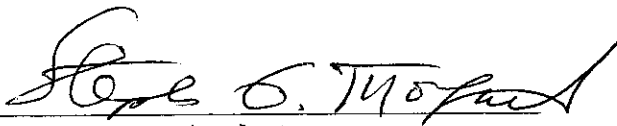

Analyst

Intox EC/IR-II: Subject Test

WAKE COUNTY BAT MOBILE UNIT 5 910

Serial Number: 008698
Test Date: 11/01/2008

Citation Number: M0000000-0
Subject's Name:
PREVENTIVE, MAINTENANCE
Subject's Date of Birth: 11/11/1911
Subject's Sex: Male
Driver's License State: XX
Driver's License Number: NONE


Analyst's Name: MORGART, STEPHEN G
Permit Number: 09372E
Effective:
12/01/2007-12/01/2009

Officer's Name: NONE, NONE
Type of Agency: FTA
Agency: DHHS
Test Type: Breath Test

Lot Number: AG723401
Exp Date: 08/21/2009


Test	g/210L	Time
DIAG	Pass	10:31pm
AIR BLK	.00	10:32pm
ACCY CHK	.07	10:32pm
AIR BLK	.00	10:33pm
SUB TEST	.00	10:34pm
AIR BLK	.00	10:35pm
SUB TEST	.00	10:36pm
AIR BLK	.00	10:37pm

Reported AC: .00 g/210L



Signature of Chemical Analyst

Court CVR

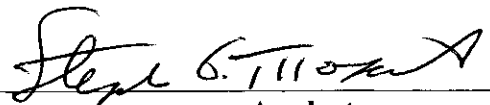

Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH

**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**


County WAKE Instrument Location BAT MOBILE UNIT #5

Instrument Serial No. 008755 W6-DH

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 15th day of November, 2008 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

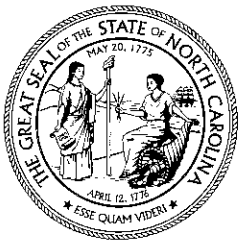

Leah C. Thayer

Signature of Certifying Official

636

Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Preventive Maintenance

WAKE COUNTY BAT MOBILE UNIT 5 910

Serial Number: 008788 Test Record Number: 122
Test Date: 11/01/2008 Test Time: 10:52pm EDT

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	10:52pm
FLO	Pass	10:52pm
FC	Pass	10:53pm

Temperature Tests

Test	Status	Time
FC1	Pass	10:53pm
SRC	Pass	10:53pm
DET	Pass	10:53pm
BAR	Pass	10:53pm
BT	Pass	10:53pm

Blank Tests

Test	Status	Time
AIR	Pass	10:53pm


Printer Tests

Test	Status	Time
PRNT	Pass	10:53pm

CRC Tests

Test	Status	Time
COMP	Pass	10:54pm
CAL	Pass	10:54pm

Preventive Maintenance
Status: Pass

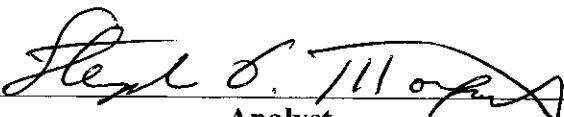

Analyst

Intox EC/IR-II: Subject Test .

WAKE COUNTY BAT MOBILE UNIT 5 910

Serial Number: 008788

Test Date: 11/01/2008

Citation Number: M0000000-0

Subject's Name:

PREVENTIVE, MAINTENANCE

Subject's Date of Birth: 11/11/1911

Subject's Sex: Male

Driver's License State: XX

Driver's License Number: NONE

Analyst's Name: MORGART, STEPHEN G

Permit Number: 09372E

Effective:

12/01/2007-12/01/2009

Officer's Name: NONE, NONE

Type of Agency: FTA

Agency: DHHS


Test Type: Breath Test

Lot Number: AG722802


Exp Date: 08/15/2009

Test	g/210L	Time
DIAG	Pass	10:36pm
AIR BLK	.00	10:37pm
ACCY CHK	.08	10:38pm
AIR BLK	.00	10:39pm
SUB TEST	.00	10:39pm
AIR BLK	.00	10:40pm
SUB TEST	.00	10:42pm
AIR BLK	.00	10:43pm

Reported AC: .00 g/210L

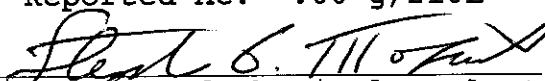
Signature of Chemical Analyst

Court CVR

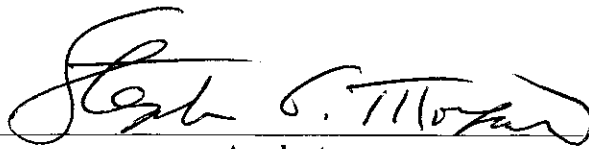

Analyst

This form is used when performing Preventive Maintenance procedures
Forensic Tests for Alcohol Branch
Department of Health and Human Services
Rev. 12/2007

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH

**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**


County Gaston Instrument Location Belmont PD

Instrument Serial No. 003733 201 Chronicle Street, Belmont
704-825-3792

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 24th day of November, 2008 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

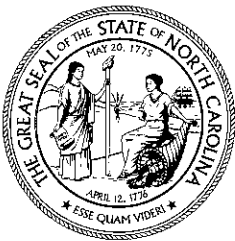

Joseph C. Hatcher
Signature of Certifying Official

650
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

GASTON COUNTY BELMONT PD 350

Serial Number: 008733
Test Date: 11/24/2008

Citation Number: M0000000-0

Subject's Name:

PREVENTIVE, MAINTENANCE

Subject's Date of Birth: 11/11/1911

Subject's Sex: Male

Driver's License State: XX

Driver's License Number: NONE

Analyst's Name: HUTCHINSON, JOSEPH E

Permit Number: 19951E

Effective:

12/01/2007-12/01/2009

Officer's Name: NONE, NONE

Type of Agency: FTA

Agency: DHHS


Test Type: Breath Test

Lot Number: AG722501

Exp Date: 08/12/2009

Test	g/210L	Time
DIAG	Pass	1:16pm
AIR BLK	.00	1:17pm
ACCY CHK	.07	1:17pm
AIR BLK	.00	1:18pm
SUB TEST	.00	1:19pm
AIR BLK	.00	1:20pm
SUB TEST	.00	1:21pm
AIR BLK	.00	1:22pm

Reported AC: .00 g/210L

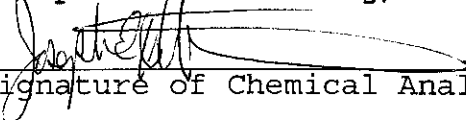

Signature of Chemical Analyst

Court CVR

Analyst

This form is used when performing Preventive Maintenance procedures
Forensic Tests for Alcohol Branch
Department of Health and Human Services
Rev. 12/2007

Intox EC/IR-II: Preventive Maintenance

GASTON COUNTY BELMONT PD 350

Serial Number: 008733 Test Record Number: 153
Test Date: 11/24/2008 Test Time: 1:12pm EST

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	1:12pm
FLO	Pass	1:12pm
FC	Pass	1:12pm

Temperature Tests

Test	Status	Time
FC1	Pass	1:13pm
SRC	Pass	1:13pm
DET	Pass	1:13pm
BAR	Pass	1:13pm
BT	Pass	1:13pm

Blank Tests

Test	Status	Time
AIR	Pass	1:13pm


Printer Tests

Test	Status	Time
PRNT	Pass	1:13pm

CRC Tests

Test	Status	Time
COMP	Pass	1:13pm
CAL	Pass	1:13pm

Preventive Maintenance
Status: Pass

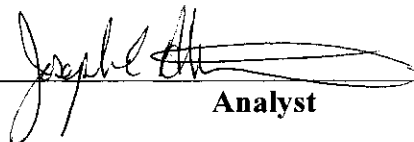

Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH


**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**

County Gaston Instrument Location Gaston County SD
Instrument Serial No. 008643 425 N. Marietta Street, Gastonia
704-869-6800

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 24th day of November, 2008 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

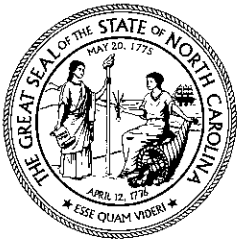

Joseph E. Hutchins
Signature of Certifying Official

654
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

GASTON COUNTY GASTON COUNTY SD 350

Serial Number: 008643
Test Date: 11/24/2008

Citation Number: M0000000-0
Subject's Name:

PREVENTIVE, MAINTENANCE

Subject's Date of Birth: 11/11/1911

Subject's Sex: Male

Driver's License State: XX

Driver's License Number: NONE

Analyst's Name: HUTCHINSON, JOSEPH E

Permit Number: 19951E

Effective:

12/01/2007-12/01/2009

Officer's Name: NONE, NONE

Type of Agency: FTA

Agency: DHHS


Test Type: Breath Test

Lot Number: AG722503

Exp Date: 08/12/2009

Test	g/210L	Time
DIAG	Pass	12:03pm
AIR BLK	.00	12:03pm
ACCY CHK	.07	12:04pm
AIR BLK	.00	12:05pm
SUB TEST	.00	12:05pm
AIR BLK	.00	12:06pm
SUB TEST	.00	12:09pm
AIR BLK	.00	12:09pm

Reported AC: .00 g/210L

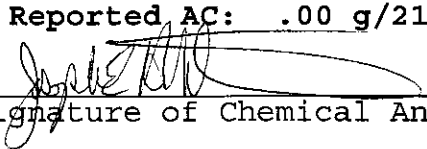

Signature of Chemical Analyst

Court CVR



Analyst

**This form is used when performing Preventive Maintenance procedures
Forensic Tests for Alcohol Branch
Department of Health and Human Services
Rev. 12/2007**

Intox EC/IR-II: Preventive Maintenance

GASTON COUNTY GASTON COUNTY SD 350

Serial Number: 008643 Test Record Number: 291
Test Date: 11/24/2008 Test Time: 11:58am EST

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	11:59am
FLO	Pass	11:59am
FC	Pass	11:59am

Temperature Tests

Test	Status	Time
FC1	Pass	11:59am
SRC	Pass	11:59am
DET	Pass	11:59am
BAR	Pass	11:59am
BT	Pass	11:59am

Blank Tests

Test	Status	Time
AIR	Pass	12:00pm


Printer Tests

Test	Status	Time
PRNT	Pass	12:00pm

CRC Tests

Test	Status	Time
COMP	Pass	12:00pm
CAL	Pass	12:00pm

Preventive Maintenance
Status: Pass

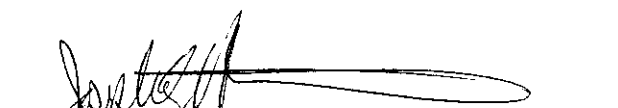

Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH

**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**


County Gaston Instrument Location Gaston County SD

Instrument Serial No. 008881 425 N. Marietta Street, Gastonia
704-869-6800

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 24th day of November, 2008 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

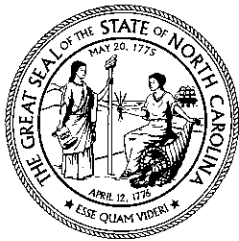

Joseph H. Hales
Signature of Certifying Official

650
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

GASTON COUNTY GASTON COUNTY SD 350

Serial Number: 008881
Test Date: 11/24/2008

Citation Number: M0000000-0
Subject's Name:
PREVENTIVE, MAINTENANCE
Subject's Date of Birth: 11/11/1911
Subject's Sex: Male
Driver's License State: XX
Driver's License Number: NONE


Analyst's Name: HUTCHINSON, JOSEPH E
Permit Number: 19951E
Effective:
12/01/2007-12/01/2009

Officer's Name: NONE, NONE
Type of Agency: FTA
Agency: DHHS
Test Type: Breath Test

Lot Number: AG722503
Exp Date: 08/12/2009


Test	g/210L	Time
DIAG	Pass	11:41am
AIR BLK	.00	11:41am
ACCY CHK	.07	11:42am
AIR BLK	.00	11:43am
SUB TEST	.00	11:43am
AIR BLK	.00	11:44am
SUB TEST	.00	11:46am
AIR BLK	.00	11:47am

Reported AC: .00 g/210L

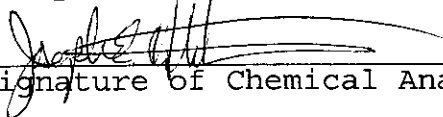

Signature of Chemical Analyst

Court CVR

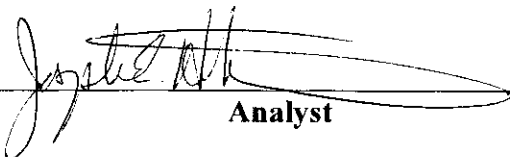

Analyst

This form is used when performing Preventive Maintenance procedures
Forensic Tests for Alcohol Branch
Department of Health and Human Services
Rev. 12/2007

Intox EC/IR-II: Preventive Maintenance

GASTON COUNTY GASTON COUNTY SD 350

Serial Number: 008881 Test Record Number: 108
Test Date: 11/24/2008 Test Time: 11:48am EST

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	11:48am
FLO	Pass	11:48am
FC	Pass	11:48am

Temperature Tests

Test	Status	Time
FC1	Pass	11:48am
SRC	Pass	11:48am
DET	Pass	11:48am
BAR	Pass	11:48am
BT	Pass	11:48am

Blank Tests

Test	Status	Time
AIR	Pass	11:49am


Printer Tests

Test	Status	Time
PRNT	Pass	11:49am

CRC Tests

Test	Status	Time
COMP	Pass	11:49am
CAL	Pass	11:49am

Preventive Maintenance
Status: Pass

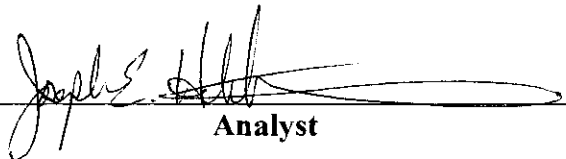
Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH


**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**

County Gaston Instrument Location Gaston County SD
Instrument Serial No. 008706 425 N. Marietta Street, Gastonia
704-869-6800

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 24th day of November, 2008 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

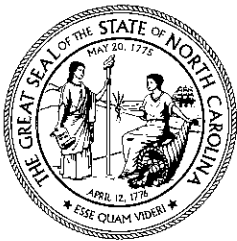

Joseph. Smith
Signature of Certifying Official

650
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

GASTON COUNTY GASTON COUNTY SD 350

Serial Number: 008706
Test Date: 11/24/2008

Citation Number: M0000000-0
Subject's Name:
PREVENTIVE, MAINTENANCE
Subject's Date of Birth: 11/11/1911
Subject's Sex: Male
Driver's License State: XX
Driver's License Number: NONE


Analyst's Name: HUTCHINSON, JOSEPH E
Permit Number: 19951E
Effective:
12/01/2007-12/01/2009

Officer's Name: NONE, NONE
Type of Agency: FTA
Agency: DHHS
Test Type: Breath Test


Lot Number: AG722502
Exp Date: 08/12/2009

Test	g/210L	Time
DIAG	Pass	12:22pm
AIR BLK	.00	12:23pm
ACCY CHK	.07	12:23pm
AIR BLK	.00	12:24pm
SUB TEST	.00	12:26pm
AIR BLK	.00	12:27pm
SUB TEST	.00	12:28pm
AIR BLK	.00	12:29pm

Reported AC: .00 g/210L

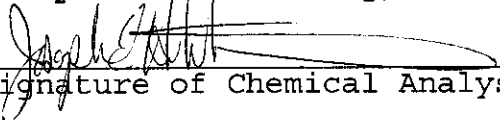
Signature of Chemical Analyst

Court CVR

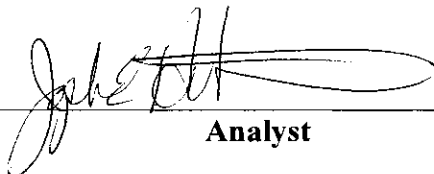
Analyst

Intox EC/IR-II: Preventive Maintenance

GASTON COUNTY GASTON COUNTY SD 350

Serial Number: 008706 Test Record Number: 366
Test Date: 11/24/2008 Test Time: 12:18pm EST

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	12:18pm
FLO	Pass	12:18pm
FC	Pass	12:19pm

Temperature Tests

Test	Status	Time
FC1	Pass	12:19pm
SRC	Pass	12:19pm
DET	Pass	12:19pm
BAR	Pass	12:19pm
BT	Pass	12:19pm

Blank Tests

Test	Status	Time
AIR	Pass	12:19pm


Printer Tests

Test	Status	Time
PRNT	Pass	12:19pm

CRC Tests

Test	Status	Time
COMP	Pass	12:20pm
CAL	Pass	12:20pm

Preventive Maintenance
Status: Pass

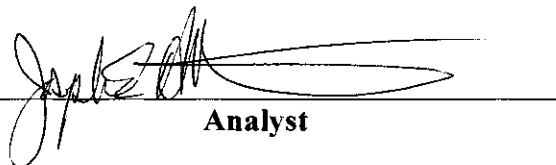

Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH

**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**


County Catawba Instrument Location Catawba County SD

Instrument Serial No. 008687 100 B Southwest Blvd, Newton
828-464-5241

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 25th day of November, 2008 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

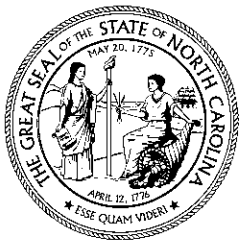

Joseph E. Hutchinson
Signature of Certifying Official

650
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

CATAWBA COUNTY CATAWBA COUNTY SD 170

Serial Number: 008687
Test Date: 11/25/2008

Citation Number: M0000000-0
Subject's Name:

PREVENTIVE, MAINTENANCE

Subject's Date of Birth: 11/11/1911

Subject's Sex: Male

Driver's License State: XX

Driver's License Number: NONE

Analyst's Name: HUTCHINSON, JOSEPH E

Permit Number: 19951E

Effective:

12/01/2007-12/01/2009

Officer's Name: NONE, NONE

Type of Agency: FTA

Agency: DHHS

Test Type: Breath Test

Lot Number: AG723302

Exp Date: 08/20/2009

Test	g/210L	Time
DIAG	Pass	10:37am
AIR BLK	.00	10:38am
ACCY CHK	.08	10:39am
AIR BLK	.00	10:40am
SUB TEST	.00	10:40am
AIR BLK	.00	10:41am
SUB TEST	.00	10:43am
AIR BLK	.00	10:44am

Reported AC: .00 g/210L



Signature of Chemical Analyst

Court CVR



Analyst

**This form is used when performing Preventive Maintenance procedures
Forensic Tests for Alcohol Branch
Department of Health and Human Services
Rev. 12/2007**

Intox EC/IR-II: Preventive Maintenance

CATAWBA COUNTY CATAWBA COUNTY SD 170

Serial Number: 008687 Test Record Number: 445
Test Date: 11/25/2008 Test Time: 10:45am EST

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	10:46am
FLO	Pass	10:46am
FC	Pass	10:46am

Temperature Tests

Test	Status	Time
FC1	Pass	10:46am
SRC	Pass	10:46am
DET	Pass	10:46am
BAR	Pass	10:46am
BT	Pass	10:46am

Blank Tests

Test	Status	Time
AIR	Pass	10:47am


Printer Tests

Test	Status	Time
PRNT	Pass	10:47am

CRC Tests

Test	Status	Time
COMP	Pass	10:47am
CAL	Pass	10:47am

Preventive Maintenance
Status: Pass

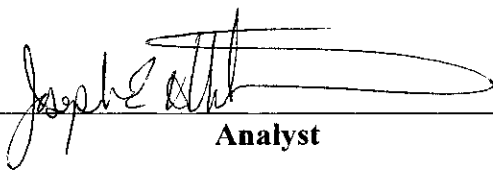
Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH

**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**


County Catawba Instrument Location Catawba County SD

Instrument Serial No. 008821 100 B Southwest Blvd, Newton
828-464-5241

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 25th day of November, 2008 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

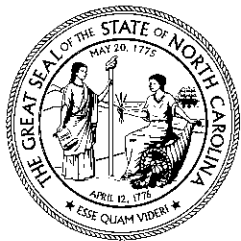

Joseph E. Hutchins
Signature of Certifying Official

650
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

CATAWBA COUNTY CATAWBA COUNTY SD 170

Serial Number: 008821
Test Date: 11/25/2008

Citation Number: M0000000-0
Subject's Name:

PREVENTIVE, MAINTENANCE

Subject's Date of Birth: 11/11/1911

Subject's Sex: Male

Driver's License State: XX

Driver's License Number: NONE

Analyst's Name: HUTCHINSON, JOSEPH E

Permit Number: 19951E

Effective:

12/01/2007-12/01/2009

Officer's Name: NONE, NONE

Type of Agency: FTA

Agency: DHHS


Test Type: Breath Test

Lot Number: AG723302


Exp Date: 08/20/2009

Test	g/210L	Time
DIAG	Pass	10:36am
AIR BLK	.00	10:37am
ACCY CHK	.07	10:38am
AIR BLK	.00	10:38am
SUB TEST	.00	10:39am
AIR BLK	.00	10:40am
SUB TEST	.00	10:41am
AIR BLK	.00	10:42am

Reported AC: .00 g/210L

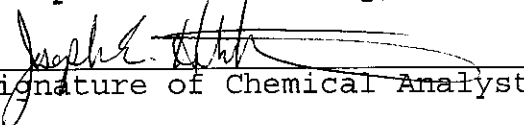
Signature of Chemical Analyst

Court CVR

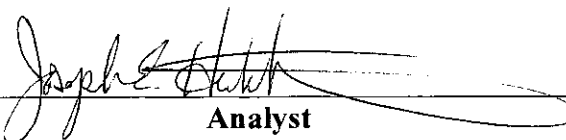
Analyst

Intox EC/IR-II: Preventive Maintenance

CATAWBA COUNTY CATAWBA COUNTY SD 170

Serial Number: 008821 Test Record Number: 224
Test Date: 11/25/2008 Test Time: 10:43am EST

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	10:44am
FLO	Pass	10:44am
FC	Pass	10:44am

Temperature Tests

Test	Status	Time
FC1	Pass	10:44am
SRC	Pass	10:44am
DET	Pass	10:44am
BAR	Pass	10:44am
BT	Pass	10:44am

Blank Tests

Test	Status	Time
AIR	Pass	10:45am


Printer Tests

Test	Status	Time
PRNT	Pass	10:45am

CRC Tests

Test	Status	Time
COMP	Pass	10:45am
CAL	Pass	10:45am

Preventive Maintenance
Status: Pass

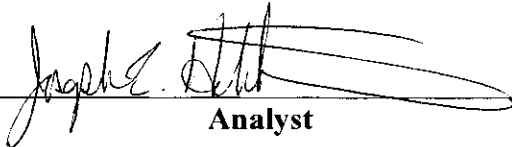
Analyst

DEPARTMENT OF HEALTH AND HUMAN SERVICES
FORENSIC TESTS FOR ALCOHOL BRANCH

**PREVENTIVE MAINTENANCE RECORD
INTOXIMETERS, MODEL INTOX EC/IR II**


County Catawba Instrument Location Hickory PD

Instrument Serial No. 008841 347 2nd Ave SW, Hickory
828-324-2060

The preventive maintenance procedures for the Intoximeters, Model Intox EC/IR II to be followed at least once every four months are:

1. Verify the ethanol gas canister displays pressure, or the alcoholic breath simulator thermometer shows 34 degrees, plus or minus .2 degree centigrade;
2. Verify instrument displays time and date;
3. Initiate breath test sequence;
4. Enter information as prompted;
5. Verify instrument accuracy;
6. When "PLEASE BLOW" appears, collect breath sample;
7. When "PLEASE BLOW" appears, collect breath sample;
8. Print test record;
9. Verify Diagnostic Program; and
10. Verify that the ethanol gas canister is being changed before expiration date, or the alcoholic breath simulator solution is being changed every four months or after 125 Alcoholic Breath Simulator tests, whichever occurs first.

I certify that on the 25th day of November, 2008 the forgoing preventive maintenance procedures were performed on the instrument indicated above, in accordance with current regulations of the N.C. Department of Health and Human Services, and the instrument is functioning properly.

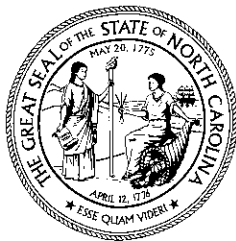

Joseph C. Hutchinson
Signature of Certifying Official

650
Certificate Number

A signed original of the preventive maintenance record shall be kept on file for at least three years.

Intox EC/IR-II: Subject Test

CATAWBA COUNTY HICKORY PD 170

Serial Number: 008841

Test Date: 11/25/2008

Citation Number: M0000000-0

Subject's Name:

PREVENTIVE, MAINTENANCE

Subject's Date of Birth: 11/11/1911

Subject's Sex: Male

Driver's License State: XX

Driver's License Number: NONE

Analyst's Name: HUTCHINSON, JOSEPH E

Permit Number: 19951E

Effective:

12/01/2007-12/01/2009

Officer's Name: NONE, NONE

Type of Agency: FTA

Agency: DHHS


Test Type: Breath Test

Lot Number: AG723302


Exp Date: 08/20/2009

Test	g/210L	Time
DIAG	Pass	12:41pm
AIR BLK	.00	12:42pm
ACCY CHK	.07	12:43pm
AIR BLK	.00	12:43pm
SUB TEST	.00	12:44pm
AIR BLK	.00	12:45pm
SUB TEST	.00	12:46pm
AIR BLK	.00	12:47pm

Reported AC: .00 g/210L

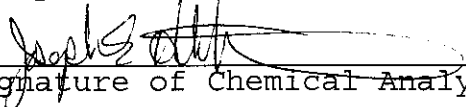
Signature of Chemical Analyst

Court CVR

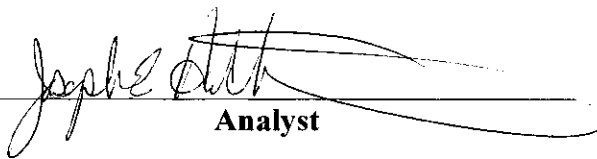
Analyst

Intox EC/IR-II: Preventive Maintenance

CATAWBA COUNTY HICKORY PD 170

Serial Number: 008841 Test Record Number: 224
Test Date: 11/25/2008 Test Time: 12:48pm EST

System Check: *Passed*

Baseline Tests

Test	Status	Time
IR	Pass	12:48pm
FLO	Pass	12:48pm
FC	Pass	12:49pm

Temperature Tests

Test	Status	Time
FC1	Pass	12:49pm
SRC	Pass	12:49pm
DET	Pass	12:49pm
BAR	Pass	12:49pm
BT	Pass	12:49pm

Blank Tests

Test	Status	Time
AIR	Pass	12:49pm


Printer Tests

Test	Status	Time
PRNT	Pass	12:49pm

CRC Tests

Test	Status	Time
COMP	Pass	12:50pm
CAL	Pass	12:50pm

Preventive Maintenance
Status: Pass

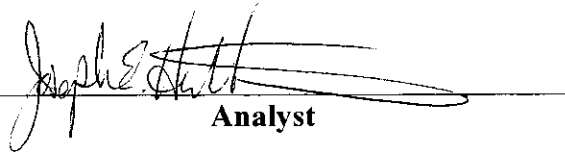
Analyst